


# PlatinumHEALTH

## Carousel Sliding Shower Chair PHB 3300(BLUE)/PHB 3400G(GRAY) Assembly and Instruction Manual


**Step 1.** Remove parts from box and take out of plastic bags.

**Step 2.** Insert the two bent legs into frame of shower chair.  
Insert until the copper tension spring buttons pop out of holes.  
Note: only two tension spring buttons will show. The other two  
tension spring buttons need to be pressed so the leg slides into  
the tube.

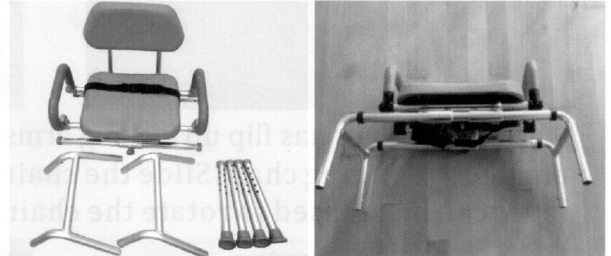


**Step 3.** Insert each frame leg into one of the four height adjustable feet.

**Note:** The rubber feet can be adjusted to the contour of the tub.

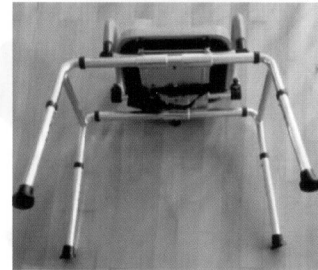


**Step 4.** Insert the backrest into the frame of the sliding chair until  
the spring tension buttons engage in the hole.

**Step 5.** Adjust legs so they will best accommodate floor heights and  
desired height for client. Legs are adjustable to the ½ inch  
by pushing the spring tension buttons.

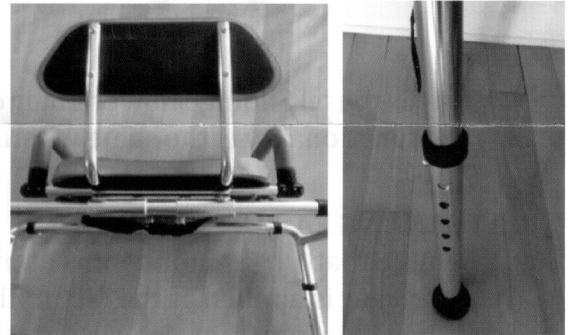


### Steps of Operation


**Step 1.** Press the grey button to engage the sliding chair. Slide the chair until it locks into place.  
To rotate chair push down on the red lever and it will allow the chair to rotate. Chair will lock  
into place every 90 degrees.

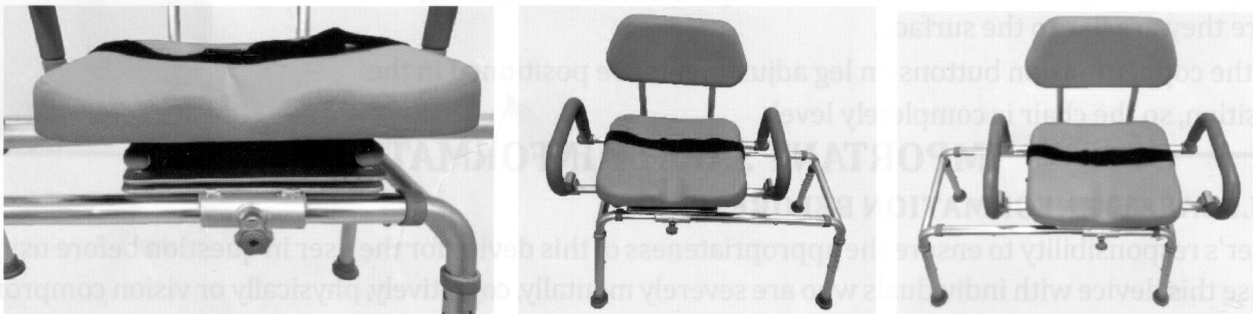

**Step 2.** Prior to showering push red lever under seat and turn the chair so the chair is facing  
the client. Once client is seated and strapped into the chair, push the red lever and have the  
client rotate into the tub one foot at a time. Once rotated push the grey button under the seat.  
This will disengage the seat so it can slide into the tub. Slide the chair until it locks into place.

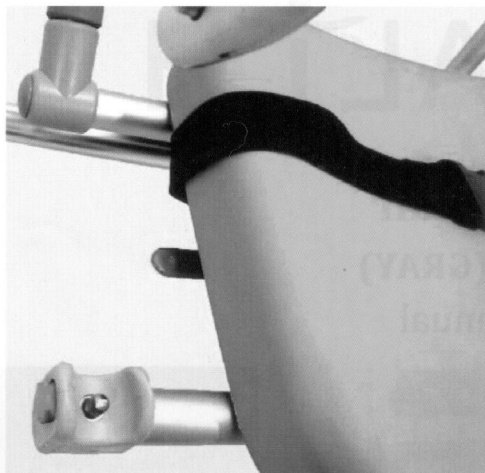

**Step 3.** The chair has flip up locking arms for easy transfer. Upon exiting the tub push the grey button to engage the sliding chair. Slide the chair until it locks. Once the chair locks into place, the red lever can be pressed to rotate the chair. Unlatch the seatbelt and exit the chair.

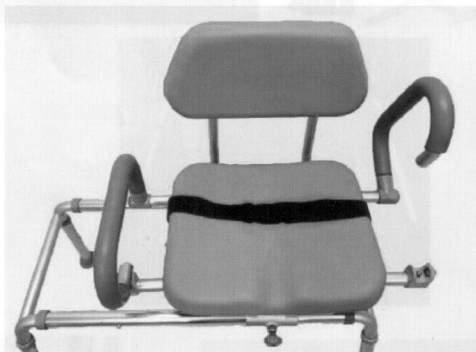

#### NOTE:

1. The Carousel Sliding Shower Chair can also be used in a shower as shown here.
2. The functionality is the same whether used in a shower or bathtub.

#### NOT SLIDING EASILY? -HOW TO FIX

Please note that a light lubricant is applied to horizontal rails that the chair slides on. This lubricant is applied at the factory before the item is shipped out. Sometimes this lubricant can evaporate if high temperatures occur during the shipping process. If your chair is not sliding easily, simply APPLY A LIGHT COATING OF WD40 LUBRICANT to both horizontal rails. This will rectify the problem.



#### Shower Chair Not level in the shower or bathtub? -HOW TO FIX

Rotate the feet to accommodate the contour or shape of the tub or shower, making sure they are flat to the surface.

Make sure the copper tension buttons on leg adjustments are positioned in the correct position, so the chair is completely level.

### IMPORTANT SAFETY INFORMATION

#### READ ALL SAFETY INFORMATION BEFORE USING

It is the user's responsibility to ensure the appropriateness of this device for the user in question before using.

- NEVER use this device with individuals who are severely mentally, cognitively, physically or vision compromised.
- ALWAYS attach safety belt and adjust to a snug fit around user's hips.
- NEVER leave a compromised user unattended.
- ALWAYS ensure that the legs are properly adjusted to ensure that the bench is level side to side and front to back.
- NEVER exceed the maximum weight capacity of 330 lbs.

NOTE: Failure to follow the above noted safety instructions can lead to severe injuries.