



ADVANCED
COMPRESSION
THERAPY
TECHNOLOGY

2023 COMPRESSION THERAPY PRODUCTS CATALOG



AIROS[®]
MEDICAL

TABLE OF CONTENTS

ABOUT

About AIROS® Medical	2
Who We Help	4
Product Overview	5

PRODUCTS

AIROS® 6/6P Sequential Compression Device	7
AIROS® 8/8P Sequential Compression Device	9
Lower Extremity Compression Garments	11
Upper Extremity Compression Garments	13
Truncal Garments	15

ADDITIONAL

Indications/Contraindications	17
Quick Reference: Devices & Garments	18
Notes	19

ABOUT AIROS® MEDICAL

AIROS Medical's team is comprised of medical device veterans and experienced clinicians specializing in compression devices and garments that effectively treat patients.

Our medical technology engineers, quality control managers, and clinical experts continually strive to develop easy-to-use devices that drive positive outcomes for patients living with chronic lymphatic and venous disorders.

We work with Durable Medical Equipment (DME) suppliers across the United States to deliver our devices and garments to the patients that need them. Our team members are always available to provide product, clinical, and billing support to our DME partners as well as the clinicians and patients who benefit from our devices and garments.





Making Compression Therapy Bearable

WHO WE HELP

Our team of medical technology designers, engineers, quality control experts, and clinical advisors collaborate to develop devices designed to drive positive patient outcomes.



Clinicians

Lymphedema therapists, vascular surgeons, wound care professionals, podiatrists, family doctors, cardiologists, and various other healthcare providers rely on the therapy our medical technology provides.



Patients

Individuals living with lymphatic and venous disorders utilize our device and garment systems every day in the comfort of their homes to reduce swelling and pain, increase mobility, and live more comfortably.



Durable Medical Equipment Suppliers

Our team includes DME veterans well-versed in billing and insurance. We pride ourselves on offering the high-quality products, consistent communication, and responsive support DMEs need to be successful.

COMPRESSION THERAPY PRODUCTS

Medical Technology & Ease of Use Meet Clinical Efficacy for Individualized Treatment

Overview

Compression therapy products, namely gradient compression devices and garment systems, have long been clinically proven to help patients suffering from lymphatic and venous disorders.

Drawing on decades of experience working with medical professionals to treat patients battling leg and arm lymphedema, venous ulcers, and other chronic conditions, AIROS' team has engineered compression devices that blend top medical technology with sound clinical direction.

Our intuitive AIROS 6/6P and AIROS 8/8P sequential compression devices feature various modes, cycle times, treatment times, and pressure ranges allowing for individualized treatment - from the most basic patient needs to complex diagnoses and prescriptions.

We are proud to offer versatile, user-friendly compression therapy products equipped with top features and functionality recommended by vascular surgeons, lymphedema therapists, physical therapists, and other healthcare professionals.

Ease of use for patients and caregivers is at the core of our product line design. Our industry experience has shown that patient compliance - and therefore everyday outcomes - is substantially enhanced when medical products are easy for patients to understand and operate.

✔ Gradient & Peristaltic Compression Therapy Modes

Pressures, cycle times, therapy times, and sequences aligned with years of therapeutic guidance from clinicians. Devices now available with both gradient sequential and peristaltic inflation sequences.

✔ Intuitive User Interface & Design

Raised buttons, large LCD screen, and versatile operating modes to promote patient compliance.

✔ Lower Extremity, Upper Extremity & Truncal Treatment

Effective therapy for the legs, arms, arm/shoulder/chest, and truncal body areas.

✔ Lymphedema Preparation Therapy (LPT) Modes

Optional secondary preparation modes to mimic manual lymphatic drainage (MLD) or tenderize fibrotic tissue prior to main therapy.

✔ Individualized Garment Fitting

42 garment sizes across 6-chamber and 8-chamber configurations allow patients' individual sizing needs to be met.

✔ Strategic Design

Carefully chosen materials for patient comfort & sturdy performance.

✔ Thoughtful Design Features

Garments include pull handles and O-ring zippers, and devices utilize a vacuum to make garment application and removal easier.

DEVICES

Reliable, Versatile Compression Therapy Devices



AIROS 6/6P

Advanced E0651 6-Chamber Device for Leg, Arm & Truncal Treatment. Now Available with Peristaltic Mode Functionality for Upper Extremity Therapy.



AIROS 8/8P

High-Performing E0652 8-Chamber Device Treating Upper & Lower Extremities Plus Truncal Area. Now Available with Peristaltic Mode Functionality for Upper Extremity Therapy.

GARMENTS

Durable, Comfortable Compression Garments



Lower Extremity Garments

For leg therapy, now available in larger sizes
13 6-chamber sizes & 11 8-chamber sizes



Upper Extremity Garments

For arm therapy and also arm/shoulder/chest area therapy
6 6-chamber sizes & 6 8-chamber sizes



Truncal Garments

Pants garments with pull handles for the truncal region
3 6-chamber sizes & 3 8-chamber sizes

AIROS 6 Sequential Compression Device

Reliable, Easy-to-Use 6-Chamber Device Ideal for the Medicare Population

The FDA-cleared AIROS® 6 Sequential Compression Device delivers compression therapy to help patients manage lymphatic and venous disorders. The device provides effective intermittent pneumatic compression (IPC) therapy for treatment of the lower and upper extremities.

The AIROS 6 provides gradient sequential compression that increases blood flow and encourages extracellular fluid clearance. The device provides cycles of compressed air at certain adjustable pressures and then inflates and deflates garments to move fluid.



Overview & Features

Regulatory Clearances: FDA 510(k) & Medicare (HCPCS Code E0651); also reimbursable through commercial insurance

Usage Tracking: Built-in timer records total usage time to measure patient compliance

Easy Operation: Large, clear LCD screen and one-button operation for starting and stopping therapy

Physical Benefits: Reduce pain and swelling, increase blood flow, move lymphatic fluid

Various Sizes of Comfortable Garments: 10 garments for lower extremity treatment & 6 garments for upper extremity treatment, including arm/shoulder/chest (NOTE: Not available for use with double extra wide leg or pants garments. These garments must be used with the AIROS 6P device)

Lymphedema Preparation Treatment (LPT) Modes: Optional Clearing Mode (mimics lymphatic drainage) or Focal Mode (fibrotic tissue relief) can be used before main Gradient Mode therapy

Powerful Vacuum: Air is quickly vacuumed from the garment after therapy is completed to assist in garment removal

Additional Features for AIROS 6P Model

- ✓ **Universal Electric:** The AIROS 6P Model features 100-240 VAC, 50/60Hz electrical rating, which allows the device to be used anywhere in the world.
- ✓ **Peristaltic Mode:** The Peristaltic Mode provides a wave-like feeling to the patient, as the pressure is not held in all chambers throughout a cycle of treatment. Some lymphedema therapists prefer this type of therapy for patients with upper extremity swelling.

The user can set the starting pressure in Chamber 1 near the hand. Pressure decreases at a fixed 5% gradient distally to proximally. Chambers 1 through 3 first inflate to their target pressures sequentially and hold. Then Chamber 1 deflates, relieving pressure near the hand and wrist, and Chamber 4 inflates to its target pressure and holds. This process continues as the 6 chambers are inflated moving up the extremity.

NOTE: The AIROS 6P must be used for the pants garments which are prescribed for truncal swelling, and some of the larger sized leg garments.



AIROS 6/6P
Instructions for
Use Guide (PDF)

AIROS 6P Sequential Compression Device

6-Chamber E0651 Model for Use with Peristaltic Mode, Truncal Treatment & Larger Leg Garments

The AIROS 6P model is based on the AIROS 6 device and provides all the same therapeutic benefits and most of the same technical characteristics. The main difference is that the 6P features an additional primary operating mode, the Peristaltic Mode, for use with the Arm and Arm plus garments (upper extremity only).



AIROS 6/6P Garment Use Reference Guide

Garment Type/Model	AIROS 6	AIROS 6P
LOWER EXTREMITY		
6-SL01, 6-SWL01, 6-SXWL01	YES	YES
6-SDXWL01	NO	YES
6-ML01, 6-MWL01, 6-MXWL01	YES	YES
6-MDXWL01	NO	YES
6-LL01, 6-LWL01	YES	YES
6-LXWL01	NO	YES
6-HL01, 6-HWL01	YES	YES
UPPER EXTREMITY		
6-MA01, 6-LA01	YES, GRADIENT MODE ONLY	YES, GRADIENT OR PERISTALTIC MODE
6-APL01, 6-APLP01, 6-APR01, 6-APRP01	YES, GRADIENT MODE ONLY	YES, GRADIENT OR PERISTALTIC MODE
TRUNCAL (PANTS)		
6-SP01, 6-MP01, 6-LP01	NO	YES

AIROS 6/6P Specifications & Information Comparison

Model	AIROS 6	AIROS 6P
Dimensions	10.2"W x 10.5" D x 4.3"H	10.2"W x 10.5" D x 4.3"H
Weight	6.4 lbs.	6.4 lbs.
Electrical	120 VAC, 60Hz, 100VA MAX (USA Only)	100-240 VAC, 50/60Hz, 100VA MAX (USA & International)
Primary Operating Modes	1 (Gradient Mode)	2 (Gradient Mode and Peristaltic Mode)
Secondary Lymphedema Preparation Treatment (LPT) Modes	2 (Focal Mode and Clearing Mode)	2 (Focal Mode and Clearing Mode)
Pressure Range Setting (Primary Operating Modes)	30-80mmHg (increments of 1mmHg)	30-80mmHg (increments of 1mmHg)
Cycle Time (Inflation) Setting	Gradient Mode and Pressure Mode: 40 or 60 seconds	Gradient Mode and Pressure Mode: 40 or 60 seconds Peristaltic Mode: No setting
Therapy Time Setting	15, 30, 45 or 60 minutes	15, 30, 45 or 60 minutes
Device Warranty	3 years	3 years
Error Indicators	3	3
Air Blocker (Unilateral Therapy)	1 included	1 included

AIROS 8 Sequential Compression Device

High-Performing 8-Chamber Device For Upper & Lower Extremity Treatment

The FDA-cleared AIROS® 8 Sequential Compression Device is an advanced compression therapy device utilized for the treatment and management of lymphedema and venous insufficiencies. The device provides effective intermittent pneumatic compression (IPC) therapy for both the upper and lower extremities.

The AIROS 8 features two main therapy modes – Gradient Mode and Pressure Mode – that provide gradient sequential compression. The device provides cycles of compressed air at certain adjustable pressures to increase blood flow and encourage extracellular fluid clearance. When utilizing the Pressure Mode, the pressure in each chamber can be individually adjusted to accommodate different therapy needs.



Overview & Features

Regulatory Clearances:

FDA 510(k) & Medicare (HCPCS Code E0652); also reimbursable through commercial insurance

Usage Tracking: Built-in timer records total usage time to measure patient compliance

Easy Operation: Large, clear LCD screen and one-button operation for starting and stopping therapy

Physical Benefits: Reduce pain and swelling, increase blood flow, move lymphatic fluid

Various Sizes of Comfortable

Garments: 8 garments for lower extremity treatment and 6 garments for upper extremity treatment, including arm/shoulder/chest (NOTE: Not available for use with double extra wide leg or pants garments. These garments must be used with the AIROS 8P device)

Lymphedema Preparation Treatment (LPT) Modes: Optional Clearing Mode (mimics lymphatic drainage) or Focal Mode (fibrotic tissue relief) can be used before main Gradient Mode or Pressure Mode therapy

Powerful Vacuum: Air is quickly vacuumed from the garment after therapy is completed to assist in garment removal

Additional Features for AIROS 8P Model

- ✓ **Universal Electric:** The AIROS 8P Model features 100-240 VAC, 50/60Hz electrical rating, which allows the device to be used anywhere in the world.
- ✓ **Peristaltic Mode:** The Peristaltic Mode provides a wave-like feeling to the patient, as the pressure is not held in all chambers throughout a cycle of treatment. Some lymphedema therapists prefer this type of therapy for patients with upper extremity swelling.

The user can set the starting pressure in Chamber 1 near the hand. Pressure decreases at a fixed 4% gradient distally to proximally. Chambers 1 through 3 first inflate to their target pressures sequentially and hold. Then Chamber 1 deflates, relieving pressure near the hand and wrist, and Chamber 4 inflates to its target pressure and holds. This process continues as the 8 chambers are inflated moving up the extremity.

NOTE: The AIROS 8P must be used for the pants garments which are prescribed for truncal swelling, and some of the larger sized leg garments.



AIROS 8/8P
Instructions for
Use Guide (PDF)

AIROS 8P Sequential Compression Device

8-Chamber Model for Use with Peristaltic Mode, Truncal Treatment & Larger Leg Garments

The AIROS 8P model is based on the AIROS 8 device and provides all the same therapeutic benefits and most technical characteristics. The main difference is that the 8P features a third primary operating mode, the Peristaltic Mode, for use with Arm and Arm plus garments (upper extremity only).



AIROS 8/8P Garment Use Reference Guide

Garment Type/Model	AIROS 8	AIROS 8P
LOWER EXTREMITY		
8-SL01, 8-SWL01, 8-SXWL01	YES	YES
8-SDXWL01	NO	YES
8-ML01, 8-MWL01, 8-MXWL01	YES	YES
8-MDXWL01	NO	YES
8-LL01, 8-LWL01	YES	YES
8-LXWL01	NO	YES
UPPER EXTREMITY		
8-MA01, 8-LA01	YES, GRADIENT MODE/PRESSURE MODE ONLY	YES, GRADIENT, PRESSURE OR PERISTALTIC MODE
8-APL01, 8-APLP01, 8-APR01, 8-APRP01	YES, GRADIENT MODE/PRESSURE MODE ONLY	YES, GRADIENT, PRESSURE OR PERISTALTIC MODE
TRUNCAL (PANTS)		
8-SP01, 8-MP01, 8-LP01	NO	YES

AIROS 8/8P Specifications & Information Comparison

Model	AIROS 8	AIROS 8P
Dimensions	10.2"W x 10.5" D x 4.3"H	10.2"W x 10.5" D x 4.3"H
Weight	6.7 lbs.	6.7 lbs.
Electrical	120 VAC, 60Hz, 100VA MAX (USA Only)	100-240 VAC, 50/60Hz, 100VA MAX (USA & International)
Primary Operating Modes	2 (Gradient Mode and Pressure Mode)	3 (Gradient Mode, Pressure Mode and Peristaltic Mode)
Secondary Lymphedema Preparation Treatment (LPT) Modes	2 (Focal Mode and Clearing Mode)	2 (Focal Mode and Clearing Mode)
Pressure Range Setting (Primary Operating Modes)	30-80mmHg (increments of 1mmHg)	30-80mmHg (increments of 1mmHg)
Cycle Time (Inflation) Setting	Gradient Mode and Pressure Mode: 40 or 60 seconds	Gradient Mode and Pressure Mode: 40 or 60 seconds Peristaltic Mode: No setting
Therapy Time Setting	15, 30, 45 or 60 minutes	15, 30, 45 or 60 minutes
Device Warranty	3 years	3 years
Error Indicators	3	3
Air Blocker (Unilateral Therapy)	1 included	1 included

Lower Extremity Garments

Durable, Comfortable Garments for Effective Therapy

AIROS Medical lower extremity compression garments are available in both 6-chamber and 8-chamber configurations for use with AIROS 6/6P and AIROS 8/8P devices.

These compression garments were designed with patient comfort, ease of use, and clinical effectiveness in mind. From durable material to a wide variety of sizes, our garments were carefully engineered to support patients with various needs.

Models span small, medium, and large inseam lengths and wide sizes are also available. There are also 2 half leg garments for use only with the AIROS 6/6P device.

Key Features

Various Sizes:

6-Chamber Garments:
11 Full Leg sizes, 2 Half Leg sizes

8-Chamber Garments:
11 Full Leg sizes

Ease of Use: Pull handles and O-ring zippers to facilitate easy application and removal

High-Quality Material:
Durable material and allows for simple cleaning

Reimbursement:

Commercial Insurance
Medicare billing codes are E0667 (Full Leg Garment) and E0669 (Half Leg Garment)



Lower Extremity Compression Garments

NOTE: Please refer to the sizing chart as some garment models are only for use with the AIROS 6P and AIROS 8P devices.





AIROS
6-Chamber
Lower
Extremity
Garment

AIROS
8-Chamber
Lower
Extremity
Garment

AIROS 6-Chamber & 8-Chamber Lower Extremity Garment Sizing Chart

Garment	Model	Inseam	Circumference (Ankle)	Circumference (Knee)	Circumference (Thigh)
Small Leg	6-SL01/8-SL01	66cm (26")	53cm (21")	64cm (25")	76cm (30")
Small Wide Leg	6-SWL01/8-SWL01	66cm (26")	61cm (24")	74cm (29")	89cm (35")
Small Extra Wide Leg	6-SXWL01/8-SXWL01	66cm (26")	71cm (28")	89cm (35")	102cm (40")
Small Double Extra Wide Leg	6-SDXWL01/8-SDXWL01	66cm (26")	81cm (32")	107cm (42")	124cm (49")
Medium Leg	6-ML01/8-ML01	81cm (32")	56cm (22")	69cm (27")	79cm (31")
Medium Wide Leg	6-MWL01/8-MWL01	81cm (32")	64cm (25")	79cm (31")	89cm (35")
Medium Extra Wide Leg	6-MXWL01/8-MXWL01	81cm (32")	74cm (29")	91cm (36")	102cm (40")
Medium Double Extra Wide Leg	6-MDXWL01/8-MDXWL01	81cm (32")	81cm (32")	109cm (43")	124cm (49")
Large Leg	6-LL01/8-LL01	91cm (36")	53cm (21")	69cm (27")	79cm (31")
Large Wide Leg	6-LWL01/8-LWL01	91cm (36")	56cm (22")	79cm (31")	102cm (40")
Large Extra Wide Leg	6-LXWL01/8-LXWL01	91cm (36")	84cm (33")	109cm (43")	127cm (50")
Half Leg	6-HL01	46cm (18")	56cm (22")	71cm (28")	-
Half Wide Leg	6-HWL01	46cm (18")	61cm (24")	76cm (30")	-

Upper Extremity Garments

Advanced, Easy-to-Use Garments for Patients with Lymphedema in the Arm, Shoulder, and Chest

AIROS Medical's upper extremity garments provide an effective treatment option for patients with lymphedema in the upper extremities, especially for post-mastectomy/breast cancer patients with lymphedema.

These compression garments were designed so that patients themselves - even elderly ones - can put on the garments and remove them without the help of a family member or caregiver. This ease-of-use increases patient compliance and drives positive outcomes.



Key Features

Various Sizes

Arm Garments: Medium and Large inseam lengths

Arm Plus Garments: Available in left and right configurations, and in two colors

Both garment types available in 6-chamber (E0651) and 8-chamber (E0652) sizes

6-chamber garments can be used with the AIROS 6 (Gradient Mode) and the AIROS 6P (Gradient or Peristaltic Mode)

8-chamber garments can be used with the AIROS 8 (Gradient Mode or Pressure Mode) and the AIROS 8P (Gradient, Pressure, or Peristaltic Mode)

High-Quality Material: Durable material and allows for simple cleaning

Ease of Use: Arm Plus garment features support padding in the back to aid in donning and removal

Arm garment simply slides on and reaches the axilla area

Reimbursement:

Commercial Insurance
Medicare billing code used for both Arm and Arm Plus garment is E0668

Quality Therapy: Patients prescribed an Arm Plus garment benefit from quality therapy immediately, so no lengthy trials are needed for billing purposes



AIROS 6-Chamber & 8-Chamber Upper Extremity Compression Garments



AIROS 6-Chamber & 8-Chamber Upper Extremity Garment Sizing Chart

Garment	Model	Inseam	Body Circumference	Color	Body Side
Arm Plus Garments					
Arm Plus Left	6-APL01/8-APL01	74cm (29")	76-137cm (30"-54")	Gray	Left
Arm Plus Pink Left	6-APLP01/8-APLP01	74cm (29")	76-137cm (30"-54")	White/Pink	Right
Arm Plus Right	6-APR01/8-APR01	74cm (29")	76-137cm (30"-54")	Gray	Left
Arm Plus Pink Right	6-APRP01/8-APRP01	74cm (29")	76-137cm (30"-54")	White/Pink	Right
Arm Garments					
Medium Arm	6-MA01/8-MA01	74cm (29")	56cm (22")	48cm (19")	33cm (13")
Large Arm	6-LA01/8-LA01	81 cm (32")	58cm (23")	51cm (20")	36cm (14")

Our Arm and Arm Plus models can be used with the AIROS 6/6P or AIROS 8/8P devices. Patients utilizing the AIROS 6P or AIROS 8P device can utilize the Arm and Arm Plus garments with the Peristaltic Mode, if the patient's doctor or lymphedema therapist believe this wave-like pressure inflation would be better than traditional gradient sequential compression therapy.

Truncal Garments

Pants Configuration Features Bilateral Leg Treatment with Compression Chambers Above the Waist to Treat Abdominal Swelling

AIROS Medical's pants garments treat the entire lower body, delivering compression therapy to the legs, abdomen, lower back, and hip area. Lymphatic drainage in these areas is critical for many patients.

These garments are made from the same durable, comfortable material as our leg, arm and Arm Plus garments, and feature pull handles for easy application and removal. The pants garments are available in both 6-chamber (used with the AIROS 6P) and 8-chamber (used with the AIROS 8P) configurations.



Key Features


Full Lower Body Therapy:

Compression chambers deliver treatment to the legs, abdomen, lower back, and hip area.

Ease of Use: Pull handles and O-ring zippers allow easy application and removal.

High-Quality Material: Durable, comfortable material that is easy to clean.

Reimbursement: Medicare and commercial insurance options.

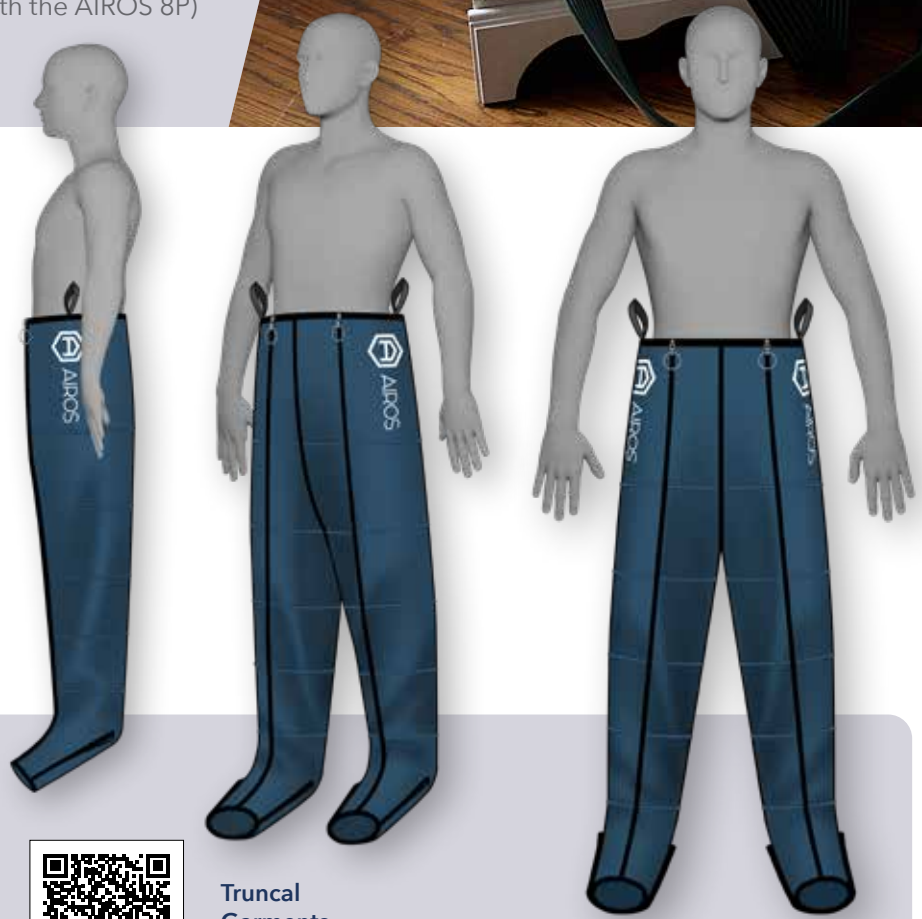

Truncal Garments



**AIROS 6-Chamber & 8-Chamber Truncal
Garment Sizing Chart**

Garment	Model	Device Compatibility	Inseam Length	Abdomen Length	Total Length	Circumference (Ankle)	Circumference (Knee)	Circumference (Thigh)	Circumference (Abdomen)
Small Pants (6-Chamber)	6-SP01	AIROS 6P	66 cm (26")	46 cm (18")	117cm (46")	53cm (21")	79 cm (31")	79 cm (31")	123 cm (48")
Medium Pants (6-Chamber)	6-MP01	AIROS 6P	81cm (32")	46 cm (18")	127cm (50")	64 cm (25")	81cm (32")	81cm (32")	125 cm (49")
Large Pants (6-Chamber)	6-LP01	AIROS 6P	91 cm (36")	46 cm (18")	142cm (56")	71 cm (28")	94 cm (37")	94 cm (37")	145 cm (57")
Small Pants (8-Chamber)	8-SP01	AIROS 8P	66 cm (26")	46 cm (18")	117cm (46")	53cm (21")	79 cm (31")	79 cm (31")	123 cm (48")
Medium Pants (8-Chamber)	8-MP01	AIROS 8P	81cm (32")	46 cm (18")	127cm (50")	64 cm (25")	81cm (32")	81cm (32")	125 cm (49")
Large Pants (8-Chamber)	8-LP01	AIROS 8P	91 cm (36")	46 cm (18")	142cm (56")	71 cm (28")	94 cm (37")	94 cm (37")	145 cm (57")

INDICATIONS FOR USE & CONTRAINDICATIONS

Below are the indications and contraindications for use for both the AIROS 6/6P and AIROS 8/8P devices, as well as associated garments.

AIROS 6/6P Sequential Compression Device

Indications for Use

The AIROS 6 Sequential Compression Device utilizes gradient pneumatic compression, which is intended for the treatment of patients with the following conditions:

- Lymphedema
- Venous stasis ulcers
- Venous insufficiency
- Peripheral edema

The device is safe for both home and hospital use.

Contraindications

Pneumatic compression **IS NOT** recommended for use in patients with the following conditions:

- Infections in the limb, including cellulitis without appropriate antibiotic coverage
- Presence of lymphangiosarcoma
- Deep Vein Thrombosis (DVT)
- Inflammatory phlebitis or episodes of pulmonary embolism
- Congestive Heart Failure (CHF)

Any local conditions in which garments would interfere, including untreated wounds, infected wounds, gangrene, recent skin graft, or dermatitis.

AIROS 8/8P Sequential Compression Device

Indications for Use

The AIROS 8 Sequential Compression Device utilizes gradient pneumatic compression, which is intended for the treatment of patients with the following conditions:

- Lymphedema
- Venous stasis ulcers
- Venous insufficiency
- Peripheral edema

The device is safe for both home and hospital use.

Contraindications

Pneumatic compression **IS NOT** recommended for use in patients with the following conditions:

- Infections in the limb, including cellulitis without appropriate antibiotic coverage
- Presence of lymphangiosarcoma
- Deep Vein Thrombosis (DVT)
- Inflammatory phlebitis or episodes of pulmonary embolism
- Congestive Heart Failure (CHF)

Any local conditions in which garments would interfere, including untreated wounds, infected wounds, gangrene, recent skin graft, or dermatitis.

ERROR INDICATORS FOR AIROS 6/6P & AIROS 8/8P

High Pressure Error Indicator: This can occur if there is a kink or other obstruction in the tubing connecting the garments, or if the connector port is blocked. The LCD screen will read **E01**, the yellow LED light will illuminate, and a buzzer will sound.

Low Pressure Error Indicator: This can occur if the garment is not properly secured to the connector port or there is a hole in the garment causing a leak. The LCD screen will read **E02**, the yellow LED light will illuminate, and a buzzer will sound.

Software System Error Indicator: This can occur if an internal software system error occurs causing the device to be inoperable. The LCD screen will read **E03**, the yellow LED light will illuminate, and a buzzer will sound.



NOTE: If you receive any of these error messages, press the **START|STOP BUTTON**, which will begin deflation of the garments. Then, contact your local authorized dealer or AIROS Medical.

QUICK REFERENCE GUIDE: DEVICES & GARMENTS

6-Chamber Garments

Garment	Model	Device Compatibility
6-SL01	Small Leg	AIROS 6/6P
6-SWL01	Small Wide Leg	AIROS 6/6P
6-SXWL01	Small Extra Wide Leg	AIROS 6/6P
6-SDXWL01	Small Double Extra Wide Leg	AIROS 6P Only
6-ML01	Medium Leg	AIROS 6/6P
6-MWL01	Medium Wide Leg	AIROS 6/6P
6-MXWL01	Medium Extra Wide Leg	AIROS 6/6P
6-MDXWL01	Medium Double Extra Wide Leg	AIROS 6P Only
6-LL01	Large Leg	AIROS 6/6P
6-LWL01	Large Wide Leg	AIROS 6/6P
6-LXWL01	Large Extra Wide Leg	AIROS 6P Only
6-HL01	Half Leg	AIROS 6/6P
6-HWL01	Half Wide Leg	AIROS 6/6P
6-SP01	Small Pants	AIROS 6P Only
6-MP01	Medium Pants	AIROS 6P Only
6-LP01	Large Pants	AIROS 6P Only
6-MA01	Medium Arm	AIROS 6/6P
6-LA01	Large Arm	AIROS 6/6P
6-APL01	Arm Plus Left	AIROS 6/6P
6-APLP01	Arm Plus Pink Left	AIROS 6/6P
6-APR01	Arm Plus Right	AIROS 6/6P
6-APRP01	Arm Plus Pink Right	AIROS 6/6P

8-Chamber Garments

Garment	Model	Device Compatibility
8-SL01	Small Leg	AIROS 8/8P
8-SWL01	Small Wide Leg	AIROS 8/8P
8-SXWL01	Small Extra Wide Leg	AIROS 8/8P
8-SDXWL01	Small Double Extra Wide Leg	AIROS 8P Only
8-ML01	Medium Leg	AIROS 8/8P
8-MWL01	Medium Wide Leg	AIROS 8/8P
8-MXWL01	Medium Extra Wide Leg	AIROS 8/8P
8-MDXWL01	Medium Double Extra Wide Leg	AIROS 8P Only
8-LL01	Large Leg	AIROS 8/8P
8-LWL01	Large Wide Leg	AIROS 8/8P
8-LXWL01	Large Extra Wide Leg	AIROS 8P Only
8-SP01	Small Pants	AIROS 8P Only
8-MP01	Medium Pants	AIROS 8P Only
8-LP01	Large Pants	AIROS 8P Only
8-MA01	Medium Arm	AIROS 8/8P
8-LA01	Large Arm	AIROS 8/8P
8-APL01	Arm Plus Left	AIROS 8/8P
8-APLP01	Arm Plus Pink Left	AIROS 8/8P
8-APR01	Arm Plus Right	AIROS 8/8P
8-APRP01	Arm Plus Pink Right	AIROS 8/8P



AIROS[®]
MEDICAL

AIROS Medical, Inc.

2501 Monroe Blvd.

Audubon, PA 19403

1-866-991-6956

info@airosmedical.com

www.airosmedical.com

PC.001 REV A