


# bili-hut™ Home Quick Reference Card

**NOTE:** Make sure to read and thoroughly understand the entire contents of the bili-hut™ Home Use Guide before operating the bili-hut™.


## Components

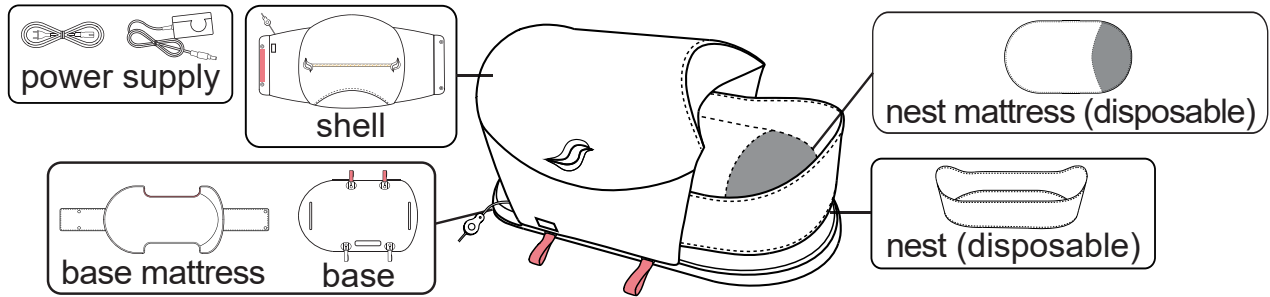


## Part 1: Assembly and Operation

**NOTE - Recommended room temperature:** select a location for treatment where the room temperature can be maintained between 70°F and 76°F (21.1°C and 24.4°C) for your baby's comfort and safety during treatment.


**a.** With the **base mattress** attached to the **base** (straps of the **base mattress** snapped closed), attach the four (4) magnetic clasps on the **shell** to the clasps on the **base** to form the **assembled hut**. The clasps should audibly “click” when in place. Gently tug upward on the **shell** to check for secure clasp attachment.

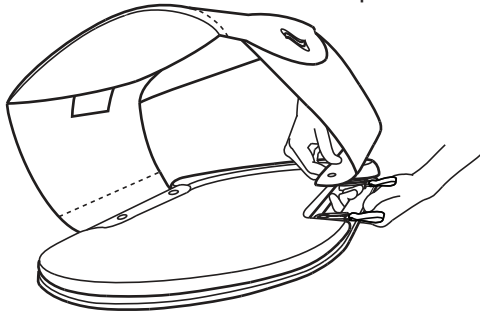

**b.** Place the **bili-hut™** on a smooth, sturdy, even surface such as a table 2 x 4 feet (0.6 m x 1.2 m) or larger.

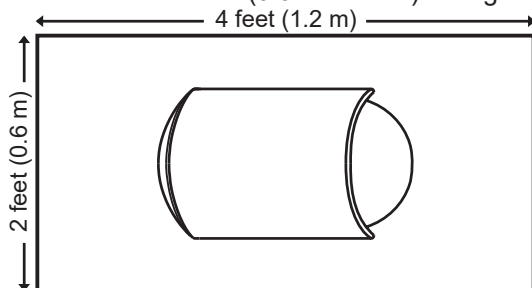


Alternatively, use a crib or other surface with barriers.

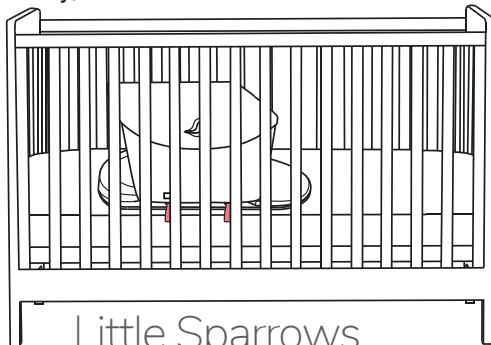

Little Sparrows  
Technologies


**c.** Place the **nest mattress** in the **nest** so that the dark section is facing up.

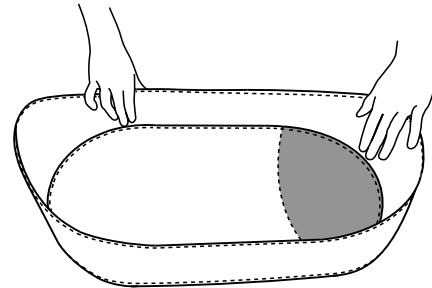


**d.** Slide the **nest** partially into the **assembled hut**. After undressing your baby except for the diaper, place eye protection over your baby's eyes. Place your baby in the **nest** so your baby's head and neck lie entirely within the dark section of the **nest mattress**.

**NOTE:** See **Section 6** of the **bili-hut™ Home Use Guide** for assistance with **disposable nest packages**.

**e.** Gently slide the **nest** further into the **assembled hut** until it touches the bumper on the **base mattress**.

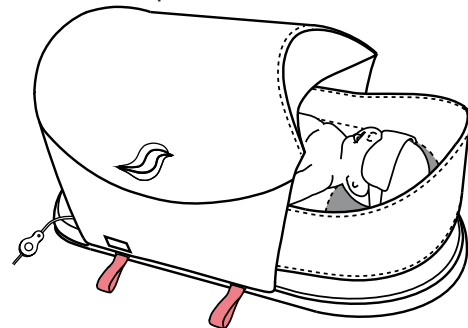

**f.** Connect the **bili-hut™** power cord to the assembled **power supply**. Plug the assembled **power supply** into a safe and active electrical outlet. Press the power switch on the power cord to turn on the **bili-hut™** and begin treating your baby.

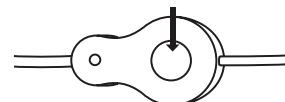


**g.** Fill out the **bili-hut™ Home Record Sheet**.

**h.** Complete prescribed phototherapy treatment.

**NOTE:** Refer to **Part 2** of this card for any treatment interruptions.

## Part 2: Removing Your Baby from the bili-hut™


- Press the power switch on the power cord to turn off the blue treatment lights before removing your baby from the **bili-hut™**.
- Gently slide the **nest** out of the **assembled hut**. Carefully lift your baby out of the **nest** and place your baby in a safe location.

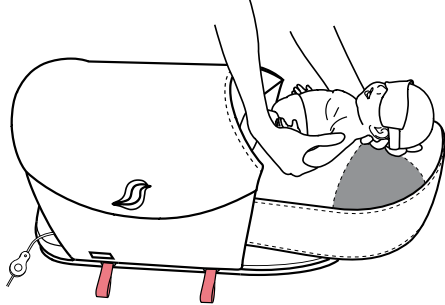

- Remove your baby's eye protection.
- Record any treatment interruptions using the **bili-hut™ Home Record Sheet**.
- To resume treatment, follow the instructions starting at **Step d** of **Part 1** on the reverse side of this card.

## Part 3: Environment High Temperature Warning

- The Environment High Temperature Warning automatically turns on when the environmental temperature surrounding the **bili-hut™** is too high.
- In this state, the blue treatment lights turn off and a yellow indicator light turns on. Reference **Section 3.4** of the **bili-hut™ Home Use Guide** for actions to take if this warning occurs.

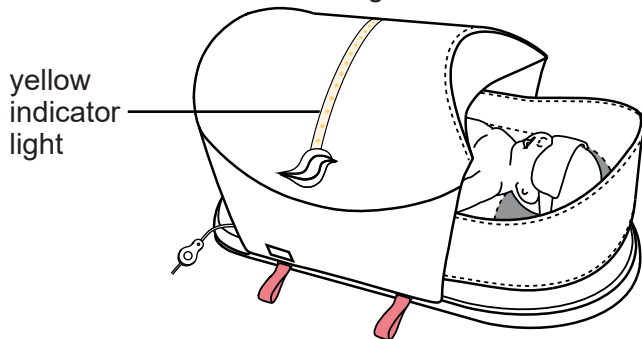

## Contact Information

Contact your healthcare provider for issues with the **bili-hut™** or health problems with your baby.

To contact Little Sparrows Technologies:

**Phone:** (781) 725-2460

**Email:** [info@littlesparrowstech.com](mailto:info@littlesparrowstech.com)

**Mailing Address:**


Little Sparrows Technologies  
288 Walnut St, Suite 300  
Newton, MA 02460

**Website:** [www.littlesparrowstech.com](http://www.littlesparrowstech.com)

**Twitter:** @TechSparrows

## Part 4: Disassembly

- Before disassembling, turn off the **bili-hut™** and remove your baby according to the instructions in **Part 2** of this card.
- Unplug the **bili-hut™**. Pull the four (4) straps to undo the clasps on the **base**.

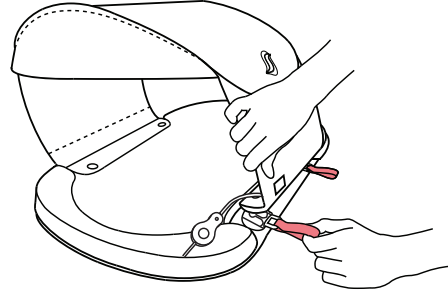

- Dispose of the **nest** and **nest mattress**; clean the **bili-hut™** according to **Part 5** of this card.








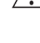
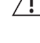
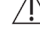
## Part 5: Cleaning

Clean the reusable components of the **bili-hut™** when soiled.

- Before cleaning, turn off the **bili-hut™**, remove your baby, and disassemble the **bili-hut™** according to **Parts 2** and **4** of this card.
- Use a disinfectant wipe or cloth dampened with disinfectant solution to clean the **bili-hut™** according to the instructions in **Section 4.1** of the **bili-hut™ Home Use Guide**.

## KEY WARNINGS (not a complete list)

**NOTE:** See the **bili-hut™ Home Use Guide** for complete safety information.

-  **WARNING:** NEVER carry your baby inside the **assembled hut** and/or inside the **nest**.
-  **WARNING:** Use eye protection to prevent damage to your baby's vision.
-  **WARNING:** Check your baby and the **bili-hut™** regularly (every 2 to 4 hours).
-  **WARNING:** DO NOT place the **bili-hut™** near anything that heats the **bili-hut™** (like a space heater). Doing so could cause the Environment High Temperature Warning to turn on.
-  **WARNING:** DO NOT place the **bili-hut™** near anything that cools the **bili-hut™** (like an air conditioner or fan).
-  **WARNING:** DO NOT place items on or in the **bili-hut™**. DO NOT spill liquids on the **bili-hut™**.
-  **WARNING:** DO NOT move the **bili-hut™** by pulling on the power cord.
-  **WARNING:** DO NOT put any creams or lotions on your baby's skin before or during treatment.
-  **WARNING:** Avoid interrupting treatment unless necessary.
-  **NOTE:** The timer only keeps track of the total lifetime use of the blue treatment lights, not your baby's treatment time.