


**Scientific Wellness
Products**

benepod

**Chronic Pain
Management**

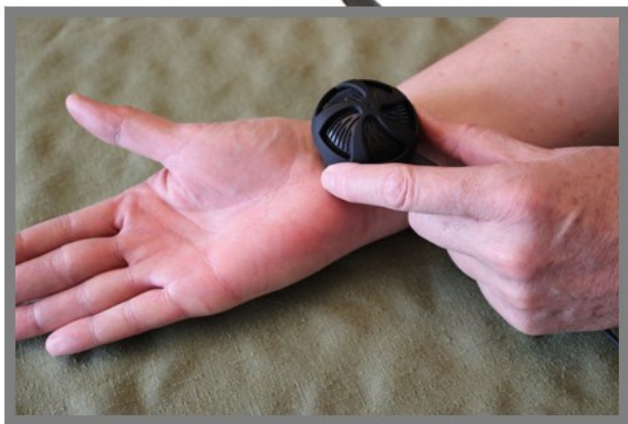

Chronic pain is more than just pain; it can interfere with your daily activities, impact your mood and prevent you from sleeping soundly. **Benepod®** is designed for people seeking relief from chronic pain including:

- ◆ **Osteoarthritis of the hands**
- ◆ **Lower back pain**
- ◆ **Chronic tendonitis conditions**
- ◆ **Chronic musculoskeletal pain**
- ◆ **Neuropathic pain**
- ◆ **Plantar fasciitis**
- ◆ **Migraines/headaches**
- ◆ **Sciatica**
- ◆ **Carpal tunnel**
- ◆ **Fibromyalgia**

Now you can regain what chronic pain has taken from you.


100% Drug Free Technology

Lightweight Slim Design



How does the Benepod® Work?

Benepod® is new type of thermal grill device which utilizes both hot and cold simultaneously to fight pain. By applying both hot and cold at the same time to a particular point on your body the **Benepod®** induces a sensation know as the Thermal Grill Illusion, targeting the same nerve receptors as drugs and competing products, in a completely noninvasive way. The sensation experienced when using **Benepod®** occurs when the body's nerve receptors can't distinguish between the hot and cold stimuli. The resulting sensation is designed to mask or even completely eliminate pain for a period of time. The patented technology used in the **Benepod®** is designed to relieve neck pain, knee pain, arthritis of the hands, headaches and various other forms of musculoskeletal related pain.

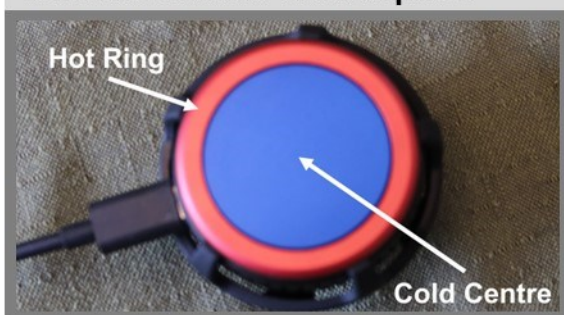

**Free Delivery in Canada & USA
(48 Lower States Only)**

Price

\$ 189.99 CDN

(Plus HST in Canada only)



Benepod® Applications

Benepod® is designed to target the same nerves as drugs, TENS units and other competing products. The Thermal Grill Illusion has been known to scientists for over a century and was developed and first studied extensively by scientists in 1896. **Benepod®** is a thermal grill device, designed to block pain by stimulating the same neuro pathways as TENS units and other similar devices.

**Patented
Technology**

**Patient
Recommended**

Advantages of Benepod® over other solutions

Cost

Benepod® requires no replacement parts, and is compatible with off the shelf accessories that many people already own. It does not require any special cleaning or sterilization, and should provide years of trouble free operation for the end user.

Safety

Benepod® operates at 5V, the same as cell phone chargers and most USB accessories, alleviating safety concerns related to products requiring a higher voltage. There are no extra wires to feed under clothing and the unit is actively cooled to prevent overheating. **Benepod®** doesn't introduce electric shocks into the body, or rely on addictive drugs to produce a result.

Ease of Use

Benepod® requires no gels, no replacement accessories, it can be used at anytime. **Benepod®** is easy to set up and use, simply plug in the power cord and wait ~30 seconds for the unit to reach its optimal operating temperature. **Benepod®** is compatible with USB accessories. Want to take **Benepod®** on the go? No problem. Because the **Benepod®** uses a standard USB-C connection. It is compatible with most of the same portable battery chargers and cell phone accessories that people already own.



Scientific Wellness Products

Division Of Cottage Camp Cabin Inc.
45 Old Trafford Drive
Hastings, Ontario, K0L 1Y0, Canada
Tel: 705-313-5460 Toll Free 800-545-2150
Email: SWP@CottageCampCabin.com
Web Site: www.ScientificWellnessProducts.com

Contact:

Gail Lebar,

Cell 705-872-6288, Email GLEbar@epna.ca

Lionel Dearing,

Cell 705-927-5555, Email

LDearling@CottageCampCabin.com