Contents

Safety Precautions	2
Warning	3
Caution	3
Contraindications	4
Indications for use	5
Read before operating the product	5
Product components	5
Product description	6
Product installation	8
Checklist before product installation	8
Product installation	8
How to use the product	9
Turning ON and OFF	9
Adjusting the product's temperature	10
After using the product	10
How to clean and store the product	11
How to clean the product	11
How to store the product	11
Trouble shooting guide	12
Check before contacting customer service	12
Product specifications	13
Product manufacturing information	14
Customer service information	16
Product warranty	17

Safety Precautions

Read and follow all instructions in this manual. Proper use of the BIO-BELT will prevent harm to the user and damage to private property.

WARNING :

Use carefully. May cause serious burns. Do not use over sensitive skin areas or in the presence of poor circulation. The unattended use of the BIO-BELT by children or incapacitated persons may be dangerous. The BIO-BELT is a medical device and should only be used by a competent adult. Elderly people and people with mobility issues need to be supervised.

- Do not alter, change or disassemble the product. Doing so may cause the product to malfunction and electric shock and fire may occur.
- Make sure to use the proper voltage for your country. Improper voltage may cause electric shock, fire or damage to the product.
- If there is a power outage, turn off the power and disconnect the plug.
- If the BIO-BELT catches fire or begins to smoke, turn off and unplug the device. Report the problem to our Customer Service Center.
- Do not allow children near the BIO-BELT unsupervised. Loss of life or injury may occur if a child wraps the power cord around their neck or trips and falls over the power cord.

CAUTION :

Property damage, device malfunction, or injury may occur if the following instructions are not followed properly. If you do not feel well when using the BIO-BELT, discontinue use and consult with a physician. Only use the BIO-BELT for heat therapy. If it is not functioning properly, discontinue use and contact Customer Service.

- Do not modify or tamper with the BIO-BELT.
- Do not jump, stand, or place heavy objects on the BIO-BELT.
- Do not place the controller on the BIO-BELT when in use. Keep the controller away from the BIO-BELT.
- Do not place the BIO-BELT on surfaces that are sensitive to temperature (i.e. leather).
- Do not cover the controller. Doing so may cause overheating of the controller.
- Do not secure the BIO-BELT with any metal products such as pins.
- Do not drop, throw or mishandle the BIO-BELT controller.
- Do not use the BIO-BELT with any ointments, or creams that generate heat. Doing so may cause burns or injury to your skin.
- Do not use the BIO-BELT with any oils without the waterproof cover. Doing so may cause damage to the BIO-BELT.
- Do not operate the BIO-BELT near objects that generate heat (i.e. ventilation from heaters, stoves) or magnetic fields (i.e. TV, magnetic objects).

- Do not use if the power cord is damaged.
- Low temperature burns can occur if using the BIO-BELT for a prolonged period of time. Use with caution when the controller is set to 104° ~ 158°F (40° ~ 70°C). Stop using the BIO-BELT if burns or blisters occur.
- Do not place the BIO-BELT directly on your skin. Place a towel in between the BIO-BELT and your skin when operating to prevent burns.
- Do not connect the BIO-BELT to any other electronic device. This may cause injury, electric shock or fire.
- Do not use on sensitive skin or on areas where circulation is poor.
- Keep the BIO-BELT away from sharp objects. Do not allow anything to puncture the folds of the BIO-BELT.

CONTRAINDICATIONS :

Do not use the BIO-BELT on infants and young children. The BIO-BELT is not indicated for the elderly and should consult with a physician prior to use.

Do not use the BIO-BELT if you have any of the following conditions:

- Person with organ transplant
- Person with high fever
- Person with cognitive disability
- Person with skin problems i.e. insensitivity to temperature change
- Person with acute disease, acute tumor, external pacemaker, the visually impaired, and any heart problems

If you have any of the following conditions, only use the BIO-BELT if directed by a physician.

- Person with immobility
- Person with internal pacemaker/defibrillator
- Person with Adrenal Suppression Disorders, Addison's Disease, Adrenal Insufficiency, Systematic Lupus Erythematosus, and Multiple Sclerosis (MS)
- Person with bypass surgery
- Person with diabetes
- Person with high blood pressure (hypertension)
- Person who is pregnant
- Person who is sick
- Person who recently received radiation/chemotherapy
- Person with surgical or silicone implants


INDICATIONS FOR USE :

For use in adults over 18 years of age. The BIO-BELT is indicated for the temporary relief of minor muscle and joint pain and stiffness; the temporary relief of joint pain associated with arthritis; the temporary relief of muscle spasms, minor sprains and strains, and minor muscular back pain; the relaxation of muscles; and the temporary increase of local circulation where applied.


Read before operating the product

Product components

※ The components may change without notice due to product improvement.

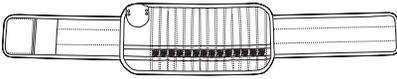BIO-BELT

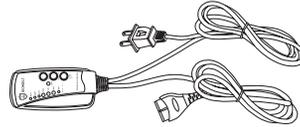
BIO-BELT Controller / Power Cable


User Manual

Product description

BIO-BELT

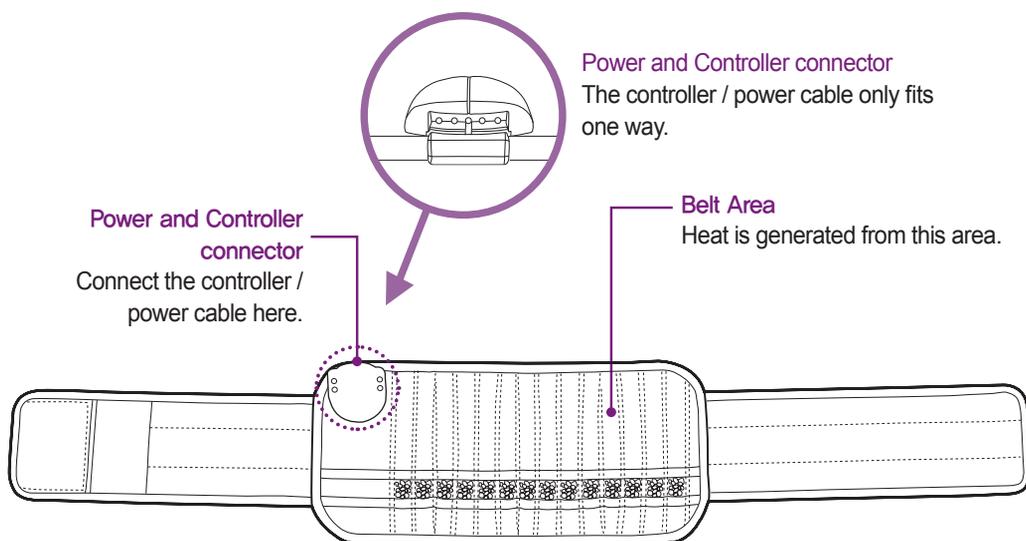Controller

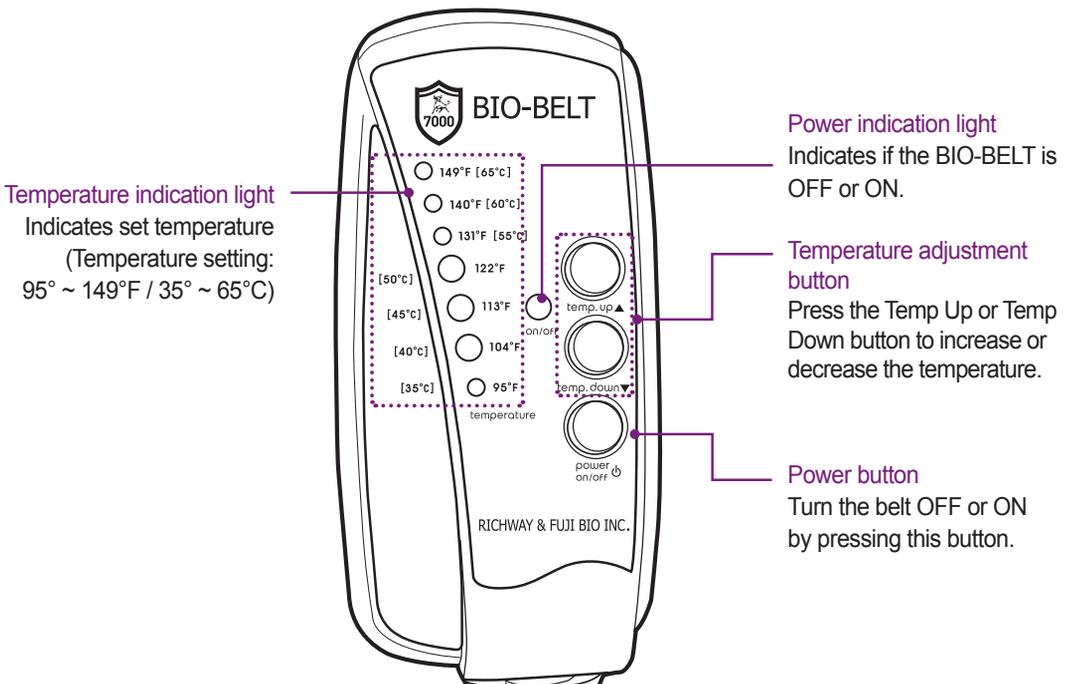Product installation

Checklist before product installation

- Do not use the product in areas with high humidity or significant dust.
- Handle the BIO-BELT with care. Any physical impact such as dropping the BIO-BELT can damage internal components.
- Check the voltage before turning on the BIO-BELT.
 - Voltage for the product: 120V
 - Electric shock or damage to the BIO-BELT can result from improper voltage.
 - The BIO-BELT and the controller must be the same voltage.
(i.e. 120V BIO-BELT + 120V Controller / 220V BIO-BELT + 220V Controller)
- Do not use the BIO-BELT if the plug, power cable, or power outlet has any damage. Always check for cuts and scrapes to the power cable.
- Do not pull on the power cord or controller when unplugging the BIO-BELT.
- Use a dedicated outlet when operating the BIO-BELT. Problems may occur when sharing a power strip with other devices.


Product installation

1. Connect the controller / power cable to the connection port of the BIO-BELT. To ensure a secure connection, there should be no space between the connector and the outlet on the mat.
The connection may be tight at first.


2. Plug the controller / power cable into an outlet. Make sure the controller is rated for the same voltage as the outlet.

NOTE: Use care when handling the controller / power cable to prevent accidental breakage. Do not drop, stand on or place heavy objects on the controller / power cable. Never pull the power cable to remove the plug from the receptacle, be sure to hold the plug.

<IMPROPER CONNECTION>
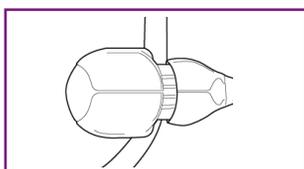
<PROPER CONNECTION>

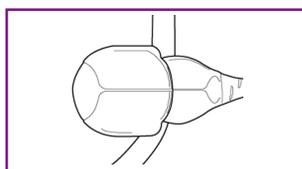How to use the product

Turning ON and OFF

Check the following to make sure the BIO-BELT is operating properly.

The BIO-BELT can be used on the legs, arms, shoulders, and waist. Its small size makes it easy to transport.

1. Turn the BIO-BELT on by pressing the Power button. The indicator light will turn on. Press the button again to turn off the BIO-BELT.
.....
2. If the light does not turn on, check for any plug/power system problems.
.....
3. Press the temperature adjustment button to see if it is working properly. Contact our Customer Service Center if the controller is not working properly.
.....
4. The BIO-BELT will automatically turn off after 2 hours of continuous use. To continue using the belt, turn it on again and set the temperature.

Caution

- Please contact our Customer Service Center if the power indication light does not turn on when the power is connected.

Adjusting the product's temperature

1. When the power is turned on, press the [Temperature Up / Down] button to set the desired temperature.
 - Range of temperature: 95 ~ 149°F / 35 ~ 65°C
 - Each time the temperature button is pressed, the temperature setting changes in increments of 9°F / 5°C.

2. When adjusting the temperature, the indicator light will change color depending on the temperature selected. Once selected, the BIO-BELT will heat up to the proper temperature. If the button of the controller is not working properly, please contact our Customer Service Center.
 - Green: 95° ~ 104°F (35° ~ 40°C)
 - Yellow: 113° ~ 131°F (45° ~ 55°C)
 - Red: 140° ~ 149°F (60° ~ 65°C)

3. Place the BIO-BELT around your legs, arms, shoulders, or waist with the gemstones faced inward (gemstones are visible through the clear window). Adjust the strap to secure the BIO-BELT. We recommend placing a towel between the BIO-BELT and your skin to prevent low temperature burns.

4. Do not place the controller on the BIO-BELT when in use. Keep the controller away from the BIO-BELT.

Warning

- Do not place the BIO-BELT directly on your skin. Doing so may cause burns or injury to your skin. Place a towel between the BIO-BELT and your skin during use.

After using the product

- After using the BIO-BELT, turn it off by pressing the Power button.
- When the power is turned off, unplug the BIO-BELT from the outlet.

How to clean and store the product

How to clean the product

1. Unplug the BIO-BELT before cleaning it. Use a moist towel to remove any dust on the surface and a soft dry fabric to clean the controller.
.....
2. We recommend wiping the clear strip on the BIO-BELT once a month with a moist towel to maintain good condition.
 - Do not leave any moisture on the BIO-BELT after cleaning the product.
 - If the BIO-BELT has any moisture or perspiration on it, use a dry towel to wipe it off.
.....
3. Do not use detergents, solvents, or flammable chemicals when cleaning the BIO-BELT and controller.

Caution

- When cleaning the BIO-BELT, do not wash with water, submerge in water, or dry clean.
 - Doing so may cause electric shock or a short circuit.
- When cleaning the BIO-BELT, make sure to turn off the power and unplug it.

How to store the product

- When not in use, turn off the power and unplug the BIO-BELT.
- Before storing the BIO-BELT, be sure to turn it off and wait until it has cooled down.
- Store the BIO-BELT lying flat (not folded) and in a well-ventilated area.
- Unplug the controller from the BIO-BELT and place it in a bag.
- Do not expose the device to water, rain or extreme humidity. It may cause electric shock or fire.

Product Storage/ Transportation

- Storage Temperature: 32 ~ 113°F / 0 ~ 45°C
- Transportation Temperature: -13 ~ 158°F / -25 ~ 70°C
- Humidity: 15 ~ 93%
- Air pressure: 700 ~ 1060 hPa

- Keep out of direct sunlight.
- Store in an area free of dust and humidity.

Trouble shooting guide

If the alarm light on the controller flashes while using the BIO-BELT, stop using it immediately. Use the trouble shooting guide to eliminate the problem. If recommended steps did not eliminate the problem and the alarm light is still flashing, stop using the BIO-BELT and contact Customer Service.

Indicator Light	Cause	Resolution method
149°F light is flashing	BIO-BELT connection is not secure	<ul style="list-style-type: none"> Check if the BIO-BELT cable is properly connected to the controller.
	Defective temperature sensor	<ul style="list-style-type: none"> Temperature sensor wire is broken. Contact our Customer Service Center.

Check before contacting customer service

If there is any malfunction during the use of the product, check the following contents before contacting our Customer Service Center.

Symptom	Resolution method
The power switch is turned on but the power indicator light is not turned on.	<ul style="list-style-type: none"> Check if the power plug is securely connected to the outlet. Check if the proper voltage is applied Connect the power plug correctly, but if it still does not properly work, contact our Customer Service Center.
The temperature is set but the BIO-BELT is not working.	<ul style="list-style-type: none"> Check if the BIO-BELT and the controller are correctly connected. Turn off the power of the controller and separate the power cable and then re-connect again. Connect the power cable correctly, but if it still does not properly work, contact our Customer Service Center.
During the use of the BIO-BELT, the power suddenly turns off.	<ul style="list-style-type: none"> Turn the power off and turn it on again.
The alarm light is on and the BIO-BELT is not working.	<ul style="list-style-type: none"> If the alarm light is on, turn off the power immediately and unplug the power cable from the power outlet and then contact our Customer Service Center.

Product specifications

Category	Contents
Product	BIO-BELT
Product size	210 x 460 mm
	1) main body of the BIO-BELT containing semi-precious stones: Length 46 cm x width 21 cm 2) Total length including the belt: 140cm (Regular) or 160cm (XL)
Weight or packing unit	1) Controller: 0.38 KG 2) Main body of the belt: 0.93KG 3) Overall weight: 1.6 KG
Maximum temperature	149°F / 65°C (The BIO-BELT is programmed to prevent overheating.)
Voltage, frequency, power consumption	120V~, 60Hz, 50W
	Operating temperature : 95 ~ 149°F / 35 ~ 65°C Variance: +/- 9°F / 5°C
Type of protection against electric shock	Class II ME equipment
Degree of protection against electric shock	Class BF applied parts - Applied parts: Top of Mat
Degree of protection against harmful ingress of water	Ordinary equipment (IPX0) - Equipment not suitable for use in presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.
Mode of operation	Continuous operation


※ The BIO-BELT will automatically turn off after 2 hours.

Class II ME equipment (single fault safe)		
	Normal Mode	Single Fault Mode
Touch Current	< 100 microA	< 500 microA
Ground Leakage Current	< 5 mA	< 10 mA


Testing in accordance with IEC 60601-1-2 indicates the following:
The BIO-BELT needs no special precautions regarding EMC and needs to be installed and put into service.
Portable and mobile RF communications equipment have no adverse effects on the BIO-BELT.

Product manufacturing information

Product label is attached to the bottom of the BIO-BELT and the controller.

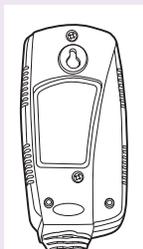Label location on the BIO-BELT

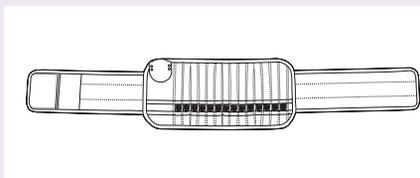Label location on the controller.

BIO-BELT product

Front

Side

Customer service information

If your BIO-BELT is not working properly call our Customer Service Center.

CUSTOMER SERVICE CENTER

North America

Richway & Fuji Bio Inc.

Tel: 1-808-589-2800

Toll-Free: 1-855-338-6410 (USA, Canada, Guam, Saipan)

Homepage: www.richwayandfujibio.com

Email: rma@richwayusa.com

Product warranty

Order information	Product name	BIO-BELT	Model name	
	Warranty period	Limited 3 years from the date of product delivery	Manufacturing number	
	Date of purchase		Branch of purchase	
Customer information	Name		Purchase price	

Richway & Fuji Bio warrants the covered product to be free from defects and malfunctions at the time sold. "Defects" as used in this warranty, is defined as imperfections, which impair or impede the function and use of the product.

How Long the Warranty Lasts

This limited warranty lasts 3 years from the date of delivery. This warranty extends to the original purchaser only. Coverage terminates if you sell or otherwise transfer ownership.

What This Warranty Covers

This warranty covers any defects or malfunctions in the electronic components, heating element, and top cover replacement. If Richway & Fuji Bio determines that such goods are defective, within the period of this warranty, Richway & Fuji Bio will repair or replace the parts free of charge for the first year. Repair or replacement for the second and third year requires a fee, which can be found on our website. Repair or replacement is the sole and exclusive remedy available to the buyer, and shall constitute the fulfillment of all Richway & Fuji Bio obligations in respect to the goods furnished under this agreement.

What This Warranty Does Not Cover

This warranty does not cover abuse, misuse, natural disaster, insufficient or excessive electrical supply, ordinary wear, fading, failure to follow instructions and improper maintenance. Richway & Fuji Bio shall not be liable for any incidental or consequential damages as a direct result of negligence or misuse. This warranty does not extend to the purchase of opened, used, repaired, repackaged, and/or resealed products, including but not limited to sale of such products on internet sites and/or sales of such products on internet auction sites and/or sales of such products by surplus or bulk resalers. Any and all warranties or guarantees shall immediately cease and terminate as to any products or parts thereof which are repaired, replaced, altered, or modified, without the prior express and written consent of Richway & Fuji Bio.

How to Get Service

To obtain service under this warranty, contact your local Customer Service Center and provide the purchaser's name, RI (Richway Invoice) number, and a detailed description of the problem you are experiencing. They will provide you with a return merchandise authorization number (RMA). This number is required and must be placed in a visible area on the front of the package. Ship your product postage paid, include a copy of the original sales slip, credit card receipt or other proof of the date of the original retail purchase. If servicing within the first two (2) months, we will send you a prepaid shipping label to cover your shipping costs. If servicing from the third month and on, the customer is responsible for shipping costs to return the product for repair or replacement. We will cover the shipping costs to return the product back to the customer.

How Do State Laws Apply

This warranty gives you specific legal rights, and you may also have other rights which vary from state to state. Any implied warranties last the duration of this warranty. Some states do not allow limitations on how long an implied warranty lasts and some states do not allow the exclusion or limitation of incidental or consequential damages, so the limitations or exclusions may not apply to you.

**THE
RICHWAY
GROUP**

RICHWAY & FUJI BIO INC.

1314 South King Street Suite 520, Honolulu, HI 96814 USA

Tel: 1-808-589-2800 / Toll-Free: 1-855-338-6410

Homepage: www.richwayandfujibio.com

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Manufacturer: Richway & Fuji Bio Inc.
Contract Manufacturer: R&L Co., Ltd.