


# Simple V3 Remote Manual




**1**  
Press the  
Power button.

**2**  
Press the  
Automatic  
Mode button.

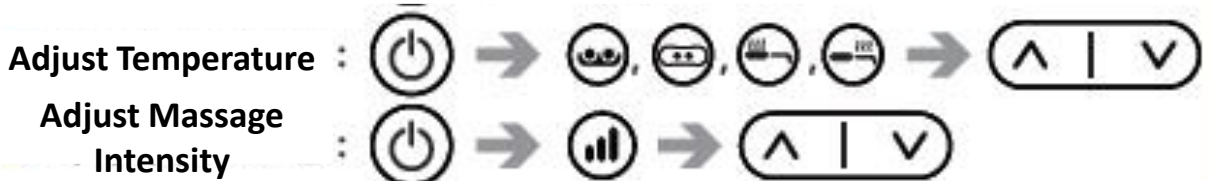


## Remote control button layout

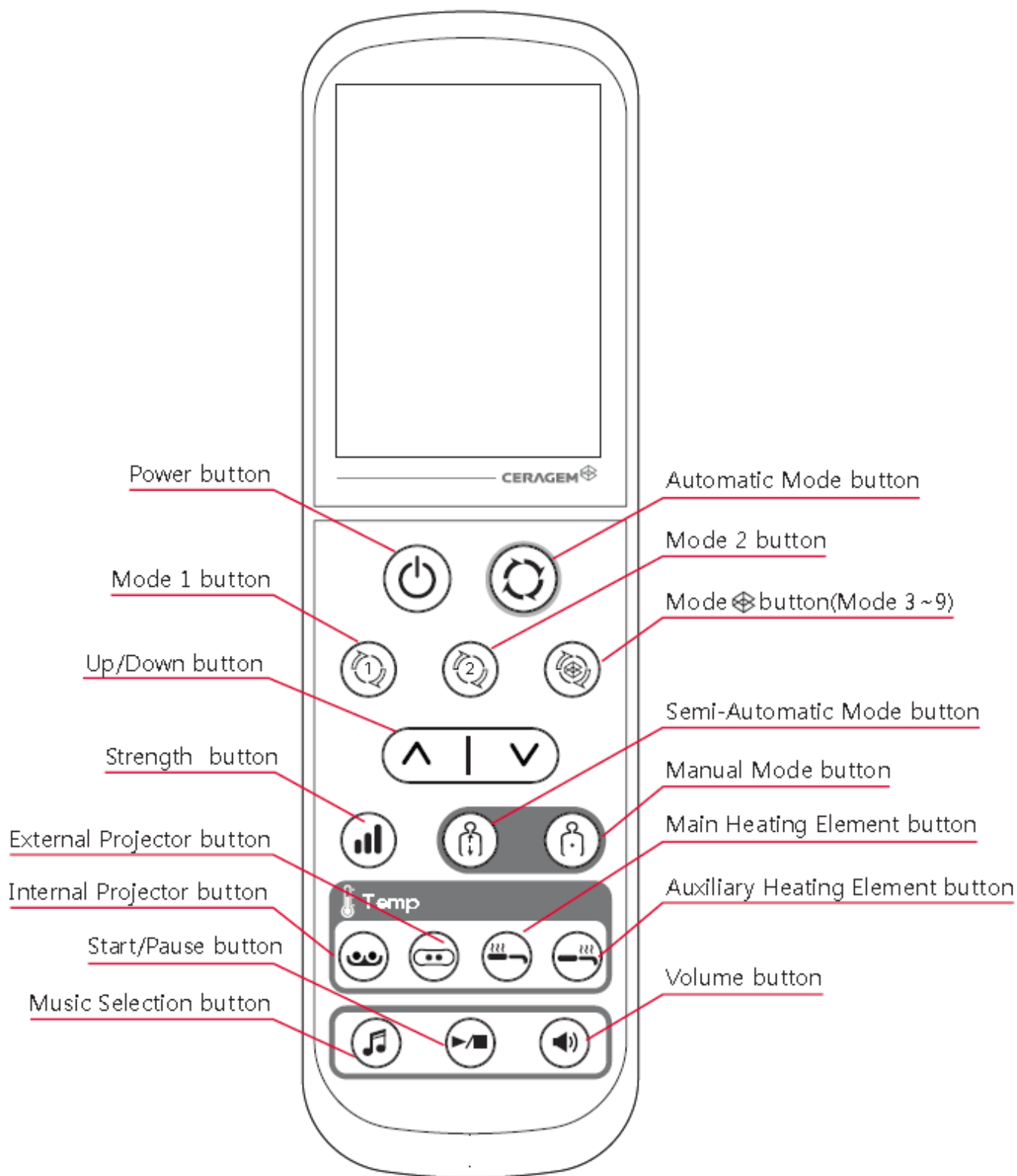



## Turning on the Remote Control


- To turn on: Press the power button.
- To turn off: Press & hold the power button for more than 1 second.

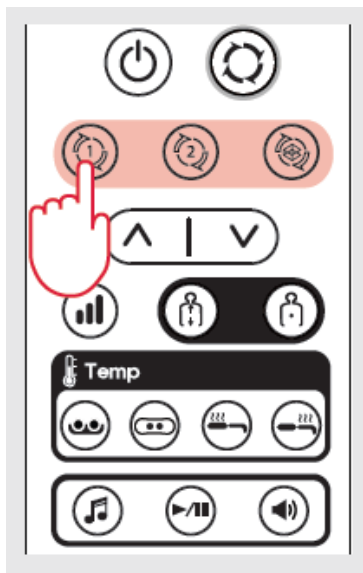

## Selecting the Modes

- Pressing the mode 1, 2, or 3-9 button twice (once for Manual mode) will activate the Thermal Massager.


	Press twice to select Mode 1
	Press twice to select Mode 2
	Long-press. (Select modes 3~9 in sequence.)

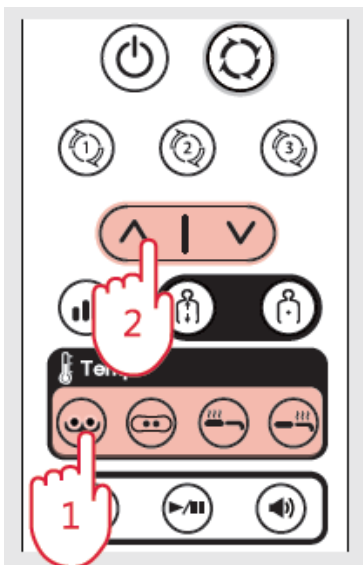

## Adjusting the Temperature


- 1) Press the temperature button.
- 2) Press the Up/Down button to set the temperature to the desired level.

Internal Projector	External Projector	Main Heating Element	Auxiliary Heating Element
OFF, 30°C ~65°C		OFF, 30°C ~60°C	

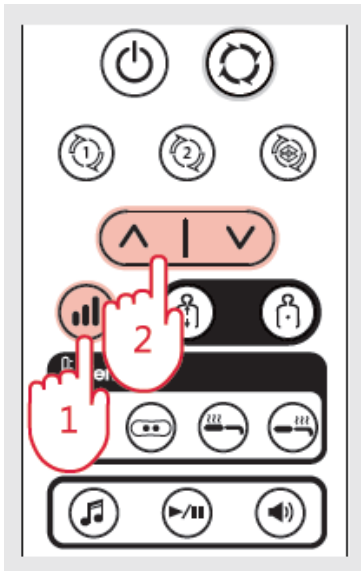




## Adjusting the Level of Intensity

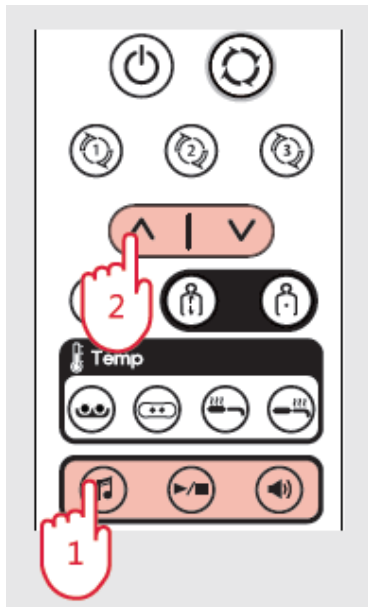- 1) Press the strength button.
- 2) Press the Up/Down button to set the level of intensity to the desired level.




## Adjusting the Sound

- 1) Press the sound button.
- 2) Press the Up/Down button to set the sound to the desired level.

		
Music Selection	Start/Pause	Volume
Ceragem Sound : 1~10 User Sound : 11~100		0~10

## Automatic Mode


(Mode Duration: Approx. 36 minutes 30 seconds)

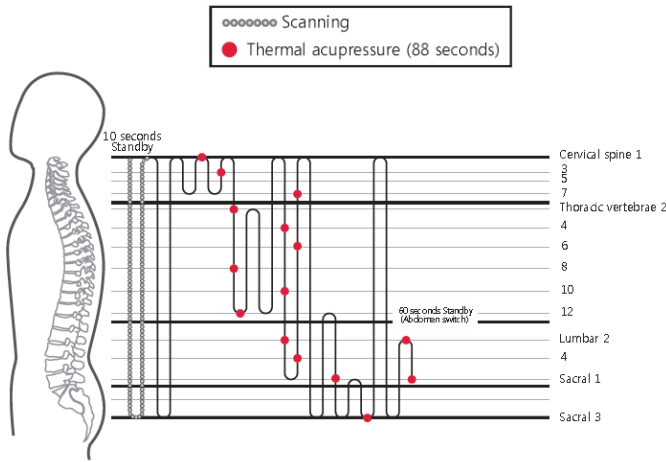

## Mode 1


(Mode Duration: Approx. 36 minutes 30 seconds)

● Thermal acupuncture (100 seconds)

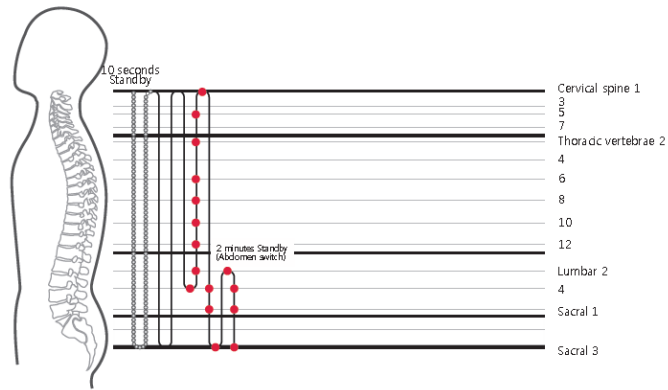


## Mode 2

(Mode Duration: Approx. 18 minutes)

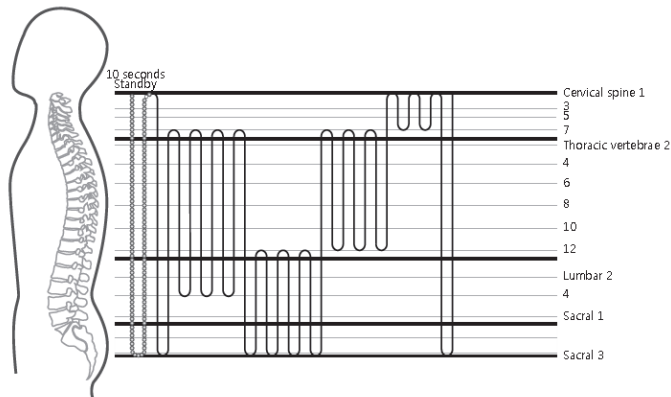


## Mode 3

(Mode Duration: Approx. 36 minutes 30 seconds)

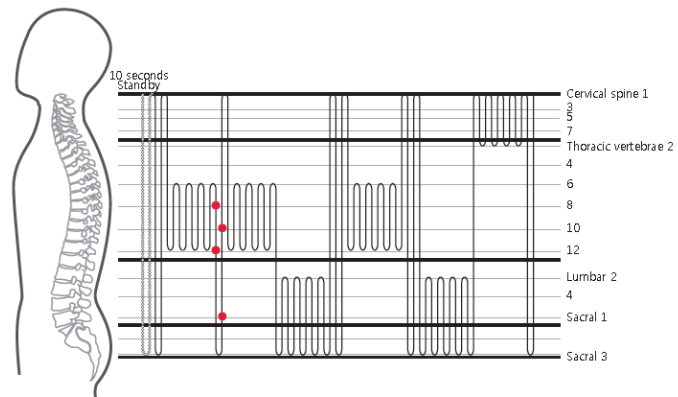


## Mode 4

(Mode Duration: Approx. 36 minutes 30 seconds)

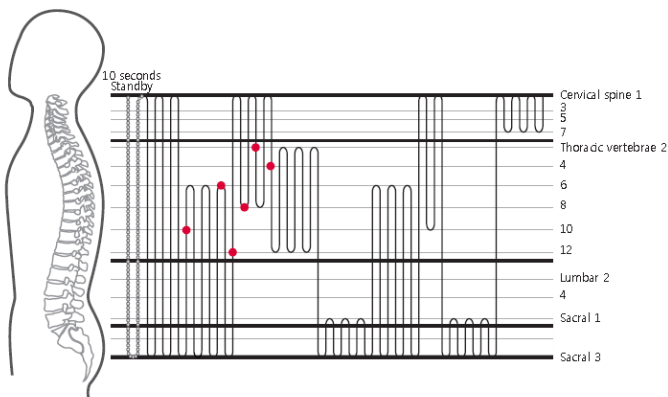


## Mode 5

(Mode Duration: Approx. 36 minutes 30 seconds)

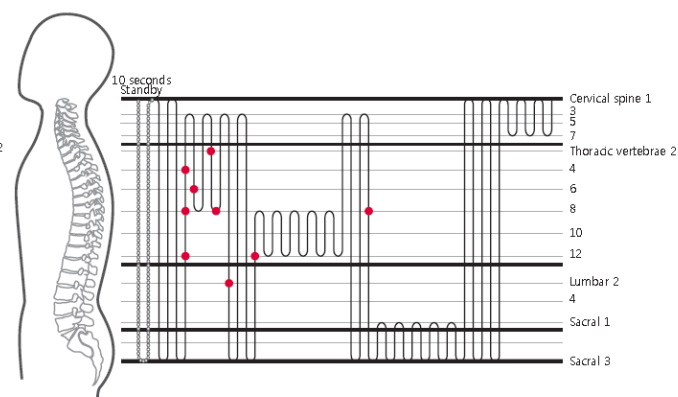


### Mode 6

(Mode Duration: Approx. 36 minutes 30 seconds)

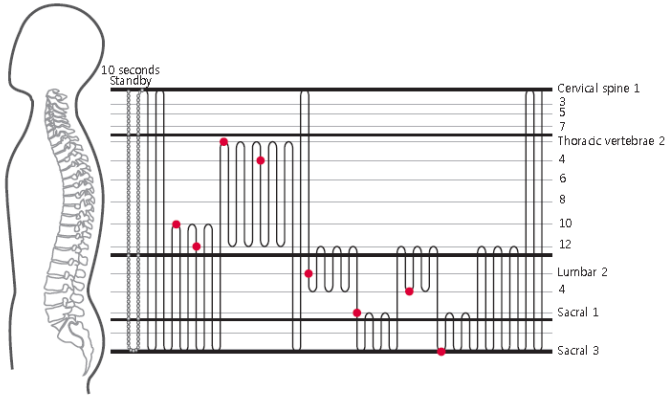


### Mode 7

(Mode Duration: Approx. 36 minutes 30 seconds)

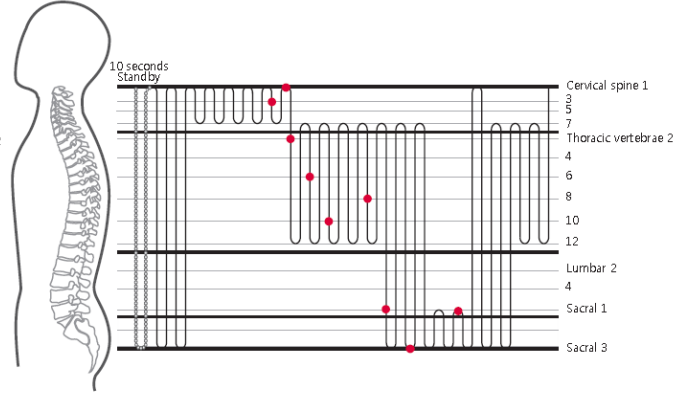


### Mode 8

(Mode Duration: Approx. 36 minutes 30 seconds)

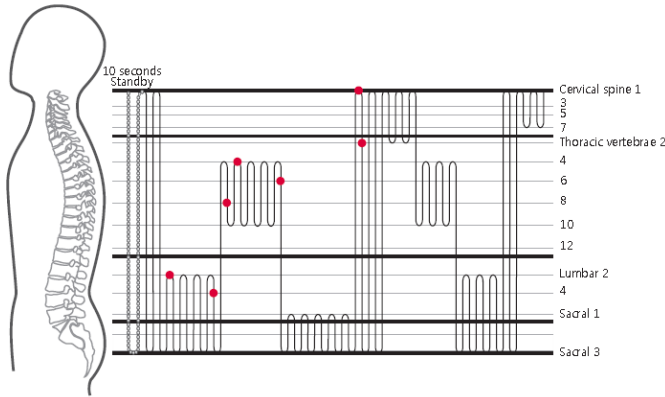

### Mode 9


(Mode Duration: Approx. 36 minutes 30 seconds)

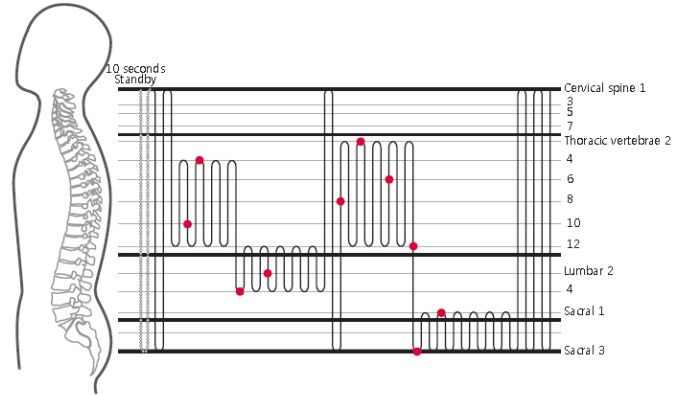


### Semi-Automatic Mode

(Mode Duration: Approx. 18 minutes)

Massaging is repeated on a specified region. This region depends on the final position of the internal projectors.

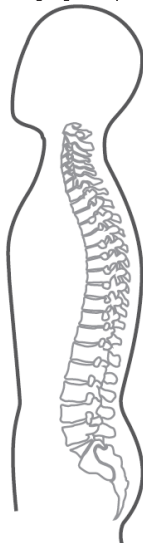

● Up&Down(3 Times) + Thermal acupressure (60 seconds)

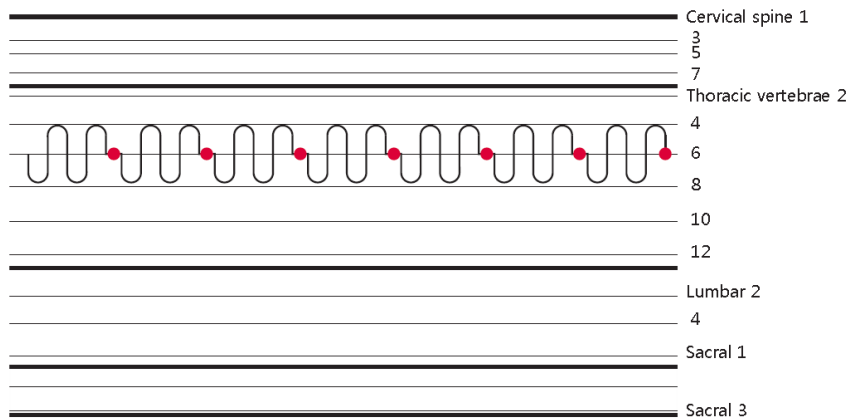

※ Users may adjust the position.