



## Microneedling System

Microneedling helps to stimulate blood flow into the scalp and promote hair growth. It can also help topical hair growth treatments to penetrate more efficiently into the scalp.

**FINE NEEDLES • ECO-FRIENDLY HANDLE • FDA & CE & ISO**




NEEDLE SIZE

0.50MM

**Compatible with Triple hair, RIZN and Plenty Natural topical hair growth treatments**

### DIRECTIONS FOR USE (ADULTS)

1. Clean/Sterilize the Triple Hair Derma Roller with alcohol and let it dry before using. Always use on a clean, towel dry scalp. Wash your hands.
2. Part hair and gently glide the Triple Hair Derma Roller back and forth 3 times. Repeat vertically, horizontally and diagonally on the targeted scalp areas.
3. Apply a Triple Hair growth spray on the targeted scalp areas for optimal results. Massage lightly until absorbed.
4. Use the Triple Hair Derma Roller two to three times per week.

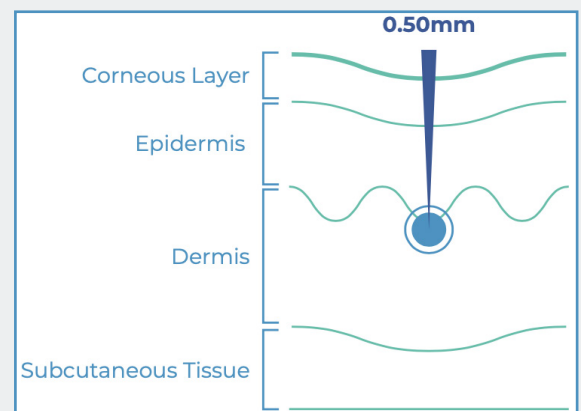

### CAUTION

Use the Derma Roller on the scalp only. Do not use the Derma Roller on an infected or irritated scalp. If there is any scalp sensitivity, stop using the Derma Roller immediately and consult a doctor or a dermatologist. If you are unsure about your scalp condition, consult a dermatologist before using the Derma Roller. Do not share the Derma Roller as it should be used exclusively by the same person. Do not use the Derma Roller through loose hair. Do not use the Derma Roller if there are bent needles. Store the Derma Roller carefully in the box provided and keep out of reach of children.

Learn more about our hair growth treatments on our website: [www.triplehair.ca](http://www.triplehair.ca)