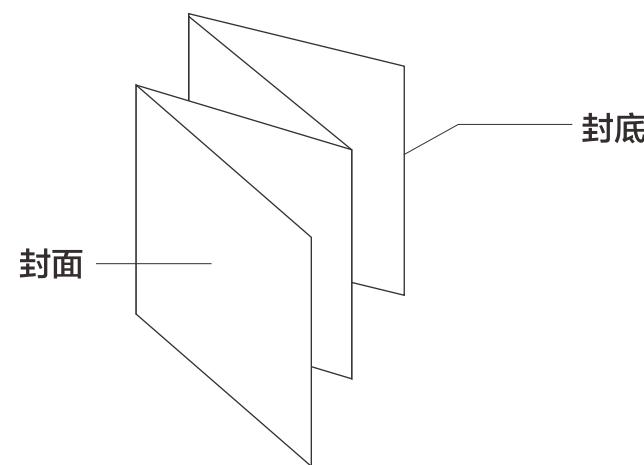



成品尺寸：85x120mm

展开尺寸：595x120mm



USING METHOD

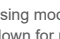Kindly note that turn Water Flosser off and charge over 3 hours before first use.

- Insert nozzle into the center of the knob at the top of the water flosser handle. You will hear a "click" once the nozzle is correctly locked into (Suggest to replace the nozzle every 3 months for better results. If the nozzle is loose or blocked by foreign objects, it should be replaced).
- Disconnect charging cable from the unit prior to filling the reservoir and product use. Lift flip top on reservoir or take the water tank, and fill the reservoir with lukewarm water.
- Close the nozzle collection cover or rotate and locked the water tank to prevent the water tank from overflowing.
- Shortly press the mode button"  " to switch the flossing mode.  
Note: The device will keep same mode with last shutdown for next working due to the memory function.

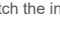
1 Standard Spray Mode: Holding the power button to spray water, flossing after 1 second, and release the button to stop spraying water.

2 Soft Mode: For the users with sensitive teeth or first time users.

3 Strong Mode: For the normal cleaning of teeth with most of peoples.


4 Pulse Mode: Alternating strong and weak, deep cleaning, and massaging gums.


5 Pressure Control: Water pressure adjustment. Low in Blue light, Medium in White light, Strong in Orange light.


- Shortly press the pressure control button"  " to switch the intensity of the water spray.
- Ensure the nozzle kept parallel to your teeth and gums. Direct the jet stream at a 90-degree angle to your gumline. Slightly close lips to avoid splashing but allow the water to flow out into your sink.
- Shortly press the power/flossing button, start flossing with white lighting. Support to switch gears and modes when working. The longest flossing time is 90 seconds then turn off automatically.
- Aim the tip toward teeth, glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.


USING METHOD

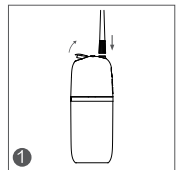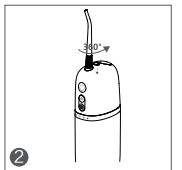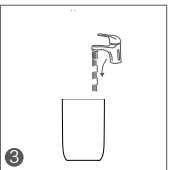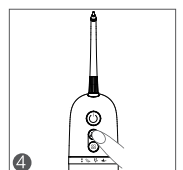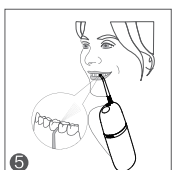
Kind Note:  
The device could be power on with Power Button/Flossing Mode Button/Pressure Control Button. Adjust any mode and water pressure to cleaning. And then start flossing after pressed the Power/Flossing Button. The device will automatically turn off after 5 seconds if no any operation.

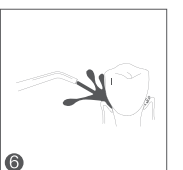













CLEANING & CARE


Do not use any corrosive or abrasive detergent as it will damage the water flosser.

**Cleaning**

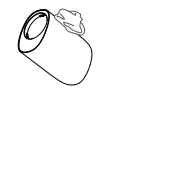- Press and hold the nozzle release button to remove the nozzle, rinsed with water. Please clean every time after used.
- The water tank can be detached, clean it with a soft cloth to wipe exterior and inside.

- Please do the clean as above steps and store in a dry environment if you plan not to use this product for a long time.





TROUDIESHOOTING

Problem	Possible Cause	Solution
Do not work.	(a) New one, bought just now. (b) Stored over 3 months.	(a) Charge the water flosser over 3 hours. (b) Contact service.
The water flosser only operates for several minutes after charged.	(a) No fully charged, low-battery. (b) At the end of battery lifetime.	(a) Charge the water flosser over 3 hours. (b) Contact service.
The water pressure is low.	(a) The nozzle is blocked. (b) The nozzle is deformed. (c) The pump is blocked.	(a) Replace the nozzle. (b) Contact service.
No flossing happened.	(a) There is no water in the reservoir. (b) There is not correct angle when using the flosser. (c) The water pipe is not esssembled well. (d) The pump is blocked.	(a) Fill water into the reservoir. (b) Direct the nozzle at 90-degree angle to gumline. (c) Esssembled the water pipe well. (d) Contact service.

DECLARATION FOR HAZARDOUS SUBSTANCES

The hazardous substances in this product, classified by component, are declared below:

Declaration for Hazardous Substances							
Part Name	Pb	Hg	Cd	Cr(VI)	PBBs	PFDEs	
Plastic Casing	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Metal	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
PCB	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Motor Components	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Battery	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	
Charging Cable	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	

**NOTE:**  
O: Indicates that the content of the hazardous substance in all homogeneous materials of the component is below the limit requirement specified in GB/T 26572.  
X: Indicates that the content of the hazardous substance in at least one homogeneous material of the part exceeds the limit requirement of GB/T 26572.

**Disposal**  
This appliance should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

WARRANTY

\* There is some delay about warranty period on the time of selling to different countries.  
\* We grant one year warranty on the water flosser no including nozzle starting on the date of purchase.  
\* We grant one year warranty on the charging cable starting on the date of purchase.

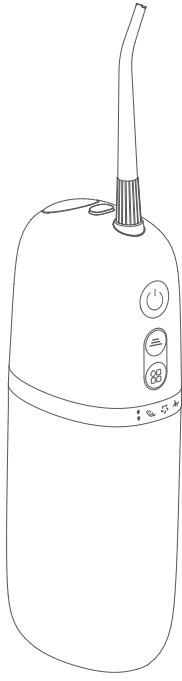
Within the warranty period, we will check faults in materials or workmanship free of charge and replace the complete device, at our discretion.

This warranty does not cover damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The flash tube is considered as wearable part of the device and is therefore not covered under the warranty.

www.imass.com.au

12

IMASS Cordless Portable Water Flosser



Model No.: Oral Care Express Pro

Ensure you have read and fully understood all these instructions. Retain this user guide for future reference.

CONTENTS

Important Safeguards ----- 2

Product Parameters ----- 4

Product Description ----- 5

Charging ----- 6

Using Method ----- 7

Cleaning & Care ----- 9

Troubleshooting -----10

Declaration of Hazardous Substances ----- 11

Warranty ----- 12

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed.

**DANGER:**

To reduce the risk of electrocution:

- Always unplug product after using.
- Do not handle plug with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord if damage before the first use and during the life of the product.

**WARNING:**

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device.
- If product is opened/disassembled for ANY reason, warranty is VOID.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any water. Contact manufacturer or your local distributor.
- Do not directly wash under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by original factory.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

IMPORTANT SAFEGUARDS

- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. They can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Do not use this product for more than 5 minutes in each two-hour period.

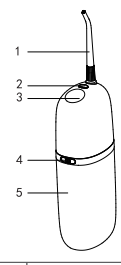
**NOTE:**

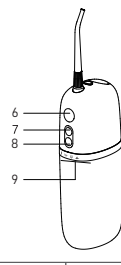
- Closely instruct and supervise children 8 years and above and individuals with special needs in the proper use of this product.
- Children should be supervised to ensure that they do not play with the appliance.
- There are no consumer serviceable parts in this appliance and it does not require maintenance.
- Fill reservoir with warm water or other dental professional recommended solution only.

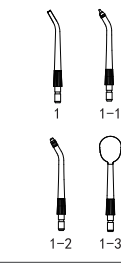
PRODUCT PARAMETERS

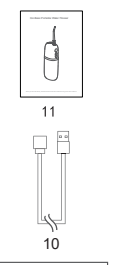
Product Name	Water Flosser
Product Size	66*49.5*166mm(Excluding nozzle height)
Reservoir Capacity	200ML
Water Pressure	About 60~130psi
Frequency	1400-1800times/min
Nozzle Diameter	0.06mm
Flossing Mode	Standard(Spot-spray), Soft, Strong, Pulse
Rated Input	5V 1A
Rated Power	5W
Battery Capacity	1000mAh
Charging Port	Type-C Charging Cable
Waterproof	IPX7
Time of Long-lasting	About 30 days(1 minute per day based on Strong mode)
Accessories	Device * 1pc Nozzle * 4pcs User Manual * 1pc Type-C Charging Cable * 1pc

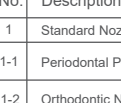
PRODUCT DESCRIPTION

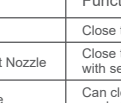


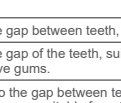


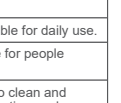


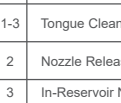


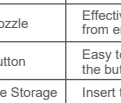


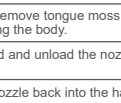


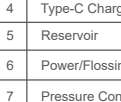





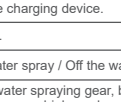












No.	Description	Function
1	Standard Nozzle	Close to the gap between teeth, suitable for daily use.
1-1	Periodontal Pocket Nozzle	Close to the gap of the teeth, suitable for people with sensitive gums.
1-2	Orthodontic Nozzle	Can close to the gap between teeth to clean and washable braces, suitable for orthodontic people.
1-3	Tongue Cleaner Nozzle	Effectively remove tongue mucus to avoid bacteria from entering the body.
2	Nozzle Release Button	Easy to load and unload the nozzle after pressing the button.
3	In-Reservoir Nozzle Storage	Insert the nozzle back into the handle.
4	Type-C Charging Port	Connect the charging device.
5	Reservoir	Store water.
6	Power/Flossing Button	Start the water spray / Off the water spray.
7	Pressure Control Button	Adjust the water spraying gear, blue low-grade, white mid-grade, orange high-grade.
8	Mode Button	Adjust the point of spray, left, strong, pulse.
9	Mode Indicator	Select the mode corresponding indicator lights up.
10	Type-C Charging Cable	Connect the product and power adapter for charging.
11	User Manual	Product instructions of guiding users.

CHARGING

**Charging The Water Flosser:**

Plug the Type-C charging into the charging port of water flosser, and USB connector to a power adapter with 5V output(No included). Suggest to charge over 3 hours if meet below situation:

(1) It is the first time to use this device.  
(2) The device is to be stored for an extended period of time (3 months or more).  
(3) The device is low in power or stops working during use.

- It is low battery when the power indicator light in orange begin blinking.
- Do not immerse both connectors of Type-C charging cable in any liquid.
- Make sure the Type-C charging cable is dry before connecting to the power supply.
- Plug the Type-C cable into the charging port of the flosser, and USB connector to 5V power adapter (Phone adapter and USB port on computer).
- It is fully charged when the power indicator light remian on.
- The water flosser could be used for 50 minutes roughly after being fully charged.
- It may shorten the lifespan of the batteries after used for a long time.

**Note:**

- Please put the metal objects such as coins or clips away from USB charging cable while charging, as this may cause the device to overheat.
- It is normal with heating when charging. If still shorun working after fully charged, it maybe at the end of the battery life. Please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.
- Use this device at the temperature between 0 °C~35 °C.