


Oral Irrigator Instruction Manual

MODEL: ORAL CARE EXPRESS

Thank you for choosing IMASS
Read thoroughly before operating
Keep for future reference

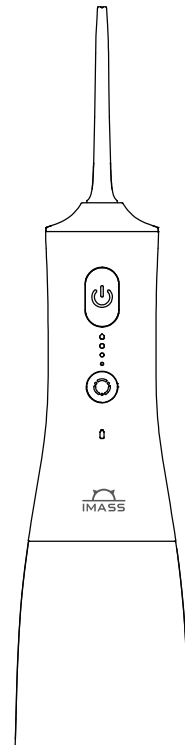

www.imass.com.au



1. Important Safeguards

All information in this section is provided to instruct you on the correct and safe operation of this device so that you can avoid injury to yourself and others as well as damage to property.

- Toddlers, infants, people who are unable to operate the appliance should not use it by themselves.
- Keep away and out of reach from children.
- Never try to alter, disassemble or repair the appliance.
- Do not place or store the appliance where it will fall or end up in a tub or sink.
- Do not reach for the appliance if it has fallen into water. Unplug immediately.
- Do not use the appliance while bathing.

2. Warnings

- This appliance can be used by children from age 8 and above; and people with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The appliance can be only be used with the provided power supply unit.
- Do not fill the tank with water warmer than 40 °C.
- If the power cord is damaged, it shall be replaced by the manufacturer, a service agent or a similarly qualified person in order to avoid any hazards.
- Keep the cord away from heated surfaces.
- Do not use the appliance without a nozzle.
- Do not submerge the appliance in water.
- Do not use for any purpose other than oral cleaning.
- Keep nozzles away from young children. The nozzle is small enough to swallow and can be a choking hazard.
- Certain types of mouthwash may damage the equipment by cracking the housing or water reservoir. Do not use mouthwash in the water reservoir.
- People who experience tooth or gum pain may not be able to use this appliance and should consult with their dentist before use.
- Do not press the nozzle firmly against your teeth or gums, it will hurt your gums if you do so.
- When using the periodontal nozzle, turn on the flosser by setting it at the lowest pressure (Soft Mode).

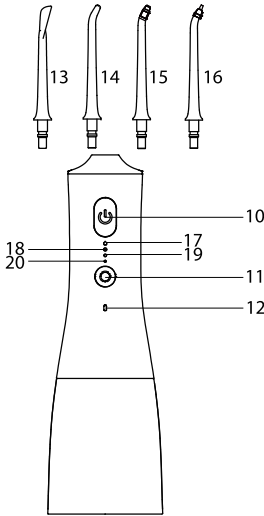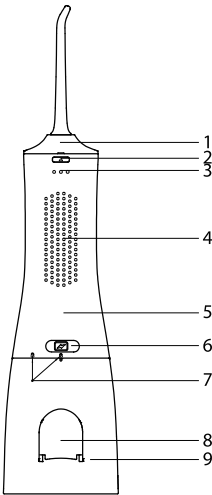


- | | |
|---|----------------------------|
| 1. Cover | 11. Mode switching button |
| 2. Nozzle release button | 12. Power indicator |
| 3. Over-blow hole | 13. Tongue Cleaner |
| 4. Anti-slip padding | 14. Standard jet tip |
| 5. Main body | 15. Orthodontic tip |
| 6. Type-C rubber plug and charging port | 16. Periodontal pocket tip |
| 7. Water tank opening and closing | 17. Pulse Mode |
| 8. Water tank door | 18. Strong Mode |
| 9. Water tank | 19. Normal Mode |
| 10. Power Switch | 20. Soft Mode |

3. How to Use

Charging the appliance (Fig. 1)

- Remove the rubber plug from the charging port and insert the charging cable.
- Attach the charging cable to the adapter and plug it into a power outlet.
- The power indicator blinks while the appliance is charging. When it is fully charged, the power indicator will stay on for 2 minutes, before turning off automatically.
- Reattach the rubber plug to seal the charging port after charging.

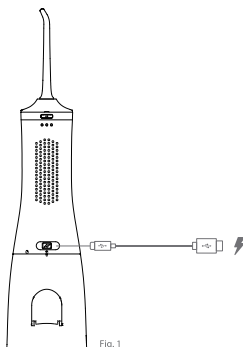



Note:

If the power indicator is blinking during or after an irrigation, that means the battery needs to be recharged. A full charge will take about 4 hours.



Using the Water Flosser

Make sure the water tank door is sealed tightly before using the water flosser to avoid water leakage.

- Attach a nozzle onto the unit. (Fig. 2)
- Fill the tank with water:
 - a. Open the water tank door to fill the tank with water that does not exceed 40°C. Close the door once it is filled. (Fig. 3)
 - b. Or, slide the water tank down from the main body and fill the tank with water that does not exceed 40°C. Reattach the water tank to the main body once it is filled. (Fig. 4)

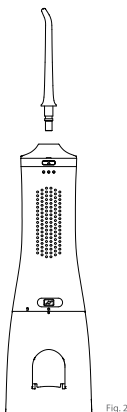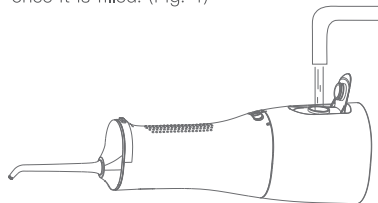

- Power on the appliance by pressing the “⏻” button.
 - Press the “⊙” button to cycle the active mode between Normal, Soft, Strong and Pulse. If the appliance has been used before, the previous mode will be retained upon powering on. (Fig. 5)
- Hold the appliance in an upright position with the nozzle tip pointing toward your teeth.
- Keep your mouth opened slightly so that the water can flow out.

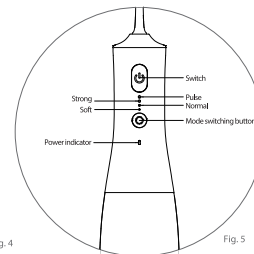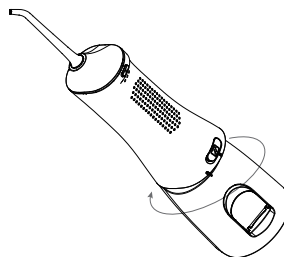

Mode Introductions

Soft: for gentle cleaning of sensitive gums

Normal: for normal cleaning


Strong: for deep cleaning, with high water pressure

Instant strong mode: To temporarily enter Strong mode, start at either Soft or normal mode. Hold down the “⊙” button to engage Strong mode. Once the button is released, the active mode will revert back to what you originated with.

Pulse: for cleaning and massaging the gums

Press “⏻” button to start irrigating. (Fig. 6)

- 1) Direct the water jet slowly at a 90° angle to the teeth and gums.
- 2) Direct the water jet slowly along the teeth.
- 3) Direct the water jet along the gum line with the nozzle tip at a 90° angle to the gum tissue.
- 4) Clean between the teeth and around any braces, crowns or bridgework thoroughly.

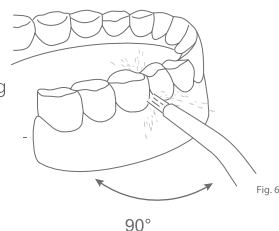

- Press the “⏻” button to power off and then press the nozzle release button to install other nozzles as needed.(Fig.7)
- Once irrigating is complete, press the “⏻” button again to turn off the appliance.

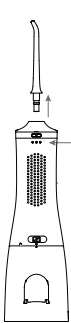


Fig.7

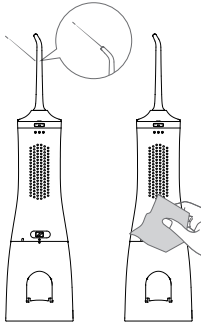
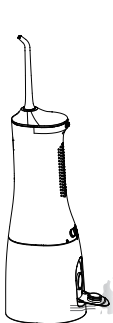



Fig.8

After using the appliance (Fig.8)

The appliance must be powered off prior to removing any excess water.

- Drain the water inside the appliance by detaching the water tank or through the water tank door.
- Reattach the water tank to the main body and close the water tank door.
- Wipe the appliance with a dry cloth.

4. Specifications

Model: **Oral Care Express**
 Battery: 1500mAh
 Pulse frequency: 1800rpm
 Noise: Max 65db
 Working temperature: 0–50°C
 Input voltage: DC5.0V  1.0A
 Waterproof rating: IPX7

5. Troubleshooting

Problem	Possible Cause	Solution
It does not work	You have just purchased the product or it has been idle for more than 3 months	Charge the unit for at least 4 hours
It only works for a few minutes even after charging	The battery has reached the end of its life cycle	Contact sales for help
	Charging time was insufficient	Charge the unit for at least 4 hours
Water pressure is very low	There is dirt stuck in the nozzle	Replace the nozzle
	The nozzle is deformed	
	There is dirt stuck in the filter	Clean the filter
Water does not come out	The tank is empty	Fill up tank
	The appliance is tilted too much during usage	Hold the appliance vertically during usage

6. Warranty

IMASS offers an extensive warranty for all products that extends through 1 YEAR of the order date. This warranty covers quality related issuesand product defects.

If you are past the 30–day return period on Amazon, contact us directly via the “Contact seller” option on your Order Details. If the product is not working properly, we can send you a replacement product. This warranty does not apply in case of abuse, mishandling, or unauthorized repair of the product.

If you have any question about this warranty or would like additional information, please contact us via the “Contact Seller” option or directly via **web-site www.imass.com.au**

Please keep the original box, packing materials and all documentation in the event that service is required.