# AVID IF<sub>2</sub><sup>TM</sup>

INTERFERENTIAL SYSTEM

## Universal Patient Protocol Protocolo Universal para Pacientes

**STOP** **DO NOT USE** if you have a pacemaker or are pregnant. For complete list of warnings and precautions, see back page.

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### ENGLISH SETUP INSTRUCTIONS

**Step 1:** Attach lead wire to Avid IF2.

**Step 2:** Attach each electrode to the lead wire—connect red electrode tips to red lead wire tips and black electrode tips to black lead wire tips.

**Step 3:** Place electrodes on skin at least 1-2 inches apart (see Electrode Placement section).

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## ENGLISH STARTING TREATMENT

**Note:** You may be asked to answer questions at the beginning and end of each treatment. Press the **+** and **-** to scroll through answers and the **✓** to select an answer.

**STEP 1:** Press the **⏸** button at the bottom of the device.

**STEP 2:** To select a preset, press the **◀** or **▶** arrow until the preset number you want is shown on the screen. Press **✓** to select.

- **Preset #1** provides 40 minutes of treatment (two 20-minute segments) for circulation and pain control. Use this preset two to three times per day, as needed throughout your rehabilitation.
- **Preset #2** provides continuous pain control. Use this preset as needed for pain control throughout the day or night.
- **Preset #3** provides 20 minutes of treatment for circulation followed by continuous treatment for pain control.

**STEP 3:** Press the **+** to start treatment by increasing the stimulation to a comfortable intensity. To reduce intensity, press the **-**. The screen will display the intensity and count down your treatment time (for continuous presets, the timer counts up). If there is an asterisk after the timer, it means that another preset will follow once the time reaches 0.

**NOTE:** The Avid IF2 has a patient safety feature that locks the screen after about 30 seconds to prevent accidental increases in intensity. If you get a **🔒** icon on the screen, press the **🔒** button, then press the **+** and **-** to adjust intensity.

**To Pause treatment:** Press the **⏸** button on the screen. To unpause, press the **▶** button

**Vectoring:** Vectoring the stimulation changes the intensity of the stimulation of each channel to “target” your pain. Vectoring is an option on the Avid IF2 when the intensity is greater than 10 mA. To vector, press the **✖** symbol in the upper right corner of the treatment screen. Then press either the **◀** or **▶** arrow to vector towards one channel (either reds or blacks). Press **✓** to return to the treatment screen.

**STEP 4:** To stop treatment and return to the main menu: Press the **🏠** symbol at the top of the screen. Then press either **●** to go to the main menu or **▶** to resume treatment.

## COMPONENT USE AND CARE

**ELECTRODES: IMPORTANT** — Read the instructions for the application and use of electrodes found on the electrode package and in the Patient Operating Manual.

- Reusable, self-adhering electrodes are included in the kit.

- Do not use alcohol to clean skin surface.
- Apply electrodes to dry, clean, unbroken skin surfaces only, checking all edges for proper adhesion.
- Do not apply skin lotion prior to application of electrodes.
- Electrodes should always be placed at least one inch apart to reduce risk of skin irritation.
- Remove electrodes before showering.
- Once electrodes are removed from skin, adhesive remover towelettes can be used to clean off any adhesive residue.
- If you experience symptoms of skin irritation while using the Avid IF2, discontinue use immediately and contact VQ OrthoCare’s Patient Care Department at **800.452.7993**, or consult your physician.
- Do not remove electrode by pulling wire; peel away from skin by grasping the edge of the electrode

**CAUTION:** Electrodes should be discarded and replaced if damaged, when proper adhesive tack can no longer be obtained, the gel has separated, or you sense a change in stimulation intensity. If in doubt about the integrity or proper function of any electrode, replace it before proceeding with treatment. **IF USING PAIN MEDICATION AND/OR OTHER PAIN MANAGEMENT DEVICES, EXERCISE CAUTION WHEN INCREASING AMPLITUDE ON THE Avid IF2.**

**WARNING:** This device is designed for use only with electrodes approved for medical electrical stimulation.

## BATTERIES AND AC ADAPTER

- The Avid IF2 can be operated with two types of power sources: an internal, rechargeable, 7.4V Lithium Ion battery or an AC adapter.
- A low-battery indicator is shown in upper left corner of screen when battery is low. Plug in the AC adapter.
- A completely drained battery will take approximately 6 hours to charge.
- The Avid IF2 can be used while the battery is charging but will take longer to charge.

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## ENGLISH APP INSTRUCTIONS

- Download the Avid IF2 app from the App Store or Google Play.
- Sign up and complete the information requested. Confirm your registration via e-mail.
- Keep the app open and it will connect to the device overnight.
- **The app will not disrupt your treatment and will not shut the device off to upload information.**
- You may receive reminders to upload your data.

**TO MANUALLY UPLOAD YOUR INFORMATION:**

- Turn on your Avid IF2 device and go to settings **⚙**. Choose Bluetooth. Make sure Bluetooth is enabled.

- Open the app on your phone. It should connect automatically. If not, press “connect” next to your device. After upload, press the two plugs at the top to disconnect **🔌**.
- For complete instructions on the app, see the patient manual. If you have any problems with the app, please contact customer service at **(800) 452-7993**.

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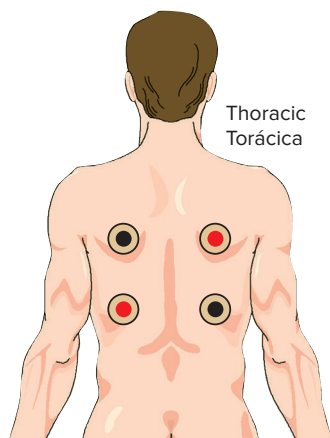
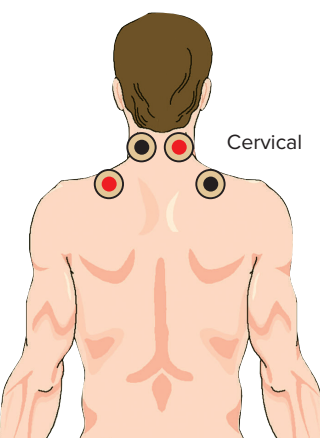
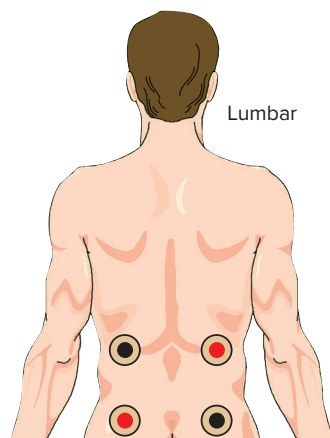
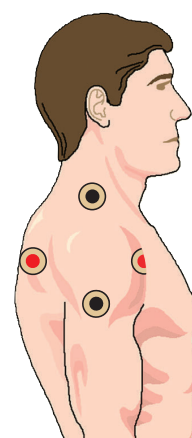
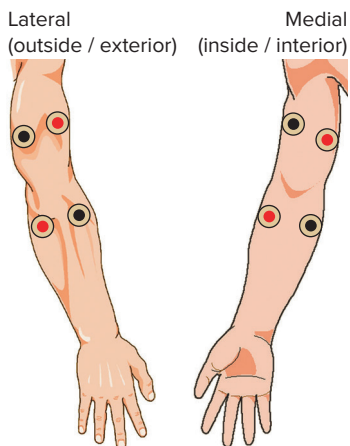
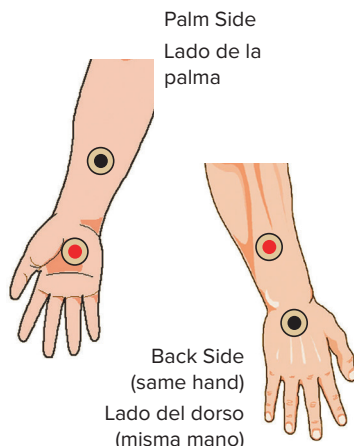
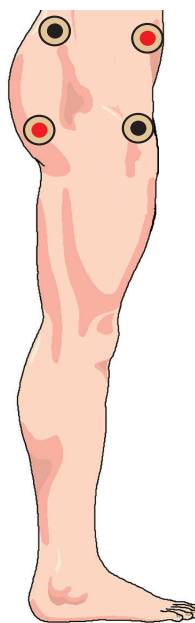
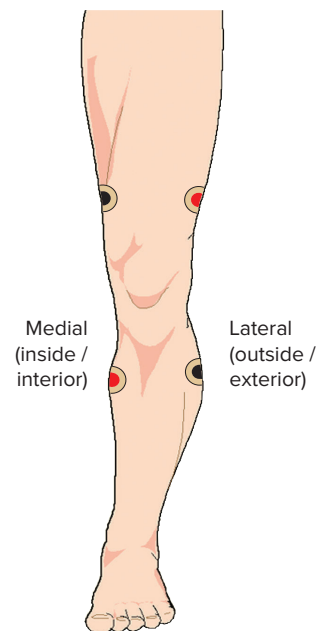
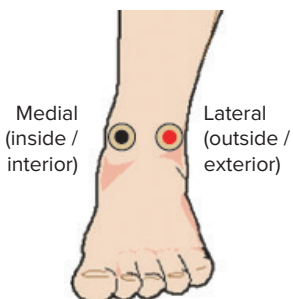
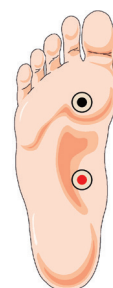
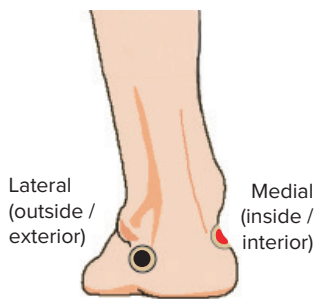
## ENGLISH

### Indications

Interferential stimulation can be used in the following applications:

- Symptomatic relief of post surgical and/or post traumatic acute pain
- Symptomatic relief of chronic intractable pain
- Relaxation of muscle spasms
- Maintaining and increasing range of motion
- Increase local blood circulation

### Contraindication

- Stimulators should not be used on patients with cardiac demand pacemakers

### Warnings

- The long-term effects of chronic electrical stimulation are unknown.
- Stimulation should not be applied over the carotid sinus nerves, particularly in patients with a known sensitivity to the carotid sinus reflex.
- Stimulation should not be applied over the eyes, neck or mouth; severe spasm of the laryngeal and pharyngeal muscles may occur and the contractions may be strong enough to close the airway or cause difficulty in breathing.
- Electrodes should not be placed on the chest in such a configuration that would allow electrical current delivered by the device to flow through the chest because it may cause a cardiac arrhythmia.
- Stimulation should not be applied transcerebrally.
- Stimulation should not be applied over swollen, infected, or inflamed areas or skin eruptions, e.g. phlebitis, thrombophlebitis, varicose veins, etc.
- Stimulation should not be applied over, or in proximity to, cancerous lesions.
- Due to the potential for causing skin burns, this device should only be used with electrodes which limit the current density to  $0.25\text{W}/\text{cm}^2$
- Stimulation should not be applied transthoracically in that the introduction of electrical current may cause cardiac arrhythmias.
- Simultaneous connection of a patient to a high frequency surgical medical electronic equipment may result in burns at the site of the stimularelectrodes and possibly damage to the Avid IF2.
- Application of electrodes near the thorax may increase the risk of cardiac fibrillation.
- Any electrodes that have current densities exceeding  $2\text{mA}/\text{cm}^2$  may require the special attention of the operator.

### Precautions

- Safety of stimulators for use during pregnancy has not been established
- Caution should be used for patients with suspected or diagnosed heart problems
- Caution should be used for patients with suspected or diagnosed epilepsy
- Caution should be used in the presence of the following:
  - When there is a tendency to hemorrhage following acute trauma or fracture;
  - Following recent surgical procedures when muscle contraction may disrupt the healing process;

- Over the menstruating or pregnant uterus;
- Over areas of the skin which lack normal sensation
- Some patients may experience skin irritation or hypersensitivity due to the electrical stimulation or electrical conductive medium; the irritation can usually be reduced by using an alternative conductive medium or alternate electrode placement
- Electrode placement and stimulation setting should be based on the guidance of the prescribing practitioner
- Keep out of the reach of children: the length of the lead wire could cause entanglement of a baby or small child, potentially causing strangulation.
- Use only with the leads and electrodes recommended for use by the manufacturer
- Portable stimulators should not be used while driving, operating machinery, or during any activity in which involuntary muscle contractions may put the user at undue risk of injury
- Do not immerse electrodes in water or trim electrodes as cut edges may affect even distribution of stimulation.
- Electronic monitoring equipment (such as ECG monitors and ECG alarms) may not operate properly when the Avid IF2 is in use.
- Wireless communications equipment such as wireless home network devices, mobile phones, cordless telephones and their base stations, walkie-talkies can affect the Avid IF2 and should be kept at least 1 (one) meter away from the equipment.

### Adverse Reactions

Skin irritation and burns beneath the electrodes have occasionally been reported with the use of stimulators.

### Caution

Federal law (USA) restricts this device to sale by, or on the order of, a practitioner licensed by the law of the State in which he/she practices to use or order the use of this device.

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