



RECOVERY — READY —

EQUIPMENT. RECOVERY. PERFORMANCE.



TRUSTED BY
HEALTHCARE
PROFESSIONALS.
PROVEN BY
RESULTS.

RECOVER WITH CONFIDENCE

Recovery Ready provides patients a rental service of the Game Ready system, for use at home after surgery, injury or intense training.



Advanced cold and compression therapy designed to support your recovery.



IMPORTANT INFORMATION



HOW TO ORDER:

It is the patients responsibility to pre book your Game Ready.
At least 1 weeks notice required prior to surgery to secure your booking.



SINGLE PATIENT USE:

Devices are provided with each booking for infection control purposes.



NOT PRE BOOKED?

If you change your mind after surgery and want a Game Ready device, please contact us and we will endeavour to allocate you a machine as soon as possible.



HIRE COSTS (INCLUDING GST)

2 WEEK HIRE incl. 1 wrap selection **\$300**

4 WEEK HIRE incl. 1 wrap selection **\$550**

DAILY EXTENSION **\$22** per day



FREE DELIVERY BRISBANE,
GOLD COAST & SURROUNDING AREAS




BRISBANE COLLECTION
AVAILABLE UPON REQUEST



SMARTER THERAPY. STRONGER RECOVERY.

- ✓ Helps manage post-operative pain and discomfort
- ✓ Assists in reducing swelling and inflammation
- ✓ Supports rehabilitation and physiotherapy programs
- ✓ Dynamic compression works together with cold therapy for effective treatment
- ✓ Anatomically designed wraps comfortably contour to the body
- ✓ Delivers consistent cooling to the treatment area

TEMPERATURE CHANGE DURING THERAPY*

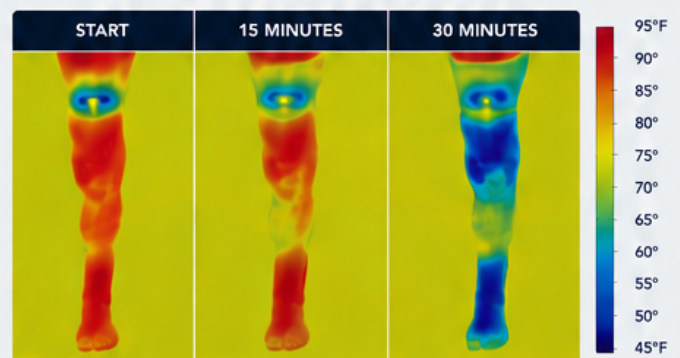

*Individual results may vary.

QUESTIONS OR READY TO BOOK?



CALL US
0428 592 572



EMAIL US
info@recoveryready.com.au



VISIT US ONLINE
www.recoveryready.com.au

EQUIPMENT. RECOVERY. PERFORMANCE.



RECOVERY
— READY —
EQUIPMENT. RECOVERY. PERFORMANCE.



TRUSTED BY
HEALTHCARE
PROFESSIONALS.
PROVEN BY
RESULTS.

HOW THE GAME READY® WORKS

IMPORTANT: The sleeve is only the outer cover of the selected wrap. The internal ATX bladder is property of Game Ready® and remains with that device after use.



1



Fill the unit with water and ice as indicated inside reservoir.

2



Connect the unit to the AC adapter, then to a grounded power outlet.

3



Insert the bladder into the patient sleeve, wrap around patient and connect hose.

CONTROL PANEL OVERVIEW



1. PRESS ON/OFF
Turn unit on or off.



**2. PRESS + BUTTON
SET TIME TO 20 MINS**
Use + or - to set time up to 20 minutes.



**3. PRESS PLAY
(ALSO PAUSE)**
Press play to start.
Press again to pause.



COMPRESSION CONTROL
Select compression level: Low, Medium or High.



TEMP CONTROL
Set target temperature.



**20 MINS ON
2 HOURS OFF
MIN 3X PER DAY**
Recommended use for best results.



ICE TIP

Game Ready® requires the unit to be filled with water and ice (see sticker levels inside reservoir). Use large blocks of ice in a Snap Lock Bags or small Tupperware Containers. Larger blocks melt slower, reducing the times needed to replace ice throughout the day.



EQUIPMENT HIRE

High quality devices professionally serviced and safety checked.



CLINICALLY TRUSTED

Used and recommended by healthcare professionals across Australia.



RECOVER AT HOME

Advanced therapy to support better recovery in the comfort of your own home.



SERVICE AREAS
Brisbane, Gold Coast
& surrounding areas



0428 592 572
Call or SMS



info@recoveryready.com.au
Email us



www.recoveryready.com.au
Visit our website

EQUIPMENT. RECOVERY. PERFORMANCE.