
OO<sup>™</sup>pik

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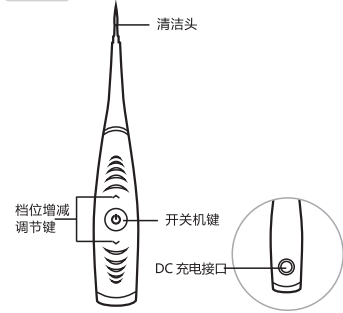
OO<sup>™</sup>pik

OOQik Calculus Remover  
User Manual

Thank you for purchasing OOQik!  
Read this manual, follow all instructions and store it with care.

Product structure

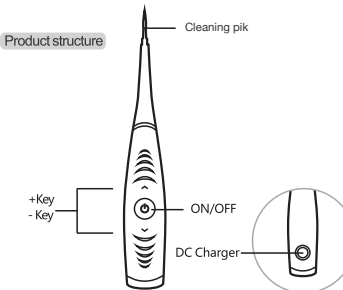OOQik is small in size and convenient to carry, allowing you to clean your teeth at home or on the go.

OOQik uses a cleaning pik just like a dental professional. It effectively removes calculus, tartar, plaque and tooth stains from your teeth.

Product features

- Medical-grade stainless steel pik for excellent cleaning
- Five adjustable modes suitable for different dental needs
- Intelligent high-frequency vibration effectively removes dental calculus, tartar, plaque and tooth stains.
- Light weight, USB chargeable and travel friendly.
- Battery is durable and long lasting.

Attention

- The cleaning pik at the top of OOQik is not replaceable.
- Do not place OOQik in direct sunlight.
- Keep OOQik indoors and between -10° c and 40° c.


- Do not wash or soak OOQik in water or any other liquids.
- OOQik is for adults only, be sure to store it properly.
- Keep out of reach of children.
- This product contains batteries. Do not place OOQik in a high-temperature environment.
- If your product has any issue, please directly contact the manufacturer or place of purchase.
- Do not use OOQik more than twice per week.
- OOQik is not a replacement for regular dental cleanings or visits to your dentist or oral healthcare professional.

Product Specifications

Rated voltage : 3.7V	Rated current : 200mA
Battery capacity : 300mAh	Standby current : ≤20uA
Product size : 183.6x30x22.3mm	Product weight : 53.6g
Charging voltage : 5V	Motor noise : ≤60dB


Usage method

Before first use, make sure to fully charge your OOQik for approximately two hours.

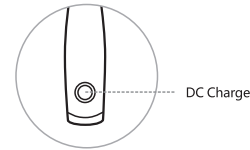Pic 1

- Press the power button to turn on OOQik.
- Adjust vibrational frequency with the upper and lower adjustment buttons. OOQik is now ready.
- The proper cleaning angle is 45° against the tooth surface (pic 2). *Tooth surface should be moist.*
- After using OOQik, wash the steel pik with hot water and wipe down the body of OOQik.

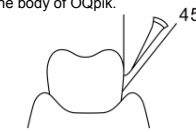Pic 2

Charger

- OOQik uses a USB charger which plugs into the charging socket on the back of your OOQik.
- OOQik should not be used while charging.
- While charging, OOQik will display a red light.
- When fully charged, OOQik will display a green light.

Pic 3

Qualification Certificate

Model: \_\_\_\_\_ P/N: \_\_\_\_\_  
LOT: \_\_\_\_\_ Ship Date: \_\_\_\_\_  
Check Date: \_\_\_\_\_ Checker: \_\_\_\_\_

Warranty Card

Name: \_\_\_\_\_  
Buy Date: \_\_\_\_\_  
Phone number: \_\_\_\_\_  
Address: \_\_\_\_\_  
Issues: \_\_\_\_\_

If you experience an issue with your OOQik, contact [hello@theOO.shop](mailto:hello@theOO.shop) and include a copy of your receipt.