

STEP 1 Wrap Cuff Firmly Around Each Calf


- Start with the middle strap ensuring that the strap is secured into place as tightly as possible.
- Do the same thing for the top strap, followed by the bottom strap.
- If cuff is not applied snugly enough, unit may give an E1 error code and will not function properly.
- Make sure the top of the unit is placed close to the knee and the unit is offset 1 inch to the right or left of the shin bone to ensure tightest fit.

STEP 2 Turn on the Device

- Press and hold the power button on each cuff for at least 3 seconds
- The device pressure level defaults to 55 mmHg
- There is a 5 second delay before inflation begins

STEP 3 Charge the Battery

- Press and hold power button for 3 seconds
- 1 red light and 3 green lights indicate a fully charged battery; will run for up to 24 hours once fully charged.
- Use the included charging cord to charge each cuff.
- Please note, the unit will not operate while charging.
- Must be charged when not in use.
- Unit takes approximately 6 hours to fully charge.

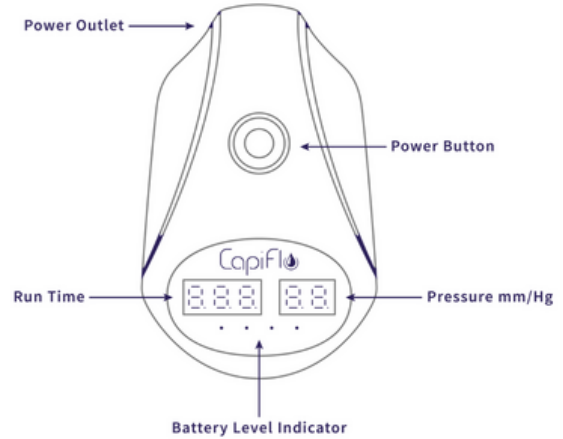

FOR MORE INFORMATION, PLEASE REFER TO THE INSTRUCTION MANUAL.

www.capiflo.com If you need assistance please contact: ☎ 813-243-2222 ✉ sales@capiflo.com

Why is an alarm sounding from the device?

An alarm is sounding because you are experiencing an E1 or E2 error code.

What does an E1 error code indicate?

An E1 error code indicates the pressure meter is not functioning properly. This can often be fixed by making sure the device is tightly secured by following detailed instructions indicated on the quick start guide.

What does an E2 error code indicate?

An E2 error code indicates there is an issue with the battery. Contact customer service to get this resolved.

How long does the battery last when fully charged?

The battery lasts up to 24 hours on a full charge.

Can the device be used while charging?

The device cannot be used while charging.

Can I use the device on one leg at a time?

Cuffs should be used as instructed by a medical professional.