


Disposable Surgical Mask

LEVEL 3

DCS004 ASTM F2100-2019




Material Type

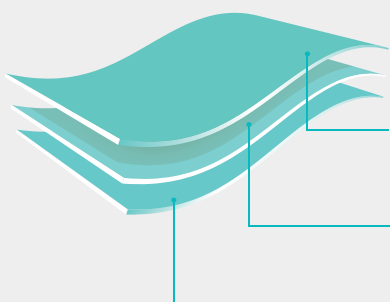

- Outer layer**
Leak proof non-woven fabric, effectively inhibit and isolate harmful particles and gases
- Middle Layer**
Fine melt-blown non-woven fabric, high density filter layer, great efficient at filtering out fine particles
- Inner Layer**
Non-woven fabric, moisture proof, breathable, non toxic, non stimulating, soft and skin friendly


Wearing directions

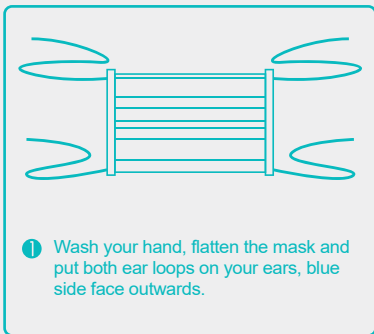

- ❶ Wash your hand, flatten the mask and put both ear loops on your ears, blue side face outwards.



- ❷ Pick up the mask by the upper ties, securely fasten the mask at the back of your head by forming a knot, press the nose clip to match your nose shape.

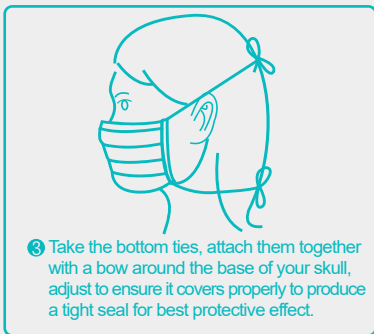

- ❸ Take the bottom ties, attach them together with a bow around the base of your skull, adjust to ensure it covers properly to produce a tight seal for best protective effect.

Disposable Surgical Mask

Indication for Use

The single use mask is intended to be worn to protect both the patient and healthcare personnel from transfer of microorganisms, body fluids and particulate material. These face masks are intended for use in infection control practices to reduce the potential exposure to blood and body fluids.
This is a single use, disposable device, provided non-sterile.

Product Information

【Product name】	Disposable Surgical Mask
【Brand】	Purism
【Size】	17.5cm*9.5cm
【Color】	Blue
【Fluid Resistance】	160mmHG
【BFE】	98%（minimum）
【PFE@0.1 micron】	98%（minimum）
【Materials】	Spun-bond non-woven fabric, melt-blown non-woven fabric, ear loops and nose clip
【Expiration date】	3 years after production
【Storage】	Store in a dry, clean environment with room temperature, avoid high temperature and exposure to flame.

Warnings

- 1.DO NOT use if the package is damaged.
- 2.Please refer to the instructions before use.
- 3.For adult use only.
- 4.The mask does not eliminate the risk of contracting any infection.
- 5.Failure to use properly and maintain the product could result in illness or even death.
- 6.Only for one-time use, please destroy and dispose according to regulations after use.
- 7.Please replace a new mask when mask is damaged, soiled, or contaminated.
- 8.Please use within the validity period.

【Important Notice】

To the extent permitted by law, DC not be liable for any loss or damage including any loss of business, loss of profits, or for any indirect, incidental or consequential loss or damage arising from reliance upon any information herein provided. Nothing in this statement will be deemed to exclude or restrict DC's liability for death or personal injury arising from its negligence.

All information presented in this packing is based on the latest data available at the time of printing. DC reserves the right to make changes at any time without prior notice.

Distributor: Purism LLC
Address: 5023 Addicks Satsuma Rd, Suite W, Houston, TX 77084
Phone: +86-400-996-5015
Website: <http://purismworld.com/>
Email: sales@purismworld.com

Chuzhou Daddy's Choice Science and Technology Co.,Ltd.
Middle of Nanjing Road, Langya District, Chuzhou City, Anhui, China 239000

MADE IN CHINA

