

RADIESSE®



Boost your
collagen
production to
work against
skin ageing
from within

RADIESSE is an injectable dermal filler used to enhance the facial area, back of the hands and diluted for smoothing lines and wrinkles of the décolleté.

All models and patients in this brochure are actual RADIESSE® patients.

MERZ AESTHETICS®



Visible results with RADIESSE



Before



After



Before



After

Patient was injected with 2.25ml of RADIESSE+ in each side of the jawline. Results immediately after treatment.



Before



After

Individual results may vary, please ask your aesthetic practitioner for a consultation.

Personalised treatment options

RADIESSE has a formulation and technology that can offer you an individualised treatment approach – tailored to your unique needs.

- Lifts for a more youthful appearance
- Provides defined facial contours and structure
- Smooths lines and wrinkles in the décolleté
- Boosts collagen and elastin in the tissues

How RADIESSE works


RADIESSE is composed of a defined range of tiny microspheres of calcium hydroxylapatite (CaHA) suspended in a gel. The CaHA microspheres activate the body's natural processes to produce new collagen and elastin. Over time, your body breaks down these tiny microspheres and the gel, leaving behind your own newly formed natural collagen.

Leaving your skin with an improvement in structure and quality depending on the treatment. For immediate and long-lasting, natural-looking results.

Immediate effect


Long-term collagen and elastin production

Immediate enhancement

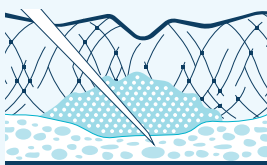


Immediate lifting and contouring effect on signs of ageing - smoothing of wrinkles and folds plus contouring.

Stimulating collagen and elastin production

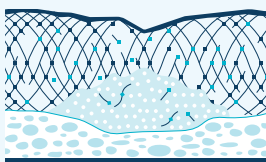

The CaHA microspheres create scaffolds that stimulate key skin cells - fibroblasts - to produce collagen and elastin.

Improving skin quality

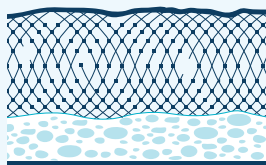

Newly formed networks of collagen and elastin provide structural support and improvement of skin quality.



Gel matrix



CaHA microspheres



Collagen fibers



Fibroblasts



Elastin



Epidermis



Dermis



Subcutis

Your questions answered

What does the treatment involve?

Following your consultation, the treatment session will typically take about 30 minutes. During the session, RADIESSE is injected into the specific area with a needle or cannula.

How soon will I see effects from my RADIESSE treatment?

Depending on the treatment performed, the effect can be visible immediately after injection due to the lifting and contouring effects. Additionally, by stimulating your body's own collagen and elastin production, the benefits of skin quality enhancement with RADIESSE will continue to develop over time.

How long do the effects of RADIESSE last?

The duration of the results depends on age, skin type, lifestyle, metabolism, and, of course, the area treated. RADIESSE results are normally particularly durable, lasting up to 12 to 18 months.

How much downtime is needed after treatment?

There is normally minimal downtime required after being treated with RADIESSE. Swelling and numbness are common following treatment and usually resolve within 7 to 10 days. You should avoid sun exposure, tanning lights, saunas, and intense facial treatments for one week after treatment.

What do I need to consider when choosing a RADIESSE practitioner?

It is important to look for a clinic or healthcare professional who is well-trained in the handling of injectables with a medical background, such as a doctor, nurse, or dentist.

Your treatment benefits

Boosts
collagen
and elastin



Shapes and
redefines facial
contours



Supports against
sagging for a lifted
appearance

Smooths lines
and wrinkles



Collagen is key in skin ageing

When we speak about skin ageing, one of the processes we are talking about is collagen loss. When we are young, collagen is organised in a dense network that keeps our skin firm, smooth and radiant. Over time, the network loosens and less collagen is produced. As a result, you begin to see the common signs of ageing – such as skin sagging and wrinkles.



Young skin



Aged appearance

Key facts about RADIESSE

Collagen boosting

Invest in your skin with this foundational treatment which acts at the source to boost key components of the skin

Well-established and favourable safety profile based on scientific and clinical experience

98% patient satisfaction rate
with results lasting up to 12-18 months

2026 New Beauty award
for Best Filler for Hand Rejuvenation



Individual results may vary, always consult a registered practitioner.

©2026 Merz Aesthetics UK Ltd – All rights reserved. MERZ AESTHETICS and RADIESSE logos are trademarks and/or registered trademarks of Merz Aesthetics UK Ltd in the United Kingdom. Registered in England No. 14506945.

Merz Aesthetics UK Ltd, Ground Floor Suite B, Breakspear Park, Breakspear Way, Hemel Hempstead, Hertfordshire, HP2 4TZ.

www.merz-aesthetics.co.uk.