

# INSTRUCTION FOR USE

## SENSIStrip Mustard 20/5 Tests

(Cat. nr. HU0030125/HU0030165)

### Lateral-flow Device for the Determination of Mustard in Food and as Cleaning Control Monitoring

Sensitivity for food matrix	2.5 ppm
Sensitivity for swabbing	0.007 µg/cm <sup>2</sup>
Sensitivity for rinse water	0.33 mg/L

#### 1. GENERAL INFORMATION

Mustard belongs to the Brassica plants. With about 30-35% the fraction of proteins in mustard seed is very high. Some of these proteins are known for being allergenic, such as Sin a 1 and Bra j 1. These proteins are predominantly heat resistant making them stable to different production processes. In addition to brown mustard (*Brassica juncea*) and black mustard (*Brassica nigra*) primarily yellow mustard (*Sinapis alba*) is used as an ingredient in many foods and food preparations. For mustard allergic persons, hidden mustard allergens in food are a critical problem. Already very low amounts of mustard can cause allergic reactions, which may lead to anaphylactic shock in severe cases. Because of this, mustard allergic persons must strictly avoid the consumption of mustard or mustard containing food. Cross-contamination, mostly in consequence of the production process, is often noticed. The sausage production process is a representative example. This explains why in many cases the existence of mustard residues in food cannot be excluded. For this reason, sensitive detection systems for mustard residues in foodstuffs are required.

The **SENSIStrip Mustard Lateral Flow Device** represents a sensitive detection system and is particularly capable of detecting mustard residues in food matrices, rinse water and swabs.

#### 2. PRINCIPLE OF THE TEST

The **SENSIStrip Mustard** test is based on the principle of immunoassay. Mustard containing sample is given into a reaction vial containing biotinylated antibody directed against mustard proteins. After 3 minutes incubation at room temperature a test strip is placed into the reaction vial. The sample migrates along the nitrocellulose membrane by

capillary forces. Along its way it releases gold nanoparticles conjugated to streptavidin. An antibody-gold complex is formed. For positive samples a red line is formed when the liquid reaches the test line area. In case of negative samples, no line is formed. In any case, above the test line area a red control line appears, indicating the validity of the test. The test is evaluated after another 5 minutes.

#### 3. PRECAUTIONS

Full compliance of the following good laboratory practices (GLP) will determine the reliability of the results:

- 1) Store the kit at 2-8°C.
- 2) Do not use the kit after its expiry date.
- 3) Prior to beginning the assay procedure, bring all samples and reagents to room temperature (20-25°C).
- 4) Extraction buffer should be mixed by gentle inversion or swirling prior to use. Do not induce foaming.
- 5) Once the assay has been started, all subsequent steps should be completed without interruption and within the recommended time limits.
- 6) Replace caps in all the reagents and samples immediately after use.
- 7) Use separate disposable consumables for each transfer of sample to the reaction vial in order to prevent cross-contamination.
- 8) Do not mix components from different batches.
- 9) Do not use reagents after the expiration date.

**NOTE:** The swab sampling device included in this kit may be supplied as sterile with a sterility expiration date printed on the device. However, this kit does not require a sterile sampling device, therefore the swab sterility expiration date does not affect the kit expiration date and can be disregarded.

#### 4. KIT CONTENTS

The kit contains components and reagents for 20 tests or 5 tests. They have to be stored at 2-8°C. Expiry data are printed on the labels of the reagent containers and the outer package.

Content	20-strip	5-strip
Test Strips, in tube with desiccant stopper	20 pcs	5 pcs

Reaction vials	20 pcs	5 pcs
Extraction tubes with caps	20 pcs	5 pcs
Extraction Buffer, 60 mL, ready-to-use.	1 pcs	1 pcs
Disposable Pipettes, 0.3 mL	21 pcs	6 pcs
Disposable Pipette, 3 mL	1 pcs	1 pcs
Disposable Spatulas	20 pcs	5 pcs
Swab Sticks	20 pcs	5 pcs
Evaluation Card	1 pcs	1 pcs
Tubes and vials racks	by kit box	by kit box
QR-Code for evaluation with RapidScan ST5 lateral flow strip reader	1 pcs	1 pcs

## 5. EQUIPMENT AND MATERIALS (NOT PROVIDED)

- 1) *RapidScan ST5* lateral flow reader for quantitative evaluation (optional)

## 6. SAMPLE PREPARATION

Due to high risk of cross-contamination all applied instruments like applicator, mortar, vials etc. have to be **cleaned thoroughly** before and after each sample. Allergen proteins adhere very strongly to different surfaces. In certain cases, they can resist a common dishwasher cleaning. To identify possible cross-contamination caused by previous extractions it is strongly recommended to note the sequence of the extractions for pattern recognition.

Chocolate and other products with high polyphenol content tend to show reduced results. To overcome this effect a special extraction additive can be ordered separately (HU0030100).

### 6.1 Solid samples / Liquid samples

- 1) Homogenize sample using appropriate methods depending on its specific nature (e.g. grind, crush, mix).
- 2) *Solid samples:* Transfer one spatula of sample to an extraction tube. Alternatively, in order to increase precision, weigh out 0.2 g of sample into an extraction tube.  
*Liquid samples:* Transfer a half spatula of sample liquid to the extraction tube. Alternatively, in order to increase precision, pipette 0.2 mL of sample into an extraction.
- 3) Add 3 mL of ready-to-use extraction buffer to the sample by using the disposable 3 mL pipette.
- 4) Close the extraction tube with a cap and shake for 1 minute.
- 5) Let the solid remain sediment. Depending on the nature of the samples this might take 1-2 minutes. Alternatively centrifuge at 2000 g or higher.

- 6) Remove cap and transfer 0.3 mL of sample supernatant into a reaction vial using a disposable 0.3 mL pipette.


### 6.2 rinse water

- 1) In case of strong acidic or basic rinse solution adjust the pH of the sample to 7 (+/- 0.5).
- 2) Transfer 0.3 mL of extraction buffer into an extraction tube using one of the disposable 0.3 mL pipettes.
- 3) Transfer 0.3 mL of rinse sample into the extraction tube using a second disposable 0.3 mL pipette
- 4) Mix the two liquids by applying the same pipet as in step 3.
- 5) Transfer 0.3 mL of mixture to a reaction vial applying the same pipet as in step 4.

### 6.3 Swabbing samples

#### DRY SURFACES

- 1) Mark out 5x5 cm area or use swab directly on (e.g. uneven) area.
- 2) Transfer 1 mL of ready-to-use extraction solution into an extraction tube by using the disposable 3 mL pipette.
- 3) Moisten a swab by dipping into the tube.
- 4) Swab marked area by using crosshatch (1. horizontally, 2. vertically, 3. diagonally) technique while rotating the tip.

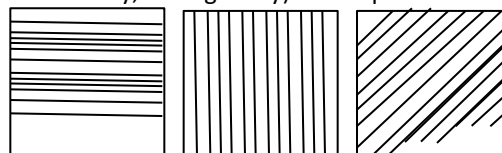

- 5) Place swab into the tube and break off the tip.
- 6) Close extraction tube with cap and shake for 1 minute to release the sample from the swab.
- 7) Remove cap and transfer 0.3 mL of sample supernatant into a reaction vial using a disposable 0.3 mL pipette.

#### WET SURFACES


Apply the same method as described for dry surfaces without prior need to moisten the swab.

## 7. ASSAY PROCEDURE

- 1) Prepare samples as described above.
- 2) After transfer of the sample to the incubation vial add cap and shake for 15 seconds. Make sure that the biotinylated antibody is completely dissolved.
- 3) Incubate for 3 minutes.
- 4) Remove cap and place one strip into the vial. For proper strip orientation make sure that the arrows on the cover foil point downwards.
- 5) Incubate for 5 minutes.
- 6) Remove the strip from the vial and evaluate immediately.

## 8. EVALUATION

SENSIS*Strip* lateral-flow devices are evaluated according to the following scheme:

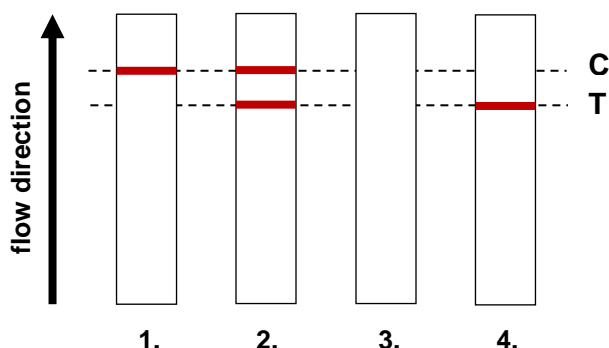

1. **Negative:** visible control (C) line, no test (T) line
2. **Positive:** visible control (C) and test (T) lines
3. **Invalid:** neither control (C) and test (T) lines visible
4. **Invalid:** no control (C) line and visible test (T) line

For a better distinguishing between negative, borderline and positive samples a colour card for evaluation is provided with the kit. The intensity of the test line has to be compared with the different increments of the colour card. Results lower than increment 3 should be treated as negative. Results according increment 3 or higher should be treated as positive. Since the increments of the colour card are ranging up to 10 a semi-quantitative evaluation is also possible. This can be improved by taking into account the results stated in the validation report of the product.

In addition, a quantitative evaluation (2.5-40 ppm) in combination with the *RapidScan ST5* lateral flow reader is possible. For further information, please contact Gold Standard Diagnostics Budapest.

## 9. PERFORMANCE

### 9.1 Sensitivity

LOD (total mustard) of the SENSIS*Strip* lateral-flow test is 2.5 ppm for food matrix, 0.33 mg/L for rinse water and 0.007 µg/cm<sup>2</sup> for swab samples applying the procedure above. The corresponding amounts of mustard protein can be calculated by anticipating a protein content of mustard of 26%.

**NOTE:** Sensitivity may vary depending on the matrix and processing of a complex food mixture. For achieving reliable results each matrix should be validated prior to routine testing.

### 9.2 Cross-reactivity

For the following foods not cross-reactivity could be

detected:

Adzuki bean	Curcuma	Peach
Almond	Dill	Peanut
Apricot	Duck	Pecan
Barley	Ewe's milk	Pepper
Bean, white	Fennel	Pine seed
Beef	Fenugreek	Pistachio
Bovine gelatin	Flaxseed	Poppy seed
Brazil nut	Garden cress	Pork
Buckwheat	Garlic	Potato
Caraway	Gliadin	Prawn, cooked
Cardamom	Goat's milk	Pumpkin seed
Carob bean	Guar gum	Radish
Carrot	Hazelnut	Rice
Cashew	Horseradish	Rye
Cayenne	Kidney bean	Sesame
Celery	Kiwi	Shrimp
Cherry	Lamb	Soy flour
Chia	Leek	Soy lecithin
Chicken	Lentil	Soy milk
Chickpea	Lupin	Split peas
Chilli	Macadamia	Sucrose
Cinnamon	Marone	Sunflower seed
Clove	Nutmeg	Thyme
Cocoa	Oats	Tomato
Coconut	Onion	Turkey
Cod	Oyster	Walnut
Corn	Paprika	Wheat
Cow's milk	Pea	White cabbage
Cumin		

The following cross-reactions were determined:

Rape seed	3.1%
-----------	------

### 9.3 High-dose-hook Effect

Reduced or absent signals can occur in case of very high concentrations. The test gives valid results up to a concentration of 5000 ppm for food samples, according 13 µg/cm<sup>2</sup> for swabs and 667 mg/L for rinse water samples.

### 9.4 Additional Performance Data

Additional data can be found in the corresponding validation report of the product, which can be inquired at Gold Standard Diagnostics Budapest.

## 10. LIABILITY

Gold Standard Diagnostics Budapest shall not be liable for any damages to the customer caused by the improper use of the kit and for any action undertaken as a consequence of results. Gold Standard Diagnostics Budapest shall not be liable for the unsafe use of the kit out of the current European safety regulation.