

EQUINE Blood Collection & Shipping Guidelines

START HERE

Equine



Blood collection should be as **STRESS-FREE** as possible

Red Top Tube
(Whole blood or Serum)
NO ADDITIVES
DO NOT USE Purple Top Tubes


Label tubes with
Tube # & Animal ID



IF NOT SENDING
RIGHT AWAY
1-2 Days Max in Fridge


1ml/cc of
Blood or Serum

Let Blood
Samples Sit for at
least 1 HOUR
Before Putting In
Fridge

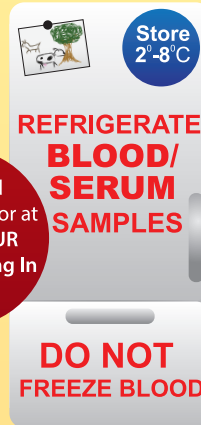

WHEN READY
TO SHIP

Bundle Samples in
Max of 10 per Set

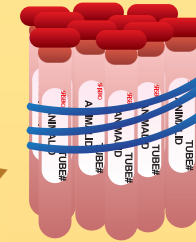

Cover Samples with
absorbent paper



Paper Towels,
Newspaper
or etc...

SHIPPING

- DO NOT send samples in an envelope. Samples must be sent in a box.
- NO Frozen Water or Dry Ice.
- Monday is the ideal time to ship samples.
- Overnight / 2-Day Shipping REQUIRED
- Email lab to provide tracking # & # of tests.
- Contact Lab for Testing Schedule (Max 2 Week Turnaround)

PAYMENT

- We accept Debit/Credit Cards, Checks, Money Orders and Cash. We do not advise sending cash by mail. Checks & money orders should be made to 'UBRL'

Strongly recommend to avoid USPS if possible.
UPS and FedEx is ideal.

REQUIRED
OVERNIGHT / 2-Day Shipping




VERY IMPORTANT:
Confirm shipping carrier,
does not label package
as 'Hazard'. It is incorrect
and it may cause shipping
delays.

Frozen Ice Pack REQUIRED For
EQUINE PREGNANCY




NO Frozen Water
or Dry Ice

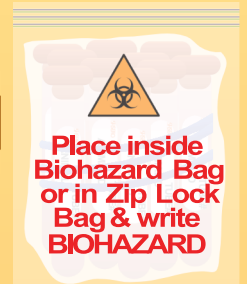

Invoice link will be emailed
once samples arrive in lab

U **UNIVERSAL**
B **BIOMEDICAL**
R **RESEARCH**
L **LABORATORY**
Livestock Diagnostics

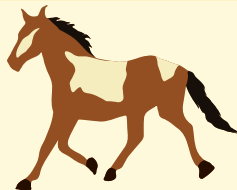

UBRL
1300 N. Fresno St.
Suite#101
Fresno, CA 93703

(559)498-0820

UBRL.ORG

Lab@UBRL.ORG

DO NOT SHIP IN
AN ENVELOPE!
(We receive broken tubes
often from samples
being sent in envelopes)

IMPORTANT:
'EXEMPT ANIMAL
SPECIMEN'
Must be written on box

Wrap in
Bubble Wrap
to protect
tubes.

Place All Documents, Payment &
ice pack inside its own zip lock bag
to avoid being damaged