

Quick Start Guide



Content

- (1) Bench Oven
- (2) Nickel-Plated Wire Shelves
- (4) Shelf Supports

Cautions

- ⚠ **Heavy Object** - Use two (2) persons when lifting or moving the oven.
- ⚠ Only use a grounded outlet rated for your unit's electrical requirement.
- ⚠ Do NOT place volatile or combustible materials inside or in close proximity to the Bench Oven.
- ⚠ Exterior surfaces may be hot to the touch during operation at high temperatures.

Before you power the unit ON

- Ensure you have removed all the packaging items from the oven before use.

Setting up the unit

- Place oven on a flat, solid surface.
- Maintain a minimum of 3 inches of separation around the oven's exterior walls and a minimum of 6 inches between the rear electrical cover and any building or vertical surface.
- Plug the oven into a grounded outlet for your unit's rated voltage. Isolate your oven to a separate, appropriately rated circuit or breaker.
- ⚠ On models 350 and 550, the 115-volt plug should not be confused with the mirror image of the 230-volt plug.

The electrical specifications are found on the label at the back of the oven.

For full manuals and guides, please visit:
<https://quincylab.com/bench-ovens-1>

For general certifications, please visit:
<https://quincylab.com/certifications>

-OR-

Scan the QR code to access the full operating manual and guides page




Contact

Quincy Lab Support:
800-482-4328 | information@quincylab.com

Analog Bench Oven


Models: 21, 31 and 51 – 250/350/550 (115V & 230V)

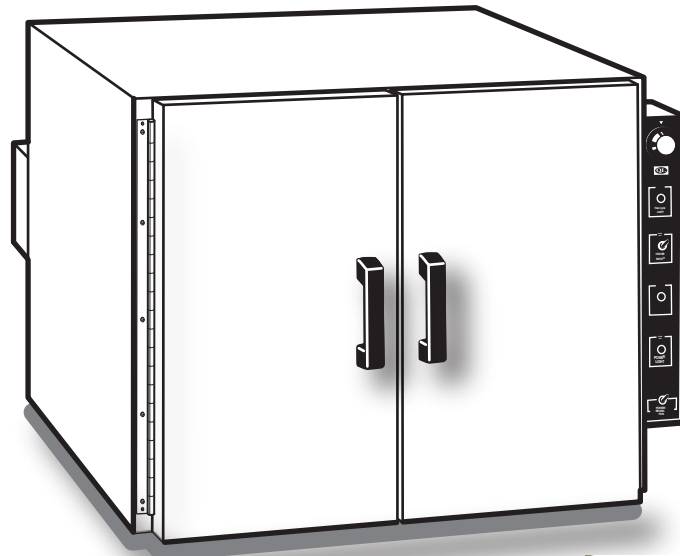

Quick start-up

- 1 Install the shelves in the oven.
See the shelf installation instructions below.
- 2 Turn the power/recirc. fan switch to the 'On' (up) position.
- 3 Turn the heater switch to the 'Heat' (up) position.
- 4 Rotate the thermostat dial to the desired temperature.

- ⚠ Once the temperature is set on the dial, the heat cycle light will illuminate until the set temperature is reached. Once reached, the heat cycle light will turn on and off with the heaters, maintaining the set temperature. Typically, the oven needs to cycle at the set temperature for at least 20 minutes to achieve equilibrium and become fully stable.

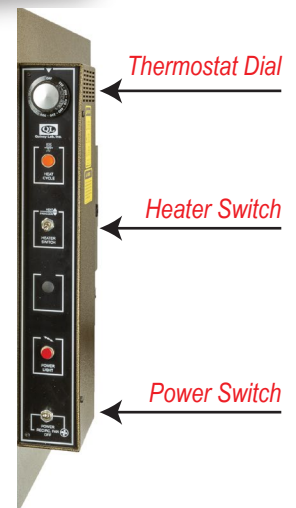

Shelf Installation

Orient the shelf support bracket in the 'down' or 'L' position, and insert the shelf bracket rivets into the corresponding keyhole supports on each inner side wall of the oven.

Place the shelf on top of the brackets.

