

# Nautilus & Medallion Stairlifts

*Stepping up your Confidence and Independence*



JunEAU Stair Lifts

(907) 275-8326

# Be Independent. Be Confident. Be Informed.

**You are not alone** - 20 million people in the U.S. have disabilities that create problems while walking, climbing, or descending stairs; making it more difficult to stay in the home they love. Homes that are so often our foundation; where good times are had and memories are created.

Compounding this challenge - according to the U.S. Department of Housing and Urban Development, 65% of single family/attached homes in the U.S. have 2 levels or more. When you couple those statistics with those from the U.S. Centers for Disease Control and Prevention - that every 11 seconds, an older adult is treated in the emergency room for a fall, modifying your home becomes a real consideration.

Don't let your staircase become an obstacle to your quality of life.



---

20 million people in the U.S. have disabilities that create problems walking, climbing, or descending stairs.

---



---

65% of single family and attached homes in the U.S. have 2 or more levels.

---



---

Every 11 seconds, an older adult is treated in the emergency room for a fall.

---

# What to Look for When Purchasing



Many stairlifts on the market have a 265 lbs - 300 lbs capacity. While that may be enough for some, there may be a need for greater range. More than one family member (or visitors) may need to use the lift. No need to strain the motor or keep guests or other family members from using the lift. Plan on what you may need 5 or 10 years from now.

If these concerns are “weighing” on you, rest assured - our 350 lb or 600 lb weight capacity will provide you peace of mind knowing your lift will handle your total current, and future needs.

## Zero Maintenance

For most people that require a stairlift, performing regular maintenance on their lift is a challenge. The Nautilus and Medallion are leading products on the market that do not require lubrication due to their nylon polymer gear rack and worm drive, featuring a full lifetime warranty on all parts. No mess, no fuss, and no expensive service contracts.



## Rail Size & Mounting

The Nautilus and Medallion rail is the best in class and sits just 3.5” - 5.5” high. The most compact lift on the market, the Nautilus folds to just 10.5” (the Medallion folds to 11”), mounting as close as 0.5” from the wall, allowing for others to easily pass by when folded.





# Safe, Quiet, Dependable, and Grease Free!

Nautilus and Medallion are a leading brand of stairlifts and are inspired by those who are mobility challenged. Jameson Medical is driven to offer mobility and accessibility products that help or enable people to go where they want, when they want, how they want while being safe ... while lasting a lifetime.

Founded in 1991, we have a proud history of growth by evolving our product lines to enrich the lives of our customers. We give people the path to independence to remain in their home, enjoy their families, and be engaged in their communities. Through innovation, we are committed to helping people of all ages overcome their limited physical mobility at home and away.

All Nautilus and Medallion stairlifts are assembled in the USA and come with a lifetime warranty on all parts, which means quality you can depend on.




# Peace of Mind Benefits Only with Nautilus and Medallion

There are many stairlifts to choose from in the marketplace. So what makes the Nautilus and Medallion different and contributes to its lower cost of ownership?

- ✔ Zero Maintenance and lower energy bills from a low friction system.
- ✔ Efficient motor that gives extended performance during power loss.
- ✔ Helical Worm Gear and Nylon Polymer Gear Rack for a smooth ride and greater stability.
- ✔ Optional Folding Rail Design.
- ✔ Compact, Grease Free, Lifetime Warranty on **ALL** Parts.

The Nautilus and Medallion provide a weight capacity of either 300, 350, or 600 lbs.

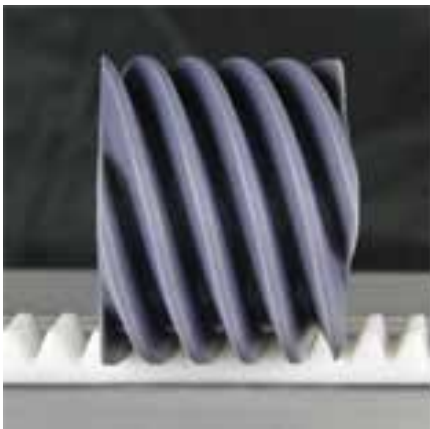

Oversized Helical Worm Gear with six points of contact for an energy efficient low friction design. That means a smooth and stable ride.



Nylon Polymer Gear Rack with steel reinforcement. Lubrication-free. That means durability and ZERO maintenance.



The Nautilus is the world's most compact stair lift, folding closer to the wall. That means others can pass by easily, leaving more room for you and fire/rescue.

# Nautilus Folding Rail Eliminates Trip Hazards



The Nautilus and Medallion Folding Rail deploys automatically to keep walkways clear. This exclusive (patented) design option eliminates the need to stop and push a button or leave the seat. And its simple mechanical design makes it reliable and minimizes system failures. The folding rail option allows the unit to be charged when parked at the top, bottom, or in the raised (hinged) position.



Nautilus and Medallion stairlifts come with two hand-held remotes so you may operate your chair without sitting on it, or a second user can call or send the chair as needed. The chair can be used as a "dumbwaiter" and items can be sent back and forth. No one needs to carry anything up and down the steps.



# The Nautilus and Medallion Advantage



**INDEPENDENT**



**CONFIDENT**



**INFORMED**

The Nautilus and Medallion stairlift gives you the freedom to go up and down the stairs with confidence. It gives Independence. It gives Peace of Mind. The Nautilus and Medallion aren't just stair lifts, they're a quality of life for those who own it and their families. Designed for a Lifetime of Dependable Use!

The Nautilus and Medallion Stairlift, made with you in mind.



*Every year, nearly 100,000 people make the decision to invest in a stair lift. Our family was one of them. We wanted the same easy access for our son that other families have without thinking of it.*

## Standard Features



Safety obstruction sensors.



Lifetime Parts warranty



Swivel seat for a safe exit




Wide footrests for stability and extra comfort.

# Measurement Guide

It really is very easy to take measurements for your new stairlift. All you will need is a good quality metal tape measure. First, take a general overview of your stairway layout and decide which side of the staircase the stairlift will be mounted on. Be aware that the stairlift rail will overhang the top landing by 7" (20cm) and meet the floor at the foot of the stairs 13" (33cm) beyond the 'nose' of the first step. Take note of any fixtures, such as window sills, hand rails, central heating radiators, doorways etc. that intrude into the stairway area - or could be blocked by the presence of the stairlift. Next, simply stand at the top of the staircase and feed the tape measure down the steps until it goes out at an angle and hits the floor. What is the measurement at your toes?

**Step 1** - Staircase Length from floor to floor.

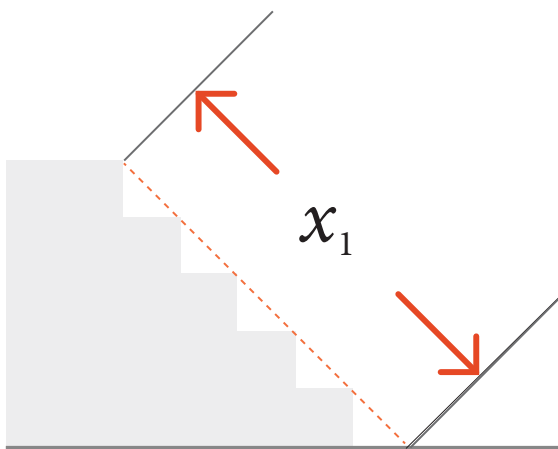

**Step 2** - If there is a door at the top – what is the distance from the face of the door to the edge of the landing when the door is closed?

**Step 3** - Standing at the bottom and looking up, which side of the staircase will your stairlift be mounted on?

**Step 4** - Is there 24 inches of clearance at the bottom of the staircase to allow the chair to fully lower to the floor?



# Nautilus and Medallion Specifications

	NAUTILUS	NAUTILUS HD	MEDALLION
WEIGHT CAPACITY	350 lbs	600 lbs	300 lbs
DRIVE SYSTEM	Nylon Polymer Helical Worm Gear	Nylon Polymer Helical Worm Gear	Nylon Polymer Helical Worm Gear
SPEED	20 FPM	14 FPM	20 FPM
SEAT	19" W x 16" H x 18" D	25" W x 16" H x 18" D	19" W x 16" H x 18" D
FOLDED WIDTH	10.5"	17"	11.0"
FOLDED WIDTH with FOLDING RAIL	13.6"	N/A	13.6"
CONTROL CIRCUIT	24V DC	24V DC	24V DC
INCLINE LIMITS	32° - 52°	32° - 52°	32° - 45°
WARRANTY	• LIFETIME on ALL Parts	• LIFETIME on ALL Parts	• LIFETIME on ALL Parts

## Other Lifts Offered by **Juneau Stair Lifts**



Curved Stair Lifts



SL350 / SL350OD



Incline Lifts

**Please feel free to call us with any questions: (907) 275-8326**

**Harmar.**  
Manufacturing Partner

 **Juneau Stair Lifts**

\*Images used for this brochure are for illustrative purposes only, please refer to local code for specific configuration requirements.