

# Quantify HS DNA Kit for Qubit measurement

Cat No: QUDNA1B002 for 500 assays

SuGenomik Quantify HS DNA Measurement Kit is ideal for quantitating DNA samples in Next Generation Sequencing applications.


This kit allows accurate dsDNA quantification in the range from 10 pg/ $\mu$ L to 100 ng/ $\mu$ L for initial sample.

The main kit component is Quantify Assay Solution. Solution includes dsDNA HS fluorescent dye,

This kit allows convenient dsDNA quantification in the broad concentration range and does not require a large

This kit is appropriate for the measurement of DNA with all versions of the Qubit fluorimeter.

All measurements are performed at room temperature and take on average 15–20 minutes (depending on the sample amount). To perform a measurement, prepare Quantify Mix, dilute DNA standards and experimental samples with it in individual tubes and perform fluorescence measurements using fluorimeter.

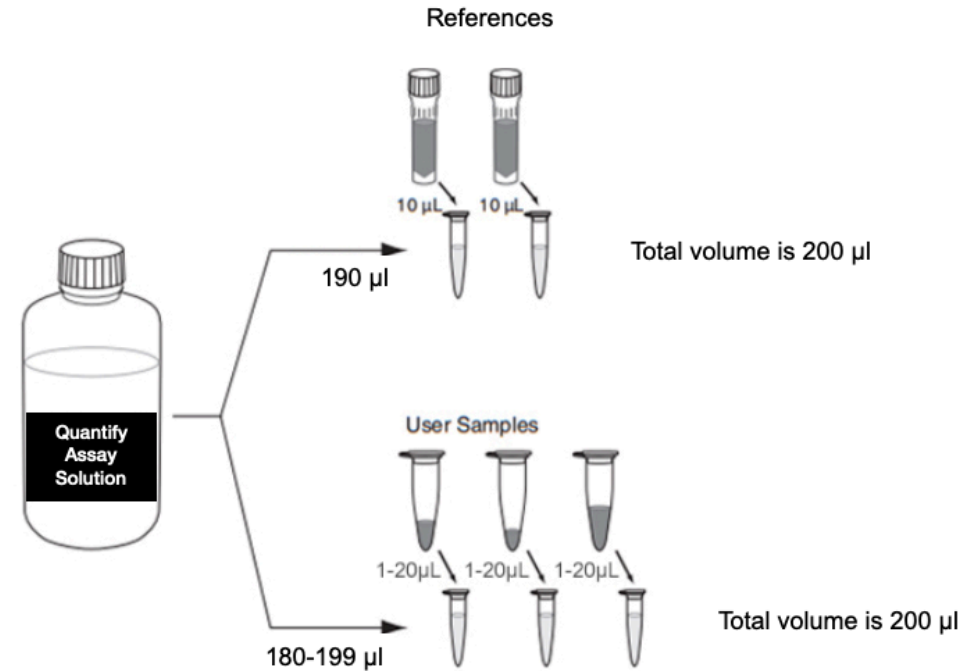

SuGenomik  
Biyoteknoloji Ltd Şti  
Address: METU  
Technopolis, Ankara.  
[www.sugenomik.com](http://www.sugenomik.com)

