







Lightweight, umbrella-fold special needs stroller for children and adults. It has innovative <code>BodyPrint</code>© seat equipped with belts that can easily be adjusted to the patient. The seat can be flat or concave depending on the user's needs and dysfunctions. Fixed <code>30°</code> backrest tilt angle (bucket) reduces the pressure on soft tissue preventing bedsores. It also gives you many other benefits.

Except for being convenient and stable, **OMBRELO** stroller is also handy and easy to maneuver.

It is ideal for transporting patients to hospitals whether for treatments or for tests. It is also good for having a long walk in the fresh air. Very simple and easy to use.

Adjustable swing-away **FlipFlap©** footplates make it easy to climb and get off the stroller. The product is available in 5 sizes, the biggest of which is designed for users who weight up to 120 kg. Have you got special requirements or look for excellent stability? Our wide range of accessories makes it easy for anyone to find the right product they are looking for.





Stroller is equipped with Light Grey rims. However there is possibility to change the rim's colour to black.

IT IS RECOMMENDED IN THE FOLLOWING CASES:

- Cerebral palsy various types and forms.
- Brain and spinal cord injuries with lower limbs paralysis and weakness, balance and motor coordination disorders with different etiology.
- Demyelinating disease (e.g. multiple sclerosis SM).
- Meningocele.
- Muscles diseases with weakness and paralysis (dystrophy, myasthenia gravis, myositis).
- Genetic disorders.
- Degeneration and neurological diseases.
- Other condition connected with muscles weakness, paralysis and musculoskeletal disorders.
- After head and spinal cord injuries and spine trauma with paralysis and paresis.

STANDARD EQUIPMENT:

Light foldable frame, folding in umbrella shape | Wheels with glass fiber, front PU swivel, rear | inflatable with 2 separate wheel locks | Multi adjustable footplates FlipFlapTM with 2 way | foot stabilizing belts | Self-tensioning upholstery BodyprintTM for perfect | positioning | Adjustable seat depth | 5 points safety belt with height adjustment | Seat - 30° fixed tilt angle | Seat to back 85° - 95° angle adjustment | Adjustable build-up headrest | Backrest | Adjustable push handles | Reflective safety marks | Comfort upholstery for size 1-4 | Flex upholstery for size 5

Our passion - your independence!

Akces-MED sp. z o.o. Jasionka 955B 36-002 JASIONKA



OMBRELO™



ADVANTAGES:



Stable structure

The device ensures adequate stability and safety of Patient.



Easy to use

Intuitive way of adjustment and easy assembly.



Lightweight construction

The product is made of light metal alloys which positively affects the low weight of the device.



Flip Flap© footrest

Special footrests system ensures a wide range of adjustment. Reclining platforms make getting on and off a device easier.



High durability

The construction ensures high durability of the device – maximum user weight is 120 kg!



Reflective elements

The product has reflective elements to ensure your safety.



24 months warranty

We offer 24 months warranty on the product. We also provide after sales service.



5-point foot straps that allow for peripheral or lower-back stabilisation.



BODY PRINT© upholstery is a unique system of belts which, when adjusted perfectly, offer optimal positioning.



Height adjustable Flip Flap® footplates with abduction and adduction angle adjustment make the stroller grow with the child.

COLORS



Blue Sea Size 1 (Comfort)



Green Grass Size 2 (Comfort)



Hot Red Size 3 (Comfort)



Brown Sugar Size 4 (Comfort)



Blue Horizon Size 5 (Flex)

You can change the colors between sizes on request. The waiting time – 30 days.




Comfort upholstery – breathable



Flex upholstery – washable and durable.





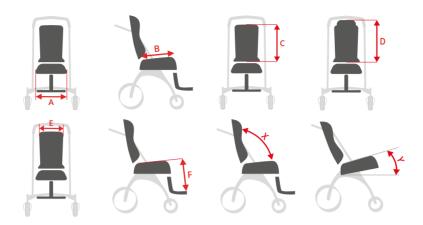

Symbol	Measurement	Unit	Size 1	Size 2	Size 3	Size 4	Size 5
А	Seat width	[cm]	23	29	34	39	43
В	Seat depth	[cm]	21-26	28-36	29-40	31-40	34-41
С	Backrest height (without hedrest)	[cm]	51	58	55	63	60
D	Backrest height (with hedrest)	[cm]	43-64	46-72	46-72	53-76	54-80
Е	Backrest width	[cm]	23	29	34	39	43
F	Footrests length	[cm]	16-64	16-64	16-64	16-64	16-64
	Footer lenght	[cm]	18	20	21	23	23
	Push handle height	[cm]	101-106	101-106	103-108	100-105	100-105
X	Seat tilt angle	[°]	90	93	91	90	90
Υ	Backrest tilt angle	[°]	31	25	26	26	28
	Max user weight	[kg]	35	45	55	75	120
	Product dimensions and weight		Size 1	Size 2	Size 3	Size 4	Size 5
	Width	[cm]	52	57	63	67	72
	Lenght	[cm]	125	125	130	137	137
	Height	[cm]	100	101	104	101	108
	Weight	[kg]	17	18	18,5	20	20,5

Note: Measurement tolerance range is +/- 1 cm / $\frac{1}{2}$ in.



Akces-MED sp. z o.o. Jasionka 955B 36-002 JASIONKA



OMBRELOTM





Safety



Led lights (2 pieces)

Front and rear LED lights for increased safety.



Soft cover for tray

Soft cover protects the user from hitting the tray. Indicated for patients with impaired limbs coordination.



Patient lift belt

Robust and ergonomic design as well as high durability facilitates daily care of the patient.



Stabilization



H vest

It acts as a trunk stabilizer. HOLD&PULL© system makes it easier to adjust the vest to the patient body. You can quickly detach the vest if necessary.



4 point vest

It acts as a trunk stabilizer. Wide range of adjustment helps to fit it to the user's needs. Equipped with HOLD&PULL© system.



6 point vest

Wide range adjustable vest that acts as a trunk stabilizer. HOLD&PULL© system makes it easier to adjust the vest to the patient body.



Thigh abduction belts

They can be easily attached to the chair by use of Velcro straps. Their use allows to keep legs in abduction and also prevents the patient from slipping off of the device.



Thigh abduction belts

They can be easily attached to the chair by use of Velcro straps. Their use allows to keep legs in abduction and also prevents the patient from slipping off of the device.



Pelvic belt

Wide belt equipped with "HOLD AND PULL" system. I provides better stabilization in the device.



Trunk stabilizing belt

It acts as an additional trunk protection which helps to keep a proper sitting position. The belt is fastened by use of Velcro straps and its size can be adjusted to the user.



Trunk belt

Universal solution applicable to many of our products. It secures the user from falling out of the device. Anatomical shape and Velcro straps allow for a comfortable fit.



Thigh strap

Helps maintain a proper sitting position. Velcro-adjustable straps as well as ergonomic design allow for a better fit for the patient.



Calf belt

It provides additional protection for lower limbs. Mounted on the front side limits the legs extension. Mounted on the rear side limits the legs bending.



Elastico headrest

Covered with thermo-active foam easily fits the shape of a head. Offers more comfortable support and better stabilisation. Attached in the place of a standard headrest.



Firm headrest made of foam

Made of poliurethane foam. Provides with basic support of the head



Akces-MED sp. z o.o. Jasionka 955B 36-002 JASIONKA



OMBRELOTM

accessories

AKCESMED



Head-neck support

Perfectly stabilizes the head. Indicated for people with limited head control. You can mount it at different height.



Lumbar support

It allows better stabilization of the lumbar back.



Abduction block

Shaped abduction block provides with correct seating position. Can be easily removed



Rigid side supports

Longer and more complete stabilization of the lower limbs indicated for patients with excessive leg abduction.



BodyMap Stabilization



BodyMap A seat cushion

Vacuum seat cushion allows for stabilisation and proper positioning of hips.



BodyMap B+ back cushion

Vacuum back cushion allows for the stabilisation of trunk. Helps maintain the correct body posture in a chair.



BodyMap D headrest

It can be easily modelled to the patient's head. It has also the self-adhesive function on one side which makes it possible to attach additional protection accesso-




Vismemo cover for BodyMap cushion Made of thermo-active foam for extra

bed sore prevention. Allows for longer and more comfortable use of the cu-



tion and absorbs moisture.

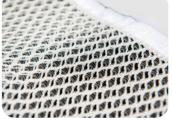

BodyMap spacer

Aerating-massaging pad to strengthen the prophylactic anti-bed sore function in patients using BodyMap. It can only be used together with the terry cover The effect of being 'suspended in the air' creates a skin-friendly microclimate.



Comfort



Plexiglas® Tray

Made of Plexiglas®. It supports the therapy and allows you to play, eat, and



Made of plastic. It supports the therapy and allows you to play, eat, and work.



Cup holder

An universal holder mounted to the frame. It enables easy and quick access to a bottle of drink.



Umbrella

Attachable to the stroller's frame, perfectly protects from the sunlight or



Akces-MED sp. z o.o. Jasionka 955B 36-002 JASIONKA







Headrest cover

It helps to keep the headrest clean. Easy to assemble/disassemble, it can be washed



Rain Cover

An excellent protection from the rain, snow, and wind either for the user and stroller



Protects against both sun and rain.

Folding canopy without side panels



Folding canopy with side panels

Side panels protect not only against sun and rain but also the wind.



Summer sleeping bag

Protects legs and trunk. This esthetic body sheet fits the height of the user covering also the footrest of the pushchair.



Winter sleeping bag

Protects legs and the torso area. A very tasteful cover of the body, adjusting to the height of the user, covering the footrest of the chair. Insulated lining provides extra protection against the cold.




Winter hood

Can be attached to the sleeping bag Protects head against frost, wind and moisture.



Backrest cushion - Elastico

It perfectly adapts to the patient's body. Made of dense shape memory foam covered with breathable and spatial fabric. It is indicated in anti-bedsore prophylaxis.

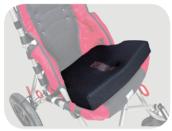

Seat cushion – Elastico

It precisely adapts to the patient's body. Made of memory foam covered with breathable and spatial fabric. It is indicated in anti-bedsore prophylaxis.



Seat cushion - profiled thighsAnatomical notch on the thighs improves sitting comfort.



Seat cushion - profiled Efficient Tilt--in-space adjustment

The use in the product of an innovative four-bar linkage, supported with the gas spring, greatly improves smooth Tilt-in--space adjustment. It perfectly adapts to the patient's body. Made of dense shape memory foam covered with breathable and spatial fabric. It is indicated in anti-bedsore prophylaxis.



Wheel Covers (4 pcs.)

Use it on the dirt tires when you want to take the products indoor or transport in the trunk.



Mosquito Net

Excellent protection against insects



Lifting Bar

It is much easier to overcome terrain obstacles. It facilitates the tilt of the trolley, curbs, steps and other vertical pavement elements.



Hand Brake

The handbrake suddenly allows stop the stroller. It also facilitates maneuvering the wheelchair. The brakes are equipped with a four-stage lock, so it functions as a parking brake





Akces-MED sp. z o.o Jasionka 955B 36-002 JASIONKA

OMBRELO™





Storage



Saddle bags (pair)

Saddle Bags attach to the side of the chair, and are perfect for storing necessary items. Made from black durable



Storage basket

Enables you to carry small objects while you stroll. It is attached to the base of the stroller.



Women' bag

This capacious bag designed for mothers will be perfect for carrying all the necessities.



Men's bag Capacious men-style bag for carrying all the necessary things.



DeLux bag

Large, durable and stylish bag useful for carrying necessary things. Made from high-quality materials.



Stroller cover 'n carry bag

Perfect solution to protect folded stro-ller against damage (e.g. during transport). Easily fits the shape and the size of a stroller.



Oxygen tank case

Oxygen tank case can be attached to the stroller by use of Velcro straps with adjustable length. Made to measure to perfectly fit the patient's oxygen tank



Build-up



Kid-sit buggypod IO

For convenient and safe transportation of a second child without having to buy a double stroller.



Buggypod LITE

Starters are very helpful if you want to take two children for a walk - in this case, you just need to mount the attachment to a larger trolley. Thanks to the special buckles, it can be easily adjusted to Ursus, Ombrelo, Hippo+



Buggypod PERLE adapter

Perle is a well-designed seat that fits into the BuggyBoard Maxi.



Buggypod Perle seat

Perle is a well-designed seat that fits into the BuggyBoard Maxi.



Rain cover for the kid-sits

An excellent protection from the rain, snow, and wind either for the user and



Extension bar

An adapter enabling to mount the Buggypod kid-sits to the strollers of other manufacturers.



Mounting set

It enables easy assembly of the Buggypod kit-sits to other strollers whose frame diameter is larger than 40mm.



Conversion kit (for The buggypod kit-sits)

Set of elements thanks to which the Buggypod kit-sit can be converted into a fully functional stroller.



Akces-MED sp. z o.o. Jasionka 955B 36-002 JASIONKA