



Your choice of nursing bed doesn't have to be a compromise

Even with the diverse demands of clinical care, staff and resident safety, and efficiency – we've managed to bring these all together in the Empresa bed.

Incorporating essential clinical features to aid safe nursing and residents' independence, the Empresa bed has people at the heart of its design.



Proven solution for falls management

Class-beating FloorBed® technology

With falls as the leading cause of death and serious injury in older people, our floor level capability provides proven effectiveness in the management of bed falls and prevention of associated injuries. Height is a crucial factor in bed selection; the impact force of a fall from a bed at 10 cm is 50% lower than at 20 cm.



Summary of product evaluation results

Description		Post FloorBed	Trend
Residents who fell	9 (of 12)	5 (of 12)	≈ 44%
Falls in bedroom (excl. bath)	15	3	≥ 80%
Notable injury and/or required diagnostic/medical intervention	2	0	≈ 100%
Falls from bed	7	0	≥ 100%

Increased safety

Where residents need extra protection, the addition of a high safety mat can help reduce the risk of falls even further, by transforming a potential fall into a simple roll.

Falls within the bedroom
were reduced by 80%
and falls from the bed
were eliminated²

80%



When moving from a supine to sitting position in bed, the body naturally lengthens. Empresa's BodyMove™ auto-regression feature


accommodates this change, so not only does it help with pressure distribution and reduction of shear for enhanced wound care support, it improves comfort too.


Prolonged mechanical loading of soft tissues covering bony prominences, such as when an individual is in bed, may lead to degeneration of skeletal muscle tissue.

Pressure redistribution and repositioning of an individual is undertaken to reduce the duration and magnitude of pressure over vulnerable areas of the body.5

Reduction in average pressure over the sacral area, when compared to a standard nursing bed4

Lower interface pressures

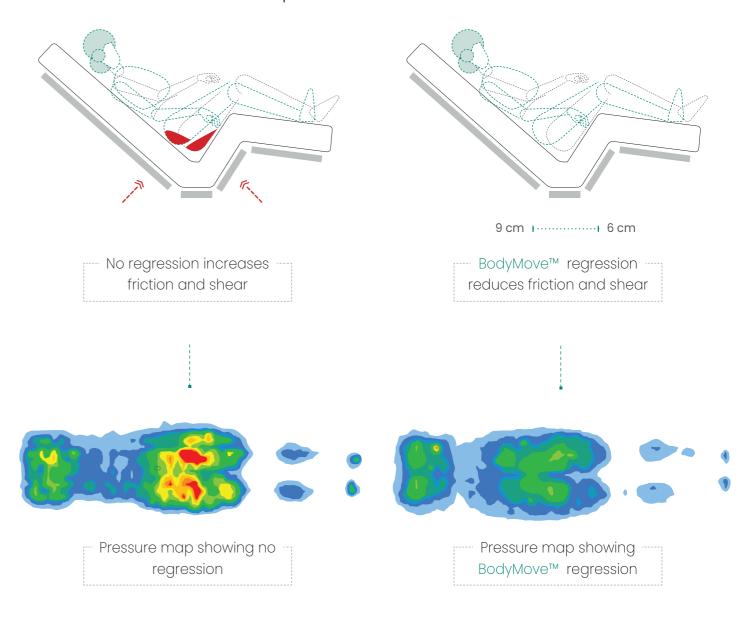

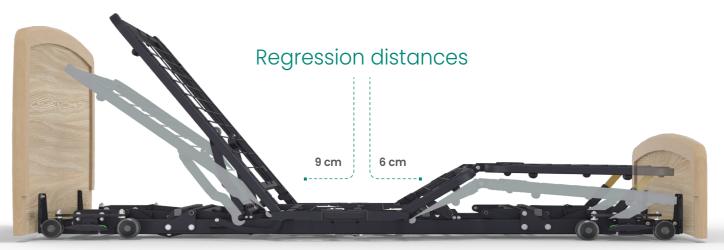



EMPRESATM FEATURING FLOORBED TECHNOLOGY

Human factors at the **centre of design**

Driven by an understanding of your requirements, a range of integrated functionality comes as standard to support caregiving and make your life easier.



Ergonomic handling

Pre-programmed "one-touch care"

Easy positioning of residents at one touch of a button, helps negate the need for lifting and saves staff time and effort. Comfortable to hold handset, accessible from anywhere around the bed.

80 cm bed height

Full height adjustment for caregiver safety and aiding the delivery of high quality care to residents.







- Comfort position --

Safety first

Lockout feature

Prevents unauthorised and inadvertent use of the handset.

Auto safety-stop

Inhibits bed from being lowered below 20 cm, without caregiver action.



Tandem castors with UniSafe™ braking system

Bed rolls effortlessly and needs only a single press of a pedal to secure each set of brakes.

Frame-mounted lifting pole -----

Ensures clear access under the bed for placement of equipment.

Quick release, -----removable bed lever

Facilitates bed mobility and transfers.






EMPRESATM FEATURING FLOORBED TECHNOLOGY

Right product, right time, right place




Raises vertically

Uses a vertical lift mechanism to raise the bed, reducing space demands and helping prevent damage to walls.

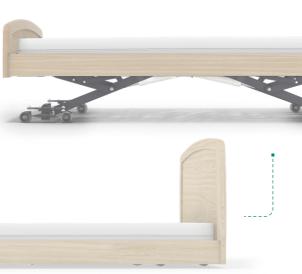

Simplified decision making

The Empresa bed enables you to standardise on bed choice, helping to streamline procurement, reduce training needs and future-proof your requirements.

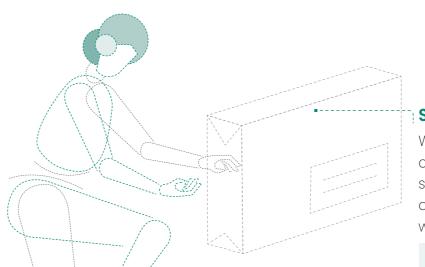

Support from start to finish

Working as part of your team, we will help capture your requirements and refine your specification. And our complete turnkey service and dedicated project management team will make implementation smooth and easy.

Quick ship wood **designs**





Quick ship wood **finishes**








Other wood finish options available.

Full range of head board **options**




Only Arc and Skandi head boards available as quick ship options

Available **accessories**

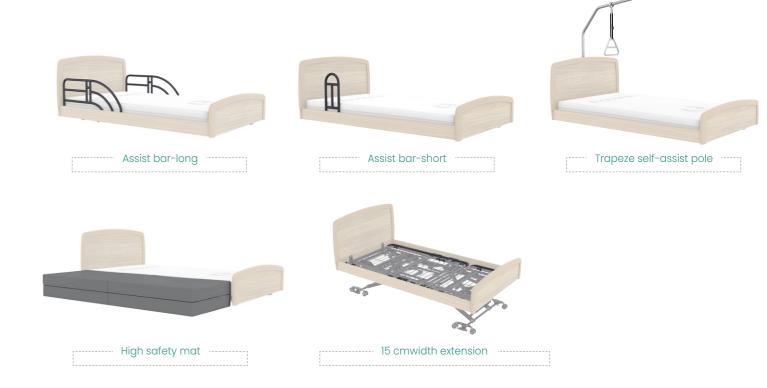

Technical specification

Warranty	Frame: 10 years Electric/parts: 3 years	
Bed weight	120 kg	
Max. patient weight	240 kg	
Safe working load	275 kg	
Liquid ingress protection	IPX6	
Trendelenburg function	12 degrees	
Legrest max. angle	36 degrees	
Backrest max. angle	70 degrees	
Auto regression	BodyMove technology	
Maximum lifting pole load	75 kg	
Castor diameter	63 mm	
8 castors	4 lockable castors; 4 castors to aid manoeuvrability	
UniSafe braking system	2 sets of brakes, each secure 4 castors	
Mattress platform height	100 mm to 800 mm	
Mattress size standard width Mattress size with width extension Mattress size with length extension	900 mm width x 2000 mm length 1050 mm width x 2000 mm length 2200 mm length	
Overall dimensions	1015 mm width x 2090 mm length	

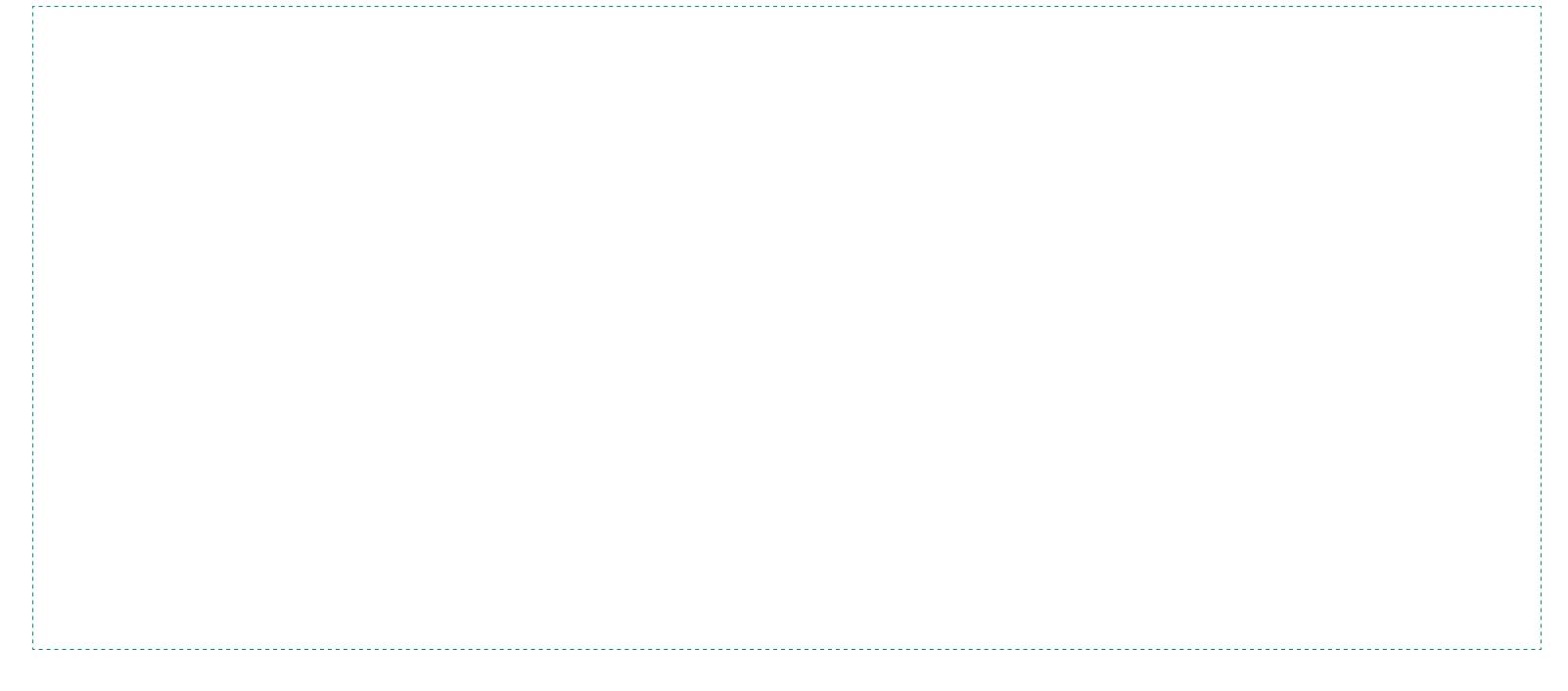

Reference

- 1. WHO (26th April 2021) Falls. Available from: https://www.who.int/news-room/fact-sheets/detail/falls [Accessed 27th May 2021]
- 2. Lutheran Home. FloorBed: resident falls analysis. Interim report. February 2019
- Zaphir J (2014) Lowest heights of floor level beds: fall velocity and impact force analysis and resulting injuries.
 Available from https://vital.net.au/sites/default/files/zaphir_study_2014_-_bed_height_impact.pdf [Accessed 5th May 2021]
- 4. Data on file. February 202
- 5. Black J et al (2007) National pressure ulcer advisory panel's updated pressure ulcer staging system.

 Dermatol. Nurs. 19(4):343–349. Available from https://link.springer.com/content/pdf/10.1007/s10439-010-0002-x.pdf. [Accessed 3rd May 2021]
- and ™ are trademarks of Accora Limited. As our policy is one of continuous improvement, we reserve the right to modify designs without prior notice.

Accora Inc.

9210 Corporate Blvd. Suite 120, Rockville, MD 20850
T: 301 560 2400 E: information@accora.care W: accora.care

Manufactured by

Accora Ltd Charter House, Barrington Road, Orwell, Cambridge SG8 5QP

T: +44 (0)1223 206 100 E: info@accora.care W: www.accora.care

