



ULISES EVO™

special stroller



ULISES EVO™ is a perfect stroller designed for daily rehabilitation of children and youth with locomotor dysfunctions. It is particularly useful for patients requiring additional stabilization.

The stroller has adjustable seat depth and backrest tilt angle which greatly facilitates its functionality.

With simple and light construction it can be easily folded for transporting purposes. Furthermore, relatively small size allows you to fit **ULISES EVO™** into an undersized car trunk.

Stable structure of the stroller, belts and supports system as well as reflective elements make the user feel safe. The backrest can be easily bend back and adjustable push handle increases the stroller maneuverability.

Rich equipment and wide range of adjustments make **ULISES EVO™** fit easily to the individual needs.



IT IS RECOMMENDED IN THE FOLLOWING CASES:

- Cerebral palsy – various types and forms.
- Brain and spinal cord injuries with lower limbs paralysis and weakness, balance and motor coordination disorders with different etiology.
- Demyelinating disease (e.g. multiple sclerosis – SM).
- Meningocele.
- Muscles diseases with weakness and paralysis (dystrophy, myasthenia gravis, myositis).
- Genetic disorders.
- Degeneration and neurological diseases.
- Other condition connected with muscles weakness, paralysis and musculoskeletal disorders.
- After head and spinal cord injuries and spine trauma with paralysis and paresis.



STANDARD EQUIPMENT:

Flex upholstery | 5 points safety belt | Adjustable footrest with 2 points foot stabilizing belts | Wheels with glass fiber, rear inflatable with 2 | separate wheel locks | Front PU swivel | Angle adjustable backrest | Adjustable push handles | BodyPrint seat © Belt system adjusting the tension of the seat upholstery and backrest | Adjustable height and tilt angle of the footrest | Adjustable depth of the seat | Headrest | Backrest extension

Our passion – your independence!

Colors may differ from the presented images. The manufacturer reserves the right to make changes in products without updating the directory. Current offer on our website.

Akces-MED sp. z o.o.
Jasionka 955B
36-002 JASIONKA

export@akces-med.com
+48 (17) 864 04 77
www.akces-med.com



ADVANTAGES:



Stable structure

The device ensures adequate stability and safety of Patient.



Easy to use

Intuitive way of adjustment and easy assembly.



Lightweight construction

The product is made of lightweight and durable materials.



24 months warranty

We offer 24 months warranty on the product. We also provide after sales service.



Five-point straps for keeping safety.



Front swivel wheels for easy maneuverability.



Inflated rear wheels provide cushioning.

High rigid lateral support can be mounted.


COLORS




Size 1
Green Grass
(Comfort)




Size 1a
Hot Red
(Comfort)

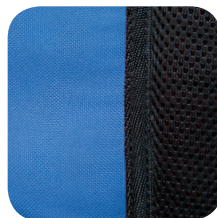

Size 2
Blue Sea
(Comfort)

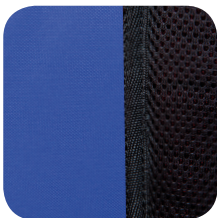


Size 2a
Blue Wave
(Comfort)

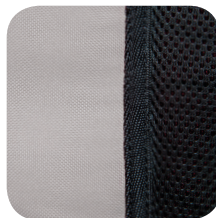

Size 3
Gray Volcano
(Comfort)



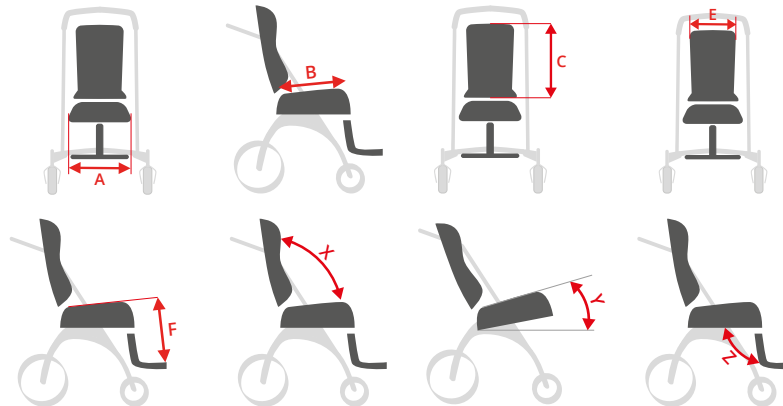

Symbol	Measurement	Unit	Size 1	Size 1a	Size 2	Size 2a	Size 3
A	Seat width	[cm]	31	31	37	37	43
B	Seat depth	[cm]	22-29	29-37	29-37	29-37	35-45
C	Backrest height	[cm]	60	68	68	68	76
E	Backrest width	[cm]	28	28	34	34	40
F	Footrests length	[cm]	19-25	24-32	24-32	30-40	30-40
	Headrest height	[cm]	40-60	43-68	43-68	51-76	51-76
	Footer length	[cm]	19	19	21	21	25
	Push handle height	[cm]	102	107	107	112	119
X	Backrest tilt angle	[°]	90-150	90-150	90-150	90-150	90-150
Y	Seat tilt angle	[°]	13	13	13	13	13
Z	Footrests tilt angle	[°]	90-180	90-180	90-180	90-180	90-180
	Max user weight	[kg]	55	55	75	75	75

Product dimensions and weight		Size 1	Size 1a	Size 2	Size 2a	Size 3
Width	[cm]	57	57	63	63	69
Folded width	[cm]	57	57	63	63	69
Length	[cm]	94	110	110	127	127
Folded length	[cm]	104	107	107	114	114
Height	[cm]	102	107	107	112	112
Folded height	[cm]	40	45	45	50	50
Weight	[kg]	18,5	19.5	20	21.5	22

Note: Measurement tolerance range is +/- 1 cm / ½ in.



Safety



Led lights (2 pieces)
Front and rear LED lights for increased safety.



Soft cover for tray
Soft cover protects the user from hitting the tray. Indicated for patients with impaired limbs coordination.




Patient lift belt
Robust and ergonomic design as well as high durability facilitates daily care of the patient.



Barrier
Makes the user safe in the stroller. Protects against falling out.



Stabilization

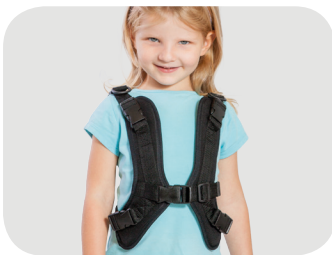

H vest
It acts as a trunk stabilizer. HOLD&PULL© system makes it easier to adjust the vest to the patient body. You can quickly detach the vest if necessary.



4 point vest
It acts as a trunk stabilizer. Wide range of adjustment helps to fit it to the user's needs. Equipped with HOLD&PULL© system.



6 point vest
Wide range adjustable vest that acts as a trunk stabilizer. HOLD&PULL© system makes it easier to adjust the vest to the patient body.



Thigh abduction belts
They can be easily attached to the chair by use of Velcro straps. Their use allows to keep legs in abduction and also prevents the patient from slipping off of the device.




Thigh abduction belts
They can be easily attached to the chair by use of Velcro straps. Their use allows to keep legs in abduction and also prevents the patient from slipping off of the device.



Pelvic belt
Wide belt equipped with "HOLD AND PULL" system. It provides better stabilization in the device.



Trunk stabilizing belt
It acts as an additional trunk protection which helps to keep a proper sitting position. The belt is fastened by use of Velcro straps and its size can be adjusted to the user.

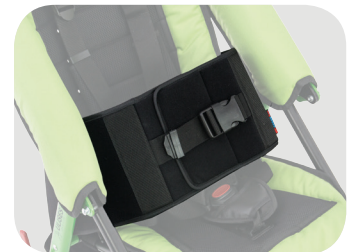

Trunk belt
Universal solution applicable to many of our products. It secures the user from falling out of the device. Anatomical shape and Velcro straps allow for a comfortable fit.



Lumbar support
It allows better stabilization of the lumbar back.



Thigh strap
Helps maintain a proper sitting position. Velcro-adjustable straps as well as ergonomic design allow for a better fit for the patient.



Calf belt
It provides additional protection for lower limbs. Mounted on the front side limits the legs extension. Mounted on the rear side limits the legs bending.



Abduction block
Shaped abduction block provides with correct seating position. Can be easily removed.



Elastic headrest

Covered with thermo-active foam easily fits the shape of a head. Offers more comfortable support and better stabilization. Attached in the place of a standard headrest.



Head-neck support

Perfectly stabilizes the head. Indicated for people with limited head control. You can mount it at different height. Head-neck support



Firm headrest

Made of poliurethane foam. Provides with basic support of the head



Seat narrowing cushions 3cm

Made of elastic foam trimmed with soft material. They are attached to the seat upholstery from both sides. Thanks to them you can narrow the seat by 6cm.



Seat narrowing cushions 5cm

Made of elastic foam trimmed with soft material. They are attached to the seat upholstery from both sides. Thanks to them you can narrow the seat by 10cm.



5-point foot straps

The ergonomic design of these criss-cross foot straps perfectly protects both foot and heel.



Low rigid lateral support with narrowing pads (3 cm).

Protects the pelvic area on the both side of the stroller.



Low rigid lateral support with narrowing pads (6 cm).

Protects the pelvic area on the both side of the stroller.



High rigid lateral support

Longer and more complete stabilization of the lower limbs indicated for patients with excessive leg abduction.



Low rigid lateral support

Protects the pelvic area on the both side OF the stroller. Made of Velcro surface to which we can fix the seat narrowing pads. We do so in case of very slim and tall patients.



Comfort



Umbrella

Attachable to the stroller's frame, perfectly protects from the sunlight or rain.



Tray

Made of plastic. It supports the therapy and allows you to play, eat, and work.



Folding Canopy

Protects against both sun and rain.



Headrest cover

It helps to keep the headrest clean. Easy to assemble/disassemble, it can be washed.



Cup holder

An universal holder mounted to the frame. It enables easy and quick access to a bottle of drink.



Rain Cover

An excellent protection from the rain, snow, and wind either for the user and stroller



Mosquito Net

Excellent protection against insects



Summer sleeping bag

Protects legs and trunk. This esthetic body sheet fits the height of the user covering also the footrest of the pushchair.



Winter sleeping bag

Protects legs and the torso area. A very tasteful cover of the body, adjusting to the height of the user, covering the footrest of the chair. Insulated lining provides extra protection against the cold.



Winter hood

Can be attached to the sleeping bag. Protects head against frost, wind and moisture.



Seat cushion - Elastico

It precisely adapts to the patient's body. Made of memory foam covered with breathable and spatial fabric. It is indicated in anti-bedsore prophylaxis.



Backrest cushion - Elastico

It perfectly adapts to the patient's body. Made of dense shape memory foam covered with breathable and spatial fabric. It is indicated in anti-bedsore prophylaxis.



Seat cushion - profiled thighs

Anatomical notch on the thighs improves sitting comfort.



Seat cushion - profiled Efficient Tilt-in-space adjustment

The use in the product of an innovative four-bar linkage, supported with the gas spring, greatly improves smooth Tilt-in-space adjustment. It perfectly adapts to the patient's body.



Wheel covers (4 pcs.)

Use it on the dirt tires when you want to take the products indoor or transport in the trunk.



Storage



Women' bag

This capacious bag designed for mothers will be perfect for carrying all the necessities.



Men's bag

Capacious men-style bag for carrying all the necessary things.



DeLux bag

Large, durable and stylish bag useful for carrying necessary things. Made from high-quality materials.



Storage basket

Enables you to carry small objects while you stroll. It is attached to the base of the stroller.



Oxygen tank case

Oxygen tank case can be attached to the stroller by use of Velcro straps with adjustable length. Made to measure to perfectly fit the patient's oxygen tank



Stroller cover 'n carry bag

Perfect solution to protect folded stroller against damage (e.g. during transport). Easily fits the shape and the size of a stroller.