

**Configura<sup>®</sup>Advance**  
by Accora

# **Comprehensive posture management**

designed for high need users

**Accora**





# Configura<sup>®</sup> Advance



Width adjustment from 350 mm to 550 mm



Dropdown armrest option – helps with transfers and sling access



Legrest elevation with 37 degrees of tilt-in-space



10 degrees of forward tilt



**BodyMove** 'Reduced shear' backrest recline



Depth adjustment from 400 mm to 550 mm



Integrated self-supply overlay



Continuous seat overlay for thigh support



Angle adjustable footplate



Height adjustable footplate



# Optional postural backrest

All cushions height and depth adjustable

Angle adjustable pillows




Optional lateral supports

ClickLock – easy-to-adjust bracket for flexible adjustment and positioning of backrest cushions

Optional profiled headrest

## Cocoon™ backrest

The **Cocoon** backrest is our latest option to fit the Advance. It combines the flexibility of the multi-adjustable pillow backrest with the support of the lateral support backrest.

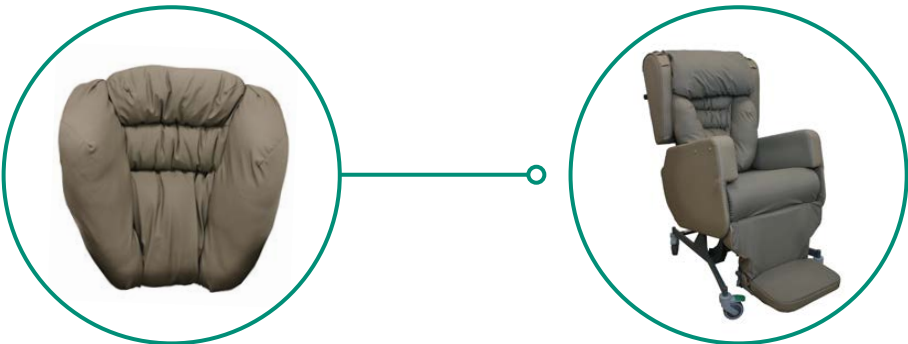

The “reduced shear backrest” mechanism on the **Configura** Advance chair complements this backrest by giving the client extra trunk support, whilst retaining comfort and pelvic stability.

## BodyMove™

Using unique **BodyMove** technology, the backrest and legrest of the Advance chair adjust in close alignment to the biomechanics of the human body. This reduces shear forces on the back and legs, improving comfort and posture.



‘Reduced shear’ recline to help maximise pelvic stability and assist in preventing tissue damage



Full-length compensation on legrest, giving superior leg elevation and support

# Benefits

The Configura Advance is a comfort chair suitable for multi-use care environments, with a high level of adjustment for size, pressure and postural requirements.



## User group

### The Configura Advance chair is appropriate for:

- Non-ambulant users
- Users with limited knee extension
- Users who need to be transported from room to room
- Moderate postural requirements
- Users up to 160 kg (25 st)
- Complex postural requirements, e.g. fixtures and contractures
- Users with stiffness or fixtures at the hips
- Stroke and neurological conditions



### Configure

Using the **Configura** Assessment Form, select the perfect chair for your service user. You can select different sizes and choose different pressure-relieving and posture management options to meet your service user's needs.



### Reissue

The **Configura** chair can be easily and cost-effectively recycled because all the contact areas can be removed for cleaning or replacement. The chair can then be re-configured for the next service user...



### Deliver

The chair is delivered within 5 working days, configured exactly to the specification chosen on your **Configura** Assessment Form. Adjustments can be made with the service user on-site to ensure maximum comfort and effectiveness.



### Adjust

If the service user's needs change over time, you can easily adjust the chair to suit. For example you can increase the specification of the pressure relief or introduce more lateral support. This can remove the need to issue a new chair.

## Specific dimensions

Description	Values
Seat height	600 mm (24")
Seat width	350 mm - 550 mm (14"-22")
Seat depth	400 mm - 550 mm (16"-22")
Footplate height adjustment	325 mm - 525 mm (13" - 21")
Armrest height	200 mm (8") as standard – adjustable with dropdown armrest selection
Overall width	690 mm (27")
Overall height – upright	1320 mm (52")
Overall length – without footplate	940 mm (37")
Overall length – with footplate	1290 mm (51")
Total length fully tilted, reclined and legs elevated	1870 mm (74")
Tilt-in-space range	7 degrees - 37 degrees
Forward tilt	10 degrees
Backrest recline	130 degrees

**Accora**

**UK:** Charter House, Barrington Road, Orwell, Cambridge SG8 5QP  
**T:** +44 (0)1223 206 100 **E:** [info@accora.care](mailto:info@accora.care) **W:** [www.accora.care](http://www.accora.care)

**IE:** 38 Main Street, Swords, Co. Dublin, Ireland K67 E0A2  
**T:** +353 (0) 1 695 0614 **E:** [info@accora.care](mailto:info@accora.care) **W:** [www.accora.care](http://www.accora.care)

\* and ™ are trademarks of Accora Limited.  
As our policy is one of continuous improvement, we reserve the right to modify designs without prior notice.