




## Next Level

A new generation of toilet aids that create possibilities

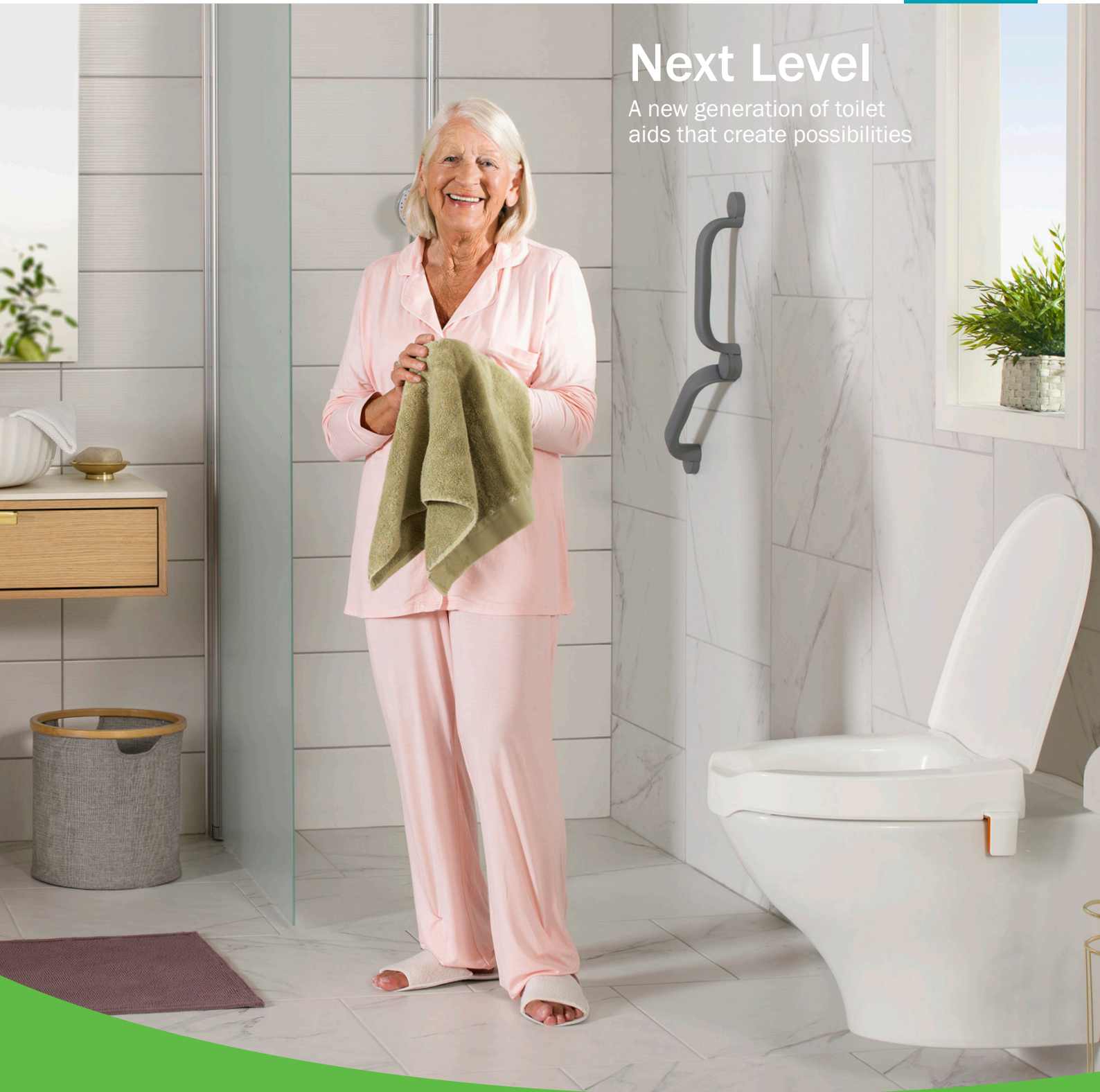

## Etac My-Loo

Raised toilet seat,  
with brackets

 **etac**<sup>®</sup>



# Etac My-Loo

## Raised toilet seat, with brackets

Together with experienced industrial designers, we have developed the new raised toilet seat with brackets – Etac My-Loo. A product that gives concepts such as hygienic, safe and appealing a new meaning. With its gently sweeping lines and excellent ergonomics, it offers extra-comfort and a design ideally suited to the modern bathroom.

### Designed for function and comfort

Etac My-Loo has been developed to blend into the bathroom environment and has a design that gives optimal comfort and accessibility. The gentle lines and the large opening make My-Loo suitable for most users. The seat raiser is offered in two heights, with or without a lid.

### Easy to keep clean and fresh

My-Loo raised toilet seat, with brackets, consists of a small number of parts made of durable materials. It is easy to wipe off with a soft cloth. Should separate cleaning be required, it is easy to remove the seat and equally simple to put it back in place.

### Safe, quick installation

It takes only a few seconds to install Etac My-Loo! You fix it in place using two specially designed brackets which are easily adapted to different toilet shapes. The extra-large friction surfaces and flexible material of the brackets provide a firm attachment to the porcelain.

Put the raised toilet seat on the toilet and push the brackets into place. Click, and the installation is finished!



Etac My-Loo raised toilet seat, with brackets	Item no.	Material
6 cm	80301520	Polypropylene and TPE
10 cm	80301521	
6 cm with lid	80301522	
10 cm with lid	80301523	
<b>Weight</b>		<b>Design</b>
6 cm: 0.8 kg	10 cm: 1.1 kg	Myra Industriell Design AB
6 cm with lid: 1.2 kg	10 cm with lid: 1.4 kg	
<b>Maximum user weight</b>		<b>Maintenance</b>
190 kg		
		Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution. The product can be safely cleaned in a cabinet washer-disinfector at 85°C for 3 minutes.



## Features



### Easy to keep clean and fresh

The clean, smooth surfaces are easy to clean with a soft cloth. Should separate cleaning be required, the seat is easy to remove and equally simple to put back in place.



### Designed for function and comfort

The soft shape of the seat gives safe support and high comfort. An extra large opening and generous recesses facilitate accessibility during personal hygiene.



### Safe, quick installation

The easily attached brackets with wide adjustability and large friction surfaces make My-Loo at all times safely fixed to, and easily mounted on, any toilet porcelain. Entirely without screws.








### Variants

My-Loo raised toilet seat, with brackets, is available in two different heights, 6 and 10 cm, with or without lid.



### Dimensions

Outer width x depth	36 cm x 40 cm
Opening width x depth	23 cm x 36 cm
Seat height	6 cm, 10 cm
Distance between brackets	29 cm - 38,5 cm

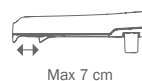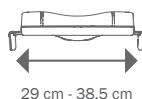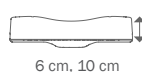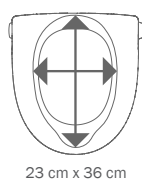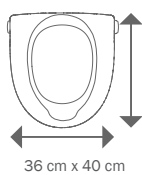



## Next level.

A new generation of toilet aids that create possibilities!

My-Loo raised toilet seat, with brackets, is part of a whole new generation of smart products from Etac. By setting the highest requirements on everything from security, comfort and flexibility to ease of installation and cleaning, we present a product range which offers entirely new and superior solutions.

With this range, we take toilet aids to a new level. 'Design drives development' has always been our guiding principle, and we are convinced that Etac My-Loo with brackets will help users and prescribers to create new possibilities.

For the latest news and continuously updated product information – visit: [www.etac.com](http://www.etac.com)

Etac R82  
Unit D4A, Coombswood Business Park East, Coombswood Way,  
Halesowen, West Midlands B62 8BH  
Tel 0121 561 2222  
Email: [enquiries@etac.uk](mailto:enquiries@etac.uk) Visit: [etac.com/uk](http://etac.com/uk)

  
Creating Possibilities