

Passive finger

Coding-7499- Silicone finger extension with rigid micro frame. 1875.00 MSRP


<u>Justification-</u> The passive finger extension by Partial Hand Solutions is designed to provide an extension for the affected finger to improve the users ability to grasp objects, preserve natural finger length and also offer a protective covering for the finger. The device consists of a custom silicone inner socket manufactured by a model of the patient's finger. A mold of the sound side finger is used to match the users size and shape to closely match the dimensions of the user. A rigid frame is applied over the silicone to create a stable and supportive device to maximize the users ability to grasp objects securely.

Casting process

Items needed.

- 1. Smooth on fast kit (Part Number) 83751 through SPS
- 2. Vaseline
- 3. Gloves
- 4. Tongue depressors
- 5. Medium size cups

Cast height for self-suspending finger. Cast both affected and non-affected side

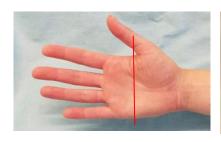





Cast should be ½ inch thick

Ship cast and work order to:

Partial Hand Solutions 86 Ladyslipper Lane Southington, CT 06489





PARTIAL HAND SOLUTIONS

86 Ladyslipper Lane Southington, CT, USA 06489 (860) 538-5532



EMERGO EUROPE
Molenstraat 15
2513 BH, The Hague
The Netherlands



Passive Finger work order

	Silicone color				Frame color		
Clinician: Contact number: Ship to address:	Otto B Black White Red Blue Green	White Red			Otto bock shade 4 Otto bock shade 9 Black White Red Blue Green Yellow		
	Affected fingers: Left Index Middle Ring Pinky Thumb						
Patient name:	Right		Middle	_	·	Thumb	

Send mold and work order to:

Partial Hand Solutions 86 Ladyslipper Lane Southington, CT 06489



Please contact Matt at Partial Hand Solutions with any questions at (860)-538-5532



PARTIAL HAND SOLUTIONS

86 Ladyslipper Lane Southington, CT, USA 06489 (860) 538-5532



EMERGO EUROPE
Molenstraat 15
2513 BH, The Hague
The Netherlands

