



SMART CAREGIVER™

Fall Prevention & Anti-Wandering



Smart Caregiver UK

info@smartcaregiver.co.uk

01455 613710

www.smartcaregiver.co.uk

Borough Court, Burrough-on-the-Hill, Leicestershire LE14 2QS

Contents Page

Corded Systems

TL-2100S—Corded Alarm System

TL-4016—Corded Silent Alarm System

TL-3100V—Dual Alarm System

Wireless Systems

TL-2100G—Wireless Alarm System

433-EC—Wireless Alarm System

433 CMU—Central Monitor Nursecall
System

Nurse Call Adapters Available

Component Products

Wireless Sensor Pads and Floor Mats

Corded Sensor Pads and Floor Mats

One-to-one Systems

TL-5102MP—Motion Sensor and Pager

TL-5102TP—Call Buttons and Pager

TL-2100S

Safety Auto-Reset Alarm Monitor



TL-2100S

Safety Auto Reset

The TL-2100S Monitor is designed for use with all Smart corded sensor pads and floor mats. Simply plug a sensor pad into the monitor, place the pad under the lower sheet in residents bed or directly on their chair. When the resident gets up and pressure is removed, the alarm will sound notifying carers that the resident is on the move. This alarm features a Safety Auto-Reset Function which will automatically reset and silence the alarm when pressed.

Optional Tamperproof Setting

When set to TamperProof mode, the monitor alert can only be reset using the Caregiver Key. Selecting this setting ensures that only the carer can silence the alarm.

Safe Always On Technology

The TL-2100S has been designed without an ON/OFF switch to ensure the monitor cannot be accidentally disarmed.

Key Features:

- Safety Auto-Reset Function: Automatically silences and resets the alarm
- Optional Setting: Caregiver Key which provides tamper proof monitoring
- Always ON Technology: No unsafe On/Off switch ensures monitor cannot be accidentally turned off
- Nurse call port, compatible with most popular nurse call systems
- Volume and tone options
- Battery operated
- AC Adapter Port

TL-4016

**Reduce Risk of Falls, Increase Carer Freedom
Use a silent monitor with wireless signal to pager**



No Alarm Noise in Room

The wireless monitor with sensor pad can be set to silent to reduce in-room alarm noise. When a resident/patient attempts to get up, a silent wireless signal is sent to the carer pager or central nurse station up to 90 meters distance (300 feet).

Flexible and Versatile

The wireless fall alarm is designed to work with Corded Bed and Chair Sensor Pads and Floor Mats. The Monitor allows carers the flexibility to monitor residents based on their unique needs.

Increase Carer Freedom

Carers can attend to other duties and receive an alert with use of a pager when someone is getting up from a bed, a chair or approaching an exit when using a floor mat.

TL-4016



Key Features:

- Monitor sends silent wireless signal to pager for noise-free carer notification
- Pager can be set on alarm or vibrate
- Call button on the front of monitor
- Nurse call port, compatible with most popular nurse call systems
- Reset button to silence the alarm
- Volume control
- Low Battery Light

TL-3100V

Recordable Dual Alarm



Unique Dual Purpose Alarm

The TL-3100V Recordable Monitor works as both a pull-string monitor or sensor pad monitor, easily transferring from bed to chair and back, eliminating the need for multiple monitoring devices. The alarm alerts when the magnetic pull-string is pulled free of the monitor or when pressure is removed from the weight sensing pad. The monitor automatically resets when the magnetic pull cord is replaced or when pressure is re-applied to the connected sensor pad, or when the reset button is pressed. Pull string is attached to residents clothing and as they lean forward to stand the string pulls the magnet from the alarm which then causes it to sound.

Personalised Recordings

Eliminate startling alarm noise by recording a personalised message to be played back when the alarm is activated. With up to 15 seconds of playback time available, record a melody or familiar voice to be heard as the alert sound, such as "Please sit down Mrs Jones, help in on the way".

Safe Always On Technology

Designed without an ON/OFF switch, the alarm is always on and cannot be turned off. Carers use the easily accessible Alarm Reset Button to silence the alarm and attend to the resident's/patients needs without the potential for accidentally leaving the alarm off.

TL-3100V



Key Features:

- Dual Functions
- Recordable Voice Option
- Magnet alarm functions
- Pad Lost, Low Battery and Monitor Status Indicator Lights, Volume Control
- Nurse call port, compatible with most popular nurse call systems
- Monitoring Options:
 - Corded Bed Pad
 - Corded Chair Pad
 - Corded Floor Mat
 - Pull String Alarm

TL-2100G

Wireless Fall Monitor Without Cords

Works with wireless bed & chair pads & floor mats



Eliminates In-Room Noise

The TL-2100G Monitor is cord free and can be mounted where you need it—away from residents - eliminating in-room alarm noise. The monitor can be wall mounted, at the bedside, in a hallway or at a nurses station up to 50 feet away (15 meters) using the wall mounting hardware provided. The TL-2100G Monitor is compatible with wireless bed and chair sensor pads, wireless motion sensors and floor mats.

Patient Safety Features

The monitor automatically lets a carer know if the signal between pad and monitor is lost or if the monitor is out of range. In addition, this monitor features an optional Tamper-proof setting which when enabled means the alarm can only be turned off by the carer using the Caregiver Key. This ensures the resident cannot accidentally turn off the monitor.

No Cords

No cords between pads and monitors, reducing tripping hazards and eliminating broken or tangled cords.

TL-2100G

Key Features:

- No alarm in the resident's/patients room
- Safety Auto-Reset Function: Automatically silences and resets the alarm
- Optional Setting: Carer with Caregiver Key can reset the monitor with the Caregiver Key only which tamper proofs the monitor
- Always ON Technology: No unsafe On/Off switch ensures monitor cannot be accidentally turned off
- Nurse call port, compatible with most popular nurse call systems
- Volume control
- Low battery alert

**Can be ordered pre-programmed to components or be programmed with instructions provided*

Wireless Paging Alarm—433-EC



**Wireless,
Rapid
Response
Alarm
Range—Up to 90
meters!**



433-EC



The Wireless Exit Paging Alarm provides affordable and effective cord free wireless monitoring with a range of up to 90 meters (300 feet), with immediate response no delay.

Easy to Use

The Alarm/Pager is simple to use, durable and light weight. **It can be paired with up to 6 wireless components** (bed & chair sensor pads, motion sensors, floor mats or call buttons) to serve as a portable carer alert.

Remove Noise from the Room

This alarm can be mounted away from the patient/resident removing in-room alarm noise and with wireless, cordfree sensor pads, you reduce entanglements and tripping hazards. And the adjustable volume and gentle chime will not startle residents and offers quieter fall prevention.


Key Features:

- Easy set up and operation
- Nurse call port, compatible with most popular nurse call systems
- Works with all wireless sensor pads, floor mats, nurse call buttons and motion sensors
- Monitor up to 6 components with one alarm
- Adjustable low, medium and high volume adjustment and gentle chime setting
- Low battery indicator
- Battery operated or optional AC power adapter

**Can be ordered pre-programmed to components or be programmed with instructions provided*

433 CMU—Central Monitor Unit

Wireless Small Facility Central Monitoring Unit
40-60 Component Capability

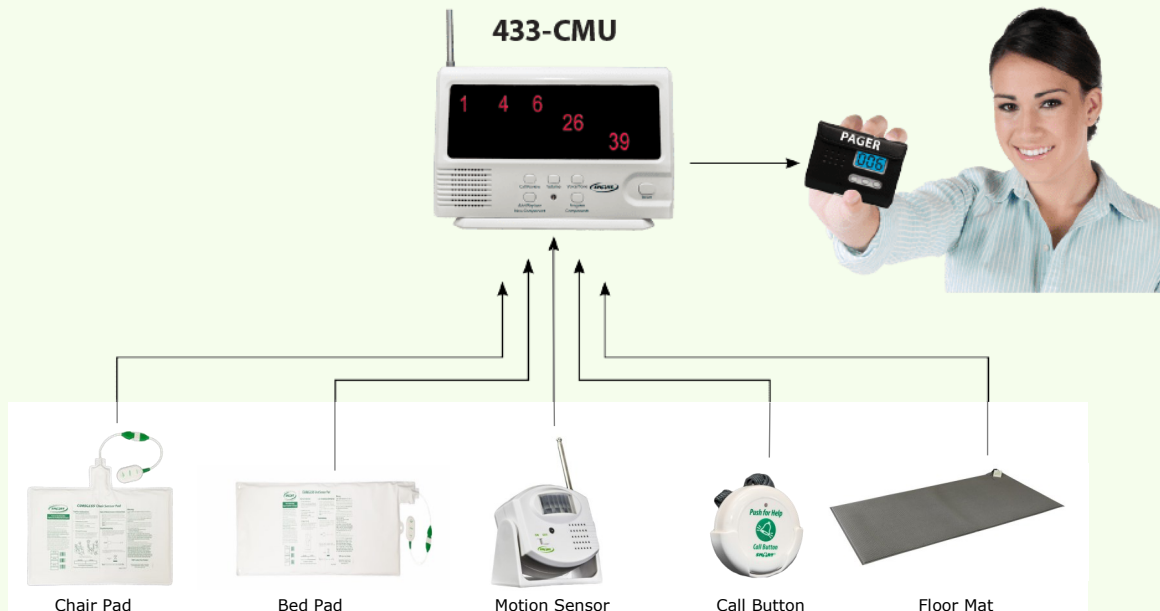

Monitor multiple residents from one central location

Eliminate In-Room Alarm Noise

The 433-CMU Central Monitoring Unit allows you to monitor multiple wireless sensor units from one central location. This eliminates bothersome in-room alarm noise and significantly reduces the tripping hazards associated with corded units.

Real Time Pager Notification

The CMU also communicates with any number of 433-Pager Units. When activated, the pager audibly alerts and the LCD screen displays which sensor unit activated the alert by number, at the same time calling out the same number. Gives carers more flexibility and freedom of movement while providing the residents/patients a quieter environment. Pager features volume, led alert, on/off switch and vibrate options with reset button to reset both CMU and Pager.

Key Features:

- Alarms audibly and visually at central location with a 300 ft range and portable caregiver pager (s) which feature an additional 150 ft range
- **40 and 60 Channel Capacity** Models - Compatible with all wireless components
- Systems may be pre-programmed prior to delivery or re-programmed using the instructions provided
- Battery operated /AC power adapter included
- Low battery indicator
- Can add additional bumper to boost range by 100 feet.

Nurse Call Systems

Nurse call systems that our products work with

We supply nurse call adapters for the following nurse call systems



Chubb

medicare
SYSTEMS



Quantec
Addressable Call System

Have a different system? Have any questions? Call our highly trained customer service team who will be able to assist you.

Cord Free Sensor Pads And Floor Mats

Wireless/Cordless Pads & Floor Mats for Wireless Exit Alarms Sizes

PTB-RI



PTB-RI—Cordless Bed Pad
10" x 30" (1 year pad)

PTB-WI



PTB-WI—Cordless Bed Pad
20" x 30" (1 year pad)

PTC-WI



PTC-WI—Cordless Chair Pad
10" x 15" (1 year pad)

PTFM-05C



PTFM-05C—Cordless Floor Mat
24" x 36" (1 year warranty)

PTFM-07C



PTFM-07C— Cordless Floor Mat
24" x 48" (1 year warranty)

Corded Pads

Corded Pads & Floor Mats for Corded Alarm System Sizes

BPB-RI



BPB-RI—Corded Bed Pad
10" x 30" (1 year pad)

PPB-RI



PPB-RI - Corded Bed Pad
10" x 30" (1 year Pad)
(With SafeTRelease Cord)

PPB-WI



**PPB-WI—Corded Bed Pad
20" x 30" (1 year pad)
(With SafeTRelease Cord)

PPC-WI



PPC-WI—Corded Chair Pad
10" x 15" (1 year pad)

FM-05



**FM-05—Corded Floor Mat
24" x 36" (1 year warranty)
(With SafeTRelease Cord)

FM-07



**FM-07— Corded Floor Mat
24" x 48" (1 year warranty)
(With SafeTRelease Cord)

**

- Unique SafeTRelease™ Cord
- No More replacing sensor pads due to broken cords (cords can be replaced if broken)
- Unique, SafeTRelease™ Cord separates instead of breaking or pulling cord away from pad
Reduces tripping, entanglement and hazards

One-to-One Motion Sensor & Paging System

Motion Activated Caregiver Paging System



TL-5102MP

Easy To Use Carer Alert

No Noise In Room

The Smart Caregiver TL-5102MP is a stand-alone motion paging system that allows you to eliminate in-room alarm noise. The Motion Sensor is placed in a doorway or by the bedside and when motion is detected, a silent wireless signal is sent to the carer pager up to 90 meters (300ft) away.

Carer Alerted By Pager

The Caregiver Pager is pocket-sized and light weight, giving carers the freedom they need.

Key Features

- Quiet non-alarming motion sensor in room
- Easy to set up
- Pocket sized pager - place in pocket, clip on belt or place on wall or table top
- Motion Sensor alarm can be placed by bed or doorway to monitor residents' activities
- Battery operated (batteries included)
 - Pager 2 x AA batteries
 - Motion Sensor 3 x AA batteries

One-to-One Call Button System

Lowest Cost Call Button Paging System



Easy To Use Two Call Button System

Increased Carer Freedom

The TL-5102TP is a stand-alone personal paging system and includes two Call Buttons and one Caregiver Pager. When pressed, the nurse call buttons send a wireless signal to the pager alerting the carer that assistance is required. This simple, streamlined system allows you to eliminate in-room alarm noise and with a range of up to 300 feet, carers have the freedom they need.

Easy Push-Button Operation

The call buttons can be worn around the neck using a lanyard or wall-mounted using the plastic cradle provided for easy access.

Versatile Pocket-Sized Pager

The pager has two audible sounds to choose from and is small enough to be carried in the pocket or clipped to a belt



TL-5102TP


Key Features

- Quiet no alarming call buttons
- Easy to set up - insert batteries and it is ready to use
- Pocket sized pager - easy to carry
- Battery operated (batteries included)

Notes



Contact us

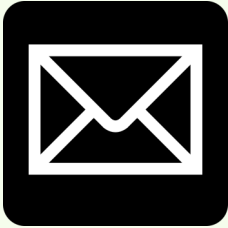

Email: info@smartcaregiver.co.uk



Phone: 01455 613710



Registered Address:

Burrough Court
Burrough-on-the-Hill
Leicestershire LE14 2QS

