



Mobility[®] in Motion

Vehicle Adaptations Enabling Independence



Life in Motion

Remarkable adaptations for everyday life.
Innovative products making car travel easy and
accessible to all, tailored to your individual needs.

autochair



BRIG-AYD CONTROLS

BraunAbility

The open road awaits. Where will it take you?



We're here to make sure you can enjoy every journey.

Whether you're reconnecting with old friends, starting adventures with new ones or riding solo, our loading and stowage adaptations give you the freedom to roam.




***"We were very pleased with
Mobility in Motion from the start.***

***The initial information, the demonstration
and the fitting have convinced us that
this was the right choice.***

***A game changer for me as getting my
small fold up scooter was getting difficult
for my husband and impossible for me!"***

5 Star Trustpilot Review

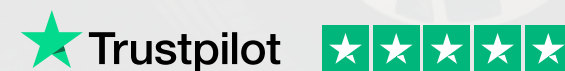

Take your mobility device further afield

Truly unlock your freedom and independence by loading your wheelchair, mobility scooter or powerchair into your car – all at the touch of a button.



autochair[®]

Autochair Smart Lifter LC
4-Way Compact Boot Hoist



Paralympian Robyn Love with the Mobility in Motion demonstrator car displaying a Boot Hoist and Wheelchair Roof Box

Loading & Stowage

4-Way Boot Hoists

4-Way Boot Hoists lift and swivel, to raise your scooter, wheelchair or powerchair up and power round to load it into your vehicle.

Fixed/folding/compact and XL models are available with lifting capacities from 40kg all the way to 200kg. We have one of the largest ranges available, which means the combination of car and scooter, or wheelchair are almost endless. We even have the Smart Lifter LC, the No.1 best selling boot hoist in the UK.

Spotlight on the Autochair Smart Lifter LC

- The LC is the No.1 Best Selling Hoist in the UK
- Capacities up to 100kg for a large range of devices
- Can lift your scooter or wheelchair into the car in less than 60 seconds
- Suitable for over 450 vehicles



2- Way Boot Hoists

The 2-Way powered lift and lower function can raise your scooter, wheelchair or powerchair up to your vehicle's boot height so you can easily swing it in.

It's the perfect adaptation for those who want a little extra help lifting their mobility device. Multiple models are available including our popular Brig-Ayd 40kg Mk2, a folding model that lays down flat in your boot freeing up space.

Spotlight on the Brig-Ayd 40kg Mk2

- Greater lifting height to fit taller devices into smaller cars
- Folds down to the boot floor so it doesn't obstruct the rear view when driving
- Ideal for smaller cars with devices up to 40kg
- Compatible with more vehicles and mobility devices



Boot Hoist Battery Pack

Independently powers the Autochair Smart Lifter range of boot hoists, eliminating the need to wire into the vehicle's electrical system - perfect for EVs.

Load and unload your mobility device up to 30% faster by speed boosting your boot hoist with a power pack. Enjoy long-lasting power & quick charging in around 2 hours. The ergonomic design makes it comfortable to hold and easy to carry. The secure mounting bracket keeps it in place on your hoist so it doesn't take up valuable boot space.

Spotlight on the Autochair Smart Lifter Power Pack

- Speed boost up to 30% faster operation
- Push-to-indicate charge level display
- Operate your hoist even if your car battery is low
- Suitable for all vehicle types - electric, hybrid, petrol, & diesel
- Award winning - Best New E-Mobility Technology 2025





Wheelchair Roof Boxes

It only takes 30 seconds to load or retrieve your folding wheelchair from the roof of your car, placing it conveniently next to your seat for an easier transfer.

Suitable for virtually all vehicles that can accept roof bars with a capacity of 75kg or more. They can be fitted for the driver or front passenger's use. Roof stowing your folding wheelchair, weighing up to 20kg, retains boot space and full seat capacity, keeping your car clean and unobstructed.

Spotlight on the BraunAbility Chair Topper

- Ideal stowage solution for independent drivers
- Retains your boot and cabin space
- Weatherproof box protects your wheelchair from the elements
- Compatible with folding wheelchairs up to 20kg
- Suitable for driver or front passenger



BraunAbility Chair Topper providing a manual folding wheelchair to the driver (sitting on a Turny Evo Swivel Seat)



Platform Lifts

Simply drive your mobility device onto the platform, step off and use the remote control to safely lift your device into the boot of your car.

A programmed platform lift allows for easy independent loading and stowing of a mobility device. With a single button it lifts and stows an unoccupied scooter or powered wheelchair inside a vehicle. The tie-down restraint system ensures your device remains secure throughout your journey.

Spotlight on The Bruno Joey Lift

- No requirement for you to physically guide or support the device while loading
- Straightforward 3-step process
- Magnetic, easy-grip, illuminated handset
- Crash tested barrier system
- Great for larger devices up to 159kg



Scooter Seats

Add functionality and practicality to your existing scooter by upgrading the seat for comfort and for ease of stowage and transport.

The Autochair Smart Scooter Seat makes your scooter easier to transport, load into smaller vehicles, and store due to its lightweight, versatile, and easily removable components. Designed to be assembled and disassembled with one hand, no part weighs more than 3kg. The Smart Scooter Seat is compatible with all Autochair Smart Lifter Boot Hoists

Spotlight on the Autochair Smart Scooter Seat

- No single part weighs more than 3kg
- Seat Rotates 360 degrees and locks every 90 degrees
- Comfortable dual-density foam padded seat and backrest for longer journeys
- Easy to remove and reassemble with just one hand
- Suitable for users up to 120kg



To the ends of the earth or just to the park.



Whether you're travelling across the country, or just down the road, Mobility in Motion has you covered.

Our range of accessible driving aids can give you confidence on the road, no matter where you want to go.





"Now that I've got the adaptations I'm able to travel longer distances by myself which is so exciting. -

I'm on cloud nine really and I can have my Independence back."

Excerpt from:

Driving with EDS Judith's Journey

Watch the full interview on
youtube.com/@mobilityinmotion





Vehicle adaptations that put you
back in the driving seat of your life.



SureDrive Next Generation Hand Controls
Mechanical Push/Pull Accelerator and Brake

Driving Aids

Paralympian Robyn Love with the new SureDrive Mechanical Hand Controls and Brig-Ayd Quick Release Steering Ball



Mechanical Hand Controls

Hand Controls move the function of the pedals to a convenient push/pull handle, behind the steering wheel, comfortably positioned for your right hand.

Enables users with partial or limited use of their legs to drive by transferring acceleration & braking to their right hand. Multiple handle styles are available to suit you. They can even include an indicator switch or secondary controls with buttons for the lights, horn and more.

Spotlight on the new SureDrive Push/Pull Hand Controls

- Adjustable handle for a comfortable driving position
- Smooth and light operation
- Passed rigorous durability, vibration & overload testing - EU CEN compliant
- Versions available with indicator control



Electronic Hand Controls

Electronic Accelerators reduce the fatigue of driving compared to a Mechanical Pull Accelerator.

They greatly reduce the force required to operate the throttle compared to a mechanical column or floor-mounted hand control. Specially designed electronic accelerators such as an Over Ring can be used independently or paired with a separate Mechanical Push Brake to suit your control requirements.

Spotlight on the BAS Over Ring Accelerator

- Smooth, light and gradual acceleration
- Moves with the steering wheel
- Park button to reduce the power of the accelerator to aid parking and town driving
- Easily removable by lifting a single pin



Floor-Mounted Hand Controls

Transfers acceleration and braking to your left hand, placing the controls in an ergonomic position beside the centre console.

A floor-mounted design leaves more leg room for the driver and is ideal for those who prefer to steer with their right hand. Like column-mounted controls, models are available with mechanical or electronic acceleration, as well as secondary controls for indicators, lights, wipers and more.

Spotlight on the BraunAbility Carospeed Classic

- Hill hold function allows you to take your hand away while stationary
- Allows your arm to rest in a comfortable position when operating
- Doesn't obstruct your knees when entering/exiting the vehicle, even via a swivel seat
- Retains full adjustability of the steering column





Steering Aids

Steering Aids allow you to continuously turn the wheel, making them ideal for steering with one hand.

Whilst they can be used independently, they are essential when hand controls are fitted to a car. Our specialist Steering Aids, like the Tetra Grip provide an alternative grip position and additional wrist support for drivers with reduced grip strength.

Spotlight on the Brig-Ayd Quick Release Steering Ball

- Quick Release Design ideal for shared vehicles
- Patented one bolt design enables it to fit closer to the spokes of the wheel
- Can be fitted wherever suits you on the steering wheel for a comfortable resting position
- Ball shape allows for a comfortable but firm grip



Remote Control Devices

Remote Controls provide drivers with access to secondary functions such as the indicators, horn, wipers and lights, in a more convenient location.

They can be integrated into a Steering Aid like a Steering Ball or a Lollipop Grip, allowing you to steer and operate secondary functions with one hand. Alternatively our Hand Control Mounted Remote Control places the secondary functions with your accelerating and braking hand.

Spotlight on Lodgeson's R200 Keypad Controls

- Various keypad configurations are available to suit you and your vehicle
- Illuminated buttons for safe operation at night
- A traditional Steering Ball shape allows for a comfortable but firm grip
- The Lollipop shape allows for an ergonomic alternative wrist position and a straight grip



Brand New fully electric VW ID. Buzz with a Brig-Ayd Quick Release Steering Ball & Left Foot Accelerator



Left Foot Accelerators

Transfers acceleration to a new left foot pedal for drivers who are left-leg dominant.

Designed for drivers unable to operate the accelerator with their right foot, an additional pedal is fitted to the left of the brake. The use of the right pedal is retained for shared vehicles. It either flips-up safely out of the way or is protected by a pedal guard, so you can easily switch between left or right foot acceleration.

Spotlight on the Brig-Ayd Twin Flip Left Foot Pedal Transfer

- Easily switch between left or right foot acceleration
- Adjustable pedal position and angle to suit the individual driver and their vehicle
- Retains the use of the right pedal for right-foot drivers
- The unused pedal stows safely out of the way allowing you to extend that leg






**Unlock a world of possibilities
and explore new destinations.**



**Mobility
in Motion**

Adapting to life's journeys

We're here for wherever you want to go with our range of accessible vehicle adaptations. So you can have the freedom to explore the world fearlessly.

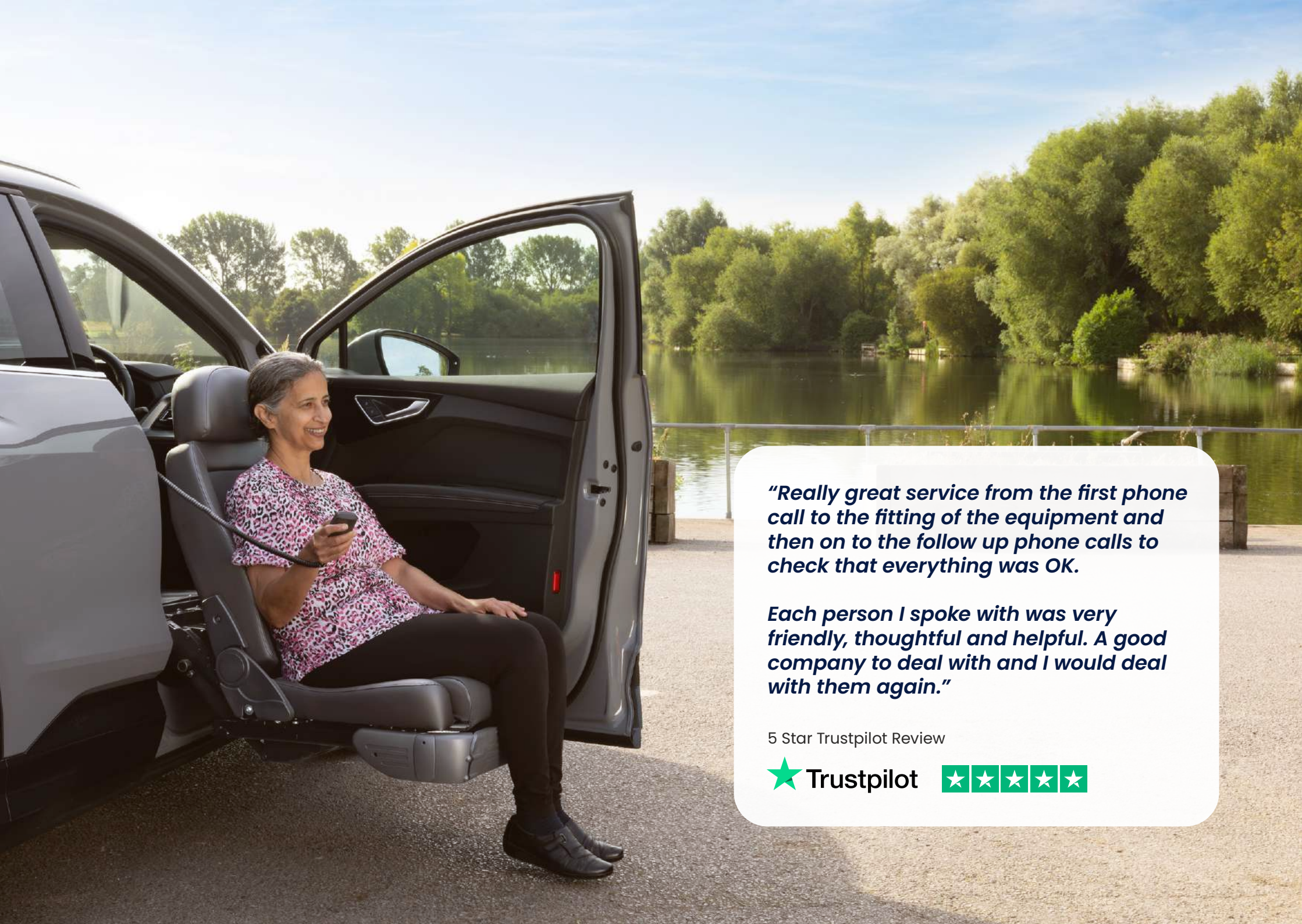

"Really great service from the first phone call to the fitting of the equipment and then on to the follow up phone calls to check that everything was OK.

Each person I spoke with was very friendly, thoughtful and helpful. A good company to deal with and I would deal with them again."

5 Star Trustpilot Review



Trustpilot



Maximise your mobility

Drivers and passengers with reduced mobility can be seated with ease by making use of our easy vehicle access solutions.



Turny® Manual + Carony®
Manual Swivel Seat & 16" Wheelchair System

Easy Vehicle Access

BraunAbility Turny Orbit Swivel Seat with additional seat accessories and postural support for a child occupant



Powered & Lowering Swivel Seats

Electronic Swivel Seats remove the physical requirement to manually rotate the seat.

The Turny Low Vehicle, provides a powered solution to rotate and extend the seat. The Turny Orbit even lowers the seat down, making it perfect for taller cars and SUVs. The Turny Evo further improves vehicle access with its programmable movement path, allowing it to be individually tailored for each user/vehicle, making a driver's side installation possible on many cars.

Spotlight on the BraunAbility Turny Evo

- Most advanced seat lift on the market today
- Fully programmable travel path to ensure maximum head & legroom
- User friendly handset with status indicator and easy to use soft touch buttons
- Ideal for a driver or passenger, including wheelchair users



Manual Swivel Seats

Swivel Seats make access to and from your car safer and easier by rotating the seat to face out of the door and extending out on sliding rails.

BraunAbility's Turny Manual swivel base is designed for those who find it difficult to enter a vehicle or make a transfer from their wheelchair. A manual lever locks the seat in place and allows for simple control of the turning mechanism. The seat can still slide forward and back in the car for comfort.



Seating Options and Accessories

Swivel bases come with a choice of seats for adults and children with various postural needs. You can even upgrade to a Wheelchair Swivel Seat.

There are also a range of seat accessories available for additional postural support depending on individual requirements. This includes; armrests, head support, backrest cushion, lumbar support, secure belts and harnesses. This will be assessed and discussed during our free home assessments.

Spotlight on the BraunAbility Turny Manual

- Manually rotate the entire seat to aid vehicle access
- Ideal for a driver or passenger, including wheelchair users
- Provides greater support than a Transfer Plate
- Seat locks in position for safe transfers
- Doesn't permanently modify your vehicle



Compact Seat
(Standard)



Caroseat Flex
(Flex)



BEV Seat
(Premium)



Kids BEV Seat
(Kids)



Options and
Accessories



Wheelchair Swivel Seat

The Carony wheelchair base removes any need for a caregiver to physically lift the user during their transfer.

This lets them stay seated in the wheelchair when moving into the car. By docking onto a swivel base the entire wheelchair seat, with the seated user, can slide into the car and swivel in where it becomes the car seat.



Spotlight on the BraunAbility Carony System

- Transfers both the occupant and the seat together
- Compatible with a range of swivel seat bases to suit your needs and your vehicle
- Options and accessories are available to tailor the seat and support your posture
- Fully crash-tested design



Transfer Plate

Make transferring easier by bridging the gap between your mobility device and the car.

Transfer Plates create a static, temporary platform on the outer edge of the seat on which to pause and pivot. Autochair's quick release Transfer Plate and BraunAbility's folding Transfer Board offer a solution for drivers and passengers in virtually any vehicle.

Spotlight on the Autochair Transfer Plate ATP

- Provides a secure platform to transfer into the vehicle
- Quickly detaches allowing for a swift transfer in and out of the vehicle
- Can be installed for passenger or driver's side use
- Reduces the transfer distance between a mobility device and the car seat



The Autochair Milford Person Lift provides a supported transfer, for adults or children, to the front or rear seats



Person Hoist

Our Milford Person Lift removes the requirement for a caregiver to physically lift the user in or out of the vehicle during their transfer.

Ideal for wheelchair users looking to get seated safely and comfortably in the original seat of a vehicle as either the driver or passenger. A great alternative to a WAV (Wheelchair Accessible Vehicle), our Person Hoist lifts you via a comfortable fitted sling that can safely accommodate children and adults alike.

Spotlight on the Autochair Milford Person Lift

- Safely lifts you via a comfortable sling
- Transfer a user in and out of the vehicle without physically lifting them
- Range of different slings to suit individual needs
- Largest weight capacity of any vehicle-based person hoist on the market, suitable for users up to 150kg



Give us a call to start your journey

Speak to our mobility experts

Call, email, or webchat with our friendly and professional UK based team, to find out how we can support your independence.

0800 288 4422

hello@mobilityinmotion.com



Your free assessment online or in person

We will guide you through a thorough online demonstration of the adaptation and answer all your questions, to find you the perfect solution.

For our more complex adaptations like our Person Hoist, one of our specialist Field Product Demonstrators will come out to you for a Free Home Assessment.

Nationwide installation

Our network of expert installers will ensure your adaptations are fitted quickly, and to an excellent standard, at your home or at one of our local install workshops. Some of our simplest adaptations can even be posted to you directly to fit yourself if you prefer.



Mobility in Motion - Information correct as of July 2025 - Product range subject to change - General



**Mobility
in Motion**



Trustpilot



Excellent