



## Regenerative Medical Technologies


#### **Speaker**

#### **Ashton Choi**

Cultivated Meat Investor / Tech Columnist
Researcher, Oxford Centre for Tissue Engineering and
Bioprocessing





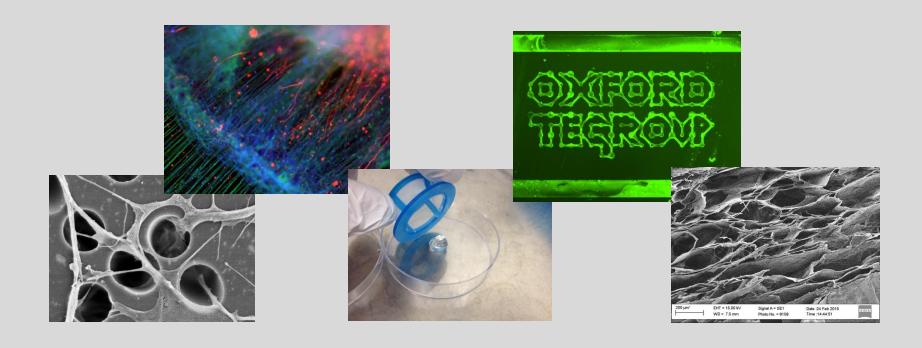

# The Oxford Centre for Tissue Engineering and Bioprocessing





# Main areas of research






- Bioreactor design
  - Large scale cell culture
  - Perfused bioreactor for cell/tissue culture
  - Hollow fibre membrane based bioreactor for generating bulky tissues
- Biomaterial selection, modification and processing
  - Cell culture substrate modification
  - Scaffolds for 3D cell culture
- 3D cell culture and characterisation
  - 3D human neural network
  - In vitro cancer model
  - Co-culture








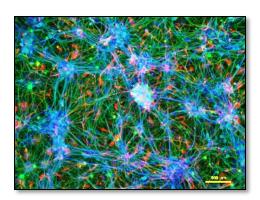

# Manufacturing process flow for cell based allogeneic regenerative medicine therapy

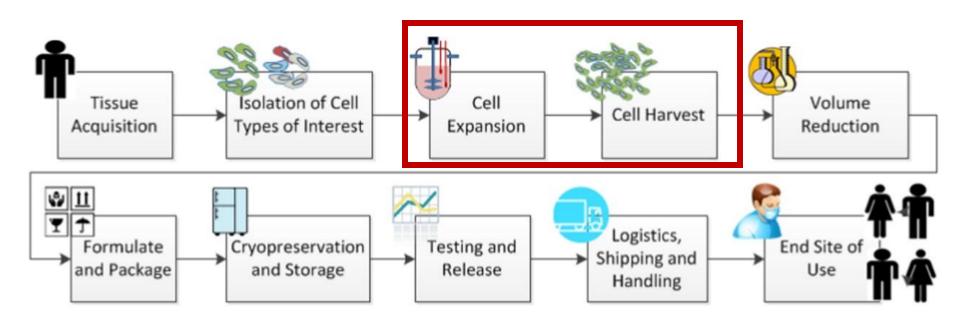

- Cell expansion step: the most expensive and time consuming step
- Cell harvesting step: not yet the optimum quality of cell harvesting

Biofabrication 7 (2015) 044103



# Bioreactor: Fluidized Bed




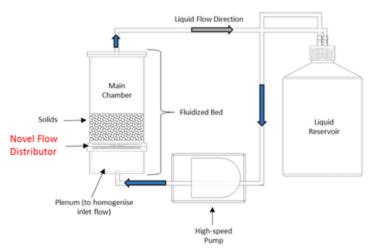


Figure: Liquid-solid fluidized bed experimental set-up



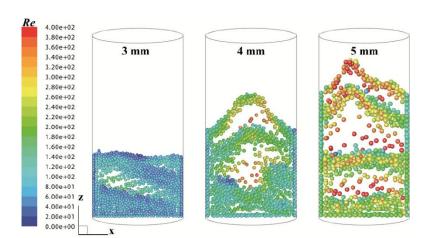

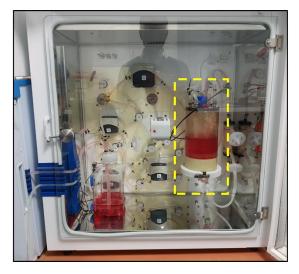


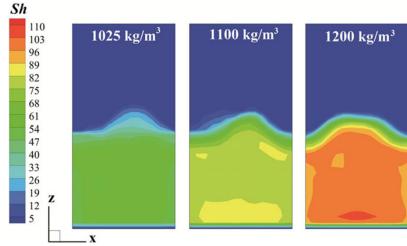

# Particle Design: Thermoresponsive Macrocarriers

A thermo-responsive polymer is able to change its properties in response to changes in temperature.


#### Poly(N-isopropyl acrylamide) – PNIPAAm or PIPAAm

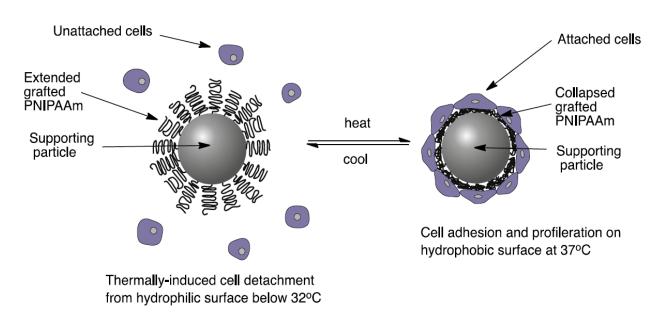



European Polymer Journal 67 (2015) 346–364



# An Automated, High-throughput, Hollow Fibre Bioreactor System for MSC-sEV Production

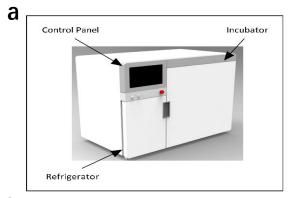

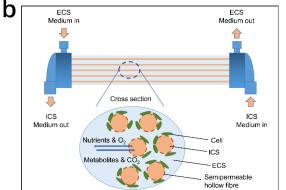

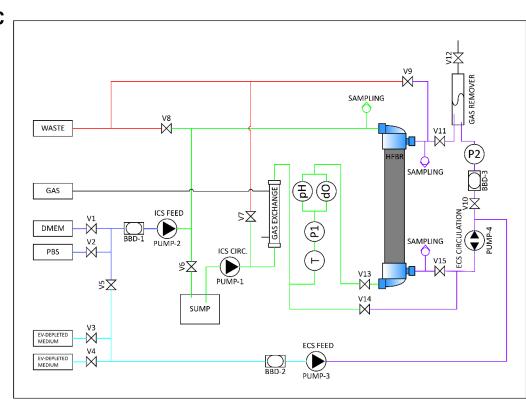

Cell Expansion System. (a) The outside appearance of the CES. (b) Schematic diagram and working principle of the HFBR cartridge. (c) Process flow diagram of the CES. BBD – bubble detector, V – valve, T – temperature sensor, P – pressure sensor, pH – pH sensor, dO – oxygen sensor.



# What if we could grow meat iust like a plant?




Lunderstand






# Ivy Farm Technologies: a spin-off company established in 2019

#### Bioreactor schematic

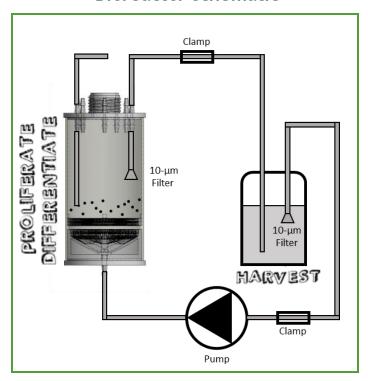


#### Exclusive license of patented thermo-responsive macrocarrier technology

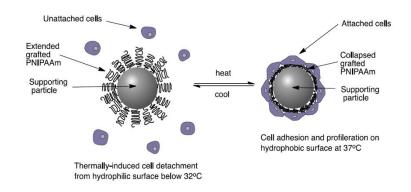

- Lower cost vs. others (fewer steps in the process, fewer additives)
- Animal-free (alternative is trypsin, an animal-derived enzyme)
- No cytotoxicity
- No adverse effects on cell morphology/physiology
- Cell phenotype remains unaltered
- Predicted regulatory-friendly methodology

#### Commercial access to advanced bioreactor design

- Angled inlets to increase axial and vertical velocity of mixture
- Encourages mixing leading to better cell quality and higher growth rate
- Defined approach to automation established

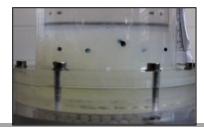





sustainable meat co.

Ivy Farm Operates the Largest Culture Meat Pilot Plant in Europe in its 18,000 Sqft Campus in Oxford, UK



Our R&D center in Oxford has a low-rent space that is capable of hosting over 120 employees, while currently we have only 50. We are thus ready to expand for the upcoming stage of development



I cannot tell the difference between traditional farmed beef and Ivy Farm cultivated beef... Its absolutely delicious" – Tom Heap

January, 29, 2024



summit

"We combine just the best bits of meat pure muscle and pure fat - to create a
product that is better for the environment." Rich Dillon

November 4, 2022

"Ivy Farm is developing a healthier alternative to traditional and existing plant-based meat products, which is a key value proposition." – Hiral Patel

January 26, 2022



"It's really good, and it's a real hot dog! Wow."- Poppy O'Toole


November 25, 2021



# First cultured burger

In 2013, this beef burger was made at the cost of eye-watering \$300k



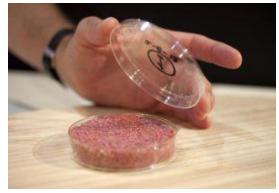



# More companies are doing it










# At much lower cost

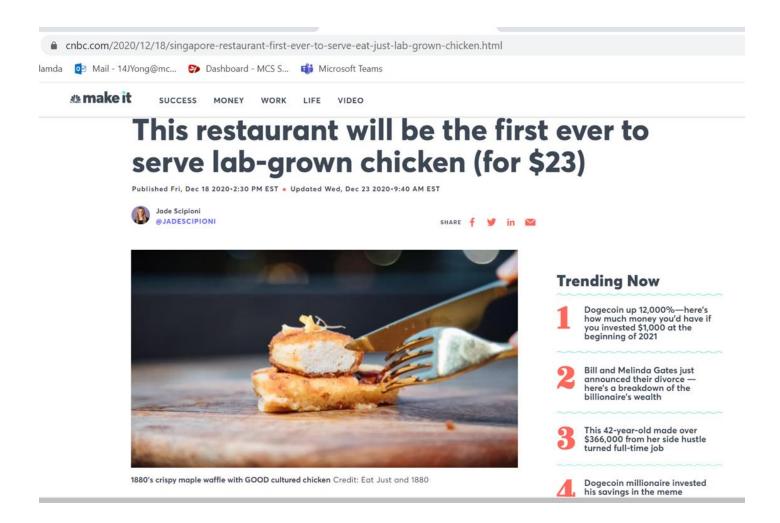



# Scale up challenge


- There are ~100 billion cells in a kg of meat
- World meat consumption: 300 billion kg/year
- If cultured meat is to supply only 10% of the world meat consumption, i.e. 30 billion kg/year, we need to produce 30 billion  $\times$  100 billion cells (3  $\times$ 10<sup>21</sup> cells /year)
- Even with high density cell culture in bioreactors, we would need at least 2  $\times$  10<sup>6</sup> m<sup>3</sup> bioreactor volume (corresponding to ~200,000 of 10,000 litre bioreactors).
- Currently, there aren't even enough amount of stainless steel materials
- Innovation is needed:
  - Bioreactor design
  - Reduce waste
  - Regulation








# Building a brain slice in a dish

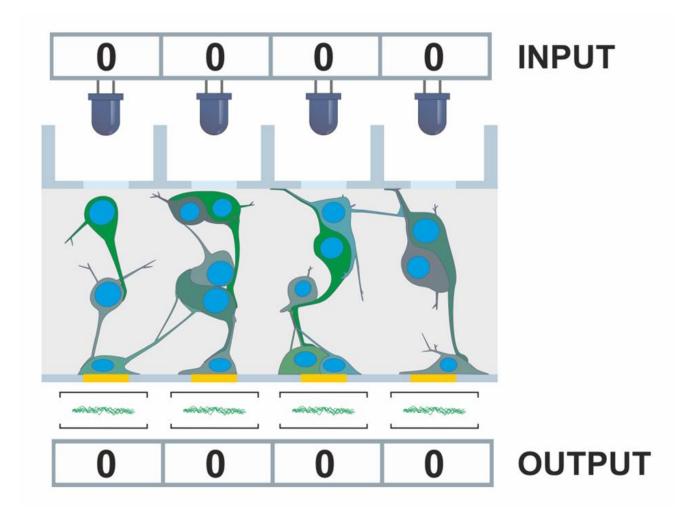

# The physical model

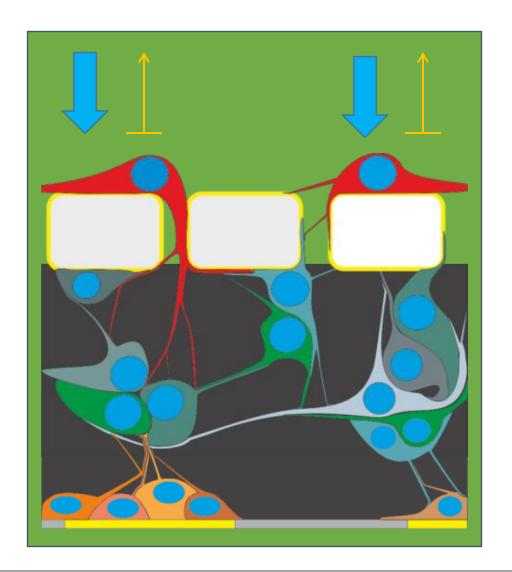

Stimuli - light


Light sensitive neurons

Porous membrane

3D Model neural network


Neurons patterned on MEA

### **Neurons like to form clusters**

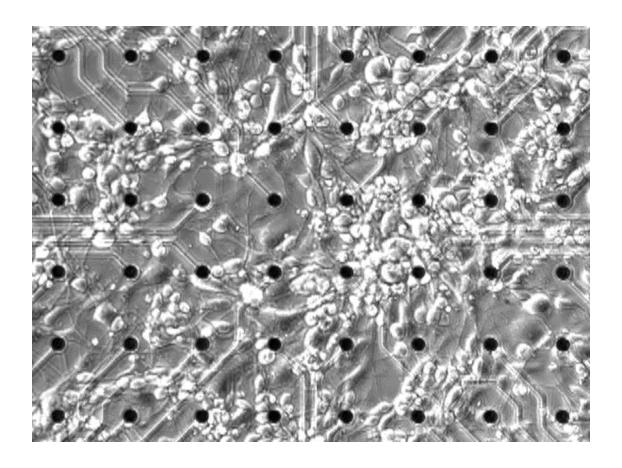

Seeding of human IPSC derived neurons on a Laminin coated MEA (720x, 4hr time-lapse)

## Neurons form networks throughout gel

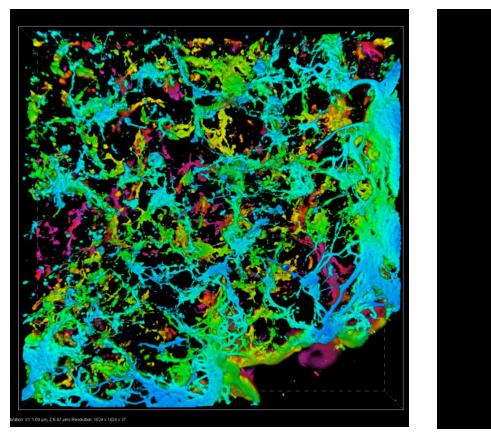



4-week culture of human iPSC derived neurons in granular HyStem hydrogel + Laminin

