
INORGANIC VENTURES



Tips for Trace Analysis by ICP-OES & ICP-MS

Mike Booth
Innovation & Applications Development Manager

What does Inorganic Ventures do?

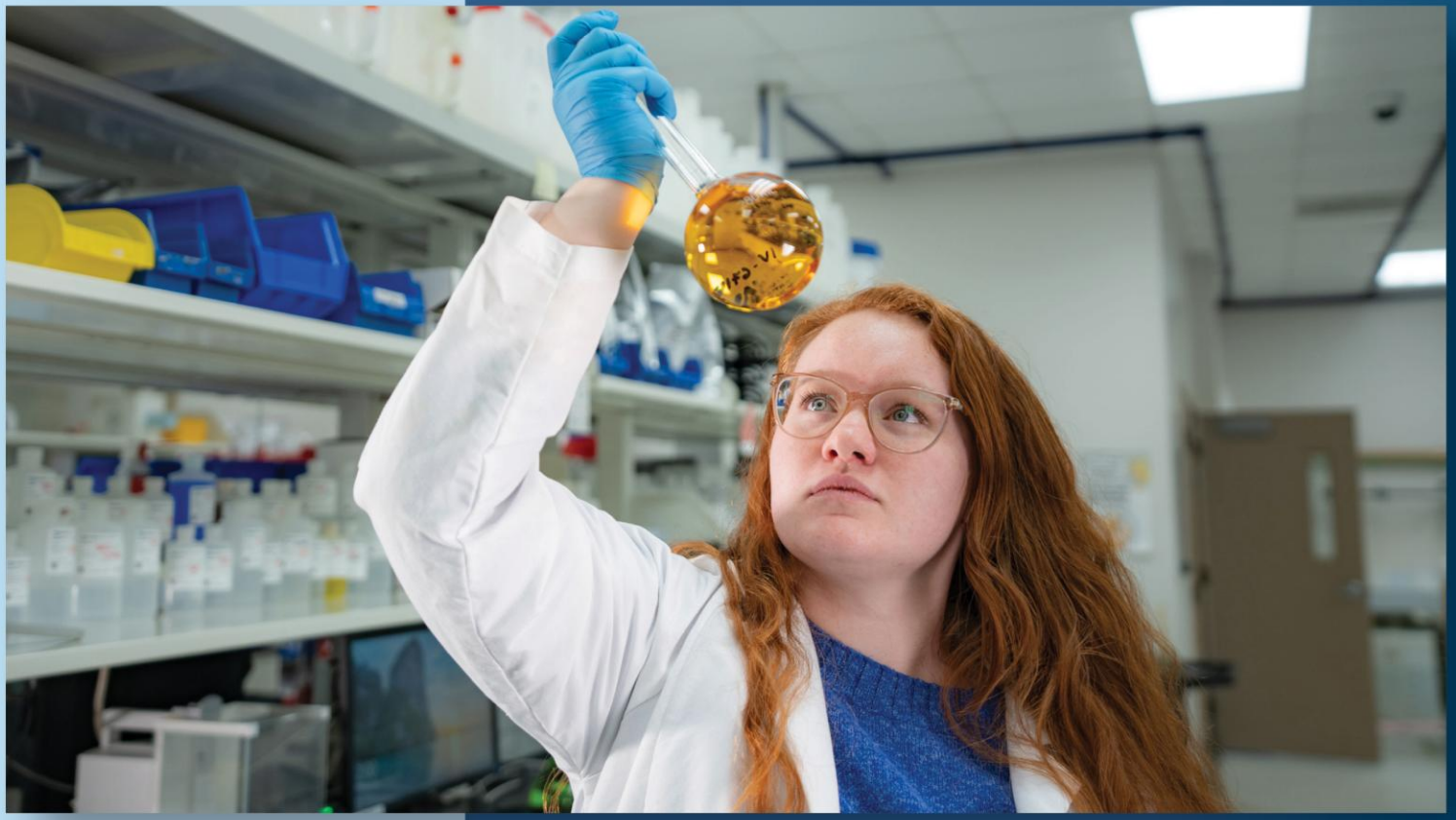

Inorganic Ventures is a manufacturer of inorganic Certified Reference Materials (CRMs) commonly referred to as a “standards”, “chemical standards” or “calibration standards”.

- Elemental Standards
  - High Purity Single Element
  - Custom Blends
- pH Standards
- Conductivity Standards
- Titrants
- Reagents



2 Day Turn-Around-Time on Stock Products!
10 Day Turn-Around-Tine on Custom Products!

01

Top Topics for Today:

- **Environmental Contamination**
  It is critical to be as contaminant free as possible for accurate trace analysis work.
- Solution Matrix
- Washout Issues
- Uncertainty of Measurements

Container Selection

- **LDPE**, HDPE, PTFE, PFA, PP, Borosilicate Glass
  - Plastic made from virgin polyethylene is critical
- Sample tubes vs. **bottles**
- **Leach containers with dilute HNO_3 (1-5% v/v)**
- **High temperature leaching is more effective (50°C)**



PP?



PP



HDPE



PTFE



Boro Glass

Container Selection

More In-Depth

Values are the ng of impurity per 125mL bottle.

*** We've seen it on occasion.**



LDPE

Element	LDPE	HDPE	PP	PTFE	Boro
Al	2	*	87		8750
Ba					5000
B					9400
Ca		28	*	25	2800
Fe		150		95	125
Mg	6	11	575	15	125
K				16	3600
Na	42	63	42	90	27500
Zn	*				125

Leaching Observations

- Un-leached bottles do contain some metallic impurities
- After high temperature leaching, impurities are significantly reduced.
- Allowing hot solution to mix within the bottles significantly reduced the amount of impurity left behind
- Chemical dishwashers deposited significant amounts of metallic impurities onto the bottle walls.
  - Chemical dishwashers contain stainless steel parts.

How Do We Leach Bottles?



- Pre-make 1-5% HNO_3
- Does not need to be highest purity acid



- Stand the bottles upright using automation



- Fill the bottles completely and cap them

How Do We Leach Bottles?



- Allow bottles to sit for at least 7 days at room temp (20°C)
  - Or 4 days at 50°C



- Empty out the bottles and rinse at least 3Xs with DI H_2O



- Clean bottle for general purpose use
  - Ex. trace analysis

Environmental Contamination

- Classed clean rooms help prevent particle contamination
- Several criteria to consider to class a clean room
  - Air changes per minute
  - Filter efficiency (ULPAs or HEPAs)
  - Ceiling type / ceiling panel material
  - Light fixture type
  - Wall and floor material
  - Air return location
  - Assessed by counting # of particles of specified size



Reagent Purity

- Clean Acids (HNO3, HCl, HF, etc.)
- Clean Bases (TEA, NH4OH, etc.)
- Clean DI Water (ASTM Type 1, 18MΩ)
- Particle Filtration (0.3µm or smaller)
- Use of reagent blanks
  - Critical for blank subtraction to determine true real-world sample impurity levels and detection limits

02

Top Topics for Today:

- Environmental Contamination
- **Solution Matrix**
  Will your sample and standards be stable for your analysis?
- Washout Issues
- Uncertainty of Measurements

Solution Matrix

- Most elements are stable using only HNO_3
- Some elements require HF for stability
  - HF stability often revolves around plastic surfaces
  - Some elements will precipitate in the presence of HF
- Some elements require HCl for stability
  - Only a few elements have issues in the presence of HCl
- Bromide/Iodide require basic matrices for stability
  - Some other elements are stable in base without HNO_3/HF

Solution Matrix – HNO_3/HF

H			Not Checked by ICP		HF Elements												He
Li	Be		HF "thieves"		Avoid HF							B	C	N	O	F	Ne
Na	Mg											Al	Si	P	S	Cl	Ar
K	Ca		Sc	Ti	V	Cr	Mn	Fe	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Br	Kr
Rb	Sr		Y	Zr	Nb	Mo	Tc	Ru	Rh	Pd	Ag	Cd	In	Sn	Sb	Te	I	Xe
Cs	Ba	*	Lu	Hf	Ta	W	Re	Os	Ir	Pt	Au	Hg	Tl	Pb	Bi	Po	At	Rn
Fr	Ra	**	Lr	Rf	Db	Sg	Bh	Hs	Mt	Ds	Rg	Cn	Nh	Fl	Mc	Lv	Ts	Og

*	La	Ce	Pr	Nd	Pm	Sm	Eu	Gd	Tb	Dy	Ho	Er	Tm	Yb
**	Ac	Th	Pa	U	Np	Pu	Am	Cm	Bk	Cf	Es	Fm	Md	No

Solution Matrix - HCl

H																	He
Li	Be											B	C	N	O	F	Ne
Na	Mg											Al	Si	P	S	Cl	Ar
K	Ca		Sc	Ti	V	Cr	Mn	Fe	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Br	Kr
Rb	Sr		Y	Zr	Nb	Mo	Tc	Ru	Rh	Pd	Ag	Cd	In	Sn	Sb	Te	I	Xe
Cs	Ba	*	Lu	Hf	Ta	W	Re	Os	Ir	Pt	Au	Hg	Tl	Pb	Bi	Po	At	Rn
Fr	Ra	**	Lr	Rf	Db	Sg	Bh	Hs	Mt	Ds	Rg	Cn	Nh	Fl	Mc	Lv	Ts	Og

*	La	Ce	Pr	Nd	Pm	Sm	Eu	Gd	Tb	Dy	Ho	Er	Tm	Yb
**	Ac	Th	Pa	U	Np	Pu	Am	Cm	Bk	Cf	Es	Fm	Md	No

Not Checked by ICP	HCl Elements
Can work w/o HF	Avoid HCl

Solution Matrix -
TEA/NH_4OH/H_2O

H																	He
Li	Be											B	C	N	O	F	Ne
Na	Mg											Al	Si	P	S	Cl	Ar
K	Ca		Sc	Ti	V	Cr	Mn	Fe	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Br	Kr
Rb	Sr		Y	Zr	Nb	Mo	Tc	Ru	Rh	Pd	Ag	Cd	In	Sn	Sb	Te	I	Xe
Cs	Ba	*	Lu	Hf	Ta	W	Re	Os	Ir	Pt	Au	Hg	Tl	Pb	Bi	Po	At	Rn
Fr	Ra	**	Lr	Rf	Db	Sg	Bh	Hs	Mt	Ds	Rg	Cn	Nh	Fl	Mc	Lv	Ts	Og

Not Checked by ICP

Basic Elements

Can work w/o HNO_3/HF, but must be basic

*	La	Ce	Pr	Nd	Pm	Sm	Eu	Gd	Tb	Dy	Ho	Er	Tm	Yb
**	Ac	Th	Pa	U	Np	Pu	Am	Cm	Bk	Cf	Es	Fm	Md	No

iNORGANIC VENTURES

Multi-Element Compatibility

(How many elements can be combined at 1,000 ppm?)

H			15% HNO_3					+ 3% HF									He
Li	Be		10% HCl									B	C	N	O	F	Ne
Na	Mg											Al	Si	P	S	Cl	Ar
K	Ca		Sc	Ti	V	Cr	Mn	Fe	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Br	Kr
Rb	Sr		Y	Zr	Nb	Mo	Tc	Ru	Rh	Pd	Ag	Cd	In	Sn	Sb	Te	I	Xe
Cs	Ba	*	Lu	Hf	Ta	W	Re	Os	Ir	Pt	Au	Hg	Tl	Pb	Bi	Po	At	Rn
Fr	Ra	**	Lr	Rf	Db	Sg	Bh	Hs	Mt	Ds	Rg	Cn	Nh	Fl	Mc	Lv	Ts	Og

*	La	Ce	Pr	Nd	Pm	Sm	Eu	Gd	Tb	Dy	Ho	Er	Tm	Yb
**	Ac	Th	Pa	U	Np	Pu	Am	Cm	Bk	Cf	Es	Fm	Md	No

Multi-Element Compatibility

(How many elements can be combined at 100 ppm?)

H			5 - 10% HNO_3 / tr. HF				+ 3% HF										He
Li	Be		10% HCl									B	C	N	O	F	Ne
Na	Mg											Al	Si	P	S	Cl	Ar
K	Ca		Sc	Ti	V	Cr	Mn	Fe	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Br	Kr
Rb	Sr		Y	Zr	Nb	Mo	Tc	Ru	Rh	Pd	Ag	Cd	In	Sn	Sb	Te	I	Xe
Cs	Ba	*	Lu	Hf	Ta	W	Re	Os	Ir	Pt	Au	Hg	Tl	Pb	Bi	Po	At	Rn
Fr	Ra	**	Lr	Rf	Db	Sg	Bh	Hs	Mt	Ds	Rg	Cn	Nh	Fl	Mc	Lv	Ts	Og

*	La	Ce	Pr	Nd	Pm	Sm	Eu	Gd	Tb	Dy	Ho	Er	Tm	Yb
**	Ac	Th	Pa	U	Np	Pu	Am	Cm	Bk	Cf	Es	Fm	Md	No

Multi-Element Compatibility

(How many elements can be combined at 10 ppm?)

H			5 - 10% HNO_3 / tr. HF														He
Li	Be		10% HCl									B	C	N	O	F	Ne
Na	Mg											Al	Si	P	S	Cl	Ar
K	Ca		Sc	Ti	V	Cr	Mn	Fe	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Br	Kr
Rb	Sr		Y	Zr	Nb	Mo	Tc	Ru	Rh	Pd	Ag	Cd	In	Sn	Sb	Te	I	Xe
Cs	Ba	*	Lu	Hf	Ta	W	Re	Os	Ir	Pt	Au	Hg	Tl	Pb	Bi	Po	At	Rn
Fr	Ra	**	Lr	Rf	Db	Sg	Bh	Hs	Mt	Ds	Rg	Cn	Nh	Fl	Mc	Lv	Ts	Og

*	La	Ce	Pr	Nd	Pm	Sm	Eu	Gd	Tb	Dy	Ho	Er	Tm	Yb
**	Ac	Th	Pa	U	Np	Pu	Am	Cm	Bk	Cf	Es	Fm	Md	No

Multi-Element Compatibility

(How many elements can be combined at 10 ppm?)

H			10% HCl / tr. HF		$(Cs,Rb)_2PtCl_6$												He
Li	Be											B	C	N	O	F	Ne
Na	Mg											Al	Si	P	S	Cl	Ar
K	Ca		Sc	Ti	V	Cr	Mn	Fe	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Br	Kr
Rb ?	Sr		Y	Zr	Nb	Mo	Tc	Ru	Rh	Pd	Ag	Cd	In	Sn	Sb	Te	I	Xe
Cs ?	Ba	*	Lu	Hf	Ta	W	Re	Os	Ir	Pt?	Au	Hg	Tl	Pb	Bi	Po	At	Rn
Fr	Ra	**	Lr	Rf	Db	Sg	Bh	Hs	Mt	Ds	Rg	Cn	Nh	Fl	Mc	Lv	Ts	Og

*	La	Ce	Pr	Nd	Pm	Sm	Eu	Gd	Tb	Dy	Ho	Er	Tm	Yb
**	Ac	Th	Pa	U	Np	Pu	Am	Cm	Bk	Cf	Es	Fm	Md	No

INORGANIC VENTURES

03

Top Topics for Today:

- Environmental Contamination
- Solution Matrix
- **Washout Issues**
  Poor washout can lead to inaccurate results. How do your blanks look?
- Uncertainty of Measurements

Causes for Washout Issues

- Solution matrix (HNO_3/HCl/HF/NH_4OH)
  - Could be from commercially available standards
  - Could be from specific sample preparation protocols
- Specific elemental affinity toward different types of plastics (Pump tubing, HF-resistant intro system)
- Spray chamber design (Single-pass, double-pass)
- Rinse protocols using the wrong acids/bases
- Allowing for the appropriate amount of time to rinse

Problem Elements

H			Not Checked by ICP				Elements with washout issues										He
Li	Be											B	C	N	O	F	Ne
Na	Mg											Al	Si	P	S	Cl	Ar
K	Ca		Sc	Ti	V	Cr	Mn	Fe	Co	Ni	Cu	Zn	Ga	Ge	As	Se	Br	Kr
Rb	Sr		Y	Zr	Nb	Mo	Tc	Ru	Rh	Pd	Ag	Cd	In	Sn	Sb	Te	I	Xe
Cs	Ba	*	Lu	Hf	Ta	W	Re	Os	Ir	Pt	Au	Hg	Tl	Pb	Bi	Po	At	Rn
Fr	Ra	**	Lr	Rf	Db	Sg	Bh	Hs	Mt	Ds	Rg	Cn	Nh	Fl	Mc	Lv	Ts	Og

*	La	Ce	Pr	Nd	Pm	Sm	Eu	Gd	Tb	Dy	Ho	Er	Tm	Yb
**	Ac	Th	Pa	U	Np	Pu	Am	Cm	Bk	Cf	Es	Fm	Md	No

Common Rinse Solutions

HNO_3

- 5-10% on an OES
- 1-2% on a MS

HCl

- 5-10% on an OES
- 1-2% on a MS

RBS-25

- 2.5% on an OES
- Not Recommended on MS due to high Sodium

H_2O

- Can be effective enough for Na, K, Ca, etc.

Rinse Solutions with HF

HNO_3

- 5-10% on an OES
- 1-2% on a MS


HF

- 0.1-2% on an OES
- 0.05-0.5% on a MS

- If using borosilicate glass nebulizer and spray chamber
  - ***Limit HF*** to a max of 0.2%
  - B and Si results will be unreliable



- If using HF Resistant nebulizer and spray chamber
  - Limit HF to a max of 3%

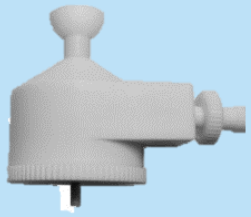

Specialty Rinse Solutions

- **NH_4OH**
  - 1-5% for OES or MS
  - Use for B, Br, I, Hg
- **HCl / Thiourea**
  - 1-10% HCl
  - 0.5% Thiourea
  - Use for Hg, Au, Os
- **HCl / Hydroxylamine·HCl**
  - 10% HCl
  - 0.5% $NH_2OH{\cdot}HCl$
  - Use for Os

Washout Study Results

Washout Times Comparison				
	Time to Washout (s)			
Rinse Solution	Au 208.209 nm	Bi 223.061 nm	Hg 184.950 nm	Os 225.585
10% HNO3	50	30	175	170
10% HCl	40	80	45	100
5% HNO3 / 5% HCl	70	170	110	120
1% NH4OH	60	60	130	60
10% HCl / 0.5% HydroxylamineHCl	50	80	50	75
ICP-TRUE-RINSE	35	25	35	30

ICP-TRUE-RINSE decreased washout time for all four elements in the study.

Troubleshooting Washout Issues

- Identify the elements that give you trouble
- Experiment with different rinse solutions
- *Keep up with routine instrument maintenance*
- Identify specific causes for washout issues
- Develop methods that include effective rinse strategies and only the elements you require


04

Top Topics for Today:

- Environmental Contamination
- Solution Matrix
- Washout Issues
- **Uncertainty of Measurements**
  What is a reasonable uncertainty for your analysis, and can you take any steps to improve it?

Standards of Comparison

- How low can you go in terms of uncertainty?
  - Review the uncertainties of your CRMs for a baseline

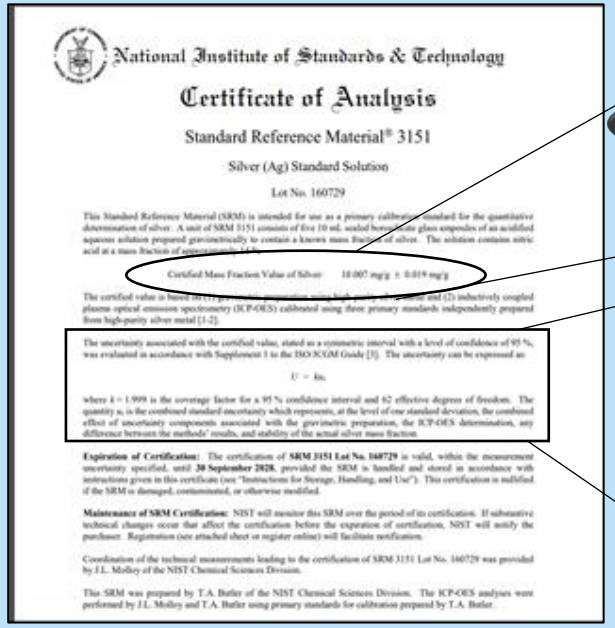

National Institute of Standards & Technology

Certificate of Analysis

Standard Reference Material® 3151

Silver (Ag) Standard Solution

Lot No. 160729

This Standard Reference Material (SRM) is intended for use as a primary calibration standard for the quantitative determination of silver. A unit of SRM 3151 consists of five 10 mL sealed borosilicate glass ampoules of an acidified aqueous solution prepared gravimetrically to contain a known mass fraction of silver. The solution contains nitric acid at a mass fraction of approximately 1.4 %.

Certified Mass Fraction Value of Silver: 10.007 mg/g ± 0.019 mg/g

The certified value is based on (1) gravimetric preparation using high-purity silver metal and (2) inductively coupled plasma optical emission spectrometry (ICP-OES) calibrated using three primary standards independently prepared from high-purity silver metal [1-2].

The uncertainty associated with the certified value, stated as a symmetric interval with a level of confidence of 95 %, was evaluated in accordance with Supplement 1 to the ISO/JCGM Guide [3]. The uncertainty can be expressed as:

$U = ku_c$

where $k = 1.999$ is the coverage factor for a 95 % confidence interval and 62 effective degrees of freedom. The quantity u_c is the combined standard uncertainty which represents, at the level of one standard deviation, the combined effect of uncertainty components associated with the gravimetric preparation, the ICP-OES determination, any difference between the methods' results, and stability of the actual silver mass fraction.

Expiration of Certification: The certification of **SRM 3151 Lot No. 160729** is valid, within the measurement uncertainty specified, until **30 September 2028**, provided the SRM is handled and stored in accordance with instructions given in this certificate (see "Instructions for Storage, Handling, and Use"). This certification is nullified if the SRM is damaged, contaminated, or otherwise modified.

Maintenance of SRM Certification: NIST will monitor this SRM over the period of its certification. If substantive technical changes occur that affect the certification before the expiration of certification, NIST will notify the purchaser. Registration (see attached sheet or register online) will facilitate notification.

Coordination of the technical measurements leading to the certification of SRM 3151 Lot No. 160729 was provided by J.L. Molloy of the NIST Chemical Sciences Division.

This SRM was prepared by T.A. Butler of the NIST Chemical Sciences Division. The ICP-OES analyses were performed by J.L. Molloy and T.A. Butler using primary standards for calibration prepared by T.A. Butler.

Certified Mass Fraction Value of Silver: 10.007 mg/g ± 0.019 mg/g

The uncertainty associated with the certified value, stated as a symmetric interval with a level of confidence of 95 %, was evaluated in accordance with Supplement 1 to the ISO/JCGM Guide [3]. The uncertainty can be expressed as:

$U = ku_c$

where $k = 1.999$ is the coverage factor for a 95 % confidence interval and 62 effective degrees of freedom. The quantity u_c is the combined standard uncertainty which represents, at the level of one standard deviation, the combined effect of uncertainty components associated with the gravimetric preparation, the ICP-OES determination, any difference between the methods' results, and stability of the actual silver mass fraction.

Sample Preparation

- Sample Preparation by weight is superior to preparation by volume.
- Example: Sample Aliquot Volume vs Weight

Preparation by Volume	Preparation by Weight		
1mL Aliquot Calibrated Electronic Pipettor	1g Aliquot 4-place Analytical Balance	1g Aliquot 3-place Analytical Balance	1g Aliquot 2-place Analytical Balance
0.3% Relative Error	0.03% Relative Error	0.2% Relative Error	1.5% Relative Error

Sample Preparation Best Practices

- Operate in units of weight/weight (Ex. µg/g) if possible

- If you must work in units of weight/volume (Ex. µg/mL) use densities to convert to weight/weight units

    1. Perform your analysis in weight/weight
    2. Convert your results back to weight/volume using your densities
    3. Remember to track the uncertainty from density measurements also

Instrumentation

Sample Intro System and Instrument Sampling Method

- This is one of the most important factors to achieve excellent results
- Double check tubing connections frequently
- This includes uptake times, stabilization delays and other settings like integration time

Matrix Matching

- Different matrices will demonstrate different nebulization efficiencies
- Density and TDS (Total Dissolved Solids) seem to have the largest impact on nebulization efficiency

Interferences

- It is important to know as much about your sample as possible
- Keep in mind contaminants that may have been introduced during sampling and/or storage
- Be mindful of interferences that may be present in your calibration curve but not present in your sample

Conclusions & Resources

- **Largest contributors to uncertainty**
  - Uncertainty of the standard of comparison
  - Measurements of small weights or volumes
  - RSDs

Resources

- GUM, Guide to the Expression of Uncertainty in Measurement
- VIM, International Vocabulary of Metrology
- EURACHEM Guide, Available from www.eurachem.org.
- Joint Committee for Guides in Metrology (JCGM)

  http://www.iso.org/sites/JCGM/JCGM-introduction.htm

Questions?

IV Ignite is the official virtual training academy of Inorganic Ventures!

TRY IT FREE

Our **Bench Boost Podcast** releases new episode every Tuesday!

TUNE IN

Visit **inorganicventures.com** to access free educational resources for chemists!

RESOURCES