OFFICE READER II

OCCUPATIONAL DESIGN










SMART

SMART ADD TECHNOLOGY ALLOWS THE EYES TO BE RELAXED WHEN VIEWING SCREEN DISPLAYS

Office Reader II incorporates Smart Add, the technology that improves the patient's visual experience when using digital devices. This technology assists with changes in focus at different working distances in a more agile and accurate way when the patient is working or reading simultaneously with different screens.

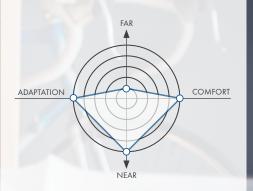

optimized areas. In addition, the soft design guarantees easy adaptation and unbeatable comfort.

4 VISION RANGES TO ADDRESS THE USER'S NEEDS.

Office Reader II offers 4 vision ranges to cover all of the patient's needs. Ask your patient what working distance is needed to focus comfortably so you can prescribe one of the 4 options.

 Office Reader II 1.3m Clear vision up to 1.3 m Office Reader II 2m Clear vision up to 2m Office Reader II 4m Clear vision up to 4m



TARGET MARKET

Wearers aged 45 and up performing near and intermediate tasks, depending on their visual range needs:

- Computer screen

- **Painting**









Smart Add

Not for driving computers

CONSUMER ADVANTAGES

Office Reader II is designed to be used as a second pair of glasses for presbyopes who perform activities at near or intermediate distances and need wider vision fields than a standard progressive lens can offer.

Softer design compared to the last

IMPROVED DESIGN

REDUCED VISUAL FATIGUE

Better ability to read text on

IMPROVED FOCUS

digital displays.

DYNAMIC VISION Easy and agile switch from the near to intermediate zone.

FASTER ADAPTATION

The adaptation to this lens practi-

NATURAL POSITION

Ergonomically comfortable

HOW TO PRESCRIBE OFFICE READER II

PROVIDE THE FOLLOWING DATA:

- Distance correction and addition power
- Range of Office Reader
- Calculation technology (Digital Ray-Path® or Surface Power®)
- Nasopupilar distance
- Frame data
- · Fitting height

All other personalization parameters are optional (pantoscopic angle, wrap angle, back vertex distance and near working distance).

LENS AVAILABILITY

Office Reader II 1.3m - Clear vision from 40 cm to 1.3 m Office Reader II 2m - Clear vision from 40 cm to 2.0 m Office Reader II 4m - Clear vision from 40 cm to 40 m Office Reader II 6m - Clear vision from 40 cm to 6.0

MINIMUM FITTING HEIGHT

14mm and 18mm

PLEASE NOTF

Available in Digital Ray-Path® or Surface Power® technologies.

Caution: Do not drive with these lenses because they do not provide distance vision.

MOUNTING INSTRUCTIONS

For proper positioning of the lens in the frame, the invisible engravings must be considered. The fitting cross should be aligned with the center of the pupil. Use the laser marks for verification.