


DELMARVA HFR High Frequency Radio Surgery



- 125 + watts of power, 4mhz in radio waves for easy, fast, and safe surgery.
- Totally isolated mono-polar configuration allowing for *Fine, Normal*, and *Coagulation* modes—with immediate *BiPolar* Switching
- Isolated circuitry as in expensive human hospital models, not ground referenced
- Isolated circuitry is compatible with all monitors
- 3 easy selectable modes of operation
 - *Cut* Fully rectified & filtered for smooth cutting and minimal hemostasis
 - *Blend* Fully rectified & not filtered for good cutting coupled with hemostasis
 - *Coag* Radio wave is modified for better coagulation controls bleeding and destroys unwanted tissue by desiccation
- Instant switching to **bipolar**, for pinpoint bleeder control or invert forceps for larger areas of coagulation
- Linear adjustable power control, assuring the correct power, whether patient is a pocket pet or equine.
- Autoclavable hand piece and cords
- Industrial grade hybrid electronics for stable, cool and safe operation
- Simple easy to use control panel
- HFR comes plug and play with all accessories and selection of electrodes.
- 3 year warranty

Frequency Chart

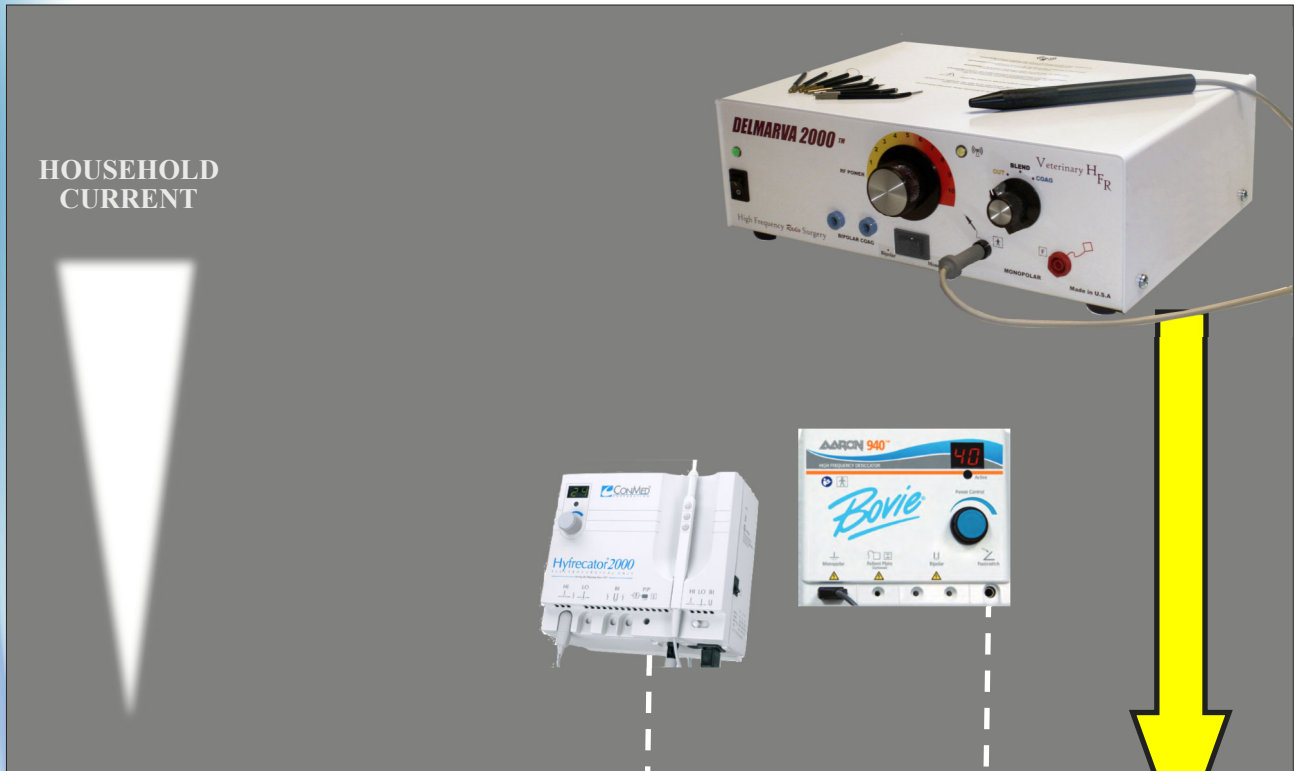

60 Hz 100 Hz 300 KHz 500 KHz 550 KHz

**4.0
MHZ**

DELMARVA 2000™

30+ Years of Distinction and
Superior Customer Service

www.delmarva2000.com

MADE IN THE USA