SWISS DOLORCLAST® MASTER

→ THE SOLUTION IN RADIAL SHOCK WAVE THERAPY





SWISS DOLORCLAST® MASTER

→ THE SOLUTION IN RADIAL SHOCK WAVE THERAPY






> Treat musculoskeletal pathologies quickly and effectively with the Radial ESWT solution.



© EMS SA FA-593-REV.A-01/EN ed. 05/2017

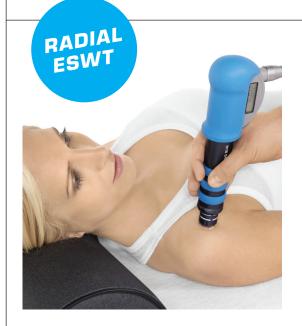

EFFECTIVE →

- > Clinically proven energy output:
 - Evo Blue® handpiece for constant energy density
 - Power⁺ handpiece for higher energy density

COMPLETE →

> Easy to master with its preset protocols, accurate manual settings, pressure ramp-up and energy density display

PATIENT-ORIENTED →

- > Inform your patients using treatment pictures, videos and method details.
- > Track clinical results by recording your patients data

PRODUCT REFERENCE	FT-217
TECHNOLOGY	Radial ESWT pneumatic with external compressor
POSITIVE ENERGY DENSITY (ED_)	Evo Blue® handpiece: 0.18 mJ/mm² - Power+ handpiece: 0.4 mJ/mm²
CAVITATION LEVEL	Increases with frequency
FREQUENCY	Up to 20Hz
CLINICAL STUDIES	26 Randomized Controlled Trials - More than 125 publications
SPECIAL FEATURES	17 preset protocols (up to 300), 3,000 patients records (VAS), pressure ramp-up, energy display, treatment pictures & videos
STANDARD CONFIGURATION	Swiss DolorClast® Master console, Evo Blue® handpiece, contact gel, cart & external compressor

EMS-MEDICAL.COM

For more information: welcome@ems-ch.com

