

POLY-LOK Polymer Ligating Clip System



Confidence, security and cost effective ligation with **POLY-LOK Polymer Clips** from **OCTO Medical**, indicated for ligating vessels and tubular structures in robotic, laparoscopic and open surgical procedures.

A dependable, CE Certified alternative with over 4 million POLY-LOK polymer clips sold worldwide in the past 5 years.

A complete range of sizes of clips and applicators, matched to the market leading Teleflex Hem-o-lok® system with the same handling, sizing and colour conventions, making it simple to switch.



3 year warranty
on reusable
POLY-LOK
Clip Applicators.

Reusable Applier Range

- Broad range of Appliers with different length and jaw configurations.
- Cost-effective Applier management programme to ensure consistent clip performance.
- 3 year warranty on reusable POLY-LOK Clip Appliers




Curved tip
Endo Appliers L, XL
Open Appliers



Right angle tip
Open Appliers



Straight tip
Endo Appliers M, ML

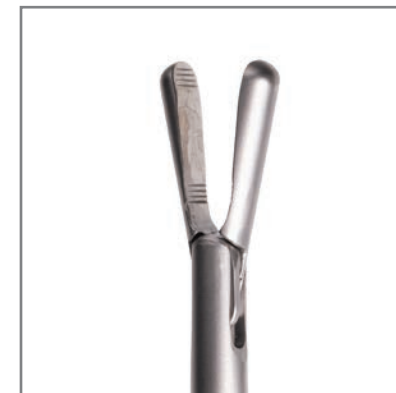

Remover
Endo Removers M/ML, L/XL
Open Removers ML/L, L/XL

Sterilisation

Prevacuum	3 times
Sterilisation temperature	132°C
Sterilisation time	1.5 min. (half cycle method)
Drying time	20 min

Sterilisation validated to ISO 17664 : 2017. Compliant with UK sterilisation standard HTM01-01.

The parameters shown were chosen to cover the widest spectrum of parameters used in reprocessing units all over the world. Other parameters with longer holding time and/or higher temperature are also validated by this procedure.

Disposable Polymer Clip Applier



Available in 3 sizes, ML, L, XL

Each device is supplied sterile complete with 2 cartridges of polymer clips.

POLY-LOK Polymer Ligating Clip System – Product Codes

POLY-LOK Clips		Endoscopic Manual Load Appliers & Removers				Open Manual Load Appliers & Removers		
		Code	Description	Configuration	Dia.	CODE	Description	Length
	M PLTX220 POLY-LOK Medium Non-Absorbable Polymer Ligating Clips 6 clips/cartridge 14 cartridges per box, sterile	PATX220E PATX220E3	POLY-LOK Endo Applier (0°) 330mm POLY-LOK Endo Applier (0°) 330mm	Single Piece Take-apart 3-piece	5mm 5mm	PATX220-20 PATX220-27 PATX220-27/70	8" Curved (9°) 11" Curved (9°) 11" Right Angle (70°)	200mm 270mm 270mm
	ML PLTX230 POLY-LOK Medium Large Non-Absorbable Polymer Ligating Clips 6 clips/cartridge 14 cartridges per box, sterile	PATX230E PATX230E3 PATX230EB PRTX230240E	POLY-LOK Endo Applier (0°) 330mm POLY-LOK Endo Applier (0°) 330mm POLY-LOK Endo Applier (0°) 450mm POLY-LOK Endo Remover ML/L 330mm	Single Piece Take-apart 3-piece Single Piece Single Piece	5mm 5mm 5mm 5mm	PATX230-20 PATX230-27 PATX230-27/70 PRTX230240-20 PRTX230240-27	8" Curved (9°) 11" Curved (9°) 11" Right Angle (70°) 8" Clip Remover ML/L 11" Clip Remover ML/L	200mm 270mm 270mm 200mm 270mm
	L PLTX240 POLY-LOK Large Non-Absorbable Polymer Ligating Clips 6 clips/cartridge 14 cartridges per box, sterile	PATX240E PATX240E3 PATX240EB PRTX230240E PRTX240250E	POLY-LOK Endo Applier (9°) 330mm POLY-LOK Endo Applier (9°) 330mm POLY-LOK Endo Applier (9°) 450mm POLY-LOK Endo Remover ML/L 330mm POLY-LOK Endo Remover L/XL 330mm	Single Piece Take-apart 3-piece Single Piece Single Piece Single Piece	10mm 10mm 10mm 5mm 10mm	PATX240-20 PATX240-27 PATX240-27/70 PRTX240250-20 PRTX240250-27	8" Curved (9°) 11" Curved (9°) 11" Right Angle (70°) 8" Clip Remover L/XL 11" Clip Remover L/XL	200mm 270mm 270mm 200mm 270mm
	XL PLTX250 POLY-LOK Extra Large Non-Absorbable Polymer Ligating Clips 6 clips/cartridge 14 cartridges per box, sterile	PATX250E PATX250E3 PATX250EB PRTX240250E	POLY-LOK Endo Applier (9°) 330mm POLY-LOK Endo Applier (9°) 330mm POLY-LOK Endo Applier (9°) 450mm POLY-LOK Endo Remover L/XL 330mm	Single Piece Take-apart 3-piece Single Piece Single Piece	10mm 10mm 10mm 10mm	PATX250-20 PATX250-27 PATX250-27/70 PRTX240250-20 PRTX240250-27	8" Curved (9°) 11" Curved (9°) 11" Right Angle (70°) 8" Clip Remover L/XL 11" Clip Remover L/XL	200mm 270mm 270mm 200mm 270mm

POLY-LOK Ligating Clips are not designed for use as a tubal occlusion device for fallopian tube contraception.
POLY-LOK Ligating Clips are contraindicated for ligating the renal artery in laparoscopic donor nephrectomies.



OCTO Medical Limited, 7 Bell Yard, London, WC2A 2JR, UK
Tel: +44 (0) 203 633 0198 | email: info@octomedical.org

A division of Winners First Trading Limited