

The health of human
The Pursuit of Yeso-med

Auto-disable Syringe (ZH-A)

Auto
-disable Syringe
Patented Product

0.5

Auto-disable Syringe (ZH-A)

Auto
-disable Syringe
Patented Product

● Main Features:

- ❶ Auto-disable feature: Syringe is auto-destroyed after injection, no extra operation is needed. No extra force is required as no obvious resistance during the activation process.
- ❷ Same operation as normal syringe, no special operation is needed.
- ❸ Air-exhaustion is allowed during the medication drawing, drawing can be continued after exhaustion.
- ❹ Fixed needle ensures the non-reusable feature of syringe and the lower medication retention level than normal syringe.

● Operation Procedures



- ❶ To avoid accidental activation, don't pull or push plunger before drawing.



- ❷ Air-exhaustion is allowed during the medication drawing, drawing can be continued after exhaustion.



- ❸ Syringe is auto-disabled by locking the metal clip once injection is completed.



Spec.	Needle spec.	Price
0.05ml		
0.1ml		
0.5ml		
1ml		

Customization is acceptable

Auto-disable Syringe (ZH-B)

Auto-disable
Disposable Syringe
Patented Product

Auto-disable Syringe (ZH-B)

Auto-disable
Disposable Syringe
Patented Product

● Main Features:

- ❶ Auto-disable feature: Syringe is auto-destroyed after injection, no extra operation is needed. No extra force is required as no obvious resistance during the activation process.
- ❷ Elastic blocking bars on the plunger can prevent accidental activation before use.
- ❸ Same operation as normal syringe, no special operation is needed.
- ❹ Fixed needle ensures the non-reusable feature of syringe and the lower medication retention level than normal syringe.

● Operation Procedures




❶ Same operation as normal syringe. To avoid accidental activation, don't push forward plunger before drawing.



❷ Push the plunger to the end of barrel, and activate the auto-disable feature.



❸ Plunger is broken as pulled out, needle and barrel are non-reusable.

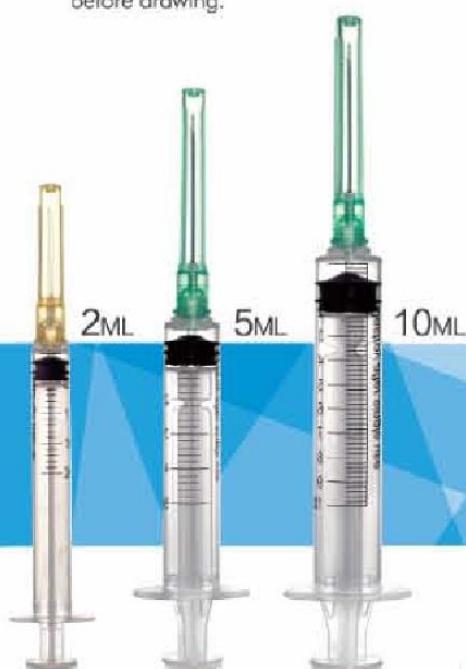

Spec.	Needle spec.	Price
2ml		
3ml		
5ml		
10ml		

Customization is acceptable

Sliding-sheath Safety Syringe (AZ-A)

**Retractable
Safety**

Disposable Syringe

Patented Product



Sliding-sheath Safety Syringe (AZ-A)

**Retractable
Safety**
Disposable Syringe
Patented Product

● Main Features:

- ❶ Same operation as normal syringe, accidental self-destruction is avoidable.
- ❷ Needle-changeable application, compatible with full range of standard needles.
- ❸ Dual-safety application
Safety before injection: needle retracted into barrel but not locked, syringe filled with medication in place, ready for injection.
Safety after injection: needle retracted into barrel and locked, put into safety box.

● Operation Procedures



- ❶ Easy operation, same as normal syringe.



- ❷ Slide sheath out and lock the needle into sheath.



- ❸ Put syringe into safety box.



Spec.	Needle spec.	Price
0.5ml		
1ml		
3ml		
5ml		
10ml		

Customization is acceptable

Safety Clip Syringe

Safety
Clip
Syringe



Safety Clip Syringe

Safety Clip Syringe

● Main Features:

- 1 Well-structured product, clip can be fixed before activating the safety feature.
- 2 Same operation as normal syringe.
- 3 Smooth movement of plunger and no extra resistance through whole process of injection.
- 4 Complete exhaustion is available before drawing.
- 5 One-hand activation of safety feature is available after injection.

● Operation Procedure

- 1 Fix the clip before drawing,
- 2 Same operation as normal syringe in drawing and injecting.
- 3 After injection, push the clip over the needle and lock it onto the clip base.



Spec.	Needle spec.	Price
0.5ml		
1ml		
3ml		
5ml		
10ml		
20ml		

Customization is acceptable



Safety needle

**Safety
needle**



- 1 Set the safety hypodermic needle on luer or luer lock cone of the standard syringe.
- 2 Pull down the protection device, pushing it on the side.
- 3 Screw the needle in order to fix it, by a rotating movement and remove the needle cap. Go on with the administration
- 4 After finishing, put up the protection device until the needle has been completely covered and "click" is heard

Spec.		Spec.		Spec.	
16G		23G		30G	
17G		24G		31G	
18G		25G		32G	
19G		26G		33G	
20G		27G			
21G		28G			
22G		29G			

Insulin syringe

Insulin syringe



Spec.	Needle spec.	Price
0.3ml		
0.5ml		
1ml		

Customization is acceptable.



0.3ml 0.5ml 1ml 0.3ml 0.5ml 1ml

Pen needle



Spec.	Needle spec.	Price
29G		
30G		
31G		
32G		

Customization is acceptable

TWO PARTS LUER SLIP SYRINGE

SPEC:2ML 3ML 5ML 10ML 20ML



THREE PARTS LUER SLIP SYRINGE

SPEC:1ML 2ML 3ML 5ML 10ML 20ML
30ML 50/60ML

THREE PARTS LUER LOCK SYRINGE

SPEC: 1ML 2ML 3ML 5ML 10ML 20ML
30ML 50/60ML



HYPODERMIC NEEDLE



SPEC: 16G 18G 19G 20G 21G 22G 23G 24G 25G 26G 27G 28G 29G 30G

- 1 YSF Series Linear Stapler
- 1 Disposable Linear Cutter Stapler
- 1 Disposable Alimentary Canal Stapler
- 1 Disposable Anorectal Stapler




WOUND PROTECTOR

PROTECTION
RETRACTOR
EXPOSURE



Mesh PerFix Plug

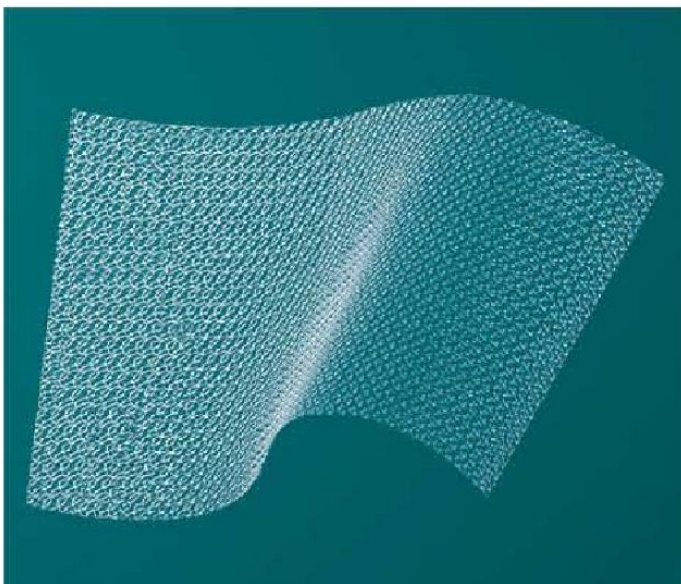



High Pressure Syringe

**High
Pressure
Syringe**





**MILLENNIUM
PRODUCTION GROUP**

Plainview, NY USA
800 280 1130
mpg-go.com
orders@mpggo.com



DENTAL SYRINGE



DENTAL SYRINGE

◀ DISPENSING TIP



Could manufacture specification as customer's requirement

Wipes – 25 pack/10 per PDQ



CLOR-O-CLEAN Disposable Disinfectant Wipes

Active Agent – Chlorine
Each pack consists of 25 wipes
Each PDQ consists of 10 packs
EPA Approved



VIR-U-CLEAN Disposal Disinfectant Wipes

Active Agent – Quaternary Ammonium
Each pack consists of 25 wipes
Each PDQ consists of 10 packs
EPA Approved



**MILLENNIUM
PRODUCTION GROUP**

orders@mpggo.com
800-280-1130

Wipes – 250 & 500

Drum of 250



Bucket of 500



Drum of 250



Bucket of 500



CLOR-O-CLEAN Disposable Disinfectant Wipes

Active Agent – Chlorine
Packed in buckets of 250 & 500
EPA Approved

VIR-U-CLEAN Disposal Disinfectant Wipes

Active Agent – Quaternary Ammonium
Packed in buckets of 250 and 500
EPA Approved



**MILLENNIUM
PRODUCTION GROUP**

orders@mpggo.com
800 – 280 – 1130

Protective Face Shields



Protective Face Shields

ALL ROUND PROTECTION - 180 degrees wide range protection.

HIGH-QUALITY MATERIAL - PET + sponge material. Anti-fog coating treatment, clear and double-sided anti-fog effect.

WIDELY USED - Use it in outdoor, kitchen and your office to protect your health. Adjustable elastic band can work for most people.

COMBINATION USE - The face shield can be worn on the glasses and can be used with disposable face masks to provide effective protection.



**MILLENNIUM
PRODUCTION GROUP**

orders@mpggo.com
800 - 280 - 1130

Safety Goggles



Safety Goggles

Fully covers eyes for full protection
Double Sided Anti-Fog, PVC Frame
Full cover with 4 valves
With elastic headband, fully adjustable
180 Degrees visual angle



**MILLENNIUM
PRODUCTION GROUP**

orders@mpggo.com
800 - 280 - 1130

Masks and Gowns



Model # MPG-0650

- Ear-Loop 3 Ply protective Mask
- Bacterial efficiency of 99%. A high quality
- Providing optimal respiratory protection
- Latex Free
- Glass free filter
- Pleated ear loop
- Level 1



Model# MSC0855

- Protective Gown with thumb hook
- Material: CPE (polyethylene)
- One size fits most
- Can be used in hospitals, laboratories
- AAMI: Level 2



Model# MPG-1100

- Protective Gown with thumb hook
- Material: CPE (polyethylene)
- One size fits most
- Can be used in hospitals, laboratories
- AAMI: Level 2

FDA
APPROVED

11-20-2020



**MILLENNIUM
PRODUCTION GROUP**

orders@mpggo.com
800 - 280 - 1130

Hand Sanitizer



8.45 oz Hand Sanitizer Liquid

- 75% Alcohol
- 1% Glycerine
- FDA Approved



4 oz Hand Sanitizer Liquid

- 80% Alcohol
- 1.45% Glycerine
- .125% Hydrogen Peroxide
- FDA Registered facility



**MILLENNIUM
PRODUCTION GROUP**

orders@mpggo.com
800 - 280 - 1130