

EnSURE[®] Touch

Smarter Monitoring. Designed for You.

EnSURE[®] Touch is an advanced monitoring system that collects, analyzes and reports data from multiple quality tests such as ATP, microorganisms, and enzymes. It features a 5" touch screen, wireless sync technology, cloud-based software and the ability to conduct a wide range of tests. EnSURE Touch is designed to adapt to your workplace, providing the data you need for audit and risk management.



Quick test: Run a test from any plan or location in as few as two touches.



Users: Manage user capabilities to differentiate personnel from managers. Enhanced administrator control.



Locations: Create testing points with customized fields and pass/fail limits on the go without needing desktop software.




Results & Reports: Easily view testing results and trend analysis.



Sync: Wirelessly synchronize data with SureTrend[®] Cloud.



MicroSnap[®]: Test and show results for organisms like Coliform, *E. coli*, *Enterobacteriaceae*, and Total Viable Count.*

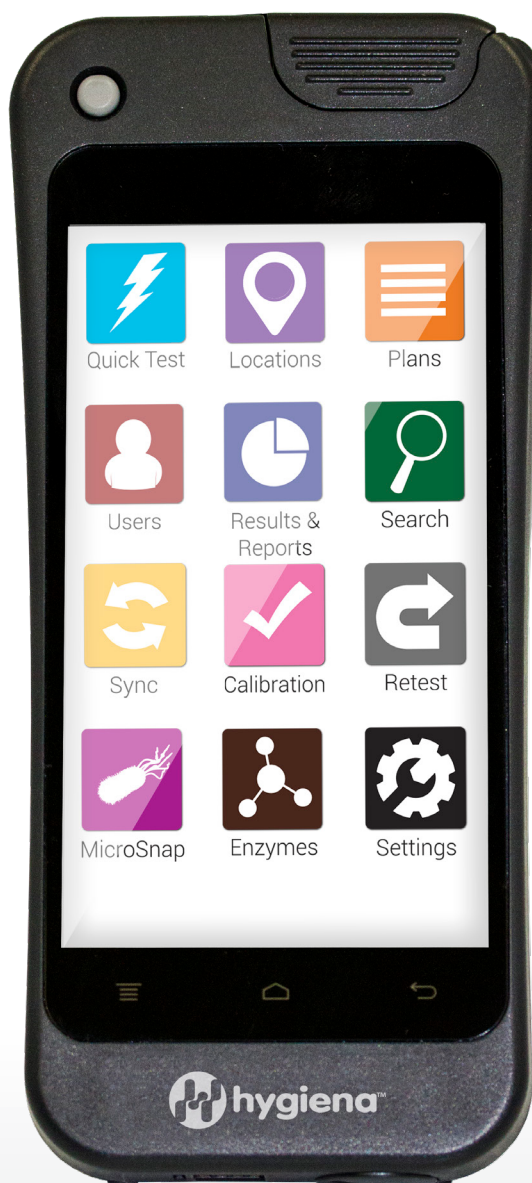

Plans: Easily group testing locations into plans and randomize test locations to ensure sampling coverage.



Search: Quickly find specific testing plans and locations on the fly, limiting time between tests.



Retest: Swipe through failed test locations and instantly run re-tests to show effectiveness of corrective action.



Settings: Modify settings for language, custom fields, tests types, simulate results from a competitor instrument, and contact support.



Calibration: Factory-calibrated and preset checks to stay GMP/ISO-compliant.



Enzymes: Test and show results for enzymes like Alkaline Phosphatase and Acid Phosphatase.*



**FOOD SAFETY SCIENTIFIC
AUSTRALIA PTY LTD**

PO Box 203, Baulkham Hills NSW 1755

T: +61 2 9674 2057 E: info@fssaust.com.au

FSSA

INNOVATIVE FOOD SAFETY SOLUTIONS

Features & Benefits

- Advanced photodiode technology and chemistry combined with unique algorithms make Hygiena's EnSURE Touch the most sensitive and durable system available
- No annual calibration requirements
- Built-in screen sharing capabilities allow for remote support directly on the device
- Wi-Fi or USB connectivity to cloud-based or desktop software
- Responsive 5" shatter-proof touch screen works even while wearing gloves
- Functions just like a smartphone and is configurable to your needs
- The ability to accurately emulate competitor instruments enables rapid system adoption
- Smart battery management features optimize device settings to maximize battery life
- Built-in tilt sensor makes sure you are running tests properly every time

SureTrend Cloud



- Fully interactive dashboard allows for quick data viewing
- Easily configure one or hundreds of instruments from a single SureTrend Account
- Monitor testing results across one or multiple facilities
- Schedule automatic reports delivered directly to your inbox
- Full online support and resource center only one click away
- Optional wireless software updates

Compatible Tests

ATP	Indicator Organism*	Process Efficiency*
UltraSnap® PERFORMANCE TESTED BY AOAC INTERNATIONAL	MicroSnap® EB	ZymoSnap® ALP
SuperSnap®	MicroSnap® <i>E. coli</i>	CROSS CHECK®
AquaSnap® FREE	MicroSnap® TOTAL	
AquaSnap® TOTAL	MicroSnap® COLIFORM	

Specifications

- Memory: 2GB internal memory and unlimited cloud-based storage
- Wi-Fi compatible
- Battery: Rechargeable lithium ion battery with USB-C charging.**
- Size: 7" x 3" x 1"
- Supports multiple languages



*Microorganism and enzyme testing will be available via future operating system updates
**Battery life varies by use and configuration