



DECORIA®

**THE PURE GENERATION OF
DERMAL FILLERS.
MADE IN SWEDEN.**

The home of Scandinavian HA


- **30+ years' experience of high-quality HA**
- **Pioneers in developing quality HA for ophthalmology, orthopedics and aesthetics**
- **Development and manufacturing in Strömstad, Sweden**
- **Strong sustainability focus with streamlined processes**

**BOHUS, Decoria® and SHAPE Technology™
(Smart Hyaluronic Acid PurE) are all
trademarks of Bohus Biotech AB.**

BOHUS
SCANDINAVIAN HA

Trädgårdsgatan 4
452 31 Strömstad, Sweden

DECORIA® PURE AND SAFE DERMAL FILLERS

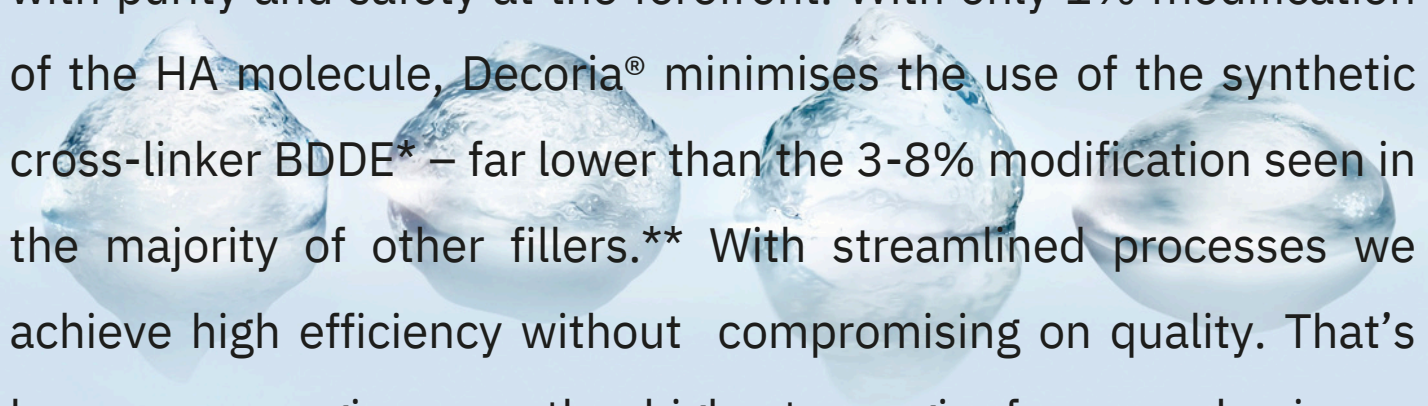


Decoria® represents a new generation of dermal fillers, designed with purity and safety at the forefront. With only 1% modification of the HA molecule, Decoria® minimises the use of the synthetic cross-linker BDDE* – far lower than the 3-8% modification seen in the majority of other fillers.** With streamlined processes we achieve high efficiency without compromising on quality. That's how we can give you the highest margin for your business. Scandinavian quality at an exceptional value.

* 1,4-butanediol diglycidyl ether.

** Wu G.T. et al., Hyaluronic Acid Basics and Rheology, Facial Plast Surg Clin North Am. 2022 Aug;30(3):301-308.

MINIMAL CROSS-LINKING MINIMISED RISK OF SIDE EFFECTS

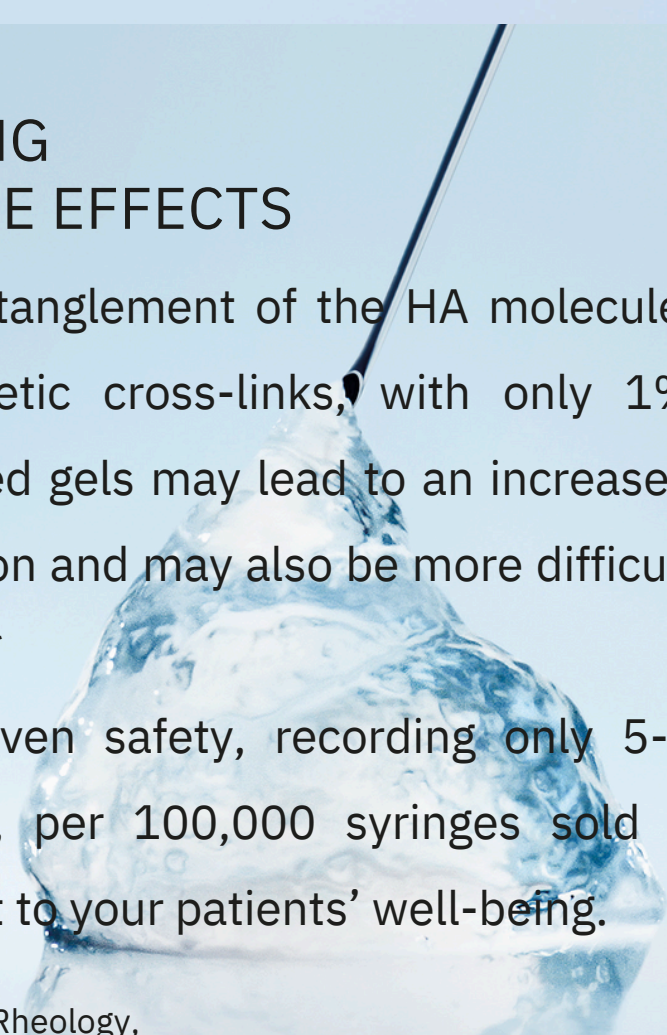

Decoria utilises the natural entanglement of the HA molecule, reducing the need for synthetic cross-links, with only 1% modification. Highly cross-linked gels may lead to an increased risk of an immunological reaction and may also be more difficult to dissolve with hyaluronidase.*


Decoria® stands out with proven safety, recording only 5-7 adverse events (non serious), per 100,000 syringes sold – demonstrating our commitment to your patients' well-being.

* Wu G.T. et al., Hyaluronic Acid Basics and Rheology, Facial Plast Surg Clin North Am. 2022 Aug;30(3):301-308.

DECORIA has a different, unique technology


Decoria has less than 1% BDDE (crosslinker) added, how is this possible?

- Decoria's unique technology uses the gel's natural tangling, making it more like tagliatelle (a tangled ball of HA strands) than spaghetti (as in uncrosslinked HA).
- Because Decoria uses the natural entanglement of the gel, it requires less than 1% BDDE (a cross-linker necessary to create a long-lasting dermal filler).
- Decoria's unique technology was developed by BOHUS based on +35 years of HA expertise.

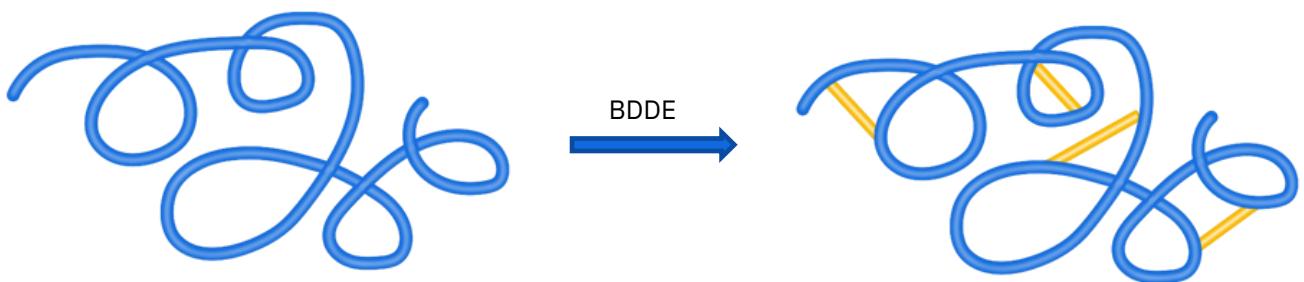


The characteristics of a filler are determined by the particle size and firmness (G') of the gel.

Firmness (G') and particle size are key factors that improve lifting capacity and better integration of the gel into soft tissue.

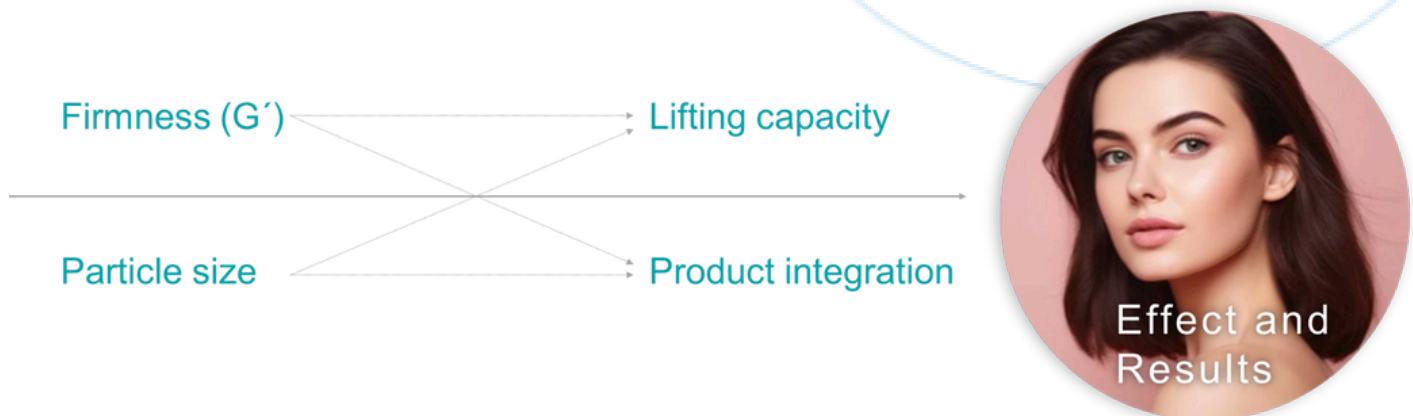

A soft gel integrates more into the tissue and is better suited for superficial injections or mobile areas (such as the lips, fine lines). A firm (hard) gel is better for deep injections and less mobile areas (such as the jawline, chin, or enhancing the cheekbones).

Decoria is designed for a complete portfolio.

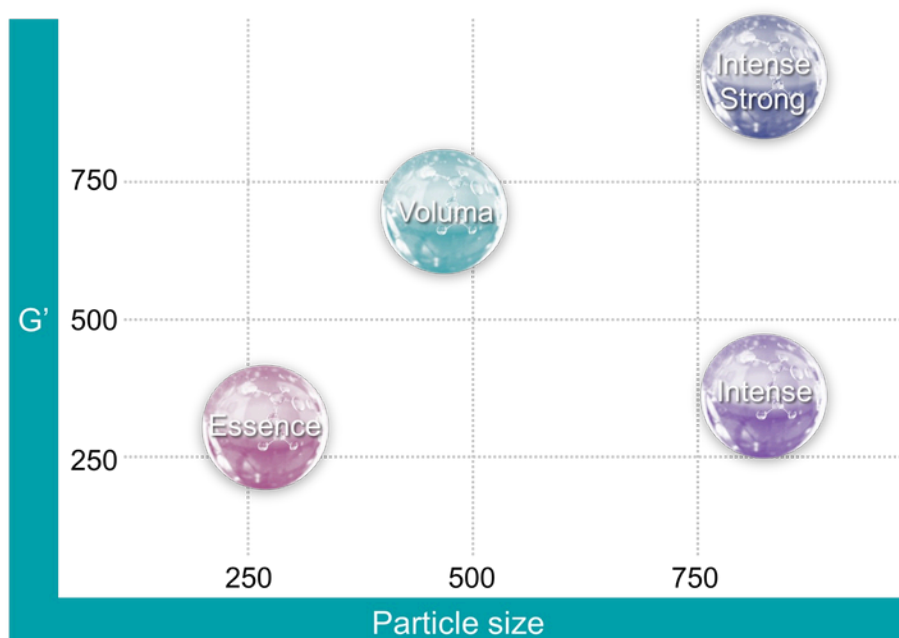

**Biocompatible
non-animal HA**

**Close to our
body's
natural HA**

**BDDE
< 1%**

SHAPE Technology™ (Global HA PrEcision)

- Helps restore youthful, natural facial features.
- Smooth and easy to inject, Helps Doctors recreate youthful, soft volume and shape lost with aging.

THE PURE DERMAL FILLER RANGE

	DECORIA® ESSENCE	DECORIA® INTENSE	DECORIA® VOLUMA	DECORIA® INTENSE STRONG
HA CONCENTRATION	20 mg/mL	20 mg/mL	20 mg/mL	20 mg/mL
COMPOSITION	Crosslinked HA	Crosslinked HA	Crosslinked HA	Crosslinked HA
TECHNOLOGY DEGREE OF MODIFICATION	SHAPE*	SHAPE	SHAPE	SHAPE
BDDE** CROSSLINKING	1%	1%	1%	1%
GEL CHARACTERISTIC	Soft voluminous gel for animated areas	Medium-firm gel for optimal volumising	Firm gel for universal support and shape	Extra-firm gel for sharp shapes and features
FOCUS INDICATIONS	Lips + Fine dynamic lines	Mid-face volume + Nasolabial folds + Higher volume lips	Cheeks + Nasolabial folds + Marionette lines	Jawline + Chin
INJECTION DEPTH	Superficial: Dermis Labial mucosa	Intermediate: Subcutaneous Labial mucosa	Intermediate to deep: Deep dermis Subcutaneous tissue to supraperiosteal administration	Deep: Deep dermis Subcutaneous tissue to supraperiosteal administration
NEEDLE	27G	27G	27G	27G



Decoria® Voluma, long lasting definition and shape

Clinical Study: Cheekbone Development

300 days. 90 patients. Maximum 5ml.

The efficacy and durability of Decoria® Voluma in cheekbone improvement were evaluated in a clinical study.


- Statistically significant improvement in zygomatic projection at day 300.
- On day 300, only a small portion of the filler was reabsorbed by 3D imaging (1.18 mL of filler at baseline, decreasing to 0.78 by the end of the study)*
- High patient satisfaction

* Torres S., Zygomatic Augmentation through Curved Cannulas with the Use of a New Stabilized Hyaluronic Acid Dermal Filler (Decoria® Voluma, BohusBioTech AB, Sweden, EC): 10 Months Follow Up, Medwin Publishers 2018 Nov. 23; Vol 3, Issue 3. All side effects were mild and resolved within 1-2 weeks.

DECORIA® ESSENCE




- Soft shape. Soft volume.
- Needle: Terumo 27G
- Hyaluronic acid concentration: 20 mg/ml

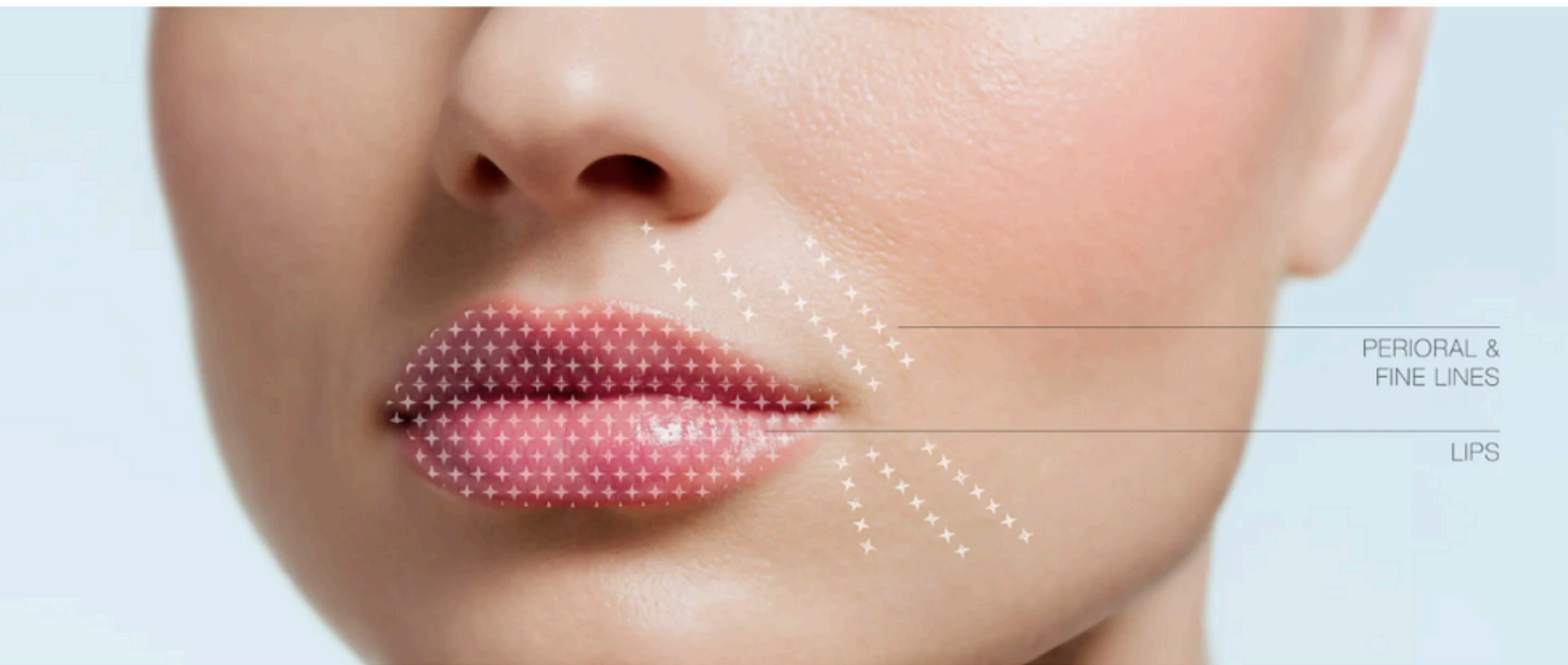

Decoria® Essence provides a natural look and feel to lips and is also suitable to treat perioral lines. Delicate adjustment. A soft voluminous gel and one of the purest lip fillers available.

Decoria® Essence is in transition to a slightly new look and feel. This new version has updated qualities specifically for lips. To be injected into the mid-to-deep dermis for correction of perioral lines and augmentation or correction of lips.

The depth of injection may vary depending on the treatment area.

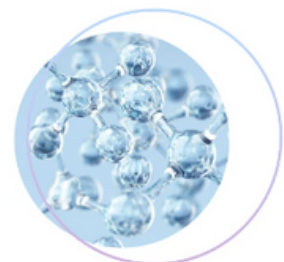

Pure, only 1% modification of the HA molecule.

DECORIA® VOLUMA




- Mid-firm gel. Universal support.
- Needle: Terumo 27G
- Hyaluronic acid concentration: 20 mg/ml



CHEEKS

NASOLABIAL FOLDS

MARIONETTE LINES

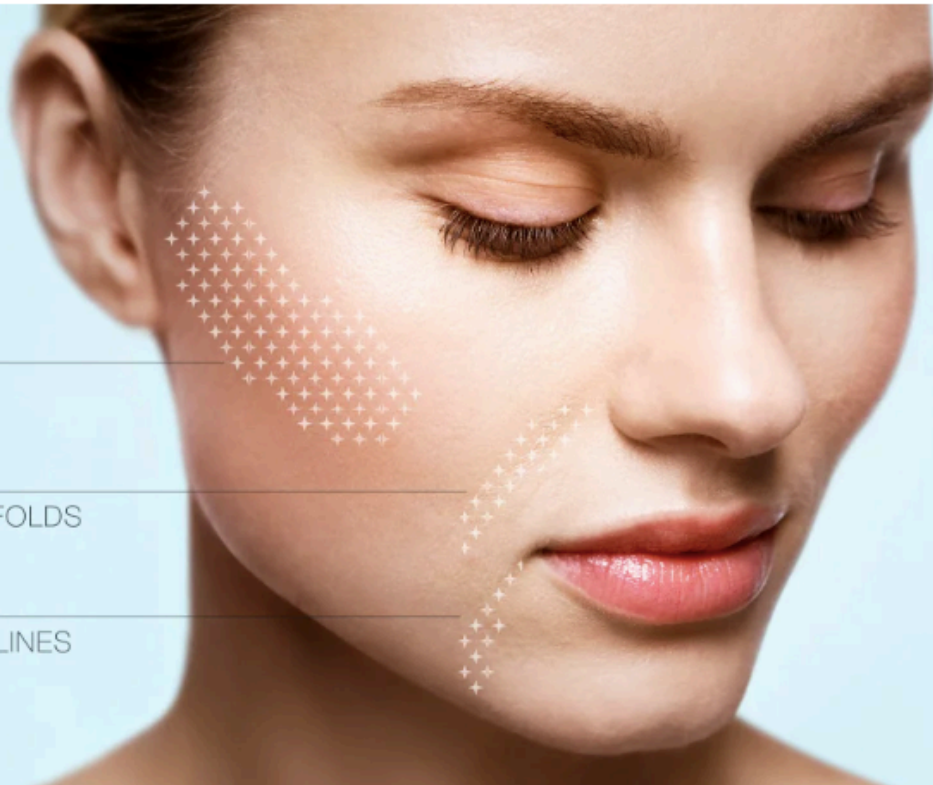

Decoria® Voluma provides universal support for cheeks, nasolabial folds and marionette lines. A firmer gel with medium lifting capacity for supportive shapes. Pure, with hyaluronic acid that is close to the body's own. Decoria® Voluma is a hyaluronic acid dermal filler to be injected into the deep dermis or the subcutaneous tissue to supraperiosteal administration. The depth of injection may vary depending on the treatment area.

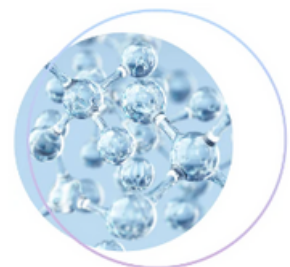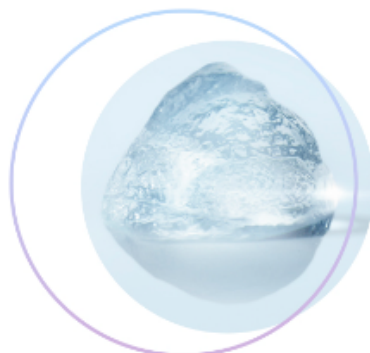

Pure, only 1% modification of the HA molecule.

DECORIA® INTENSE





- Soft shape. More volume.
- Needle: Terumo 27G
- Hyaluronic acid concentration: 20 mg/ml



CHEEKS, MID-FACE

Decoria® Intense provides a soft voluminous look to the cheeks. A medium to firm gel with large particles for optimal volumising. Soft shapes. Pure, with hyaluronic acid that is close to the body's own.. Decoria® Intense is a hyaluronic acid dermal filler to be injected into the deep layer of the dermis and/or the surface layer of the subcutis. The depth of injection may vary depending on the treatment area.

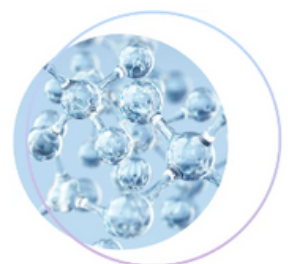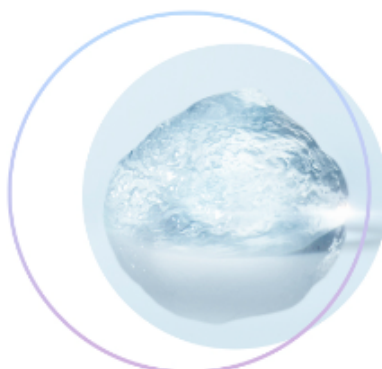

Pure, only 1% modification of the HA molecule.

DECORIA® INTENSE STRONG





- Firm gel. Strong shapes.
- Needle: Terumo 27G
- Hyaluronic acid concentration: 20 mg/ml



Decoria® Intense Strong is ideal when enhancing facial features such as the jawline and chin, to create a strong, defined profile. An extra firm gel with strong lifting capacity for a pure confident look.

Decoria® Intense Strong is a hyaluronic acid dermal filler to be injected into the deep dermis or the subcutaneous tissue to supraperiosteal administration.

The depth of injection may vary depending on the treatment area.

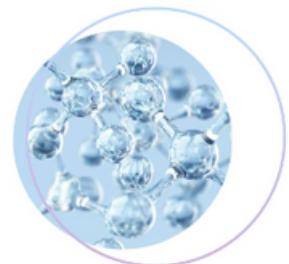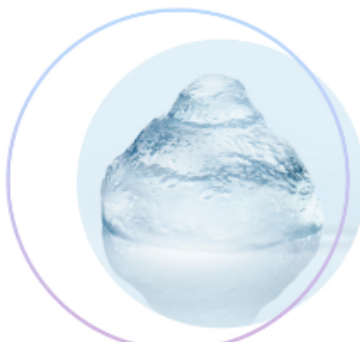

Pure, only 1% modification of the HA molecule.

DECORIA®



BOHUS
SCANDINAVIAN^{HA}


BIODERMALIST
• MEDİKAL VE ESTETİK ÜRÜNLERİ •

Türkiye yetkili
distribütörüdür.

