


Finally, a buttress that enables seamless tissue manipulation¹

Experience the difference with ECHELON ENDOPATH™ Staple Line Reinforcement.

ECHELON ENDOPATH™ Staple Line Reinforcement enables a seamless buttress experience from start to finish.¹ The simple click-and-go applicator and proprietary buttress attachment material (BAM) are designed for consistent loading every time.¹ The applicator allows you to experience precise control during attachment while the BAM provides seamless tissue manipulation and release.¹ All so you can transect as you intend.

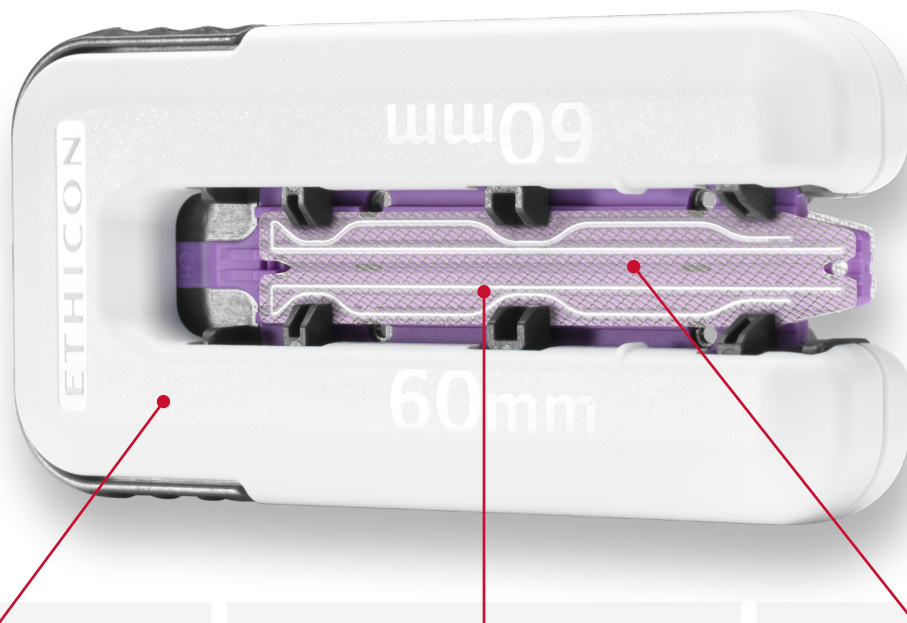

Seamless application

- **Simple click-and-go applicator** with no anvil or cartridge side for fast, easy buttress application.
- **Designed for OR efficiency.** 4 steps to attach compared to 12 steps to attach with GORE® SEAMGUARD® Staple Line Reinforcement and 16 steps with PERI-STRIPS DRY with VERITAS Collagen Matrix Staple Line Reinforcement.²

Seamless tissue manipulation

- Proprietary buttress attachment material (BAM) affixes buttress across the **entire anvil and cartridge face** for no slipping, twisting, sliding or bunching.¹
- **Delivers exceptional coverage.** 100% of ECHELON ENDOPATH™ Staple Line Reinforcement applications covered all staples with buttress vs. 28% for GORE® SEAMGUARD® Staple Line Reinforcement³ and 0% for Endo GIA™ Reinforced Reload with Tri-Staple™ Technology.⁴

Seamless performance

- Reinforcement made from the **same trusted materials used in VICRYL® and PDS® sutures.**
- **Exhibited the same hemostasis** as GORE® SEAMGUARD® Staple Line Reinforcement with 45% less material left behind.⁵

¹ Based on coverage after manipulation assessment on porcine tissue. ² Ricketts C, Pollack E. Evaluation of setup by nurses of a novel click and go staple line reinforcement system. *Glob Surg.* 2020;6. doi:10.15761/GOS.1000215. ³ In a benchtop test, ECHELON ENDOPATH™ Reinforcement covered all staple legs in (40/40) applications compared to GORE® SEAMGUARD® (11/40), p<0.05. ⁴ In a benchtop test, ECHELON ENDOPATH™ Reinforcement fully covered staple legs in (15/15) applications compared to Endo GIA™ Reinforced Reload with Tri-Staple™ Technology (0/15), p<0.05. ⁵ ECHELON ENDOPATH™ Staple Line Reinforcement vs. GORE® SEAMGUARD® Staple Line Reinforcement. Hemostasis: evaluated via gross observations of firings requiring intervention in 1.38mm to 2.32mm thick porcine ileum and jejunum tissues, evaluated via a 5-point Likert scale resulting in statistically no difference in Likert scores. Material: based on average total surface area measurements (26.79cm² vs. 49.21cm²).

Ordering information

ECHELON ENDOPATH™ Staple Line Reinforcement

Code	Description	Qty per sales unit
ECH60R	ECHELON ENDOPATH™ Staple Line Reinforcement	6



ECHELON™ 3000 Stapler

Code	Description	Jaw length	Shaft length	Qty per sales unit
ECH60C	ECHELON™ 3000 60mm Compact	60mm	280mm	3
ECH60S	ECHELON™ 3000 60mm Standard	60mm	340mm	3
ECH60L	ECHELON™ 3000 60mm Long	60mm	440mm	3



ECHELON™+ Stapler

Code	Description	Jaw length	Shaft length	Qty per sales unit
PCEE60A	ECHELON™+ 60mm Compact	60mm	280mm	3
PSEE60A	ECHELON™+ 60mm Standard	60mm	340mm	3
PLEE60A	ECHELON™+ 60mm Long	60mm	440mm	3



ECHELON ENDOPATH™ Reloads

Code	Color	Description	Open staple height	Closed staple height	Qty per sales unit
GST60W	White	White reload intended for use in thin tissue	2.6mm	1.0mm	12
GST60B	Blue	Blue reload intended for use in regular tissue	3.6mm	1.5mm	12
GST60D	Gold	Gold reload intended for use in regular/thick tissue	3.8mm	1.8mm	12
GST60G	Green	Green reload intended for use in thick tissue	4.1mm	2.0mm	12
GST60T	Black	Black reload intended for use in very thick tissue	4.2mm	2.3mm	12



How to order

All purchase orders are made to Johnson & Johnson Health Care Systems, Inc. (JJHCS):

- <https://us.jjcustomerconnect.com> or 1-866-565-4283
- Electronic Data Interchange (EDI) Helpline: 1-800-262-2888

Customer support

For product use assistance, clinical guidelines, service and repair, emergency assistance, copy of a 510(k) clearance letter, or complaints, please contact our Customer Support Center by calling 1-877-ETHICON (384-4266).

See it come to life with augmented reality:



Scan the QR code with your iOS 13 iPhone or iPad camera to visit the website. Then tap the augmented reality (AR) icon to start the experience.

Learn more at ethicon.com/Seamless