

Technical Specification

Ventilation Settings

Ventilation Modes

- VCV (A/C), PCV (A/C), PRVC (optional), PSV (optional), STANDBY
- SIMV (VCV)+PSV, SIMV (PCV)+PSV, SIMV (PRVC)+PSV
- SPONT/CPAP+PSV
- BIVENT/APRV +PSV (optional)
- NIV/CPAP, NIV-T, NIV-S/T

Enhancements

- Apnea ventilation, Pressure and Flow trigger, Automatic Tube Compensation (ATC), Smart suction
- Manual breath, Insp/Exphold, Screen freeze, Nebulization, Lung recruitment

Parameters

• Tidal volume (VT)	20-2000ml
• Respiratory rate (RR)	1 to 80 bpm
• Inspiratory time (Ti)	0.2 to 9 s (adult), 0.2 to 5 s (pediatric)
• Inspiratory flow (Flow)	0 to 100 L/min (pediatric), 0 to 180L/min (adult)
• Inspiratory pressure (P _{insp})	5 to 70 mbar (or cmH ₂ O)
• Inspiratory pressure limit (P _{max})	80 mbar (or cmH ₂ O)
• PEEP	0 to 35 mbar (or cmH ₂ O)
• Tslope	0 to 2 s
• O ₂ concentration (FiO ₂)	21 to 100 Vol%
• Trigger sensitivity	0.5 to 20 L/min (Flow trigger), -20 to 0 mbar (or cmH ₂ O) (Pressure trigger)
• I/E ratio	1/10 to 4/1
• Apnea alarm time	10-60 seconds

Monitoring

- Pressure values P_{plat}, P_{peak}, P_{mean}, P_{min}, PEEP
- Volume/Flow values VT_i, VT_E, MV, MV_e, MV_{spont}
- Time values f_{total}, f_{spont}, I:E
- Inspiratory O₂ concentration (FiO₂), End-expiratory CO₂ concentration (etCO₂)
- Compliance (dynamic & static), Resistance (R), MVleak, RSBI, WOB, I:E, V_{daw}, PEEP_i
- Pressure-Volume loop, Pressure-Flow loop, Flow-Volume loop

Alarms

Expiratory minute volume (MV) High/Low, Airway pressure (Paw) High/Low, VT_e Low, PEEP High/Low, Insp. O₂ concentration (FiO₂) High/Low, End-expiratory CO₂ concentration (etCO₂) High/Low, f_{spont} High, Apnea alarm, Disconnection, Flow sensor error, Gas supply, Electrical supply & battery failure, Exhalation obstruction, Apnea backup for low frequency alarm

Physical Specifications

• Dimensions (WxDxH)	375mm x 395mm x 430mm
• Weight	15kg (33.1lbs)
• Screen	12.1" TFT color touch screen

Remark: Above configurations include standard and option. Please check price with your Aeonmed sales representative.

AEOMED
Reliable Quality Thoughtful Service

HQ: No.4 Hangfeng Road, Fengtai Science Park, Fengtai District, Beijing, China (100070)
Science Park: No. 10, Chaobai Street, Yanjiao Development Zone, Sanhe City,
Hebei Province, China (065201)
TEL: +86-10-5841 1198 FAX: +86-10-6371 8989
Http://www.aeonmed.com E-mail: marketing@aeonmed.com

Information contained in this leaflet is correct at the date of publication. The policy of Beijing Aeonmed Co., Ltd. is one of the continued improvement to its products. Because of this policy Beijing Aeonmed Co., Ltd. reserves the right to make any change, which may affect the information in this leaflet without giving prior notices.
AVT70-1706



An Optimal Combination of Invasive and Noninvasive Ventilator

VG70 Ventilator

CE 0123

AEOMED
Reliable Quality Thoughtful Service

Superior Mobile ICU ventilator

- Comprehensive ICU ventilator including BIVENT and PRVC
- Compact, big capacity battery, no air compressor, intra-hospital mobility
- Flexible device configuration: equipped on a trolley, bed or ceiling pendant

Cost Effective Solution

- Unique metal-based, autoclavable, heated exhalation valve
- Built-in flow sensor, non-consumable design
- Upgradeable ventilation system software, with an available USB port




An Optimal Combination of Invasive and Noninvasive Ventilator


As noninvasive ventilation is used increasingly in a wide range of clinical situations, we offer a dual solution. VG70 combines the advantages of a flexible noninvasive ventilator with a full-featured invasive ventilator for the ICU.



360° Visible Alarm Lamp

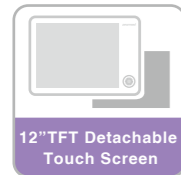

12" TFT Detachable Touch Screen

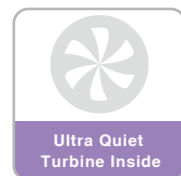

Ultra Quiet Turbine Inside



Integrated Power Supply Solution



Built-in Battery, With Extended Backup Option



Optimal patient-ventilator synchrony, increase patient comfort

- **The Unique Leak Compensation System** - Keep precise control on the tidal volume of each breath delivered to the patient by adjusting compensation dosage automatically
- **Advanced Trigger Technique** - Enhance sensitivity, avoid spurious triggering

Auto-detect and Adjust Leak Compensation

Automatically Adapt to Patient's Breathing Pattern

Multi-parameter Monitoring

Safe Ventilation Through Whole Treatment Phase

Initial Treatment Phase

- Noninvasive ventilation mode associated with decreased intubation rates, shortened patient stays, improved patient comfort, and a reduced risk of cross infection
- Preset patient's height and IBW. Reduce clinician's workload

Stable Condition Phase

- PRVC and BIVENT employ lung-protective strategies, delivering intelligent ventilation
- Comprehensive lung mechanics monitoring include compliance, airway resistance, PEEP_i and time constant
- Three waveforms & three loops with user-friendly display provide a continuous monitoring of the patient's condition

Weaning Phase

- Various ventilation modes enhance the weaning process
- The unique trigger and leakage compensation system safeguards each and every patient breath resulting in smooth and comfortable breathing, avoiding extra workload on the patient and promoting recovery
- RSBI and WOB provide accurate reference for weaning

Rehab Phase

- Data export port provides connection to hospital monitors and Patient Data Management Systems
- Provides pressure support for the patient when spontaneous breathing is present

