


# Collection & Shipping Guidelines for Blood Samples

**START HERE**

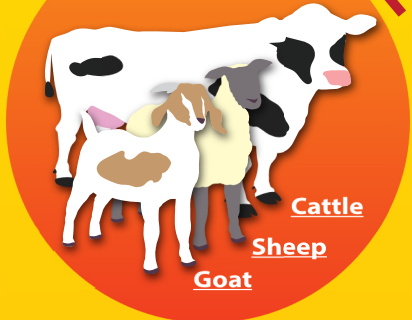

Cattle  
Sheep  
Goat


ANIMAL ID and TUBE # REQUIRED on each Red Top Tube



**IF NOT SENDING RIGHT AWAY**  
1-2 Days Max in Fridge


1ml/cc of Blood  
**OR**  
1ml/cc of Serum

Let Blood Samples Sit for at least 1 HOUR Before Putting In Fridge

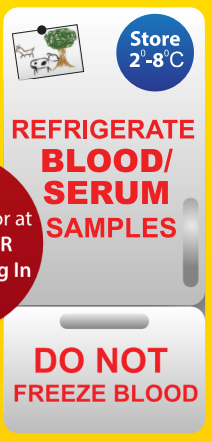

**WHEN READY TO SHIP**

Bundle Samples in Max of 10 per Set

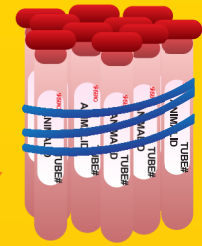

Cover Samples with absorbent paper




Paper Towels, Newspaper or etc...

**IF SENDING RIGHT AWAY**

**Frozen Ice Pack Required For**  
CL, CAE/OPPV, Johne's & Q-Fever



**NO Frozen Water or Dry Ice**

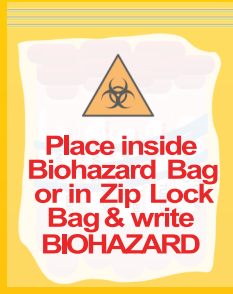

Place inside Biohazard Bag or in Zip Lock Bag & write BIOHAZARD

3-5 Day shipping for CL, CAE/OPPV, Johne's & Q-Fever.  
Pregnancy (BioPRYN) can be sent standard shipping.



**SHIPPING**

- DO NOT send samples in an envelope. Samples must be sent in a box.
- NO Frozen Water or Dry Ice.
- Friday is the ideal time to ship samples.
- 3-5 day shipping for all samples except (BioPRYN) pregnancy. Pregnancy can be sent at room temperature and standard shipping.
- Email lab to provide tracking # & # of tests.
- BE AWARE OF TEST DAYS to avoid delays.

**PAYMENT**

- We accept Debit/Credit Cards, CHECKS, MONEY ORDERS and CASH. We do not advise sending cash by mail. Checks & money orders should be made to 'UBRL'

**bioPRYN®**  
Blood Pregnancy  
(No ice pack required if only sending pregnancy samples)



**IMPORTANT: 'EXEMPT ANIMAL SPECIMEN'**  
Must be written on box

Wrap in Bubble Wrap to protect tubes.


**DO NOT SHIP IN AN ENVELOPE!**  
(We receive broken tubes often from samples being sent in envelopes)

Place All Documents, Payment & ice pack inside its own zip lock bag to avoid being damaged



Invoice link will be emailed once samples arrive in lab

**U** UNIVERSAL  
**B** BIOMEDICAL  
**R** RESEARCH  
**L** LABORATORY  
Livestock Diagnostics

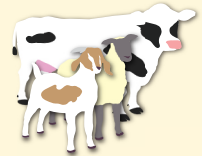

**UBRL**  
1300 N. Fresno St.  
Suite#220  
Fresno, CA 93703

(559)498-0820

**UBRL.ORG** Lab@UBRL.ORG