



High Performance Portable Suction

plus b - emergency

The V7 plus b Emergency Aspirator is rugged, light and field proven. With a superior flow-rate and pressure, it is lighter and more compact than comparable units. The pressure indicator gauge is positioned to provide excellent visibility for the operator and the pressure adjustment knob is positioned so it minimizes the chance of it getting caught or damaged while in use.



WITH EMS TRANSPORT BAG

FEATURES :

- Oil-free, no maintenance piston pump.
- Vacuum manometer with easy to turn control knob from 0 - 630mmHg
- Water resistant ON/OFF switch protector and LED battery
- Light reflecting indicators make it easy to see in low light conditions
- Weight: 9.25lb
- Dimensions: 10" x 15" x 6"
- Packing dimensions: 11 1/2"x 18" x 8 1/4"
- Optional 10G Rated Self-charging Truck/Car mount




EMS 12V Wall Mount





High Performance Portable Suction plus b - emergency

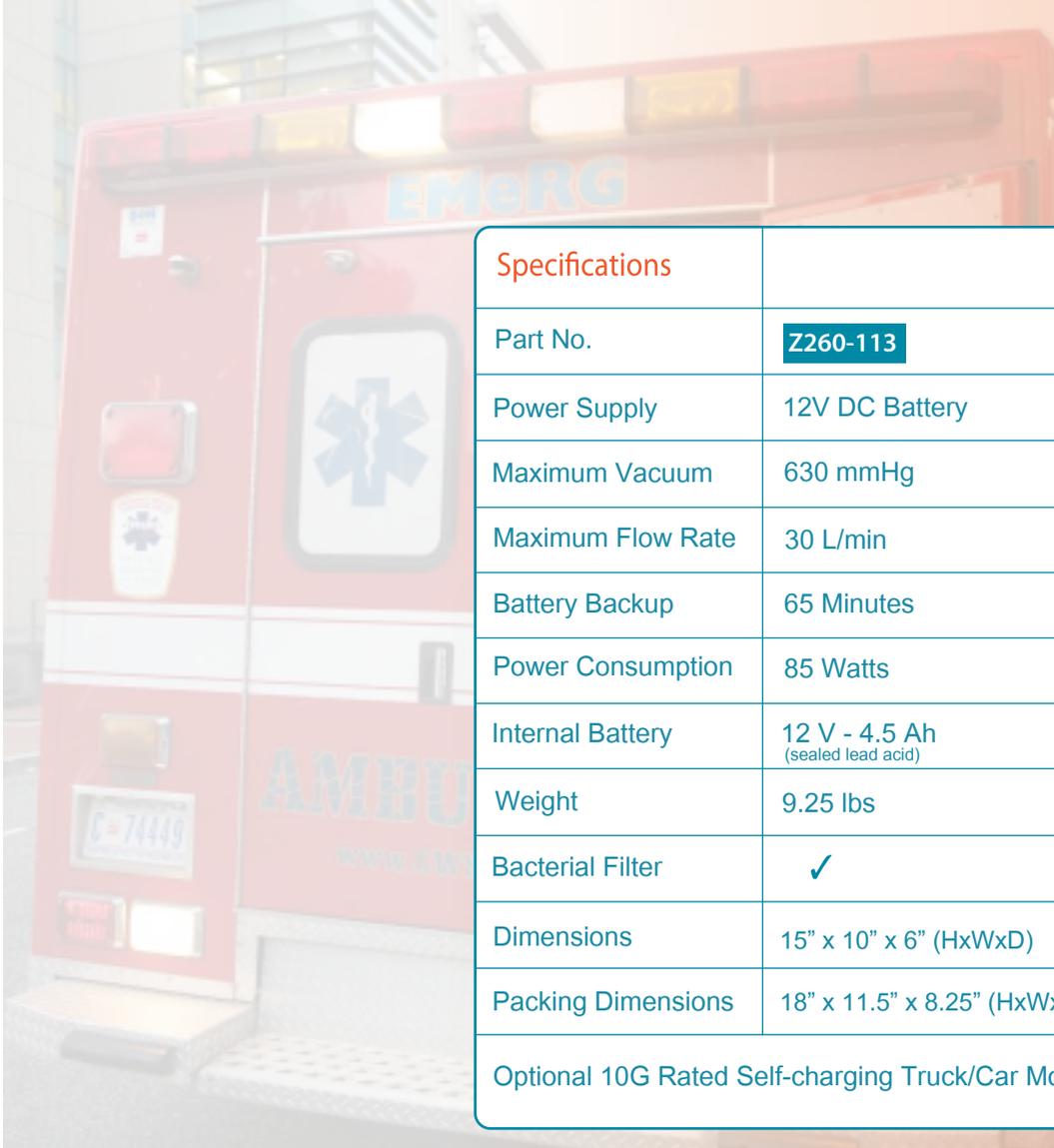

Specifications	
Part No.	Z260-113
Power Supply	12V DC Battery
Maximum Vacuum	630 mmHg
Maximum Flow Rate	30 L/min
Battery Backup	65 Minutes
Power Consumption	85 Watts
Internal Battery	12 V - 4.5 Ah (sealed lead acid)
Weight	9.25 lbs
Bacterial Filter	✓
Dimensions	15" x 10" x 6" (HxWxD)
Packing Dimensions	18" x 11.5" x 8.25" (HxWxD)
Optional 10G Rated Self-charging Truck/Car Mount	

Accessories and Parts	
Z260-302	10G Wall Bracket
Z260-310	Mobile Trolley Cart
Z260-320	Battery Charger AC/DC
Z260-319	Charging Cable 12/24V
Z260-342	Bacteria Filter
Z260-318	PC1 (1L) Reusable Collection Jar
Z260-111	Translucent Hose
Z260-103	EMS Transport Bag

For additional canister options, please contact your Vital Sales Associate.

Limited 2 year warranty. Excludes batteries, replacement tubing and filters.

