


miyo[®]
liquid ceramic

miyo[®]
PINK[®]
liquid ceramic

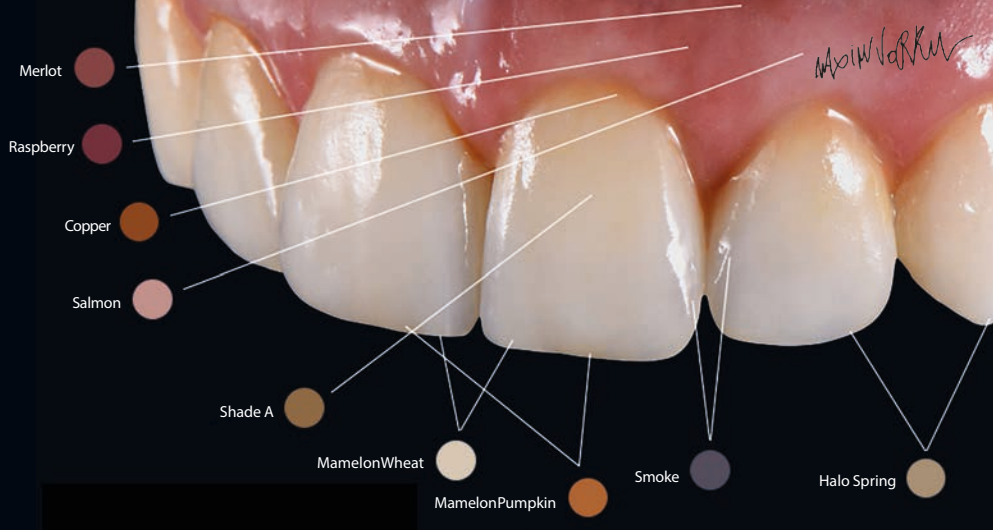

Informacije o proizvodu

srpski

Slojevita estetika - monolitni rad

Jensen trenutno nudi dva MiYO* Liquid Ceramic System seta, prijavljena za patent: MiYO Liquid Ceramic sistem za zube i MiYO pink Liquid Ceramic Gingival sistem za tkivo. Oba sistema predstavljaju jedinstven sistem za završnu obradu koji čine paste opalescentnih i fluorescentnih materijala za bojenje, strukturiranje i glaziranje.

Sa MiYO možete lako i brzo postići visoke estetske rezultate, uporedive sa slojevitim restauracijama za isto vreme koje bi bilo potrebno ako samo bojite.

MiYO Boja se sastoji od različitih tipova samo-glazirajućih boja sa različitim nivoima translucencije, od kojih je svaka jedinstveno formulisana da reprodukuje prirodnu boju i strukturu zuba i tkiva, kao što su incizalna translucencija, mameloni, linije mikro pukotina i gingivalno tkivo.

Preko MiYO boje stavljamo MiYO strukturu. Sa strukturom MiYO možete stvoriti dubinu, vitalnost i teksturu kao kod prirodne gledi i tkiva, u neverovatnim debljinama od 0.1mm do 0.2mm. Zahvaljujući visokoj stabilnosti individualne površinske strukture mogu lako biti integrisane u pastu sa četkicom.

Niska temperatura pečenja ne menja karakteristike površine, tako da se pouzdani i estetski rezultati postižu samo sa jednim ili dva pečenja.

* Prijavljena za patent

MiYO PREDNOSTI:

- + Pouzdani rezultati svaki put
- + Jednostavno rukovanje materijalima u pasti
- + Ekonomičnost i ušteda vremena kroz brze radne procese
- + MiYO boja sa savršeno podešenom fluorescencijom i neprozirnošću za kontrast i dubinu
- + MiYO Structure - jedinstvena struktura u pasti
- + Prirodni indeks prelamanja MiYO Structure masa za strukture je vidljiv na monolitnim krunicama debljine sloja već od 0,1 mm
- + Visoka stabilnost boja
- + Minimalno i redukovano skupljanje
- + Detaljni estetski rezultati, jer se sa MiYO masama ima potpuna kontrola nad oblikom, površinom i strukturom pre ciklusa pečenja

MiYO JE DIZAJNIRAN ZA:

- Monolitni cirkonijum
- Cirkonijum dioksid
- Litijum disilikat
- Cirkon-keramiku (npr. InSync Zr)
- Metal-keramiku (npr. InSync MC)
- Press na cirkon

MiYO PINK JE DIZAJNIRAN ZA:

- Monolitni cirkonijum
- Cirkonijum dioksid
- Litijum disilikat
- Cirkon-keramiku (npr. InSync Zr)
- Metal-keramiku (npr. InSync MC)




1

MIYOCOLOR



Različiti oblik, boja i efekat počevši od jedne baze

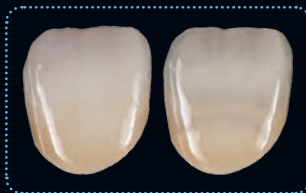

Od visoke do niske vrednosti sa
MiYO Color Storm i Smoke



Od A3 do A1 sa MiYO Color Lumín

2

MIYOSTRUCTURE



Strukturiranje površine


MiYO Structure pre
pečenja

MiYO Structure
posle pečenja




1

MIYO PINK COLOR

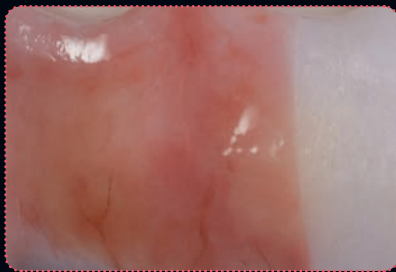


Neograničene mogućnosti
za detaljan dizajn.

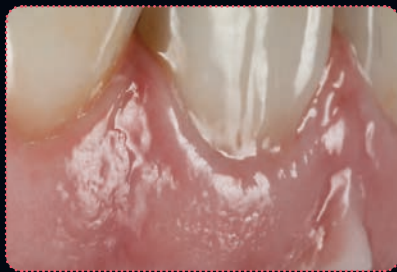

Dizajn prirodnog tkiva sa
MiYO pink Structure.

2

MIYO PINK STRUCTURE

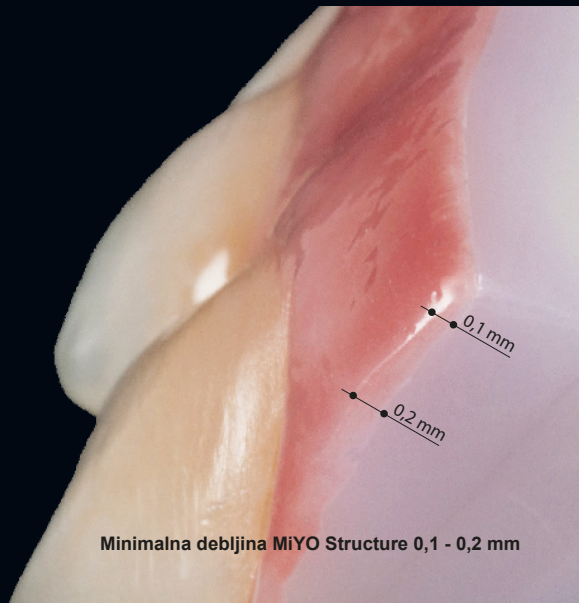

Minimalna debljina MiYO Structure 0,1 - 0,2 mm



KOMPONENTE SISTEMA MiYO:

TRANSLUCENT	COLOR	STRUCTURE
 Trans Shade A	 Snow	 Structure Window
 Trans Shade B	 Linen	 Structure Ghost
 Trans Shade C	 Fissure	 Structure Enamel
 Trans Shade D	MAMELON	 Structure Ice
 Trans Straw	 Mamelon Wheat	 Structure Blush
 Trans Sage	 Mamelon Coral	 HT Structure Window
 Trans Lotus	 Mamelon Pumpkin	 HT Structure Enamel
 Trans Clementine	HALO	
 Trans Storm	 Halo Spring	 Glaze Paste - fluor -
 Trans Cobalt	 Halo Autumn	 InSync Liquid - one for all -
 Trans Lumin		
 Trans Lumin plus		
 Trans Smoke		
 Trans Slate		
 Trans Sunflower		
 Trans Garnet		

KOMPONENTE SISTEMA MiYO PINK:

TRANSLUCENT	COLOR	STRUCTURE
 Trans Raspberry	 Flamingo	 Structure Orchid
 Trans Copper	 Crimson	 Structure Rouge
 Trans Midnight	 Plum	 Structure Frost
 Trans Carnation	 Merlot	
 Trans Garnet	 Sorbet	 Glaze Paste - no fluor - InSync Liquid - one for all -
	 Salmon	
	 Sable	
	 Thistle	
	 Venule	

TEHNIČKI PODACI:

	CTE (25-500 °C) [10 ⁻⁶ K ⁻¹]
MiYO Color	7,5 ± 0,5
MiYO Structure	7,4 ± 0,5
Glaze Paste	7,3 ± 0,5





Vaš lokalni zastupnik:

Mikodental Dental Depot
Stojana Navakovica 19
15000 Šabac
+381 (0) 15345123
info@mikodental.com
www.mikodental.rs



Prodavač MiLO-a u Evropi

Proizvođač:

Chemichl AG
Landstrasse 114
9490 Vaduz, Liechtenstein
info@chemichl.com



Prodaja i tehnička podrška u Evropi:

Jensen GmbH
Gustav-Werner-Straße 1
72555 Metzingen, Germany

Tel. +49 7123 922 60
info@jensendental.de
www.jensendental.de

...Više o MiYO

