


INJECTABLE SYRINGES




An Injectable Syringe is comprised of a barrel for containing the product; a plunger, which is used to evacuate the contents; and a cap, which secures itself on the nozzle of the barrel.

Injectable Syringe Barrels range in size from 2 ML to 100 ML and are manufactured using polypropylene material. These translucent polypropylene barrels are often decorated with dosage marks which helps with an accurate injection. The Injectable Syringe Barrels are available with luer lock or luer slip nozzles allowing for needle attachment commonly used for injecting medicants.

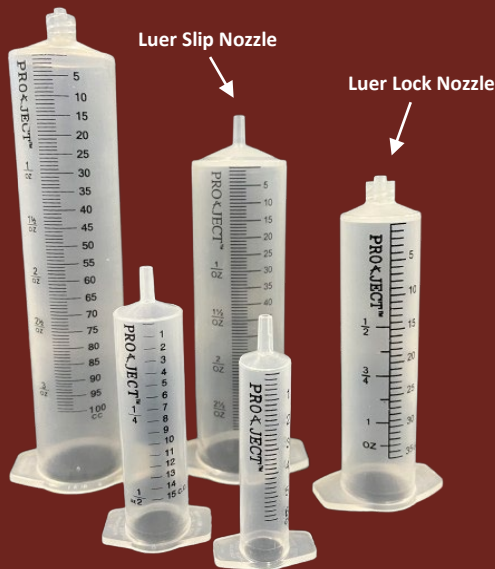

Plungers can be either one piece, meaning the seal is molded onto the end of the plunger; or two piece, where a seal is installed on the end of plunger.

We specialize in syringe decoration including hot stamping and laser printing on plungers to represent dosage marks or customer specific graphics. Natural polypropylene barrels used for injection or dispensing are decorated with graphics using silk screening, pad printing or offset printing techniques. White or natural polyethylene barrels can be offset printed or pressure sensitive labels can be applied.

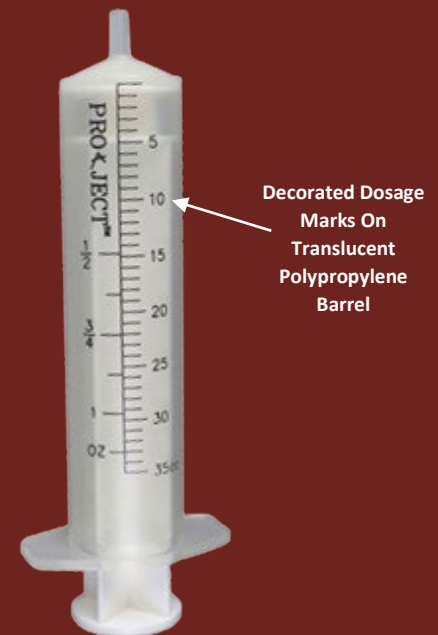

For 40 years, Genesis Industries, Inc. has been manufacturing precision products that exceed our customers' expectations. We work directly with R & D Managers, Production Managers and Engineers to invent and design high quality products.

INJECTABLE SYRINGES



2 ML INJECTABLE SYRINGES

Available Nozzle:
Luer Slip (B13)

Standard Material:
Polypropylene

Available Cap:
Luer Slip (C13)




6 ML INJECTABLE SYRINGES

Available Nozzles:
Luer Slip (B13), Luer Lock (B14)

Standard Material:
Polypropylene

Available Caps:
Luer Slip (C13), Luer Lock (C14)

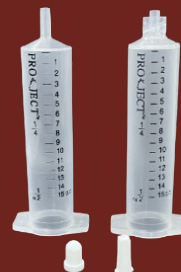

15 ML INJECTABLE SYRINGES

Available Nozzles:
Luer Slip (B13), Luer Lock (B14)

Standard Material:
Polypropylene

Available Caps:
Luer Slip (C13), Luer Lock (C14)



35 ML INJECTABLE SYRINGES

Available Nozzles:
Luer Slip (B13), Luer Lock (B14)

Standard Material:
Polypropylene

Available Caps:
Luer Slip (C13), Luer Lock (C14)



60 ML INJECTABLE SYRINGES

Available Nozzles:
Luer Slip (B13), Luer Lock (B14)

Standard Material:
Polypropylene

Available Caps:
Luer Slip (C13), Luer Lock (C14)



80 ML INJECTABLE SYRINGES

Available Nozzles:
Luer Slip (B13), Luer Lock (B14)

Standard Material:
Polypropylene

Available Caps:
Luer Slip (C13), Luer Lock (C14)



100 ML INJECTABLE SYRINGES

Available Nozzles:
Luer Slip (B13), Luer Lock (B14)

Standard Material:
Polypropylene

Available Caps:
Luer Slip (C13), Luer Lock (C14)

Genesis Industries, Inc. – USA
1.800.826.3301
Sales@GenInd.com

Genesis Manufacturing Ltd. – United Kingdom
+44 (0) 1225 791 184
CustomerServices@GenesisIndustries.co.uk

Genesis Industries (Aust) Pty. Ltd.
+61 (0) 2 6372 6622
Sales@GenesisIndustries.com.au

Genesis Industries (New Zealand) Ltd.
+64 (0) 9 444 2115
NZ.Sales@GenesisIndustries.co.nz