




ORAL & TOPICAL SYRINGES

An Oral & Topical Syringe is comprised of a barrel for containing the product; a plunger, which is used to evacuate the contents; and a cap, which secures itself on the nozzle of the barrel.

Oral & Topical Syringe Barrels range in size from 1 ML to 120 ML and are manufactured using High Density Polyethylene (HDPE) or Low Density Polyethylene (LDPE) material.

Oral & Topical Syringe Barrels are available with a variety of nozzles which are used for oral or topical products and dispense single or multiple doses. Nozzle options include wide mouth, serrated tip, short oral, long nose, short nose, extended nose, large cannula, paste tip, broad nose and enema tip.

Oral & Topical Syringe Plungers can be either one piece, meaning the seal is molded onto the end of the plunger; two piece, where a seal is installed on the end of the plunger; or three piece, meaning a dosage ring is included with the plunger. Plungers can be either straight ribbed or notch ribbed accepting metered dosing rings. Dosage Rings are available in Multi-Dose™ and Accu-Dose™. The Multi-Dose™ ring is rotated on the plunger to the prescribed dosage mark on the plunger rib, while the Accu-Dose™ ring is notched to reflect one rotation of the ring for a pre-determined dose. Seals include flat faced, tapered face and EZ Glide™.

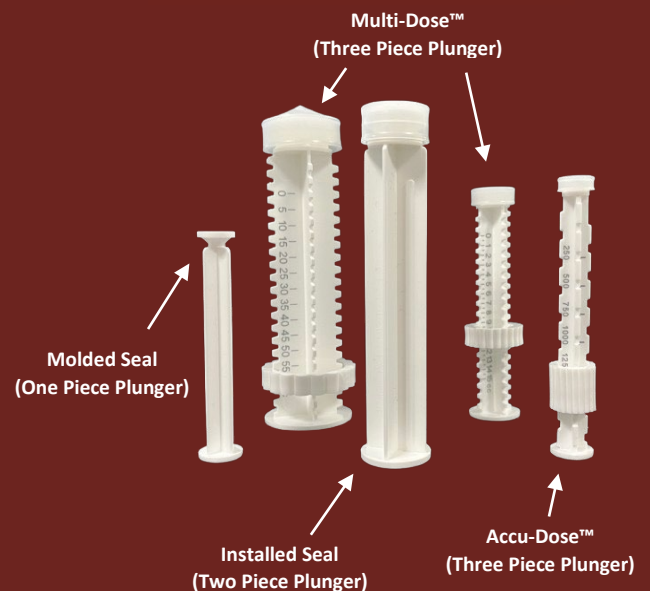

Genesis Industries offers highly accurate air-powered filling machines made with stainless steel construction. These filling machines are ideal for filling a variety of containers and syringes manufactured by Genesis. The SVS is ideal for large production runs and includes a 12 gallon hopper and built-in scale. The SVS can fill from fractional ounces up to 157 ounces. The SVS-J is ideal for smaller production runs and can fill from fractional ounces up to 8.5 ounces. Both filling machines can be operated using auto cycle mode or single foot pedal control.

Over 40 years of experience in precision manufacturing allows Genesis Industries to consistently deliver exceptional products and solutions. We take pride in our ability to produce quality products that exceed our customers' expectations.



1 ML BARRELS

Available Nozzle:
Wide Mouth (B02)

Available Cap:
Wide Mouth (C01)



4 ML BARRELS

Available Nozzle:
Enema Tip (B20)

Available Cap:
Enema Tip (C20)



6 ML BARRELS

Available Nozzles:
Short Oral (B11),
Enema Tip (B20)

Available Cap:
Enema Tip (C20)



7 ML BARRELS

Available Nozzles:
Short Nose (B22),
1.8" Extended Nose (B52),
2" Extended Nose (B53),
Extended Nose (B72),
Sure-Grip™ Short Nose (B125)

Available Caps:
Long Nose (C12), Sure-Grip™ 7 CC (C125-07)



10 ML BARRELS

Available Nozzles:
Wide Mouth (B02),
Paste Tip (B15)

Available Caps:
Wide Mouth (C02),
Paste Tip (C15)



15 ML BARRELS

Available Nozzles:
Wide Mouth (B02),
Short Oral (B11), Long Nose (B12),
Paste Tip (B15),
Short Paste Tip (B151),
Threaded Paste Tip (B156),
Enema Tip (B20)

Available Caps:
Long Nose (C12), Paste Tip (C15)



32 ML BARRELS

Available Nozzles:
Serrated Tip (B09),
Long Nose (B12),
Paste Tip (B15),
Large Cannula (B21)

Available Caps:
Long Nose (C12),
Paste Tip (C15),
Luer Slip (C13)



36 ML BARRELS

Available Nozzles:
Long Nose (B12), Paste Tip (B15),
Sure-Grip™ Short Paste Tip (B155),
Short Broad Nose (B281),
Short Nose (B32)

Available Caps:
Long Nose (C12), Paste Tip (C15),
Sure-Grip™ Short Oral (C155-36),
Broad Nose (C28300SS)



60 ML BARRELS

Available Nozzles:
Long Nose (B12), Paste Tip (B15),
Broad Nose (B28)

Available Caps:
Long Nose (C12), Paste Tip (C15),
Broad Nose (C28300SS)



80 ML BARRELS

Available Nozzles:
Long Nose (B12), Paste Tip (B15),
Broad Nose (B28)

Available Caps:
Long Nose (C12), Paste Tip (C15),
Broad Nose (C28300SS)



100 ML BARRELS

Available Nozzles:
Long Nose (B12), Paste Tip (B15),
Broad Nose (B28)

Available Caps:
Long Nose (C12), Paste Tip (C15),
Broad Nose (C28300SS)



120 ML BARRELS

Available Nozzle:
Broad Nose (B28)

Available Cap:
Broad Nose (C28300SS)