

# Silicone Foam Dressing

Gentle silicone contact layer foam dressing for use on moderately to heavily exuding wounds.

Happy   
Healthcare



**Silicone Foam Dressing consists of a gentle silicone wound contact layer, a super absorbent pad and a vapour permeable and waterproof film.**

It is ideal for wound exudate management. The multi-layered construction facilitates exudate management to provide the optimal moist wound environment which promotes wound healing, minimising trauma and reducing pain during dressing changes.

No residue is left after removal.

## Protective Backing

- Smooth and flexible
- Water proof & bacteria barrier

## Super Absorbent Core

- Exceptional absorption capability
- Helps minimise maceration risk

## Non-woven Layer


- Exudate laterally wicks to enhance absorption

## Absorbent PU Foam

- Great fluid handling capacity
- Offers additional comfort and cushioning protection


## Gentle Silicone Wound Contact Layer

- Gentle silicone and perforation technology
- Distributes wound exudate fast and evenly

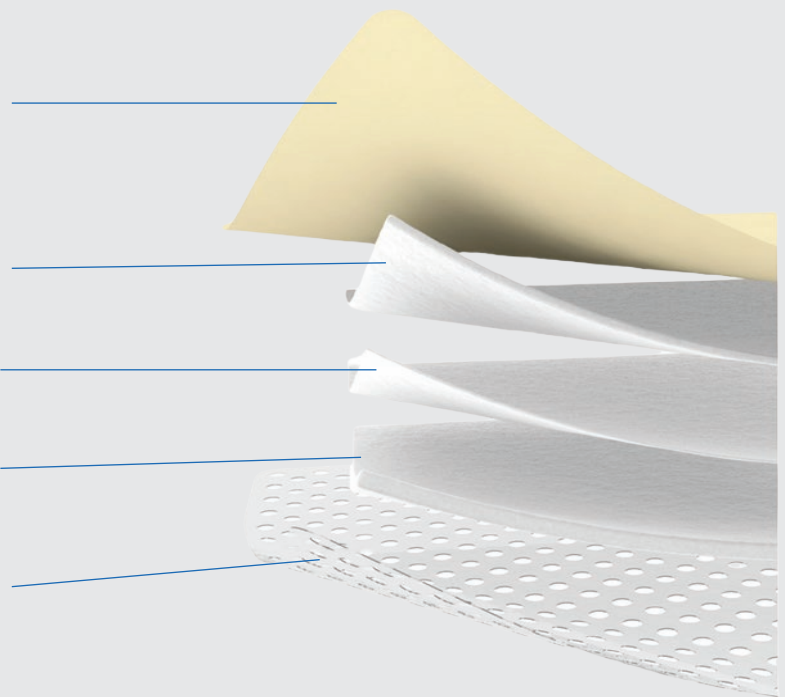

## Benefits:

- Core silicone gel technology provides a gentle silicone contact layer for fragile skin.
- Excellent wound exudate absorption, reduces the risk of maceration and leakage.
- Minimises pain on removal; no second dressing needed.
- Various shapes adapt to different requirements.
- Maintains its integrity during use

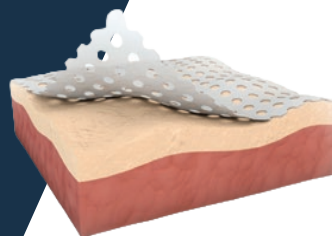

## Soft & Tacky

Able to create multiple points of contact between the uneven surface of the skin and wound which creates a tight seal, preventing exudate leakage and maceration.



Various shapes and sizes are available to accommodate different wound types.