


Veti-Screen NICOTINE

In-House Urinary Drug Screen for Nicotine

Veterinary Adjusted Drug Screen for Nicotine

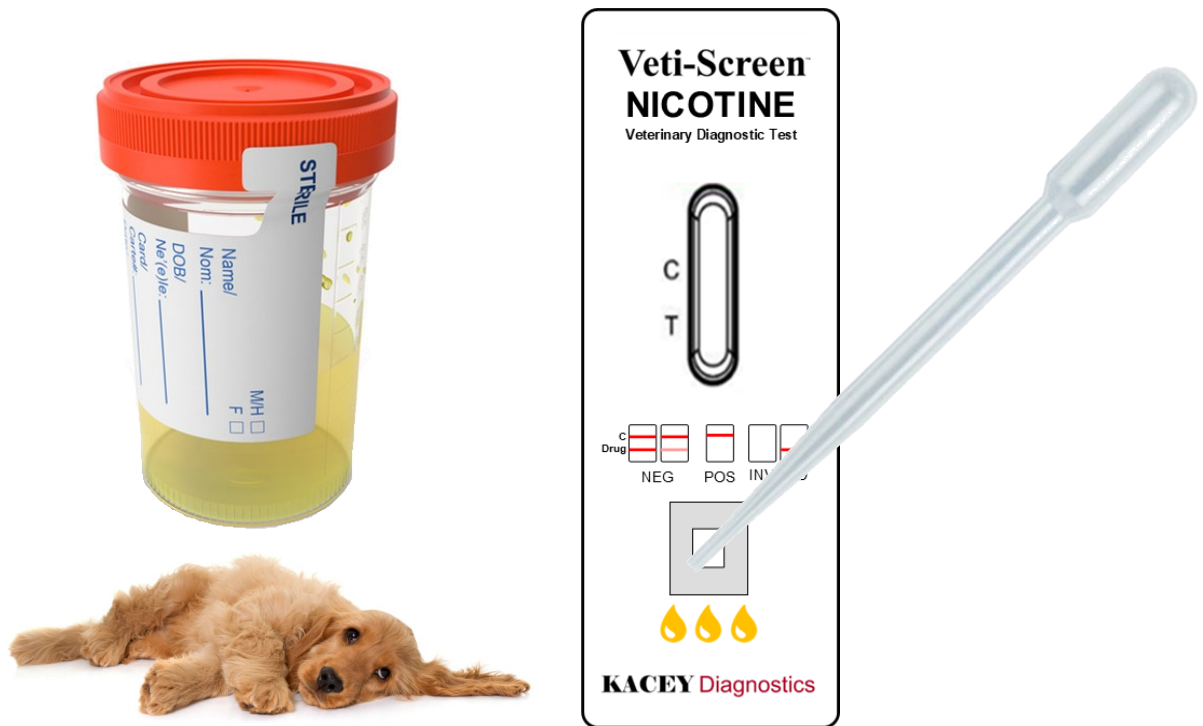

- Adjusted to the lowest threshold in the Veterinary market
- Rapid urine test utilizing monoclonal antibodies to detect levels
- Simply place a drop of urine in the well, read in 5 minutes
 - 14 months expiration, individually packaged
- Available in 5 (#30342) & 10 (#30343) cassette packs




"Kacey Products Providing Better Animal Care"

Kacey Vet-Screen Animal Nicotine test compared to the Standard Human Nicotine Test

There is a tremendous difference between the Kacey Nicotine tests and standard common human Nicotine tests. Standard Nicotine Human tests have a much higher LOD (level of detection) for human levels of detection than that of the Kacey Nicotine test for animals.

The Kacey Vet-Screen Nicotine test was developed for animals and not humans therefore, the Kacey Nicotine tests should never be used to determine Nicotine values in humans. However, the Kacey DOA tests have a lower LOD which detects clinical levels which would otherwise appear as a **“False Negative”** in human tests.

The difference is shown here in the Level Of the Detection (LOD) chart for Nicotine :

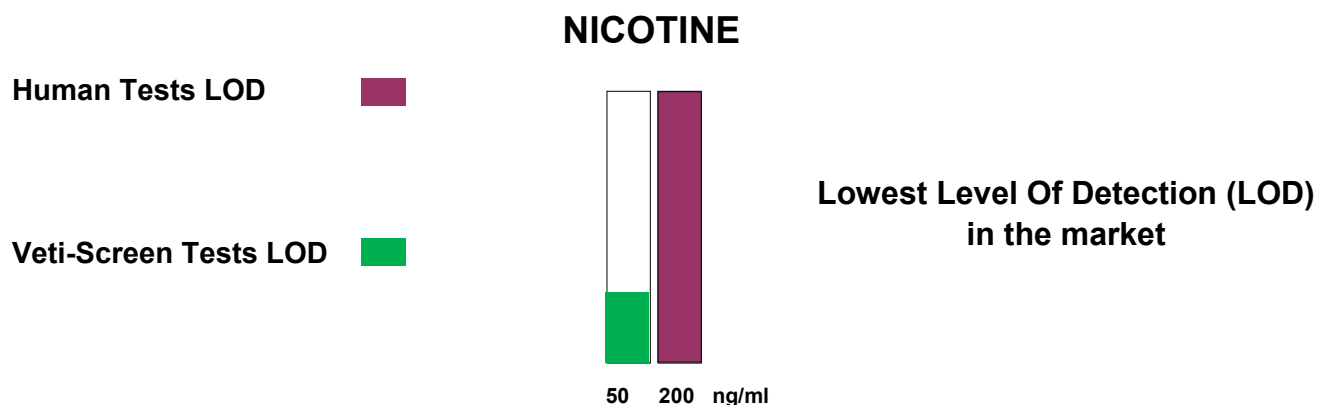

All measurements for the Kacey Nicotine test are in nanograms / mL (ng/mL).

In all cases the Kacey Nicotine detection level is four (4x) times lower at the LOD of 50ng/ml versus the standard human cutoff values of 200 ng/ml for Nicotine tests available. In almost if not all cases the human Nicotine tests would indicate a **FALSE NEGATIVE**, while the Kacey Nicotine test would indicate a **POSITIVE** as part of the Veterinary diagnostic Nicotine screen test.

Note: The Kacey Nicotine tests are tested only on verified Animal Urine for QC (Quality Control) and QA (Quality Assurance) testing. All controls are tested for their performance characteristics of precision, reproducibility, and accuracy are performed on



1854 A Hendersonville Rd, Asheville, NC 28803

Phone: 1(828) 685-3569 Fax: 1(828) 685-3571

Email: sales@kaceyinc.org

www.KaceyDiagnostics.com