



# NIGHTHAWK DENTAL



**Manufacturer**                      **Implants**                      **Ncm**                      **Driver Tip**                      **ASC**                      **MU**                      **Driver Tip**                      **Ncm**



Special Driver for ASC

<b>BIOHORIZONS</b>	All Implants	30 Ncm	0.050 - 1.25 hex	YES 15°	YES	same	15 Ncm
--------------------	--------------	--------	------------------	---------	-----	------	--------

<b>BlueSkyBio</b>	Bio/One Stage = Straumann Tissue Level Compatible	30 Ncm	0.048 hex	N/A	N/A		
	Bio/Trilobe = NobelReplace Compatible		0.050 hex		YES		20 Ncm
	Bio/Internal Hex = Implant Direct / Zimmer Compatible						
	Bio/Conus 12 = Astra OsseoSpeed Compatible	30 Ncm	.065 Star	YES 15°	YES	0.048 Hex	
	Bio/Forte =		Unigrip				
	Bio/Quattro = Straumann bone level Compatible						
Bio/Max = NobelActive Compatible							

<b>camlog</b>	Camlog & Conelog	20 Ncm	0.050 - 1.25 hex	YES 15°	YES	0.050 hex	15 Ncm
---------------	------------------	--------	------------------	---------	-----	-----------	--------

<b>Dentsply Sirona Implants</b>	Ankylos c/x - c/ = 3.5, 4.5, 5.5, 7.0 (non-hex)	15 Ncm	1.0mm	NO	N/A		
	Ankylos c/x - /x = 3.5, 4.5, 5.5, 7.0 (hex)			YES 15°			
	OsseoSpeed/ OsseoSpeed TX = 3.0 - Yellow	20 Ncm	0.050 hex	YES 15°	N/A		
	OsseoSpeed/ OsseoSpeed TX = 3.5, 4.0 - Aqua			YES 20°			
	OsseoSpeed/ OsseoSpeed TX = 4.5, 5.0 -Lilac	YES 25°					
	OsseoSpeed TX Profile = 4.5 , 5.0	NO					
	OsseoSpeed EV = 3.0 Green	25 Ncm	0.050 hex	YES 25°	YES	same	15 Ncm
	OsseoSpeed EV = 3.6 Purple, 4.2 Yellow, 4.8 Blue, 5.4 Brown				N/A		
	OsseoSpeed Profile EV = 4.2 Yellow, 4.8 Blue						
Xive S = 3.0, 3.4, 3.8, 4.5, 5.5	24 Ncm	1.22 hex	YES 24°				

<b>DentiumUSA</b>	Super Line & Implantium	25-30 Ncm	1.28 hex				
-------------------	-------------------------	-----------	----------	--	--	--	--

<b>HIOSSEN IMPLANT</b>	Hiossen/Osstem HG mini = 3.5	20 Ncm	1.2 hex				
	Hiossen/Osstem HG Standard = 4.0, 4.5, 5.0	30 Ncm					

<b>ids integrated dental systems</b>	Reflect Recover = Compatible with Nobel Replace Select	35 Ncm	1.2 hex	N/A	YES	same	15 Ncm
	Reflect Certus = Compatible with Biomet 3i						
	Reflect Rapid = Compatible with Novel Active						
	Reflect Aspire = Compatible with Astra Tech OsseoSpeed TX						
	Reflect Tapered Screw = Compatible with Zimmer TSV						

<b>IMPLANT DIRECT</b>	Legacy = 3.0 Blue	20-25 Ncm	0.050 hex	YES 25°	YES	0.050 hex	20 Ncm
	Legacy = 3.5 Green, 4.5 Purple, 5.7 Yellow						
	InterActive/Swish Active = 3.0 Purple, 3.4 Yellow						
	Replant/Trilobe = 3.5 Purple, 4.3 Yellow, 5.0 Blue, 6.0 Green						
	Swish = 4.8 Brown, 6.5 Blue						
	30 Ncm	N/A					

<b>KEYSTONE DENTAL GROUP</b>	Prima	30 Ncm	0.047 □ Tapered	NO	YES	0.047 □	20 Ncm
	Genesis						
	Tilobemaxx						
	Paltop Internal = NP Blue	25 Ncm	0.050 hex	YES 25°	YES	0.050 hex	30 Ncm
	Paltop Internal = SP Silver, WP Purple	35 Ncm					

<b>MEGAGEN IMPLANTS</b>	AnyRidge = 3.5 Magenta, 4.0-5.5 Blue, 6.0-8.0 yellow	30-35 Ncm	1.2 hex	NO	YES	same	15 Ncm
	AnyOne Internal	35 Ncm			N/A		
	Rescue	45 Ncm					
	Mini	15 Ncm					

<b>MIS</b>	C1 = Conical Connection = NP Yellow, RP Purple, WP Green	30 Ncm	0.050 hex	YES 20°	YES	0.050 hex	25 Ncm
	V3 = Conical Connection = NP Yellow, RP Green						
	M4 = Internal Hexagon = NP Yellow, RP Purple, WP Green						
	Seven = Internal Hexagon = NP Yellow, RP Purple, WP Green						



**Manufacturer**                      **Implants**                      **Ncm**                      **Driver Tip**                      **ASC**                      **MU**                      **Driver Tip**                      **Ncm**



Special Driver for ASC



GM Line	20 Ncm	Neodent driver	YES 25°	yes	same	10 Ncm
CM Line			N/A			
Zygomatic Implant						



3.25 Implant	20 Ncm	Neoss Driver	YES 25°	N/A	same	20 Ncm
3.5 - 6.5 Implants	32 Ncm		YES			



NobelActive = 3.0	15 Ncm	Unigrip	N/A	N/A	N/A	N/A
NobelReplace = NP, RP	35 Ncm		YES 25°	YES	same	15 Ncm
NobelActive = NP, RP, WP			N/A			
Branemark = NP, RP, WP			YES			
Angled Multi-unit Abutments are 15 Ncm into implant	15 Ncm	Special Driver for abutment				
Straight Multi-unit Abutments are 35 Ncm into implant	35 Ncm	ASC Driver	YES 25°			
Procrabutment or crowns		mini ASC Driver	YES 25°	YES	mimi ASC	15 Ncm
Multi-unit Abutments ASC option for Procera bridges	15 Ncm	Special Driver	N/A		N/A	
NobelPearl = RP, WP						



Tissue Level = 3.3 NNC, 3.3 RN, 4.8 WN	35 Ncm	Flat Star	YES 25°	N/A		
BL/BLT = 3.3 NC Yellow, 4.1 RC Purple,				YES	same	15 Ncm
BLX = RB purple, WB Brown						
Zygomatic				N/A		
Pure Ceramic	15 Ncm		N/A			




Thommen = 3.5 Yellow	10 Ncm	4-lobe	YES 25°	YES	same	15 Ncm
Thommen = 4.0 Green, 4.5 Blue, 5.0 Gray, 6.0 Purple	25 Ncm					



Eztetic implant = 2.9 Blue	15 Ncm	0.050 hex	N/A			
Tapered Screw-Vent Implants = 3.5 Green,	30 Ncm		YES 25°			15 Ncm
Certain Internan Connection	20 Ncm	Special Driver	YES	0.050 hex	10 Ncm	
Straight Multi-unit Abutments		0.050 hex				
External Hex Connection		Special Driver				
Straight Multi-unit Abutments						

The screw will be captive inside the solution - Do NOT remove

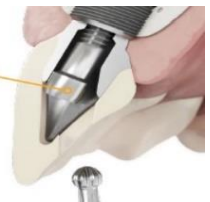

Clones parts are available for most implant. Torque values may differ from OEM parts



0.048 Hex  
Tapered  
Short: HNT  
Long: HNL



0.050 Hex  
Tapered  
Short: HTT  
Long: HTLT



Torx  
Short: PROC-S  
Long: PROC-L  
NobelBiocare



Torx  
Short: ITS  
Long: ITL  
I/VStraumann



0.017  
Slot Driver  
Short: SN  
Long: SNL