


ABOUT US

Spring Health Products, Inc. is a recognized industry leader in the manufacture of diamond burs, LED curing lights, hand-pieces and laboratory rotary and is committed to research and the continued advancement of materials science to bring innovative dental products to the market. Our products are made in the USA and have received the highest clinical evaluations. We have been producing premium quality diamonds for more than 30 years! Our factory is located in Norristown, Pennsylvania. Backed by Spring Health Products 100% satisfaction guarantee.



Spring LED Curing Lights

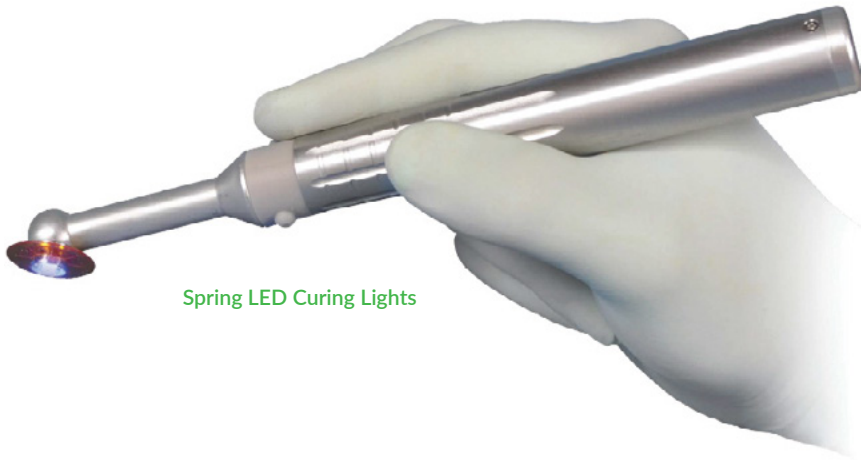

Handpieces



Lab Rotary

Spring Health Products, Inc.

705 General Washington Avenue
Norristown, PA 19403 USA
Phone: 610.630.9171
Toll-Free: 1.800.800.1680
Fax: 610.630.9730
customer@springhealthproducts.com
www.springhealthproducts.com

Request a **FREE** sample of your favorite diamond bur shape and grip by emailing or calling us.

SAMPLE
DIAMOND
ENCLOSED

SPRING

HEALTH PRODUCTS



MADE IN USA



Celebrating
30 Years!

AWARD WINNING
SINGLE-PATIENT USE DIAMONDS

"A quality diamond that has the speed and durability to compete with any diamond on the market."

Andrew Deak, MD, DMD - Elyria, OH

"The burs cut like butter... even after the first use."

Joseph Lichter, DDS - Brooklyn, NY

"Felt like I was using an Electric Handpiece when I really wasn't!"



Spring's Single-Patient diamond burs provide a fresh, fast-cutting bur every time while improving infection control, expediting procedures, reducing patient trauma and decreasing handpiece wear and tear.

A sharper tool every time!

Cuts like a new multi-use bur at a fraction of the cost!

Diamonds stay sharp with no decrease in cutting efficient.

Always a fast, smooth, vibration free cut!

Premium Diamonds at Economical prices!

FASTEST & SHARPEST CUT! LESS CHAIRTIME, LESS TOOTH TRAUMA

Spring's Advanced Diamond Bonding Technology secures diamonds to the shank with maximum staying power thus ensuring that there is NO loss of diamonds. The diamond partials are ingrained at the optimum depth for the most cutting surface possible and no bald spots. Bald spots or over-plating on the bur reduce speed causing heat and risks trauma to the tooth. Spring's long experience in diamond plating has made them excel at diamond covering the tips of the bur - the hardest part to cover.

CONCENTRIC SHANK NO WALK-OUT

Concentric bur Shanks that are micro-precisely calibrated with the strictest tolerances on shank diameter in the industry for a secure fit and made with premium 420 heat treated stainless steel ensures less turbine wear giving a smooth, fast cut with no chatter or vibration.

SPRING DIAMOND BENEFITS

- Highest quality all natural diamonds, blocky in shape, tougher and longer lasting.
- Uniform diamond partial size for better cutting, smoother preparations and less vibration.
- True Super Coarse Grit - a claim most competitors cannot state. Plating 180µ is not easy. The bur cuts much fast with larger grits.
- Precision tip coverage the Spring Diamonds always cut crisp clear margins.
- Made in USA- Manufactured in America with pride and adhering to tough FDA rules and regulations.

TURBO DIAMONDS

Spring turbo diamonds can reduce chair time by eliminating the use of an additional instrument. The spiral design cuts a smooth axial wall as a single instrument - diamond-free cooling zones that allow the cooling water to flow. This bur cuts significantly cooler, making it suitable for use in the bulk reduction required for crown preps and little or no finishing is required.



Heat Dissipates



GRIT LEGEND

- Xtra Fine 30µ=XF.....
- Fine 100µ=F.....
- Medium 130µ=M.....
- Coarse 145µ=C.....
- Super Coarse 180µ=SC.....