


# Assembly and Setup Instructions

## CL-30ZR-60



Step 1. Install The Head Part And Face-plate.....	1
Step 2. Install The Rear-view Mirrors.....	2
Step 3. Install The Storage Box With Armrest.....	3
Additional crutch holder installation.....	4
Additional seat belt installation.....	5

# Step 1

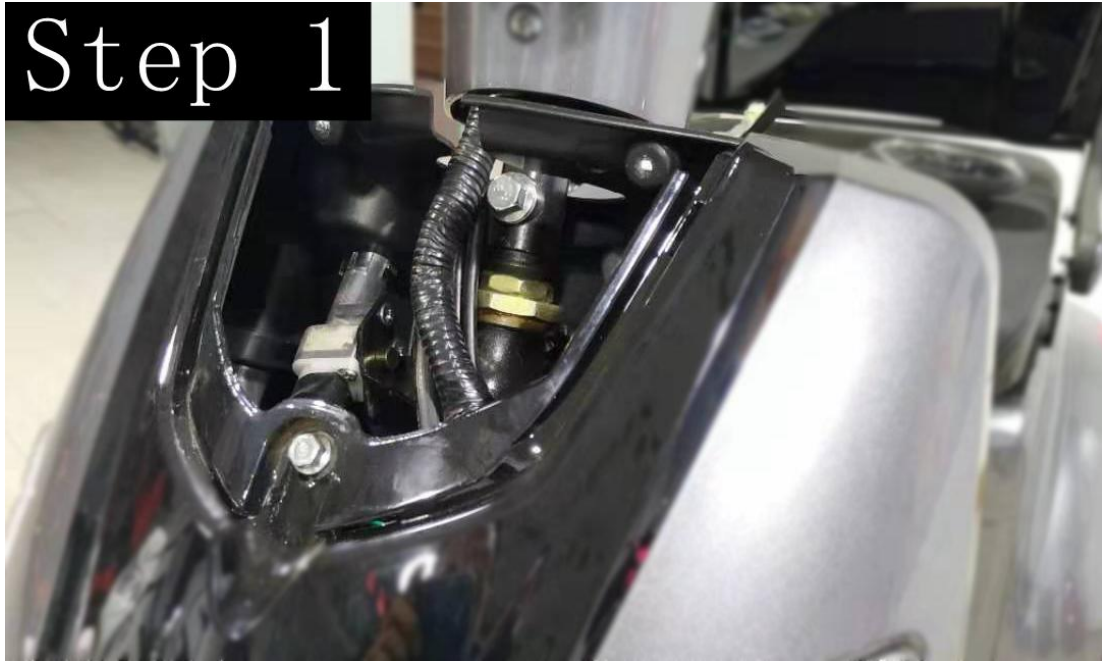

1. Install the head part with screw 10\*45 and latch segment as above pictures. Make sure the screw is tight. (Note: No need nut here, because the nut is fixed on the frame)



2. Put the face-plate in the right position and tap lightly, fix it to the vehicle with self-tapping screw 4.2\*16 on both left and right sides.

## Step 2




1. Check the letter on the rear-view mirror, “L” means left and “R” means right.
2. Turn clockwise to secure the right rear-view mirror to the vehicle, tighten the nut to prevent it from loosening.
3. Install the left rear-view mirror.

## Step 3



1. Install the armrest to slidaway with hexagon socket screw 8\*20 and self-locking nut M8, position as above picture.
2. Install the decoration covers.

## Additional crutch holder installation

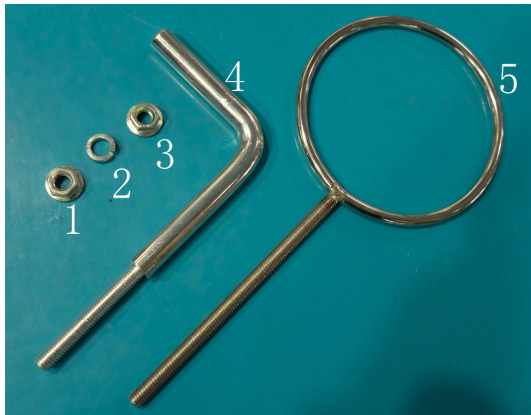

1. Install nut 3 (M8) , spacer 2 and part 4 to position A.
2. Install nut 1 (M8) and part 5 to position B.

# Additional seat belt installation



1. Install the seat belt to slidaway with hexagon socket screw 8\*25 and self-locking nut M8, position as above picture.
2. Install the decoration covers.