



Instruments for Surgery



TABLE OF CONTENTS

About ACE Medical Company	2
Why Buy Our Bipolar Forceps	2
Why Buy Our MIS Bipolars	2
Why Buy our Xclusive Bipolars	3
Why Buy Our Orbitaris Laparoscopic Bipolar	3
Why Buy Our Reusable Cables	3
Bipolar Details	4
Bipolars	5
Dissecting Ring	5
Xclusive Bipolar Forceps	6
MIS Bipolars	7
Isabella Port Tumor Forceps	8
Orbitaris	9
Orbitaris Electrodes	10
Powergrip 3.0	11
Powergrip Electrodes	12
Trays	13
Cables	14
Contact Information and Terms of Sale.....	15

ABOUT ACE MEDICAL COMPANY

ACE Medical Company provides our customers with the highest quality surgical instruments that are as reliable and precise as the care they provide their patients. Our instruments are hand finished by skilled craftsmen, who use the latest technology and finest raw materials.

Beyond that, we understand that each surgical suite is unique and has its own set of needs, and so every effort is made to accommodate those to create a long-term relationship built on trust and efficiency. With dedication to cutting-edge science and technology combined with an understanding of our customers' requirements, ACE Medical Company strives to provide an unparalleled service specifically tailored to each individual customer.

After browsing our catalog, we are confident you will appreciate our product offering. Please call our office at **800.822.0314** or visit our website at **www.acemedicalco.com** to learn why we are preferred in most advanced surgical suites.

WHY BUY OUR BIPOLAR FORCEPS

During bipolar coagulation, a build-up of tissue on the forcep tips is a dangerous potential effect. Using the ACE Bipolars nearly removes this risk completely, and saves the difficult and time-consuming task of cleaning the forceps in between procedures. This also reduces possible interruptions during surgery that would be necessary if cleaning had to take place.

The non-stick effect used with ACE & MIS Bipolars is achieved by using a unique solid silver alloy tip which has excellent thermal conductivity. Over time the natural oxidizing effect will blacken the tips, however this does not impact their non-stick properties or performance whatsoever as it remains permanent despite repeated cleaning and sterilization cycles for the life of the instruments.

Our forceps, the first of their kind, have a unique ergonomic grip, are weighted, and balanced excellently to serve as extensions of your surgical hands. Sturdy enough to grasp tissue, yet precise enough to navigate and visualize the most cavernous anatomy.

WHY BUY OUR MIS BIPOLARS

Working with Dr. Daniel Prevedello to develop forceps specifically meant for minimally invasive neurosurgery, we focused on two objectives: reducing the profile of the tines and changing the angle of the forcep. To reduce the profile, we ensured that when closed, the tines would completely converge. This allowed narrower anatomy to be navigated with improved visualization. We then altered the angle of the tines to give surgeons better control and move their hands away from the line of sight through a microscope. Last, we opted for two angles of the tips, 15° and 45°, three lengths and solid silver non-stick tips so that these bipolars could be used in port mis procedures without issue. With these modifications in hand, our ACE MIS bipolars became a top choice whenever neurosurgeons prefer a less invasive approach.

WHY BUY OUR XCLUSIVE BIPOLARS

The Claris Xclusive Bipolar Forceps are an essential forcep for delicate procedures. Made of high-quality materials, these forceps contain two micro silver tips on our thinnest tines ever manufactured. They are perfect for narrow anatomy that requires exceptional performance.



WHY OUR ORBITARIS LAPAROSCOPIC BIPOLAR

Do you need an ergonomic and adjustable handle for laparoscopic electrocautery? Look no further than the Orbitaris. These forceps have been precision-crafted to offer superior control when performing procedures such as hysterectomies and gynecologic procedures. The handle, which is designed to twist can always be brought into an ergonomic position while reaching tough anatomy. Browse our extensive line of electrodes made to find the perfect addition to your surgical armamentarium. You can be confident that this device will provide long lasting functionality, performance, and reliability.


WHY OUR REUSABLE CABLES

Looking for a healthy and sustainable alternative to the wasteful single-use products? Please consider our reusable bipolar cables, designed with Best-in-Class Manufacturing for up to 200 uses. They are optimally created to fit any generator, thereby reducing hospitals' and surgeons' environmental footprint. In offering an economical and care-free solution, these products combine sustainability with convenience. Experience the ultimate combination of affordability and efficiency today!


BIPOLAR DETAILS

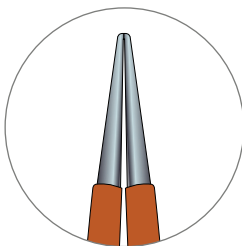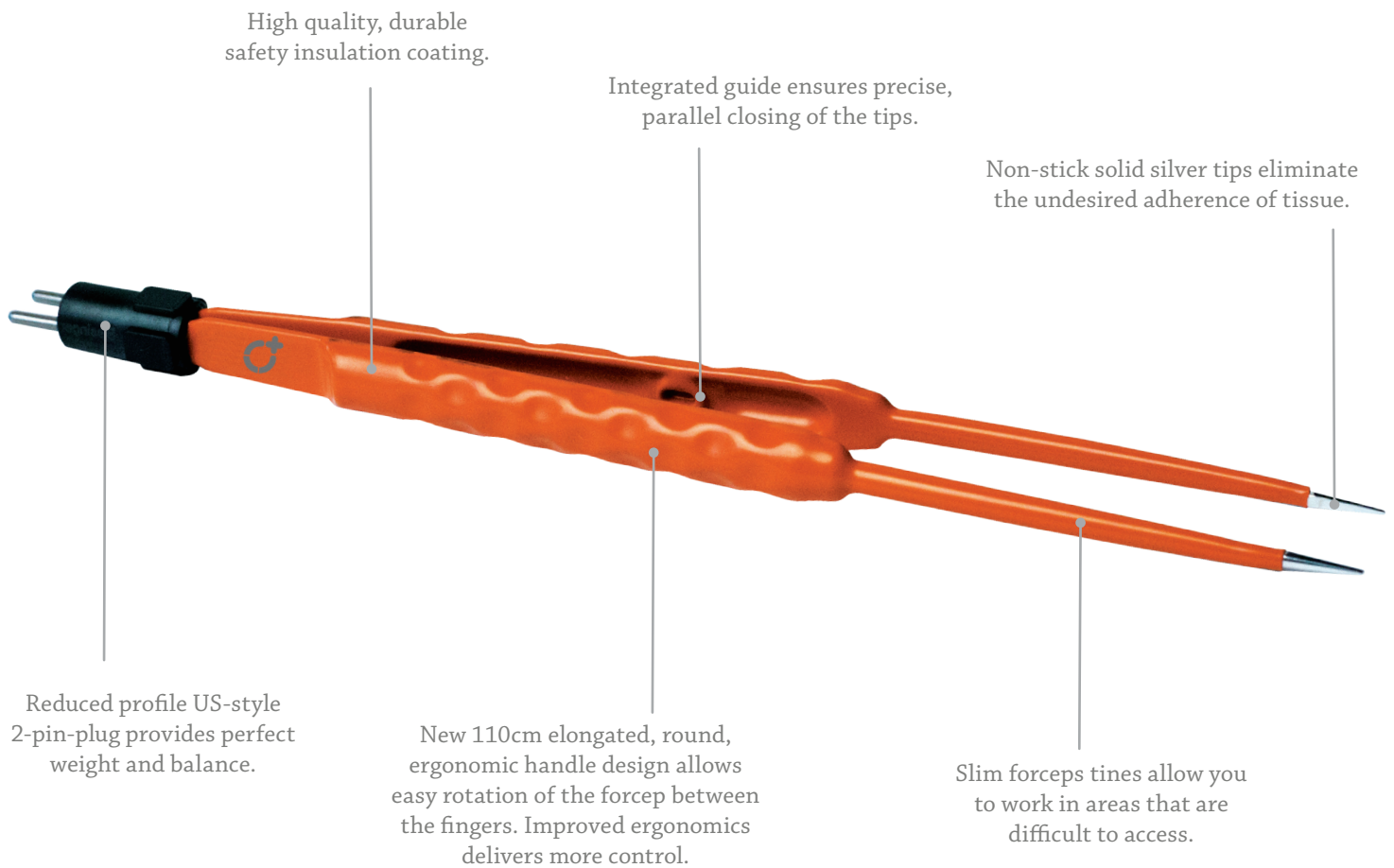


0.3 mm fine tip for precision surgery and dissection

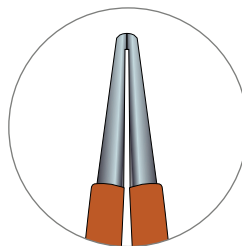


0.6 mm fine tip for precision surgery, dissection and coagulation power

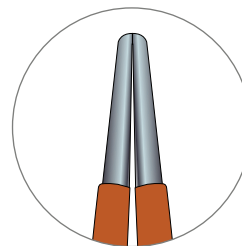

1.0 mm blunt tip for grasping and coagulation of specific tissue

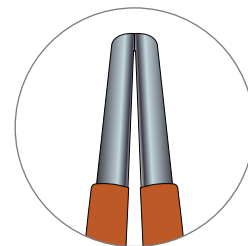

2.0 mm blunt tip for grasping and coagulation

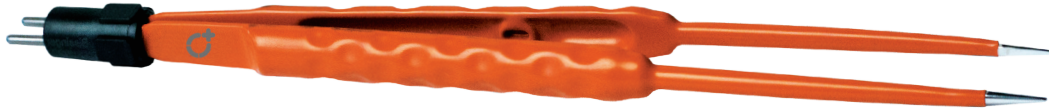

STRAIGHT

(total length/working length)

tip	(140/30) mm	(165/55) mm
0.3 mm	81460100	81460108
0.6 mm	81460102	81460110
1.0 mm	81460104	81460112



BAYONET

(total length/working length)

tip	(165/55) mm	(185/75) mm	(210/100) mm	(235/125) mm	(260/150) mm
0.3 mm	81460132	81460140	81460148	81460156	
0.6 mm	81460134	81460142	81460150	81460158	
1.0 mm	81460136	81460144	81460152	81460160	81460190
2.0 mm	81460138	81460146	81460154	81460162	81460191



BAYONET ANGLED UPWARDS

(total length/working length)

tip	(210/100) mm	(235/125) mm
0.6 mm	81460181	81460182

DISSECTING RING


#80636000



Add dissecting force to any bipolar with a reusable silicone ring sold in a pack of 10





XCLUSIVE

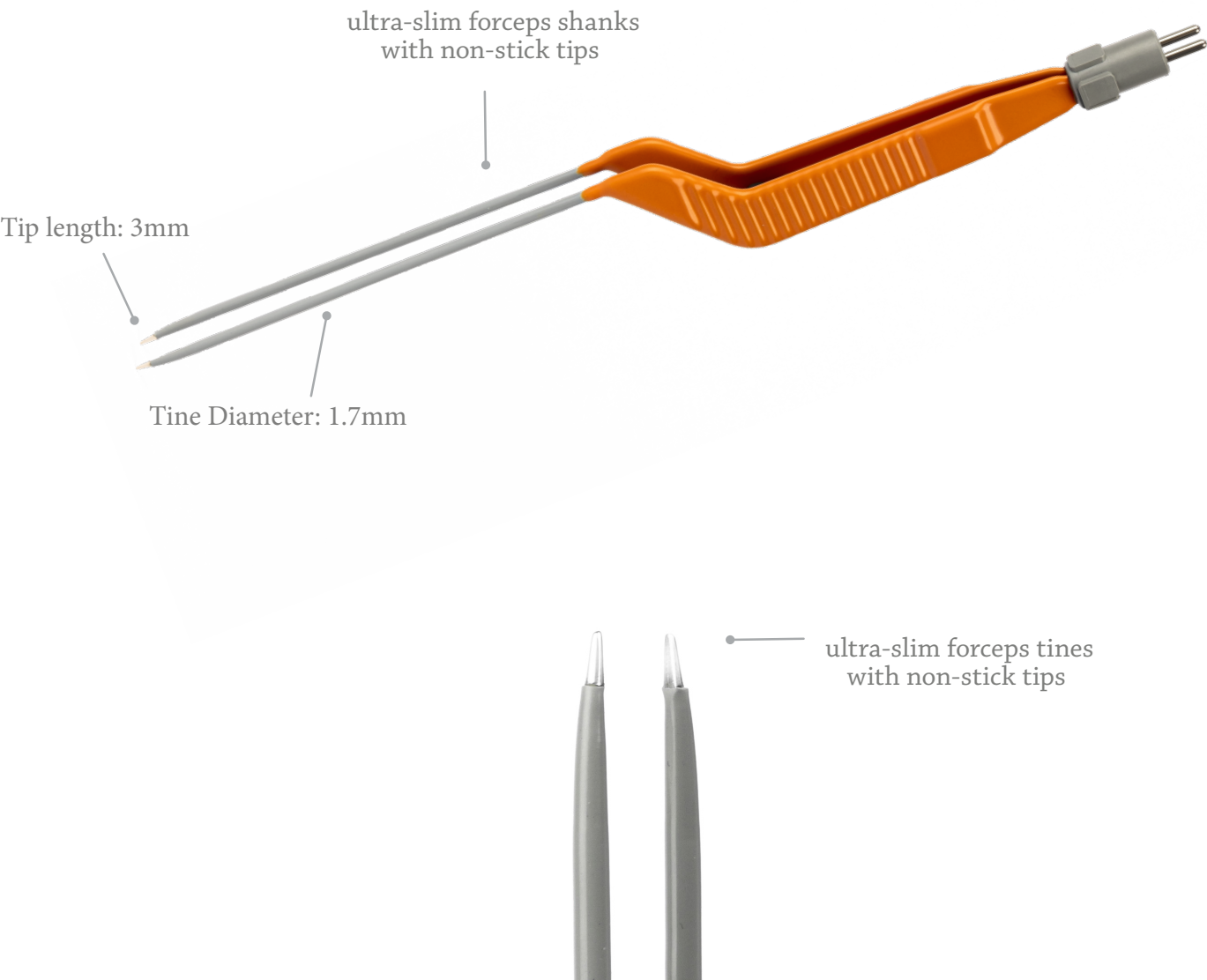

2-PIN PLUG


tip	200 mm	220 mm	240 mm
0.3 x 3 mm	80620397	80622348	80624337
0.6 x 3 mm	80620398	80622349	80624342

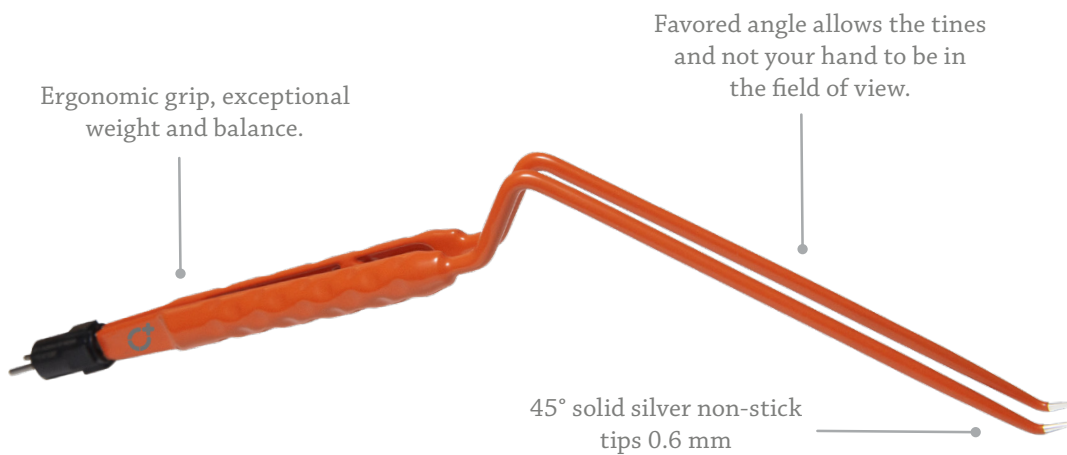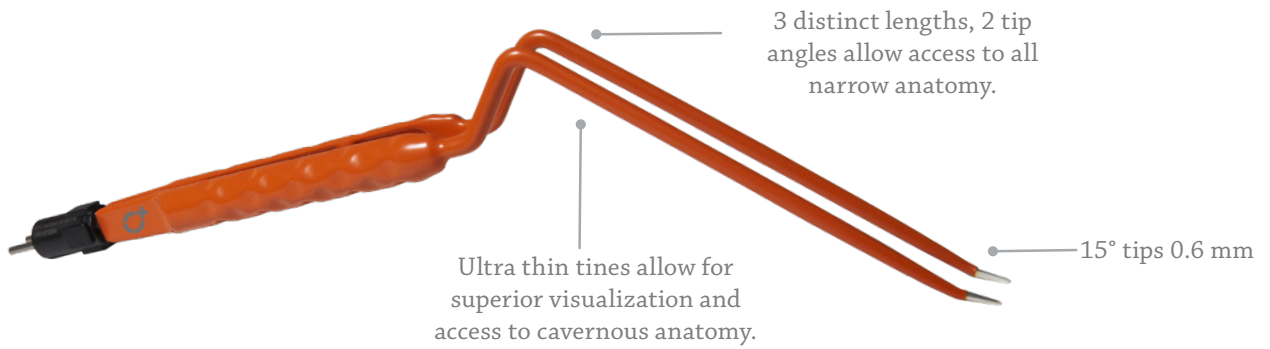


A forcep to accomodate minimally invasive Surgery

MIS BIPOLARS

degree	75 mm	90 mm	105 mm
15°	81470010	81470020	81470030
45°	81470015	81470025	81470035

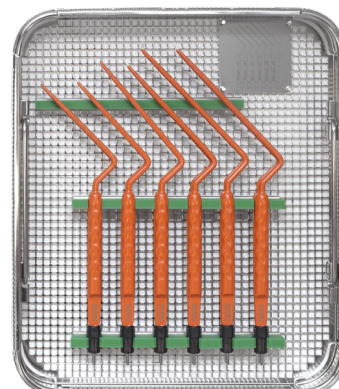


ISABELLA PORT TUMOR FORCEPS

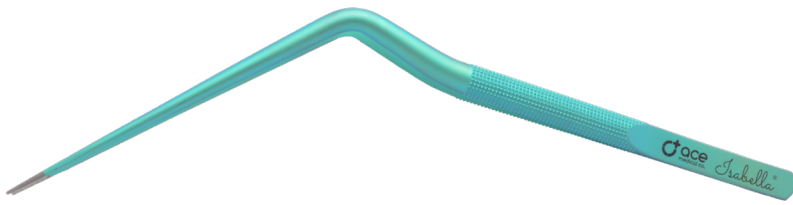


FORCEPS

AM-513-01

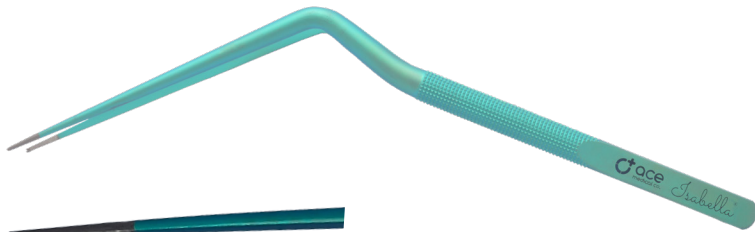

INTERLOCKING TEETH 1x2

AM-513-02

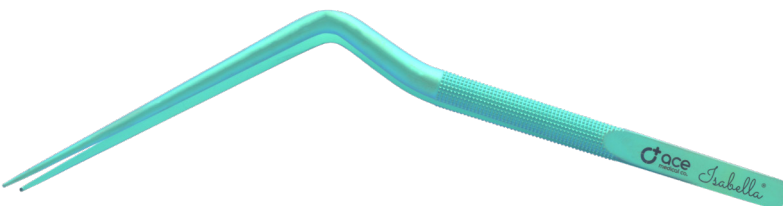


CUP FORCEPS

AM-513-03



Sold in a set of 3

ORBITARIS

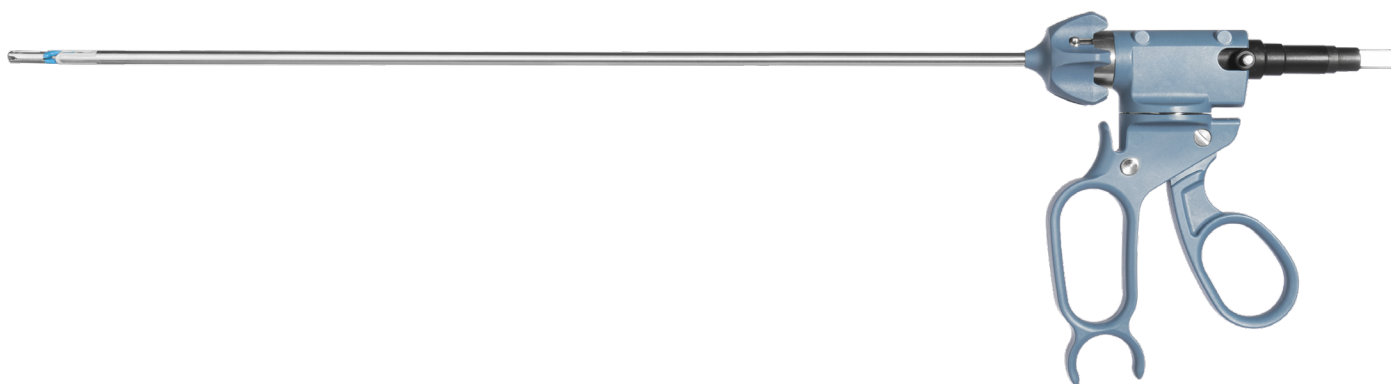

HANDLE

#82460002



SHAFT

#82400034



ELECTRODES

#82408010



ELECTRODES



Shaft Length: 340mm

Shaft Width: 5mm

LARGE GRASPING FORCEPS

#82408010



FORCEPS, CURVED

#82408011



GRASPING FORCEPS, CURVED

#82408015



SCISSORS, CURVED

#82408018



MARYLAND DISSECTOR, CURVED

#82408019



eCUT


#82408039



sideview of the eCUT



POWERGRIP 3.0

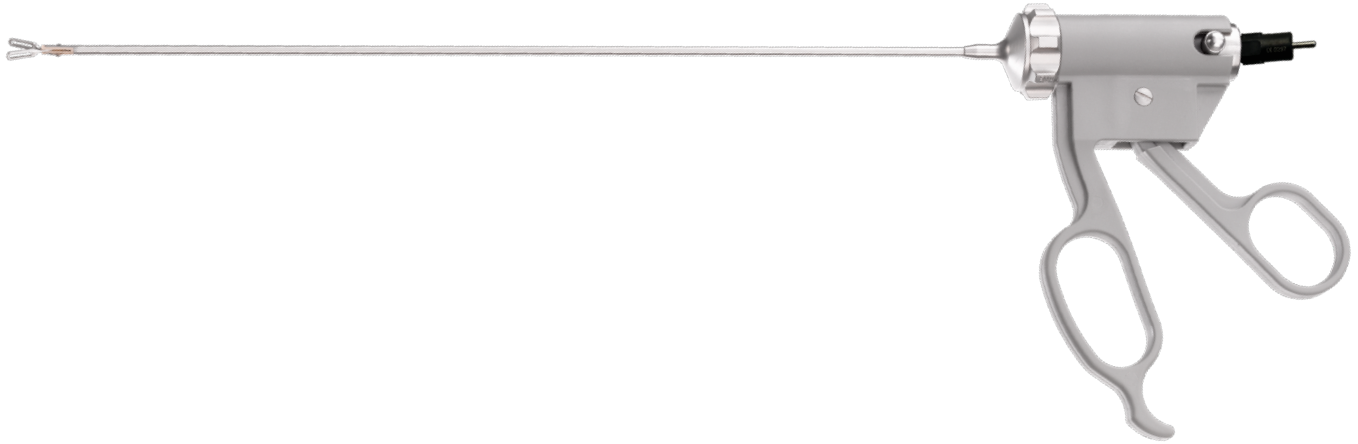

HANDLE

#82460152



SHAFT

#82460204



ELECTRODES

#82460152



ELECTRODES



KELLY CLAMP

#82460386



MARYLAND DISSECTOR, CURVED

#82460302



GRASPING FORCEPS, FENESTRATED

#82460300



SCISSORS, CURVED

#82460324




TRAYS

INSTRUMENT TRAY


Width	240 mm
Length	340 mm
Height	52 mm
Material	Stainless Steel
Part #	BT-10

The tray can store up to 10 instruments.

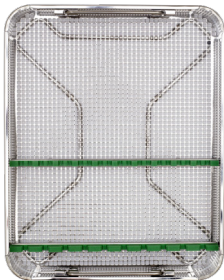


Width	70 mm
Length	290 mm
Height	60 mm
Material	Stainless Steel
Part #	MIS-BT

The tray can store up to 6 instruments.

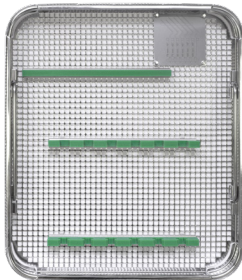

Width	70 mm
Length	290 mm
Height	60 mm
Material	Stainless Steel
Part #	BT-2

The tray can store up to 2 instruments.

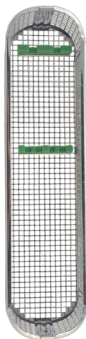

Width	340 mm
Length	240 mm
Height	50 mm
Material	Stainless Steel
Part #	MSS-1

The tray can store up to 2 handles, 2 shafts, and 10 electrodes.

CABLES

CONNECTION CABLES (RE USABLE - UP TO 200 USES)

Plug, generator side	Length	2-Pin	2-Pin Angled
2-Banana	5m / 16.4ft.	80100027	80100129
2-Banana	3m / 9.8ft.	80100025	80100128
Valleylab, Lamidey, EMC	5m / 16.4ft.	80100055	80100154
Valleylab, Lamidey, EMC	3m / 9.8ft.	80100053	
Bipolar Scissors (2-Banana)*	3m / 9.8ft.	80100081	
Bipolar Scissors (2-Banana)*	5m / 16.4ft.	80100085	

*For use with single shaft bipolar

GENERATOR CONNECTION



2-Banana



Valleylab

BIPOLAR CONNECTION



2 Pin Angled



Standard 2 Pin



Single Shaft Bipolar Cable

CONTACT INFORMATION AND TERMS OF SALE

PURCHASE ORDERS MAY BE SENT DIRECTLY TO:

ACE MEDICAL COMPANY
2230 Park Ave., Suite 202
Cincinnati, Ohio 45206

Fax: (513) 321-0787
Customer Service: (800) 822-0314
info@acemedicalco.com

PRICES

Prices are listed in the Ace Medical Company current price list and are subject to change without notice. Freight is FOB destination, prepaid and added to the invoice.

WARRANTY

All instruments are warranted by the manufacturer against defects in workmanship and materials. If the products are not used in accordance with the intended purpose, altered, repaired or modified by persons not authorized by Ace Medical Company, the warranty will be void. Ace Medical Company's sole obligation and customer's exclusive remedy for breach of any warranty, shall be at Ace Medical Company's option, to repair or replace the product. There are no other expressed or implied warranties, including any warranty of merchantability or fitness for a particular purpose.

REPAIR PROGRAM

Ace Medical Company offers a full-service, warranted repair program. Check with customer service for details.

CUSTOM INSTRUMENTS

Ace Medical Company will discuss the feasibility of a custom instrument. We are happy to quote any customization that is deemed possible to manufacture.

RETURN GOODS POLICY

1. All returns must be authorized by Ace Medical Company
2. Returns must be in resalable condition within 30 days of the invoice date. A 25% reprocessing fee will be charged unless the return is due to Ace Medical error.
3. All returns are for credit only.
4. All returns must have an RGA number issued by the Ace Medical Company customer service department.

NON-RESALABLE MERCHANDISE INCLUDES THE FOLLOWING:

Products not purchased from Ace Medical, products that have been sterilized and/or used, products attempted to return 30 days beyond the invoice date, sterile products opened or removed from their original container, custom made or discounted products, and products not in the original packaging.



Visit Ace Medical Company:
www.acemedicalco.com

2230 Park Ave., Suite 202
Cincinnati, Ohio 45206

Phone: (800) 822-0314
Fax: (513) 321-0787
info@acemedicalco.com